



IMAGING SYSTEM ERROR MESSAGE GUIDE

Version 3.0

March 2007

Department of Veterans Affairs
System Design and Development
VistA Imaging

Property of the US Government

No permission to copy or redistribute the software described in this document is given. Use of unreleased versions of this software requires the user to execute a written test agreement with the VistA Imaging Development Office.

This is a controlled document. No changes to this document may be made without the express written consent of VistA Imaging Development Office.

While every effort has been made to assure the accuracy of the information provided, this document may include technical inaccuracies and/or typographical errors. Changes are periodically made to the information herein and incorporated into new editions of this document.

Product names mentioned in this document may be trademarks or registered trademarks of their respective companies, and are hereby acknowledged.

VistA Imaging Software Design & Development
Department of Veterans Affairs
1335 East-West Highway 3rd Floor
Silver Spring, MD 20910
(301) 734-0100

Revision History

10 Nov 2003	Changed title page/footer dates.
13 Nov 2003	Edited document to remove Reference Source errors.
20 Nov 2003	Corrected references to C-Store.
20 Mar 2004	Patch 11 updates
14 Oct 2004	Patch 30 updates
1 Mar 2005	Patch 51 updates . E. deMoel, K. Harding.
6 Mar 2007	Patch 69 updates for conversion to Caché. E. deMoel, A. McFarren. Affected sections: 1.1 – 1.2.1, 1.5.1.2, 1.7.2, 1.7.3, 1.8.3, 1.11, and 1.12.2. Also added sections 1.16.6.20-29; these sections describe previously undocumented error messages introduced in Patch 26.

Preface

The purpose of this guide is to provide information about error messages as they relate to the Veterans Health Information Systems and Technology Architecture (VistA) Imaging V. 3.0 package (i.e., files, routines, and configuration that comprise the VistA Imaging System).

Food and Drug Administration Policy

The Food and Drug Administration classifies this software as a medical device. As such, it may not be changed in any way. Modifications to this software may result in an adulterated medical device under 21CFR820, the use of which is considered to be a violation of US Federal Statutes.

VA Policy

Those components of a national package (routines, data dictionaries, options, protocols, GUI components, etc.) that implement a controlled procedure, contain a controlled or strictly defined interface or report data to a database external to the local facility, must not be altered except by the Office of Information (OI) Technical Services (TS) staff. A controlled procedure is one that implements requirements that are mandated or governed by law or VA (Department of Veterans Affairs) directive or is subject to governing financial management standards of the Federal Government and VA or that is regulated by oversight groups such as the JCAHO or FDA. A controlled or strictly defined interface is one that adheres to a specific industry standard, will adversely affect a package and/or render the package inoperable if modified or deleted. For national software that is subject to FDA oversight, only the holder of the premarketing clearance (510(k)) is allowed to modify code for the medical device. The holder of a premarketing clearance is restricted to specifically designated TS staff that are located at the registered manufacturing site and operating in the designated production environment. Modifying FDA regulated software under any other conditions is a severe violation of the Code of Federal Regulations. Local, that is field-based, developers are prohibited from modifying national software that is certified by the FDA.

Caution: Federal law restricts this device to use by or on the order of either a licensed practitioner or persons lawfully engaged in the manufacture or distribution of the product.

Table of Contents

Chapter 1	DICOM Error Messages	1
1.1	License Limits.....	1
1.1.1	Limit on Number of C-Store Processes	1
1.1.2	Limit on Number of Telnet Sessions	2
1.2	Improper Start-Up Sequence.....	3
1.2.1	Start C-Store Program before MUMPS Server is Running	3
1.2.2	Start C-Store Program while it is already Running	3
1.2.3	Transmit Image before C-Store Program is Running	3
1.2.4	Too many Instruments for one Image Gateway.....	4
1.2.5	Images are acquired, but are not moved to the File Server.....	4
1.3	Menu Options not Available.....	5
1.3.1	System Configuration Settings	5
1.3.2	End-User Privileges	6
1.3.3	Ad-Hoc Operational Issues	6
1.3.4	Logins Disabled	8
1.4	Privileges and Permissions	8
1.5	Disk Full.....	9
1.5.1	Full MUMPS Databases	9
1.5.1.1	Data Stored in Wrong Database.....	9
1.5.1.2	Purge MUMPS Error Log.....	9
1.5.1.3	Purge Old Modality Worklist Entries	12
1.5.1.4	Purge Old HL7 Transaction Global Nodes.....	12
1.5.1.5	Purge Old Audit Records	12
1.5.1.6	Purge DICOM Image Processing Error Logs	12
1.5.2	Full Physical Media	12
1.5.2.1	Purge Old DICOM Message Files	12
1.5.2.2	Purge DICOM Image Files	13
1.5.2.3	Other Possible Issues	13
1.6	Extremely Slow Performance	14
1.7	Network Connectivity Problems.....	15
1.7.1	Only one VistA Imaging DICOM Gateway (that used to work) is Affected ..	15
1.7.2	Several Processors are Affected.....	16
1.7.3	One or more windows show an error message	17
1.7.4	PACS isn't Sending an Image Complete Message	18
1.8	Other Configuration Issues	18
1.8.1	Loading Master Files	18
1.8.2	Incomplete Master Files.....	19
1.8.3	Multiple Instruments are Configured to Use the Same Port Number	19
1.8.4	Digitized Images Appear Small on the Diagnostic Workstation	19
1.8.5	Digitized Images Appear Dark	19
1.8.6	Modality Worklist Find Failure	20
1.9	Failover Procedure	20

1.10	Error Messages from VistA Imaging DICOM Gateway Executable Programs	20
1.10.1	Error Messages from MAG_CSTORE.EXE	20
1.10.2	Error Messages from MAG_DCMTOTGA.exe	22
1.10.3	Error Messages from MAG_ABSTRTGA.exe.....	23
1.10.4	Error Messages from MAG_DCM_COPY.exe	23
1.10.5	Error Messages from MAG_TGATODCM.exe	24
1.11	Run Time Error Messages	24
1.12	DICOM error messages	32
1.12.1	Error message: \$ZE=<STKOV>BTXT+2^DIALOG:::6:6:	32
1.12.2	Error message “Port in Use” or “Port Unavailable”	32
1.13	Signing on to the GUI.....	32
1.14	Dumping Files that Cannot be Processed by the Gateway Application	33
1.15	Microsoft Windows Error Code.....	34
1.16	Run Time Errors Reported by the DICOM Gateway	35
1.16.1	DICOM Association Errors	35
1.16.1.1	Unexpected “DATA” PDU.....	35
1.16.1.2	Peer application entity Requests “Abort Association”	35
1.16.1.3	DICOM Association – Invalid Protocol Identifier.....	35
1.16.1.4	Opcode for Unsupported Operation.....	36
1.16.1.5	DICOM Association – Invalid DICOM Standard UID	36
1.16.1.6	Invalid Context Identifier.....	36
1.16.1.7	Invalid Sub-Type	37
1.16.1.8	Illegal Unique Identifier (UID).....	37
1.16.2	Storage Provider Error – Unexpected Response.....	37
1.16.3	Errors Encountered in Reading/Writing DICOM Messages	38
1.16.3.1	Unsupported Value Representation (Input)	38
1.16.3.2	Value Representation Mismatch	38
1.16.3.3	Required Type 1 Data Item Missing.....	38
1.16.3.4	Unsupported Value Representation (Output).....	39
1.16.4	Modality Worklist – Invalid SOP Class (C-Find).....	39
1.16.5	HL7 Error Messages	39
1.16.5.1	Bad HL7 Message Header	40
1.16.5.2	Incomplete HL7 Message	40
1.16.5.3	Bad Fillers Order Number	40
1.16.6	DICOM Image Processing Errors During Input	41
1.16.6.1	Undefined Imaging Service	41
1.16.6.2	File 2005 Corruption Problem	41
1.16.6.3	Radiology Report File Corruption Problem.....	42
1.16.6.4	Group Entry Number Problem.....	42
1.16.6.5	Image Entry Number Problem	43
1.16.6.6	Imaging Patient Mismatch Problem.....	43
1.16.6.7	Radiology Patient Mismatch Problem	44
1.16.6.8	No Radiology Case Number	44
1.16.6.9	Radiology Case Not in ^RADPT	44
1.16.6.10	Cannot Create Group for Old Radiology Images	45
1.16.6.11	Cannot Create Group for New Radiology Images	45

1.16.6.12	Cannot Find Image Group Pointer in Radiology Report	45
1.16.6.13	Cannot Create Image File Entry in VistA Database	45
1.16.6.14	Cannot Create Subdirectory to Store Image File	46
1.16.6.15	Cannot Write Image File.....	46
1.16.6.16	TIU Note Patient/Report Mismatch	47
1.16.6.17	TIU External Data Link Association Error	47
1.16.6.18	Cannot Find Table with Additional Data Items	48
1.16.6.19	DICOM Element has Too Many Values.....	48
1.16.6.20	Cannot Create Targa™ File	49
1.16.6.21	Cannot Create Abstract File.....	49
1.16.6.22	No Image File Created	49
1.16.6.23	Error while Creating Image File	50
1.16.6.24	Could not Copy Image File.....	50
1.16.6.25	Could not Create Image File	51
1.16.6.26	Image File Created, but File Size not Greater than 0.....	51
1.16.6.27	Image File Created, but Size Differs from Expected File Size.....	52
1.16.6.28	Image File Created, but File Size Differs from Expected File Size....	52
1.16.6.29	Could not Delete Image File	53
1.16.6.30	Cannot find Table with Additional Data Items.....	53
1.16.6.31	Image Patient Mismatch (parent file 2006.5839)	54
1.16.6.32	Image and TIU Note Associated with different Notes.....	54
1.16.6.33	Wrong Group Object Type	54
1.16.6.34	Ignore Image File	55
1.16.6.35	Unknown Instrument Mnemonic	55
1.16.6.36	Rename Failed	55
1.16.6.37	Image File Creation Error	55
1.16.6.38	Image File Subdirectory Creation Error	56
1.16.6.39	New Files Might Overwrite Existing Ones (Single).....	56
1.16.6.40	New Files Might Overwrite Existing Ones (Multiple)	57
1.16.6.41	Cannot Delete File	57
1.16.6.42	Unknown Status (Patient Safety RPC)	58
1.16.6.43	Rollback Error (number of returned elements).....	58
1.16.6.44	Rollback error (returncode).....	58
1.16.6.45	Problem with Temporary File (1)	59
1.16.6.46	Problem with Temporary File (2)	59
1.16.6.47	Problem with Temporary File (3)	59
1.16.6.48	Image Copy Error	60
1.16.6.49	Image Write Failed	60
1.16.6.50	Targa™ Abstract Creation Error	60
1.16.6.51	Multiframe Image File Creation Error (number of returned elements)61	
1.16.6.52	Multiframe Image File Creation Error (return code not “processed”)61	
1.16.6.53	Multiframe Image File Creation Error (return code not “store”).....61	
1.16.6.54	Multiframe Image File Creation Error (error in command).....62	
1.16.6.55	Targa™ File Creation Error	62
1.16.6.56	Could not Open Image File for Write	63
1.16.6.57	Data Transfer Dictionary Missing	63

1.16.6.58	Invalid Multiplicity	63
1.16.6.59	Image Entry Number Problem	64
1.16.6.60	Routing Not Actively Used.....	64
1.16.6.61	Cannot Create Abstract.....	65
1.16.7	DICOM Image Processing Errors During Output	65
1.16.7.1	Bad Targa™ File	66
1.16.7.2	Cannot Establish Connection with MAG_VISTA_SEND_IMAGE	66
1.16.7.3	Connection to Destination Storage Provider Failed.....	66
1.16.7.4	Could Not Establish Association with Destination SCP.....	67
1.16.7.5	DICOM Network Error.....	67
1.16.7.6	Image Transmission Acknowledgement Failed.....	67
1.16.7.7	Image Transmission Failed.....	68
1.16.7.8	Cannot Acknowledge Establishment of Association.....	68
1.16.7.9	Cannot Acknowledge Release of Association	68
1.16.8	Error Messages Report by DICOM Message Queuing Software	68
1.16.8.1	Error in Directory Lookup of FIFO Queue Subdirectory	69
1.16.8.2	Cannot Create FIFO Queue Subdirectory	69
1.16.8.3	Invalid Data in Queue Pointer File	70
1.16.8.4	Cannot Create a FIFO Queue Text File	70
1.16.8.5	Unrecognized Command Code in DICOM Message.....	71
1.16.8.6	Unsupported Command in DICOM Message.....	71
1.16.9	TCP/IP Communications Errors	71
1.16.9.1	Cannot Connect to TCP/IP Socket (Set-Up Error)	71
1.16.9.2	Cannot Find DICOM Message File to Send.....	72
1.16.9.3	Invalid First 8-bytes in DICOM Message.....	72
1.16.9.4	Length-to-End Missing in DICOM Message.....	72
1.16.9.5	SOP Class UID Missing in DICOM Message	73
1.16.9.6	Invalid SOP Class UID in DICOM Message.....	73
1.16.9.7	Unknown Presentation Context ID	73
1.16.9.8	Unexpected PDU (expected 04 = P-Data)	74
1.17	Errors that may occur in Remote Procedure Calls.....	74
1.17.1	MAG DICOM AUDIT COUNT (COUNT^MAGDRPC7)	74
1.17.2	MAG DICOM FIND LOCATION (FINDLOC^MAGDRPC8)	74
1.17.3	MAG DICOM GET BASIC IMAGE (IMAGE^MAGDRPC2)	74
1.17.4	MAG DICOM GET NEXT QUEUE ENTRY (NEXTIMG^MAGDRPC4)	75
1.17.5	MAG DICOM GET PLACE (GETPLACE^MAGDRPC8)	75
1.17.6	MAG DICOM GET SERVICE INFO (SERVICE^MAGDRPC2)	75
1.17.7	MAG DICOM HL7 POINTER ACTION (HL7PTR^MAGDRPC8).....	75
1.17.8	MAG DICOM IMAGE PROCESSING (ENTRY^MAGDIR8)	76
1.17.9	MAG DICOM INCORRECT IMAGE CT (CORRECT^MAGDRPC8).....	76
1.17.10	MAG DICOM LOOKUP STUDY (LOOKUP^MAGDRPC4).....	76
1.17.11	MAG DICOM NETWORK STATUS (ONOFFLINE^MAGDRPC5)	77
1.17.12	MAG DICOM PACS CUTOFF DATE (CUTOFF^MAGDRPC1)	77
1.17.13	MAG DICOM PACS MINIMUM SPACE (MINSPACE^MAGDRPC1).....	77
1.17.14	MAG DICOM QUEUE IMAGE (QUEUE^MAGDRPC3)	77

1.17.15	MAG DICOM QUEUE INIT (INIT^MAGDRPC4)	78
1.17.16	MAG DICOM ROUTE EVAL LOG (EVALLOG^MAGDRPC6).....	78
1.17.17	MAG DICOM ROUTE EVAL START (START^MAGDRPC5).....	78
1.17.18	MAG DICOM ROUTE NEXT FILE (XMIT^MAGDRPC5).....	79
1.17.19	MAG DICOM ROUTE PURGE DONE (PURGDONE^MAGDRPC6).....	79
1.17.20	MAG DICOM ROUTE REMOVE OBSO (REMOBSO^MAGDRPC6).....	79
1.17.21	MAG DICOM ROUTE REQUEUE (REQUEUE^MAGDRPC6).....	79
1.17.22	MAG DICOM ROUTE TRANSACT STS (TRANSTS^MAGDRTIM)	80
1.17.23	MAG DICOM ROUTE VALID DEST (VALDEST^MAGDRPC1)	80
1.17.24	MAG DICOM SET PACS PARAMS (SETPACS^MAGDRPC8)	80
1.17.25	MAG DICOM TEXT PROCESSING (ENTRY^MAGDHRS1).....	80
1.17.26	MAG DICOM UPDATE GATEWAY NAME (UPDTGW^MAGDRPC8)	81
1.17.27	MAG DICOM UPDATE SCU LIST (UPDTAPP^MAGDRPC8).....	81
1.17.28	MAG DICOM WORKSTATION VERSION (STATION^MAGDRPC1)	81

Table of Contents

Chapter 1 DICOM Error Messages

This chapter contains a number of hints for resolving problems.

1.1 License Limits

The MUMPS database server used by the VistA Imaging DICOM Gateway requires a large enough license to allow all users, background tasks, and network connections to operate simultaneously. Of course, the larger the license, the higher the fee. In the VistA Imaging DICOM Gateway Installation Guide, a minimum license size is specified for MSM as **16 users** and **8 network connections**, and for Caché as **16 users**. This size system will allow for a sufficient number of job partitions. When the MUMPS engine on a VistA Imaging DICOM Gateway has a smaller license, it may not function adequately, and the error messages shown in this section may appear.

An error message may appear when an attempt is made to start an additional task on the database server, while all license slots are being utilized. There are a number of situations where such a message can show up.

- C-Store connection
- Telnet connection

The discussion below provides more details about each of these.


As the actual usage of the database server increases, the first action should be to rethink the allocation of tasks over the various VistA Imaging DICOM Gateway processors. Only when additional tasks are needed, and all of the available processors have been “maxed out”, should an additional processor be acquired, or the license count on at least one processor be increased.

Contact the vendor of the MUMPS Database server (InterSystems Corporation) or a designated reseller to acquire licenses, or to increase license limits.

1.1.1 Limit on Number of C-Store Processes

For each C-Store processor, a background task will be started on the VistA Imaging DICOM Gateway. The license quota for number of users needs to be high enough to support all C-Store processes that are assigned to a single Image Gateway.


When the quota for “users” is not high enough to start all C-Store processes, the start up log of the failing MAG_CSTORE.EXE task will contain this error message:

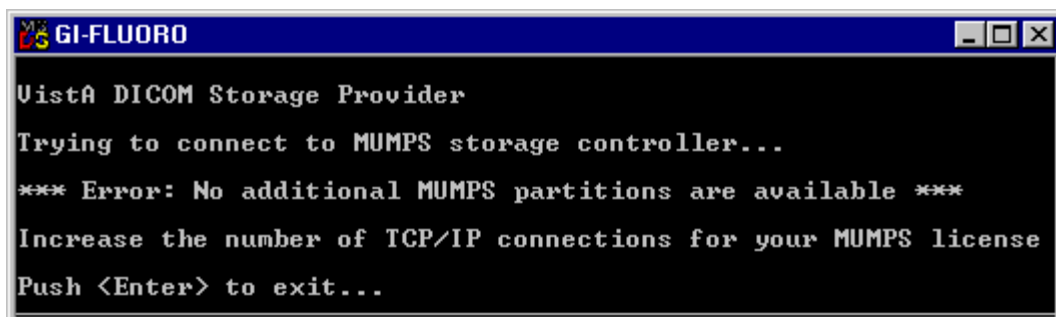

If this error message occurs, first try to re-distribute the various C-Store providers over the VistA DICOM Image Gateways. If all Image Gateways are maxed out, consider obtaining a license with a higher “user” quota.

1.1.2 Limit on Number of Telnet Sessions

Each telnet session counts as one user. When the “user” quota is not high enough to accommodate all needed interactive sessions, the attempt to launch a telnet session will produce this error message:

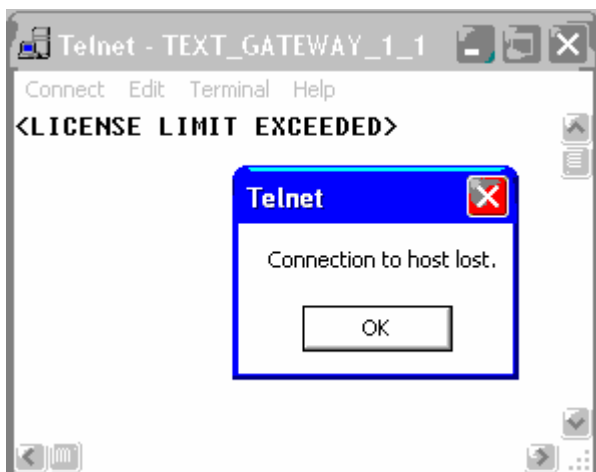


When this error occurs, the local system manager should use the Caché Control Panel (right-click on the Caché cube, select **Control-Panel**) to run a System Status, and terminate any non-essential telnet sessions. (Select the display that shows all running processes, right-click on redundant processes and select **Terminate**.)

When all running telnet sessions are essential, consider obtaining a license with a higher user quota.

1.2 Improper Start-Up Sequence

1.2.1 Start C-Store Program before MUMPS Server is Running

When an attempt is made to start a C-Store program (MAG_CSTORE.EXE) before the MUMPS (MSM or Caché) Server has been started, an error message like the following will appear:

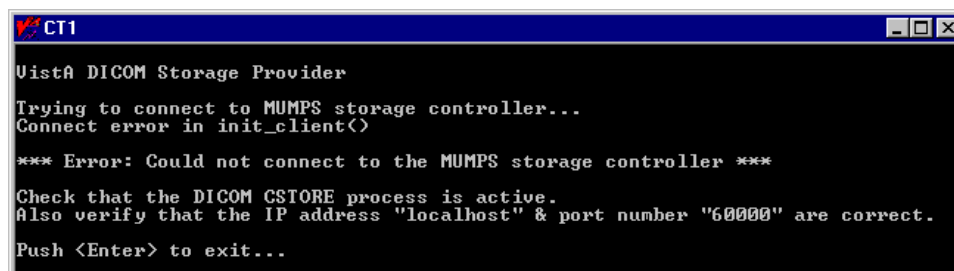

```

CT1
UistA DICOM Storage Provider
Trying to connect to MUMPS storage controller...
Connect error in init_client()


*** Error: Could not connect to the MUMPS storage controller ***

Check that the DICOM CSTORE process is active.
Also verify that the IP address "localhost" & port number "60000" are correct.
Push <Enter> to exit...
  
```

To resolve this problem, first start the MUMPS Server and then start the C-Store program.

1.2.2 Start C-Store Program while it is already Running

When an attempt is made to start an instance of a C-Store program, while a copy of that program for that instrument is already running, the following error message will be reported:

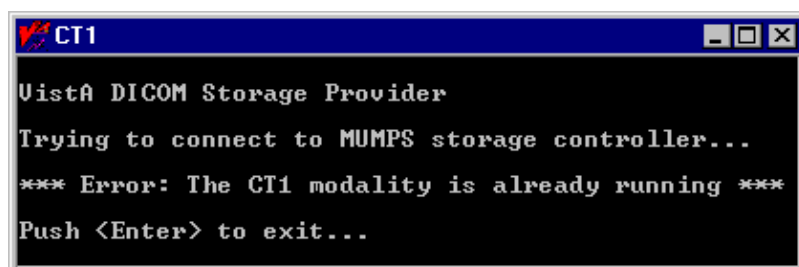

```

CT1
UistA DICOM Storage Provider
Trying to connect to MUMPS storage controller...


*** Error: The CT1 modality is already running ***

Push <Enter> to exit...
  
```

To resolve this problem, press <Enter> to terminate the second instance of the program.

1.2.3 Transmit Image before C-Store Program is Running

When an attempt is made to store an image while the C-Store program for the instrument is not running, an error message is likely to appear. The actual error message is instrument-specific. The following is a sample message produced by the Mallinckrodt Send_Image program.

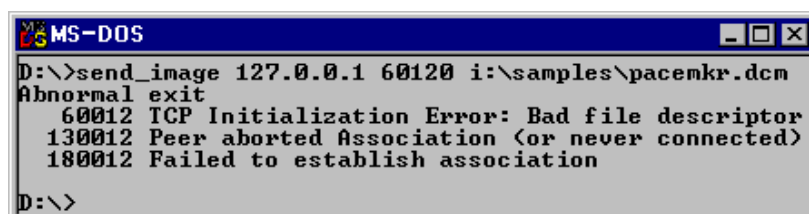

```


MS-DOS
D:\>send_image 127.0.0.1 60120 i:\samples\pacemkr.dcm
Abnormal exit
  60012 TCP Initialization Error: Bad file descriptor
 130012 Peer aborted Association (or never connected)
 180012 Failed to establish association

D:\>
  
```

To resolve this problem, first start the C-Store program by double-clicking on its icon, and then transmit the image.

1.2.4 Too many Instruments for one Image Gateway

Each instrument that transmits images to an Image Gateway uses one MUMPS background partition. When too many instruments are assigned to one single Image Gateway exceeding the maximum partition count, the following error message will appear:

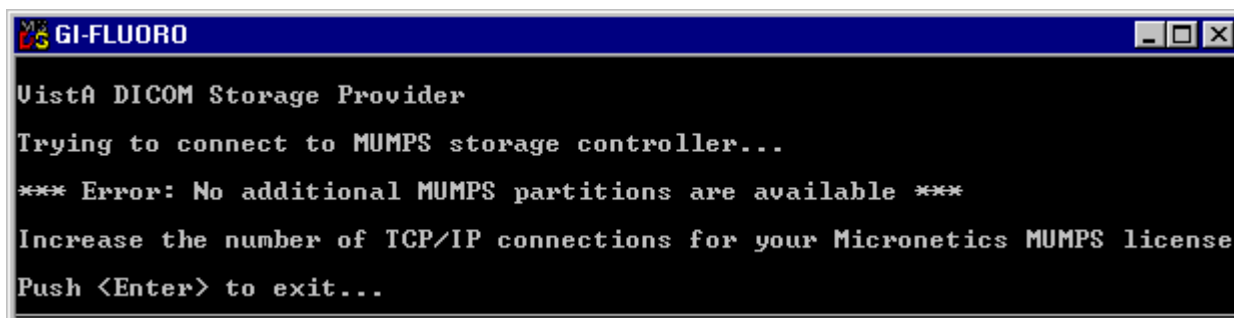

When this happens, first try to re-distribute the various instruments over the available Image Gateways. If that isn't possible, consider acquiring a license for more "users" for the processor.

1.2.5 Images are acquired, but are not moved to the File Server

When images are being acquired successfully, but are not being transported to the Windows File Server, this is most likely due to the fact that the "Process DICOM Images" task is not running. Check to make sure that it is running.

If "Process DICOM Images" is running, and images are still not moving to the Windows File Server, one of the following six problems might be happening:

- The patient/study information in the image header is incorrect and the association cannot be properly made.
- The network connection from the DICOM Image Gateway to the Windows File Server may be down.
- There may not be an entry in the Modality.DIC for the manufacturer-model-modality defined in the image header.
- There may be a run-away job on the VistA Imaging DICOM Gateway, either in MUMPS or elsewhere, which is preventing the image-processing task from having the necessary CPU resources. (Use Windows[™] Task Manager and ^%SS to further identify this problem.)
- There may be a problem with a run-away job on the VistA main HIS preventing the Kernel Broker from having the necessary CPU resources.
- There may be other problems.

1.3 Menu Options not Available

The availability of the menu options on the DICOM Gateway is affected by a number of issues:

- System configuration settings
- End-user privileges
- Ad-hoc operational issues

The very first time that the menu is displayed on a newly installed system, only those menu options that relate to set-up and configuration will be available. Menu options that are available are shown in bold face; menu options that are not available are shown in a fainter font and are enclosed in parentheses:

```
VistA DICOM Gateway Menu --
```

```
1  (Text Gateway)
2  (Image Gateway)
3  (Routing Gateway)
4  System Maintenance
5  Quit
```

```
OPTION:
```

1.3.1 System Configuration Settings

Using menu option **Update Gateway Configuration Parameters** (menu path 4-2-2), a number of parameters can be (re)defined. These parameters determine how the DICOM Gateway can be used. Three of these parameters determine which main menu groups will be available.

When these parameters are set as indicated in the sample dialog below, only the menu options for the “Text Gateway” will be available.

Sample dialog:

```
Will this computer be a DICOM Image Gateway? NO <Enter>
Will this computer be a DICOM Text Gateway? YES <Enter>
Will this computer be a Routing Processor? NO <Enter>
```

Sample menu:

```
VistA DICOM Gateway Menu -- Test System
```

```
1  Text Gateway
2  (Image Gateway)
3  (Routing Gateway)
4  System Maintenance
5  Quit
```

```
OPTION:
```

1.3.2 End-User Privileges

A user can be logged in on a DICOM Gateway as either a “fully privileged” user or as a “view only” user. When a user is “fully privileged”, typically all menu options will be available. When a user has “view only” access, all menu options that can be used to modify data will be made unavailable. For instance, a “view only” user will see the following menu on the Text Gateway:

Vista DICOM Text Gateway -- Test System

- 1 (Start Processing Text Messages from HIS)
- 2 (Send DICOM Text Messages to Commercial PACS or Broker)
- 3 Display Text Gateway Statistics
- 4 Display Modality Worklist Statistics
- 5 Modality Worklist Query
- 6 Display a HL7 Message
- 7 Display a DICOM Message
- 8 (Modify the HL7 Message Pointer)
- 9 Generate a Daily Summary Report
- 10 (Purge Old Modality Worklist Entries)
- 11 (Purge Old DICOM Message Files)
- 12 (Purge Old HL7 Transaction Global Nodes)
- 13 (Purge Old Audit Records)

OPTION:

1.3.3 Ad-Hoc Operational Issues

Certain menu options can only be executed after certain configuration steps have been completed. For instance, there are menu options that can send e-mail messages when errors occur. Those menu options will not be available until an e-mail address has been configured (using menu option **Update Gateway Configuration Parameters** (menu path 4-2-2). When menu options are made unavailable because no e-mail address has been set up yet, the following message will appear at the bottom of the menu:

Vista DICOM Text Gateway -- Test System

- 1 (Start Processing Text Messages from HIS)
- 2 Send DICOM Text Messages to Commercial PACS or Broker
- 3 Display Text Gateway Statistics
- 4 Display Modality Worklist Statistics
- 5 Modality Worklist Query
- 6 Display a HL7 Message
- 7 Display a DICOM Message
- 8 Modify the HL7 Message Pointer
- 9 Generate a Daily Summary Report
- 10 Purge Old Modality Worklist Entries
- 11 Purge Old DICOM Message Files
- 12 Purge Old HL7 Transaction Global Nodes
- 13 Purge Old Audit Records

You must first use option 4-2-2 to specify E-mail information.

OPTION:

Other menu options require that the location of the DICOM Gateway must be properly configured. The location-information is also set up using menu option **Update Gateway Configuration Parameters** (menu path 4-2-2). When the location-information is missing, the following error message will appear:

VistA DICOM Text Gateway -- Test System

- 1 (Start Processing Text Messages from HIS)
- 2 (Send DICOM Text Messages to Commercial PACS or Broker)
- 3 (Display Text Gateway Statistics)
- 4 (Display Modality Worklist Statistics)
- 5 Modality Worklist Query
- 6 Display a HL7 Message
- 7 Display a DICOM Message
- 8 Modify the HL7 Message Pointer
- 9 (Generate a Daily Summary Report)
- 10 Purge Old Modality Worklist Entries
- 11 Purge Old DICOM Message Files
- 12 Purge Old HL7 Transaction Global Nodes
- 13 (Purge Old Audit Records)

You must first use option 4-2-2 to obtain Location information.

OPTION:

Finally, the menu program can be started in two different ways. Normally, the menu program is started after an end-user logs in using an access and verify code that is valid on the VistA system. When proper credentials have been provided, the user can use the various menu options and perform actions that may require that Remote Procedures are executed on the VistA system.

The other possibility is that an end-user logs in using credentials that are used for local maintenance on the DICOM Gateway only. In this case, the DICOM Gateway will not be able to call any Remote Procedures on the VistA, but menu options will be available to perform any local maintenance that may be necessary. When a user is logged in using “local maintenance” credentials, a warning message is displayed at the bottom of the menu:

VistA DICOM Gateway Menu -- Test System

- 1 Text Gateway
- 2 Image Gateway
- 3 Routing Gateway
- 4 System Maintenance
- 5 Quit

Currently not authorized to invoke RPCs on VistA.

OPTION:

1.3.4 Logins Disabled

Usually, end-users are permitted to login into their VistA system. However, there are times when this permission is temporarily withdrawn. This may happen in a number of cases:

- automatically, when there are too many users on the VistA system
- intentionally, when maintenance on the VistA system is scheduled
- casually, for the duration of a single KIDS installation

When any of these situations occurs, logins will be temporarily disabled, and all attempts to execute Remote Procedure Calls will fail.

In such a case, a DICOM Gateway will report a failure to connect to VistA:

```
Could not connect to M-to-M Broker
(0,SIGNON - - - Signons not currently allowed on this processor.)
```

When a Remote Procedure was called for an interactive application option, the software will return to the menu, and the end-user can restart that application at a later time.

When the RPC was called from a non-interactive application, the software will wait for 5 minutes, and then attempt to reconnect to VistA.

```
*** Error: <ECODE>CALLRPC+14^MAGM2VCU:::6:27:
(,M13<NORPC> 0SIGNON - - - Signons not currently allowed on this processor.,) ***

Attempting to reconnect...

Cannot re-connect to VistA at this time.
0,SIGNON - - - Signons not currently allowed on this processor.

*** Wed 15:49 A time-out has occurred on the M2M Broker ***
*** (communication with the VistA System). ***
*** This system will wait 5 minutes and try to link again ***

Attempting to reconnect...
... reconnected
```

1.4 Privileges and Permissions

By the nature of the information that is stored and manipulated in a hospital environment, many aspects of the VistA Imaging DICOM Gateway are subject to access restrictions. It is important to keep track of the various user-accounts that have to be created and utilized to ensure proper access for the various parts of the system.

When any part of the software reports an access restriction type of error, the most common causes are that...

- The current user is logged in using a user-account that does not have the privileges needed to perform the action at hand
- A password has been changed, and the current user has not yet updated the registered copy of that password to the current value.

It is also important to be aware that...

- Some access restrictions apply to the current user
- Some access restrictions apply to the file at hand
- Some access restrictions apply to the disk-share where a file is stored
- Different access privileges may apply for reading an existing file and for creating a new file

Typical error messages that may appear are “Image file subdirectory creation error”, “Cannot create the image subdirectory: -5”, “Cannot find file xxxx” and “Access Denied”.

1.5 Disk Full

In the course of doing their work, the various processes create temporary data files which are later removed. On occasion, it may be necessary to manually remove data files that have been inadvertently left behind, perhaps by a run-away task (for example, one that is trying to read an incorrectly constructed DICOM image file). The sections below describe several utility functions that may be of help in removing this old computer litter.

There are two major groups of problems that result in “disk full” error messages:

- Situations where all space within a MUMPS database has been used up.
- Situations where all space on a physical medium has been used up.

1.5.1 Full MUMPS Databases

1.5.1.1 Data Stored in Wrong Database

The VistA Imaging DICOM Gateway locally *only* stores and maintains data in *its own* database (volume group DCM). Any error messages that indicate issues with storage in other databases (e.g. MGR or IMx) would indicate that data has been stored in a wrong place. When this happens, the spurious data should be removed from the database where it should not have been stored. (Of course, the VistA Imaging DICOM Gateway remotely stores data in ^MAG(2005) and other translated globals.)


1.5.1.2 Purge MUMPS Error Log

When a database “fills up”, **DO NOT** increase the size of the database container file. Under normal circumstances, 15MB is more than enough to hold the DICOM database. As a safety margin, the database is sized at 50MB. When the database becomes “full”, this indicates an error by itself. First, the root-cause of this error should be determined, then the extraneous data should be removed, and after this normal processing can continue with plenty of space in the 50MB database.


The most likely cause for a “database full” situation is that a burst of similar errors caused the error log to overflow the database.

When this type of error occurs, execute the following steps (preferably, call customer support for help in executing these steps):

1. Start the Caché Control Panel (right-click on the Caché cube and select **Control Panel**).
2. In the window for the Caché Control Panel, click on Processes.
3. In the right-hand pane, a list of all active Caché processes will be displayed. Identify the active telnet sessions and terminate all of them (right-click on the line for the process and select **Terminate** from the pop-up menu).
4. In the window for the Caché Control Panel, expand the tree for Local Databases, then right-click on the folder for the DICOM database, and select **Integrity Check**.

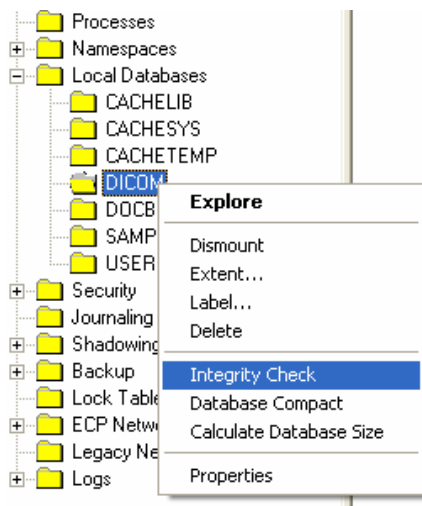

5. A new window will appear; a log of the integrity check will be compiled in this window. If all is well, the log will end with the text **“No Errors were found in this directory”**:

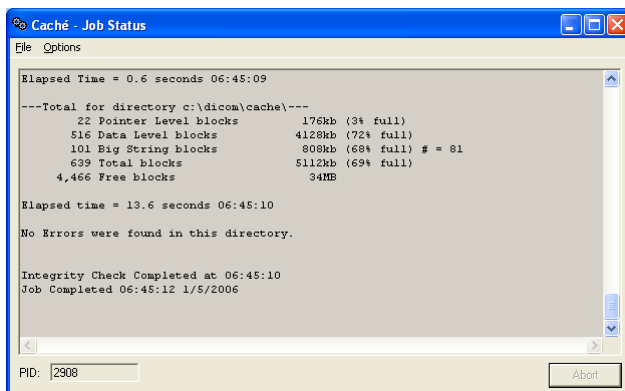


6. If the database is corrupted, first resolve all issues that are reported, using ^REPAIR. (Learning to use ^REPAIR requires a multi-day class from the database vendor. Most of the information learned in this class is “under non-disclosure restrictions”. Call customer support if a database repair is needed.)
7. Special cases are global variables ^TMP and ^MAGDMLOG. When either of these global variables has become corrupted, it’s not worth the time to try and fix it. Use

^REPAIR to remove the global variable from the Global Directory.

8. Once it is established that it is safe to use the database, proceed with the next steps, otherwise continue repairing the database.
9. Check the error-log: in the Caché Control Panel expand the tree for **Logs**, and drill down to **Application Errors** → **DICOM** → current date. If the error log indicates that there is a persistently recurring software problem, be sure to save a copy of one of those error reports for the Development Team.
10. Once the important error-logs are saved, erase all error logs.
11. Remove temporary data:
 - a. Open a console window (right-click on the Caché cube and select **Terminal**)
 - b. Issue the commands:

```
> zn "DICOM" <Enter>
> Kill ^TMP,^MAGDMLLOG <Enter>
> Halt <Enter>
```

12. Check how much space is now available (in the window for the Caché Control Panel, left-click on the text **Local Databases** in the left-hand panel). A survey of the statuses of all local databases will appear in the right-hand panel. Look at the line for the DICOM database.

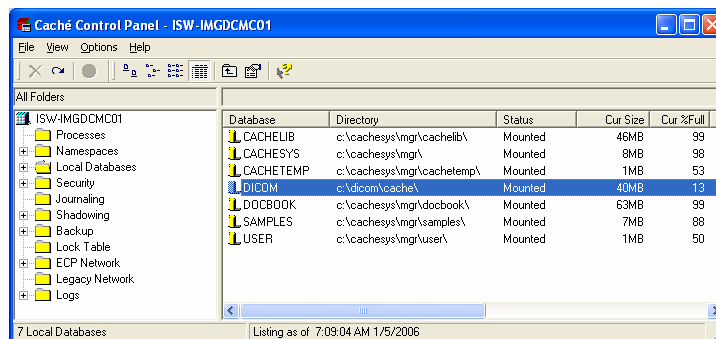

If the amount of free space in the DICOM Database is more than 15%, the system can be put back into operation. If necessary, shut down Caché and restart it, so that all background processes will be running again.

Once the system is back up and running, execute routine purge functions.

Purging the files and HL7 pointers to 30 days should clear sufficient disk space. (Options 10, 11, and 12 on the menu)

The U-V/W-X queue pointers should be checked.

MUMPS error logs reside in the MUMPS databases on each VistA Imaging DICOM Gateways in global variable ^ERROR.

Old entries in the MUMPS error log can be purged using a menu-option.

1.5.1.3 Purge Old Modality Worklist Entries

The Modality Worklist database resides in the MUMPS database on the DICOM Text Gateway in global variable ^MAGDWLST.

Old entries in the Modality Worklist database can be purged using a menu-option.

1.5.1.4 Purge Old HL7 Transaction Global Nodes

The HL7 transactions are stored in the main hospital VistA MUMPS database in global variable ^MAGDHL7.

Old nodes in that global variable can be purged using a menu-option.

1.5.1.5 Purge Old Audit Records

Audit records that pertain to the whole system are stored in the main hospital VistA MUMPS database in global variable ^MAGDAUDT.

Old audit records can be purged using a menu-option. This allows the global space to be recovered for other use.

1.5.1.6 Purge DICOM Image Processing Error Logs

There are four DICOM log files that hold images that fail to be properly processed:

- DICOM Failed Images – wrong patient/study id – ^MAGD(2006.575)
- DICOM Undefined Modalities – unknown modalities - ^MAGDINPT(2006.5713)
- DICOM Fixed Image – fixed image files - ^MAGDINPT(2006.5712)
- DICOM Error Log – holds deleted images – ^MAGD(2006.599)

Use of the DICOM Failed Image Correction procedure (see Chapter 7 of the VistA DICOM Gateway User Manual) automatically deletes entries from the first three. No purging is currently available for the DICOM Error Log, however.

1.5.2 Full Physical Media

1.5.2.1 Purge Old DICOM Message Files

On the DICOM Text Gateway, the DICOM message files reside in directories with names like **D:\DICOM\DATA n** , where n is an integer number, typically 1 and usually not higher than 2.

Old DICOM message files are purged automatically by the message generating application. They can also be purged using a menu-option.

The entire subdirectory tree can be purged by running **INIT_DICOM.BAT** in the **DATAn** directory from the command line. (This has to run on a “quiet system” – be sure to stop all MUMPS processes before doing this.)

1.5.2.2 Purge DICOM Image Files

On the DICOM Image Gateway, the DICOM image files are stored in the **D:\DICOM\IMAGE_IN** and **D:\DICOM\IMAGE_OUT** directories for incoming and outgoing images respectively.

Normally, all image deletion from these two directories is under automatic program control. Incoming images are deleted from the **D:\DICOM\IMAGE_IN** directory as they are processed. Incoming images that cannot be processed remain in the directory until they are resolved using the DICOM Failed Image Correction application. Incomplete image files will remain in this directory for at least an hour before they are renamed by appending “_incomplete” to the filename. These files (*.dcm_incomplete) will remain in the **D:\DICOM\IMAGE_IN** directory for further research by site personnel and will require manual deletion.

Outgoing images are automatically deleted from the **D:\DICOM\IMAGE_OUT** directory as they are sent. However, they may remain if there is a problem with the transmission. Old *.dcm files in this directory can be manually deleted.

In either case, on an as-needed basis, use MS Explorer to delete the old image files.

1.5.2.3 Other Possible Issues

In addition to the possibility that a disk is physically full, there are a number of other reasons why an error message like “Insufficient disk space to store images. Waiting 10 minutes.” might appear:

- Images are not processed due to lost connection with server
- Images are not processed due to incorrect information in headers.
- Images are not processed due to failure of the “Process Image” process.

Such error messages are usually displayed in the telnet window for the session labeled **PACS Request Image Transfer (2_2)**.

Make sure that the session labeled **Process DICOM Images (2_3)** is running.

It is very likely that the directory **x:\DICOM\Image_In** contains files for images that require corrections. See the window for the session labeled **Image Status (2_5)** for details. The latter window will also contain information on how badly backed up the system may be (number of files to be processed). Disk space will become available as images are processed. Just keep an eye on the status window to make sure the image count is *decreasing*.

Some image files may be as large as 7MB. When a site needs to store a large number of these files, it is imperative to ensure that sufficient disk space is available.

It may be possible to use the DICOM FIX utility to correct images with demographic information that does not match the HIS information.

The VistA Imaging DICOM Gateway uses a site parameter that establishes the percentage of disk space that should remain free. If this percentage is too high, the system will start producing error messages about insufficient disk space when there still is a lot of disk space. To ensure that this parameter is set to a reasonable value, use VA-FileMan (Enter/Edit, file 2006.1):

```
PCT FREE SPACE DICOM MSGS: 10//
```

Typically, this value should not be larger than 15%.

1.6 Extremely Slow Performance

When a normally operating system begins to run very slowly, the first question to ask is “what changed?”

If a specific system is extremely slow and others are fine, then launch Windows™ Task Manager to determine what is running on the system. Is there a run-away task? (The problem might be as simple as a user who inadvertently clicked too many times on an icon and launched multiple tasks!)

The problem might be with the main Hospital Information System. There could be a “run-away” job on it or a situation where the MUMPS-to-MUMPS Broker connection is not functioning properly.

A slow DICOM Image Gateway may be related to the Windows™ file server problem. Is there adequate free disk space on the server for new images?

The problem could also be with the network. Check whether the PC network interface is properly set up, especially the speed and duplex settings. The settings should be the same on the VistA Imaging DICOM Gateway as on the hubs and routers. In a wide-area network situation, check that the WAN is operating properly.

If you suspect that the amount of work for the system might be exceeding the available resources, check the following possible causes:

- If bad performance mainly occurs during peak times, there may be a network bandwidth issue. All devices used for VistA Imaging must be placed on a switched rather than a shared connection. This will reduce delays due to contention and collisions on the device’s Ethernet segment. If the device is already on switched Ethernet, one could try upgrading the network for a higher throughput capacity.
- There could be a bottleneck on the Imaging Storage Server. This could be the case when the Storage Server has a shared Ethernet connection or is connected via a 10 Mbps or a 100 Mbps Ethernet. In that case, upgrading to 1000 Mbps (Gigabit Ethernet) may resolve the

problem.

- There could be an IP routing bottleneck. If the device is on a different TCP/IP subnet than the Imaging Storage Server, the router might be the bottleneck. It is recommended that all VistA Imaging devices be on the same subnet, so that there will be no need for any routing. The VistA Imaging DICOM Gateways and Background Processors must be on the same subnet as the Imaging Storage Server. It may be possible to “collapse” multiple subnets into a single one by migrating to a wider subnet mask.
- If the site uses VLANs, routing bottlenecks may not be completely eliminated. Double-check the definitions of the VLANs and the assignment of processors to subnets. Review the routing set-up.
- Check the Network Interface Card (NIC) and the network switch or hub to make sure that none of the components is using “automatic sensing” of speed. If one component is set to “auto” or “auto-detect”, and another component is expecting to use a specific speed, time may be lost in unnecessary negotiations. This feature can sometimes fail.

1.7 Network Connectivity Problems

Connectivity problems can exhibit themselves in a number of ways:

- Users are unable to login. They see messages like “Domain Controller Unavailable” when they attempt to login into the system.
- The window “**Network Neighborhood**” does not show any other computers
- There are error messages while booting up the computer, or when attempting to start a program.
- There are error messages related to TCP/IP, NetBIOS name resolution or NIC hardware in the system’s “**Event Viewer**”.
- There may be errors in configuration information.

If any of the above situations occurs, the following steps usually lead to a quick assessment of the actual problem.

1.7.1 Only one VistA Imaging DICOM Gateway (that used to work) is Affected

Check first whether the problem is localized in one computer, or shared by multiple ones. If the problem occurs only in one computer, check the following:

1. Check that all cables are properly connected.
2. Check that the information in the appropriate configuration file is correct (Instrument.DIC and PortList.DIC contain port numbers, SCU_List.DIC contains IP addresses and port numbers and AE Titles).
3. Try to “ping” the affected computer from itself (first try “**ping 127.0.0.1**” and then try “**ping assigned_IP_address**”). If either of these attempts fails, the IP settings are mis-configured.

4. Check the computer's IP configuration (utility program `ipconfig.exe`). Make sure that the computer has a usable IP address, subnet mask and gateway address. Attempt to reach a server or network resource that is known to be operational. Use the "ping" utility to see if network traffic is being routed correctly to the ailing system. Use `tracert.exe` to further verify the network pathway.
5. Try to "ping" the default gateway from the affected computer. If this fails, either the address for the default gateway is not entered, or it is wrong.
6. Make sure that no other system has the VistA Imaging DICOM Gateway's IP address. This can be easily checked by disconnecting the ailing system, and then trying to "ping" its IP address from another location. If there is a response, then obviously another computer is using the same IP address.
7. Make sure that no other MUMPS system has the same 3-character system name.
8. Make sure that Microsoft's DHCP is not being used, check that the affected computer has correct WINS and DNS information. Make sure that IP addresses of any incoming clients are in the Gateway's HOSTS table. This will prevent excessive name lookup.
9. Double-check the values for the subnet mask and the default gateway. If either of these is wrong, the processor may not be able to "see outside" its subnet. As a result, Image Gateways and Workstations may be unable to access either the modalities or the VistA Hospital Information System.
10. Double-check that the TCP/IP protocol is configured properly on the system in question by looking at the network properties. Confirm that TCP/IP is in the list of protocols and verify that it is configured with an IP address and an address for a default gateway. The default gateway must have an address that is on the same subnet as the system itself.
11. Use a utility program like `netperf.exe` (combined with `netserver.exe`) to measure the effective speed of the network connection.

1.7.2 Several Processors are Affected

When multiple computers are having problems, check the following:

1. Look for problems with the network or with the router. To see if network traffic is being routed correctly, use the "ping" utility on an affected computer to attempt to reach a server or network resource that is known to be operational.
2. Check the server(s). Is the domain controller up and running? Is it reachable from the network? Can the server be reached via the network? Are the network interface cards and protocols configured properly?
3. If it is possible to login into the Windows[™] Domain, but not into VistA, check whether VistA is running and "reachable" via the network.

4. On the server, use “**Server Manager**” to check “**Services**” to make sure that the appropriate services are running.
5. Check that the information in the appropriate configuration file is correct (Instrument.DIC and PortList.DIC contain port numbers, SCU_List.DIC contains IP addresses and port numbers and AE Titles).
6. Try to “ping” the server from any of the affected processors. If this attempt fails, try to “**ping 127.0.0.1**” on the server to make sure that the server’s IP address and settings are correct. From the server, try to “ping” its default gateway. From one of the affected stations, try to “ping” the server’s gateway.

7. Use the utility program `tracert.exe` to trace the access route to the server. E.g.

```
D:\>ipconfig <Enter>

Windows NT IP Configuration

Ethernet adapter DC21X41:

    IP Address. . . . . : 11.22.33.43
    Subnet Mask . . . . . : 255.255.255.0
    Default Gateway . . . . . : 11.22.33.45

D:\>tracert 11.22.33.45 <Enter>

Tracing route to wciofo-reglwan.eth100.net.va.gov [11.22.33.45]
over a maximum of 30 hops:

    1    <10 ms    <10 ms    <10 ms    wciofo-reglwan.eth100.net.va.gov [11.22.33.45]

Trace complete.

D:\>
```

8. If `tracert` is not successful, there is most likely a routing problem. Links that are probably down are indicated by asterisks in the report from `tracert`.
9. Try using **Windows™ Explorer** to see if it is possible to browse to the server. Users with Administrator privileges on the server have access to the “administrative shares”. They could attempt to map a drive to one of these “shares”. If such attempts fail, it is likely that the WINS/DNS setup is incorrect.

1.7.3 One or more windows show an error message

One error message that may show up on occasion is “MUMPS length network_read() error: -1”.

There are several possible causes for this error message:

- Network error
- Caché console is closed or crashed
- Disk is Full

When this error message appears, follow the steps below to find the root-cause of the problem:

- Verify connectivity by PINGing in both directions
- Verify connectivity by executing a DICOM echo
- Check that the Caché Server is still up and running (Caché cube is blue).
- Verify that there is sufficient disk space on the disk that contains the directories x:\DICOM\Data1 and x:\DICOM\Image_In
- Check that there is sufficient space in the MUMPS database for normal processing

1.7.4 PACS isn't Sending an Image Complete Message

This situation typically comes about when the TCP link with the PACS is broken.

Before attempting to restart any processes, examine the log files to ascertain that the exam complete message is being received from the PACS. The correct C-Move request should be issued by VistA, but response may not be returned from the PACS. Restarting the gateway process for the PACS (**1-2-x**) usually restores normal functionality

1.8 Other Configuration Issues

1.8.1 Loading Master Files

After (re)loading the Master Files, the error message “nnn error(s) encountered while building master files” may appear. When this message appears, it indicates that some error(s) have been found in one or more of the Master Files. The site-modifiable master files are:

- Instrument.DIC
- Modality.DIC
- PortList.DIC
- SCU_List.DIC
- WorkList.DIC

Check the contents of these files to repair any reported issues.

Most commonly, the “Institution Name” is not set in WorkList.DIC. The default value that is specified in the distributed file is: “<add local name here>”. It should be either the facility number or the name, e.g.:

```
AE_Title|523| ... or  
AE_Title|BOSTON, MA| ...
```

See the VistA Imaging DICOM Gateway Installation Guide for more information about properly formatted specifications in Master Files.

1.8.2 Incomplete Master Files

It may happen that not all modalities that are transmitting images to an image gateway are described in the master file named `Modality.DIC`. When this happens, the error message “The following images have undefined modalities” will appear.

In such a case, use the information about the modality that is presented in the error message to add the missing modalities to the master file.

Note: Don’t forget to re-load the master file.

1.8.3 Multiple Instruments are Configured to Use the Same Port Number

The VistA Imaging DICOM Gateways should be configured such that each instrument has its own unique port number assigned for communication. When multiple instruments share a port number, only the one that happens to connect first will be able to communicate with its gateway. Any other instrument attempting to use the same port number will cause an error message “Cannot listen to DICOM port because it is already being used”.

When this happens, verify that there are no conflicts in the master file `Instrument.DIC`. Reload the master file after it has been corrected.

When the configuration has been corrected, restart the VistA Imaging DICOM gateway. Just closing and re-starting the modality windows will not have the effect that any ports are released. The Caché Server must be restarted to release all ports.

1.8.4 Digitized Images Appear Small on the Diagnostic Workstation

The master file `Modality.dic` contains a parameter that provides a number of settings that influence the conversion from the original file format to TARGA™ format. See the VistA Imaging DICOM Gateway Installation Guide for a complete description of these parameters.

Pay special attention to the (sub)parameters **R1**, **R2**, **R4** and **R8** that may be used to create a “small” thumbnail version of the image (also known as “abstract” or “icon”). When one of these parameters is applied to the conversion of the “full” image, a severe loss of resolution may be the result.

1.8.5 Digitized Images Appear Dark

The master file `Modality.dic` contains a parameter that provides a number of settings that influence the conversion from the original file format to TARGA™ format. See the VistA Imaging DICOM Gateway Installation Guide for a complete description of these parameters.

Pay special attention to the (sub)parameter **Bnnn**. When the number of bits specified by this parameter is incorrect, the resulting images may appear with a severely reduced number of gray-values, typically showing extremely dark colors.

1.8.6 Modality Worklist Find Failure

The registration of a modality's AE Title in the master file `Worklist.DIC` must be an exact match with the string that the modality transmits. This value is case-sensitive.

After the master file `Worklist.DIC` has been corrected, verify that the values in the master files `Instrument.DIC` and `Modality.DIC` are also correct, and be sure to re-load any modified master files into the database.

1.9 Failover Procedure

It is possible to allow a gateway processor to take over the tasks of another, in case one of the gateways should fail. To allow for a smooth “fail-over”, take the following steps:

1. Add the TCP/IP address of the *failed* gateway to the substitute gateway processor. This will be the *second* IP address on the substitute gateway's Network Interface Card (NIC).
2. On the substitute gateway, open the folder that contains the instruments for the failing gateway and click on the modality icons within that folder to begin listening for the modalities on their respective ports. The substitute gateway can continue to work with its “regular” icons for the modalities that already were assigned to that processor.

Nothing will have to change on the modality configuration. It will be completely “transparent” to the modality. There might be a small performance degradation, however, and the screen will certainly be more busy.

1.10 Error Messages from VistA Imaging DICOM Gateway Executable Programs

1.10.1 Error Messages from MAG_CSTORE.EXE

The program `MAG_CSTORE.EXE` maintains a TCP socket connection with an image acquisition device on one side and with a MUMPS Server on the other side. It processes bi-directional data streams that are encoded according the DICOM standard and causes image files to be stored on its host-PC. It is called with the following parameters:

```
Mag_cstore <IP address> <port> <modality> [<debug level>]
```

<u>Parameter</u>	<u>Meaning</u>
IP address	For the VistA Imaging DICOM Gateway, this value is always “ localhost ”.
Port	For the VistA Imaging DICOM Gateway, the port number for the on which the MUMPS Server listens, which is always 60000 .
Modality	An abbreviation for the name of the instrument that is sending the images, typically a code like “ CR1 ”, “ CT3 ”, ...
Debug Level	A code that indicates the verbosity of debug messages: 0 =none, 1 =some, 2 =lots

The program MAG_CSTORE.EXE can produce the following error codes:

Error: Could not connect to the MUMPS storage controller
 Error in reading PDU header -- connection dropped
 Error in reading dataset -- connection dropped
 Error in writing echo response -- connection dropped
 Error in reading control dataset -- connection dropped
 The WinSock initialization has failed
 Error writing control PDU to disk
 Error writing control PDU to DICOM socket
 Error reading PDV header
 Error reading PDV data, *nnn* bytes input
 Disk read error in net_write_data()
 Unknown disk read error in net_write_data()
 Error writing PDU hdr to socket
 Error writing PDV hdr to socket
 Error writing DATA PDU to socket
 select() error in network_read()
 Error in network_read(), total read: *nnn*, errno=*nnn*
 MS NT Socket Error #*nnn*
 select() error in network_write()
 Error in network_write(), Error Number: *nnn*
 Socket initialization error in init_client()
 Connect error in init_client()
 ioctlsocket() (non-block) error in init_client()
 ioctl() (non-block) error in init_client()
 setsockopt(SO_LINGER) error in init_client()
 socket() error in init_server()
 gethostid() error in init_server()
 bind() error in init_server
 Could not listen on the DICOM port because it is already in use
 listen() error in init_server()
 ioctlsocket() (non-block) error in init_server()
 ioctl() (non-block) error in init_server()
 select() error in init_server()
 accept() error in init_server()
 setsockopt(SO_LINGER) error in init_server()
 MUMPS length network_read() error: *nnn*
 MUMPS message network_read() error: *nnn*
 Disk Read error
 Unknown read error
 Write error: number_written= *nnn* count= *nnn*
 Unknown write error
 Error: Incorrect Message *xxx* It should have been *xxx*
 Error: No additional MUMPS partitions are available
 Error Reading PDU Header

1.10.2 Error Messages from MAG_DCMTOTGA.exe

The program MAG_DCMTOTGA.exe reads files that are encoded according to the DICOM standard and produces files in TARGA™ format. It is called with the following parameters:

```
MAG_DCMTOTGA {input filename} {output filename} {parameter(s)}
```

<u>Parameter</u>	<u>Meaning</u>
A <i>nnn</i>	Add <i>nnn</i> to each pixel (before min/max check)
B <i>nnn</i> †	Number of bits in a pixel (stored in the TARGA™ header)
C <i>nnn</i>	Ceiling (maximum) pixel value (any value > <i>nnn</i> becomes <i>nnn</i>)
F <i>nnn</i>	Floor (minimum) pixel value (any value < <i>nnn</i> becomes <i>nnn</i>)
I	Invert each pixel
O <i>nnn</i> †	Byte offset in DICOM file to image
R1	Reduce the size of the image file a factor of two (a two-byte 8-bit pixel is stored in one byte, loss-less)
R2	Reduce the size of the image file a factor of two (a two-byte 16-bit pixel is stored as one byte, usually not loss-less)
R4	Reduce the size of the image file a factor of four (a square of four 16-bit pixels is averaged and stored in two bytes (16 bits) as one pixel, definitely not loss-less)
R8	Reduce the size of the image file a factor of eight (a square of four 16-bit pixels is averaged and stored in one byte as one pixel, definitely not loss-less)
R16	Reduce the size of the image file a factor of sixteen (a square of four by four pixels is averaged and stored in as one pixel, definitely not loss-less)
R32	Reduce the size of the image file a factor of thirty-two (a square of sixteen 16-bit pixels is averaged and stored in one byte as one pixel, definitely not loss-less)
S <i>nnn</i>	Subtract <i>nnn</i> from each pixel (unsigned arithmetic -- before add)
X <i>nnn</i> †	X-dimension of image (horizontal width or number of columns)
Y <i>nnn</i> †	Y-dimension of image (vertical height or number of rows)
† - Required	<i>nnn</i> represents an unsigned integer number

The environment variable MAG_DCMTOTGA_VERBOSE may be used to override the default settings:

MAG_DCMTOTGA_VERBOSE=1	Selects brief verbose mode
MAG_DCMTOTGA_VERBOSE=2	Selects total verbose mode
MAG_DCMTOTGA_VERBOSE=3	Turns on density histogram

The program MAG_DCMTOTGA.exe can produce any of the following error codes:

- 1 End of Input File
- 2 Read Error detected by `ferror()`
- 3 End of Output File
- 4 Write Error detected by `ferror()`
- 5 Maximum number of columns exceeded
- 999 Unknown Error

1.10.3 Error Messages from MAG_ABSTRTGA.exe

The program MAG_ABSTRTGA.exe reads files that are encoded in Targa™ format and produces reduced resolution (thumbnail, abstract) copies of these files. It is called with the following parameters:

```
MAG_ABSTRTGA.c {input filename} {output filename} [/s:x,y] [/8]
```

<u>Parameter</u>	<u>Meaning</u>
/s:x,y (Optional)	Modify the maximum X and Y dimensions.
/8 (Optional)	Force all black & white (monochrome) pixels to be 8-bits

The following environment variables can override the default settings:

ABSTR_VERBOSE=1	Select verbose mode
ABSTR_MAX_X=nnn	Maximum X dimension is <i>nnn</i>
ABSTR_MAX_Y=nnn	Maximum Y dimension is <i>nnn</i>

The program MAG_ABSTRTGA.exe can produce any of the following error codes:

- 1 End of Input File
- 2 Read Error detected by `ferror()`
- 3 End of Output File
- 4 Write Error detected by `ferror()`
- 5 Maximum number of columns exceeded
- 999 Unknown Error

1.10.4 Error Messages from MAG_DCM_COPY.exe

The program MAG_DCM_COPY.exe is typically called when the DICOM header is being updated either to change from Explicit VR to Implicit VR or when the values of the attributes are being changed. It is called with the following parameters:

```
Mag_dcm_copy {input file} {output file} {input offset} {number of bytes to copy}
```

<u>Parameter</u>	<u>Meaning</u>
Offset	The byte address where the copy program starts reading
Number of Bytes	The number of bytes to be copied to the output file. (Any information to be copied is appended to the output file.)

The program MAG_DCM_COPY.exe can produce any of the following error codes:

- 1 End of Input File
- 2 Read Error detected by `ferror()`
- 3 End of Output File
- 4 Write Error detected by `ferror()`
- 999 Unknown Error

1.10.5 Error Messages from MAG_TGATODCM.exe

The program MAG_DCMTOTGA.exe reads files in TARGA™ format and produces files that are encoded according to the DICOM standard. It is with the following parameters:

```
Mag_tgatodcm {input filename} {output filename} {number of bytes
to copy}
```

The program MAG_TGATODCM.exe can produce any of the following error codes:

- 1 End of Input File
- 2 Read Error detected by `ferror()`
- 3 End of Output File
- 4 Write Error detected by `ferror()`
- 999 Unknown Error

1.11 Run Time Error Messages

The system error log contains a list of any errors that have occurred within the Caché Server. It is good practice to review this log on a regular basis and to report important errors.

Some entries in the error log are the result of the normal operation of the system and can be ignored.

The Caché system may issue the following error messages:

Caché Error	Description
<BLOCKNUMBER>	A reference has been made to a block outside the range of the database file.
<CANNOT GET THIS PROPERTY>	There has been an attempt to get a property of a class for which getting this property is invalid.
<CANNOT SET THIS PROPERTY>	There has been an attempt to set a property of a class for which setting this property is invalid.
<CLASS COMPILING>	There has been an attempt to instantiate a class or invoke a class method of a class which is currently being recompiled on the local system.
<CLASS DESCRIPTOR>	There has been an attempt to run a routine which is actually a class descriptor.
<CLASS DOES NOT EXIST>	A reference has been made to a nonexistent class.
<CLASS EDITED>	There has been an attempt to use an object hosted on the local system whose class has been recompiled from a remote system since the object was created.
<CLASS RECOMPILED>	There has been an attempt to use an object hosted on the local system whose class has been recompiled on the local system since the object was created.
<CLASS TOO BIG TO LOAD>	A class cannot be used because its class descriptor is too large to fit into a routine buffer.
<CLASS TOO BIG TO SAVE>	A class cannot be created because its class descriptor is too large to fit into a routine buffer.

Caché Error	Description
<CLIENT-SERVER MISMATCH>	A network request cannot be processed due to incompatibility between the client and server.
<CLUSTERFAIL>	A cluster member has failed during global buffer lock processing.
<COLLATEMISMATCH>	Subscript level mapping failed due to misconfigured collation type.
<COLLATION NOT SUPPORTED>	A reference has been made to a global whose collation type is not supported on the current system.
<COMMAND>	A command has been used improperly in this context, such as an argumentless GoTo in a routine.
<COMMITFAIL>	Received during a COMMIT when Caché receives an error while processing a TCommit. This error means that Caché is not sure whether one or more remote machines actually processed the commit.
<CONFLICTING BLOCK NUMBERS>	There has been an attempt to reserve a block that was already reserved.
<CORRUPT OBJECT>	An internal object system error occurred. Contact InterSystems Worldwide Response Center if this error occurs.
<CORRUPT VOLUME SET>	The volume set is corrupted. Usually, this means the label on the volume set is wrong. Use the ^LABEL utility to correct it.
<CP NOT STARTED>	One of the major processes required for proper operation of the system failed to start. This is potentially a very serious system error; notify your system manager.
<DATABASE MAP LABEL>	There is an invalid label in a database map block.
<DATABASE>	Caché has detected degradation in this database (this is potentially a very serious system error; notify your system manager).
<DIRECTORY>	There is no such directory on the target system, no Caché database, the Caché database is not mounted, or the database is locked by another configuration.
<DISCONNECT>	A TCP disconnect has been detected while a long-duration request is being processed.
<DISKHARD>	Caché has encountered an uncorrectable disk hardware error (this may also be the result of a database problem; notify your system manager).
<DIVIDE>	There has been an attempt to divide by zero.
<DOMAINSPACERETRY>	Repeated attempts to contact the domain space master have failed.
<DSCON>	There has been an attempt to read from a disconnected terminal.
<DYNAMIC LIBRARY LOAD>	An error has occurred during an attempt to load a dynamic library via callout. See cconsole.log for additional information.
<ECODETRAP>	A user-generated software trap was generated by setting the \$ECODE system variable to a non-null string value.

Caché Error	Description
<EDITED>	There has been an attempt either (1) to return to an edited routine or edited line, (2) to edit or save a routine was being edited, or (3) to issue a Use command on a file for the first time after opening it and the routine issuing the Use command has been edited.
<ENDOFFILE>	There has been an attempt to read past the end-of-file marker of a sequential file.
<FILEFULL>	Caché attempted to allocate a disk block for more global data or routine storage, but the attempt failed because the Caché database is full and could not be expanded.
<FRAMESTACK>	The routine has too many nested calls to Do, For, Xecute, New, or user-written functions.
<FUNCTION>	The specified function does not exist or is being used improperly.
<GARBAGE COLLECTOR FAILED>	The one of the processes that reclaims space in the database has failed. This is potentially a very serious system error; notify your system manager.
<HALTED>	An internal error message.
<ILLEGAL VALUE>	There has been an attempt to use a negative value where one is not allowed, such as, for \$X or \$Y.
<INSUFFICIENT CLASS MEMORY>	A class cannot be used because Caché has run out of shared memory.
<INTERNAL OBJECT ERROR>	An internal object system error. Contact Intersystem's Worldwide Response Center if this error occurs.
<INTERRUPT>	A user has interrupted the routine. (In many implementations, the user has pressed CTRL-C.)
<INVALID ARGUMENT>	There is an invalid argument prototype in the zfentry specification of a callout function.
<INVALID BIT STRING>	The bit string used in a bit string operation is not valid.
<INVALID CLASS>	There has been an attempt to use a class that has been corrupted. Recompile the class and try again.
<INVALID OREF>	No object with the specified OREF is currently in memory.
<INVALID TYPE>	An OREF has been used where not allowed.
<Java Exception>	An exception occurred during a call into the Java runtime environment.
<Java VM not loaded>	No Java Virtual Machine is available.
<LABELREDEF>	A routine has a duplicate label within it. Labels must be unique within the routine.
<LANGUAGE MISMATCH>	While compiling and inserting code into an existing routine, the current language mode differs from that of the routine.
<LICENSE LIMIT EXCEEDED>	There has been an attempt to exceed the number of processes allowed on your system under the current license.
<LICENSE SERVER UNAVAILABLE>	The license server is unreachable at the moment. Check your network.
<LIST>	An improperly-formed list has been used.
<LOCKLOST>	Some locks once owned by this job have been reset.
<LOGIN INHIBITED>	The system is initializing. No users are permitted to begin work.
<MAGTAPE>	A magnetic tape operation encountered an error. Check \$ZA.
<MAXARRAY>	There are too many subscripts at this level.

Caché Error	Description
<MAXNUMBER>	An arithmetic operation has produced a number larger than the implementation allows.
<MAXSCOPE>	There has been an attempt to issue more than 31 levels of New commands.
<MAXSTRING>	There has been an attempt to specify or create a data string longer than the implementation allows (32,767 characters).
<METHOD DOES NOT EXIST>	The method does not exist in the specified class or the class of the specified object.
<MNEMONICSPACE>	There has been an attempt to use control mnemonics for a device with no associated mnemonic space.
<NAKED>	There has been an attempt to use a naked global reference when the naked state was undefined.
<NAME>	There is invalid syntax in a name.
<NAMEADD>	There has been an overflow of device name table, resulting from the Open command.
<NAMESPACE>	The specified namespace is undefined or not active.
<NETFORMAT>	There has been an error in a network message. The remote system found fault with the format of a request. Call your support center to resolve this serious error.
<NETGLOREF>	There has been an error in a network message. The remote system found fault with the format of a request. Call your support center to resolve this serious error. (This can also occur if a remote global is a string-collated or 7-bit encoded global on a DSM system.)
<NETJOBMAX>	Another high-speed networking process cannot be added. This is usually due to an insufficient number of global buffers.
<NETLOCK>	A Caché ObjectScript Lock command has been attempted to a remote computer whose remote system index is greater than 31. To correct, redefine your network configuration to include fewer than 32 remote computers.
<NETRETRY>	An operation failed at the network level in a way that could be immediately retried.
<NETSRVFAIL>	During a transaction COMMIT or a Set, Kill, ZKill command, a client system has detected that one of the servers involved in the transaction has restarted while the transaction was open. Caché
<NETVERSION>	The client and server systems are running different DDP versions, which cannot accept each other's message format.
<NETWORK DATA UPDATE FAILED - BLOCKNUMBER>	Updates sent over the network were lost because the remote system attempted to refer to a block that is outside the bounds of the database; notify your system manager.
<NETWORK DATA UPDATE FAILED - CLIENT-SERVER MISMATCH>	Updates sent over the network were lost because a network request could not be processed due to incompatibility between the client and server.
<NETWORK DATA UPDATE FAILED - CLUSTERFAILED>	Updates sent over the network were lost because a cluster member failed during global buffer lock processing.

Caché Error	Description
<NETWORK DATA UPDATE FAILED - DATABASE>	Updates sent over the network were lost because Caché on the server has detected degradation in this database. This is potentially a very serious system error; notify your system manager.
<NETWORK DATA UPDATE FAILED - DIRECTORY>	Updates sent over the network were lost because the referenced directory is not on the remote system.
<NETWORK DATA UPDATE FAILED - DISKHARD>	Updates sent over the network were lost because Caché on the server has encountered an uncorrectable disk hardware error. This may also be the result of a database problem; notify your system manager.
<NETWORK DATA UPDATE FAILED - FILEFULL>	In networked operation, a <FILEFULL> error occurred.
<NETWORK DATA UPDATE FAILED - MAXSTRING>	There has been an attempt to specify or create a data string longer than the implementation allows (32,767 characters).
<NETWORK DATA UPDATE FAILED - NETFORMAT>	Updates sent over the network were lost because the remote system found fault with the format of a request. Call your support center to resolve this serious error.
<NETWORK DATA UPDATE FAILED - NETGLOREF>	Updates sent over the network were lost because the remote system found fault with the format of a request. Call your support center to resolve this serious error.
<NETWORK DATA UPDATE FAILED - NETVERSION>	The client and server systems are running different ECP versions, which cannot accept each other's message format.
<NETWORK DATA UPDATE FAILED - PROTECT>	In networked operation, a <PROTECT> error occurred.
<NETWORK DATA UPDATE FAILED - STRINGSTACK>	In networked operation, a <STRINGSTACK> error occurred.
<NETWORK DATA UPDATE FAILED - STRMISMATCH>	Updates sent over the network were lost because there has been an internal error handling big strings over the network.
<NETWORK DATA UPDATE FAILED - SUBSCRIPT>	In networked operation, a <SUBSCRIPT> error occurred.
<NETWORK DATA UPDATE FAILED - SYSTEM>	Updates sent over the network were lost because a <SYSTEM> error occurred on the server. Either there has been an attempt to do something not allowed by the operating system or ECP. Or there is an error condition in Caché, in which case you should notify your support center with as much information as possible.
<NETWORK DATA UPDATE FAILED - WIDECHAR>	In networked operation, a <WIDECHAR> error occurred.
<NETWORK DATA UPDATE FAILED>	Updates sent over the network were lost. The reasons for the loss are undetermined. Call your support center to resolve this serious error.

Caché Error	Description
<NETWORK>	Typically, one of the following has occurred: the network timeout has expired; the local port has gone down; the node being accessed is down; or the remote server connection is disabled.
<NLS TABLE>	There has been an attempt to perform NLS translation using data that is not proper for the conversion table.
<NO CURRENT OBJECT>	There is no current object.
<NO MAILBOX>	A resource needed for inter-process communication is unavailable.
<NO SOURCE>	A source line is missing from a routine in the routine source global.
<NODEV>	There has been an attempt to a write-only device, or write to a read-only device, with interjob communication.
<NOJOB>	There has been an attempt to specify an incorrect process number in a View command, or an error occurred in a Job command.
<NOLINE>	There has been an attempt to refer to a nonexistent routine line.
<NOROUTINE>	There has been an attempt to refer to a nonexistent routine.
<NOSYS>	There has been an attempt to make an extended or implicit reference to a remote system that is not reachable in the current network configuration. The remote system is not in the tables, or is a DSM system that has not advertised itself.
<NOT PRIMARY VOLUME>	Volume sequence is not 1; the volume label disagrees with the function of the volume.
<NOTOPEN>	The device cannot be opened, or there has been an attempt to use an unopened device.
<NULL VALUE>	A null string appears where one is not allowed.
<OUT OF \$ZF HEAP SPACE>	The \$ZF heap lacks the necessary available space to support one of the input or output parameters being passed between Caché and the external program invoked via the \$ZF function.
<PARAMETER>	The number of parameters passed to a labeled line by a user-written function reference or a Do command exceeded the number of formal parameters declared for the labeled line.
<PRIVATE METHOD>	There has been an attempt to invoke a private and, therefore, unavailable method.
<PRIVATE PROPERTY>	There has been an attempt to access a private and, therefore, unavailable property.
<PROPERTY DOES NOT EXIST>	The property is not part of the class of the specified object.
<PROTECT>	There has been an attempt to do something with a global (Read, Write, or Kill) for which there was no authorization; or there has been an attempt to use a View command which modifies memory, \$View or \$ZUtil(49); or there has been an attempt to use a nonexistent directory, possibly with extended global syntax, or some other protection violation occurred.
<RANGE>	A bit position is out of allowable range.
<READ>	The record cannot be read.

Caché Error	Description
<RECOMPILE>	A routine has been compiled under a different version of Caché or an Intersystem's legacy product. It cannot be loaded onto this system with ^%RIMF, which transfers object code. Transfer it as source code (using ^%RO and ^%RI, or the Caché Explorer) and then recompile it.
<REMOTE CLASS EDITED>	There has been an attempt to use an object hosted on a remote system whose class has been recompiled from a remote system since the object was created.
<REMOTE CLASS RECOMPILED>	There has been an attempt to use an object hosted on a remote system whose class has been recompiled from the local system since the object was created.
<ROLLFAIL>	Caché has encountered an error processing a call to TRollBack. This error means that Caché is not sure whether one or more remote machines actually processed the rollback.
<ROUTINELOAD>	An error occurred in loading a routine. This is database degradation; notify the system manager.
<SELECT>	A \$Select function contains no true condition.
<SLMSPAN>	There has been an attempt to kill a global across a subscript level mapping boundary.
<STACK>	The argument stack is out of room or contains an incorrect type.
<STORE>	The partition space for the process is exhausted. The process needs a block larger than the largest contiguous block of available space in the partition. Even though \$S may show available space, it is fragmented. Reduce the number or size of the local variables. Utilities may not be available until some variables are deleted in order to free up space. Issuing ZLOAD for too large a routine from a tape or sequential file can also cause this error.
<STRINGSTACK>	An expression is too long, there are too many expressions in an argument for a single command, or an expression contains many very long strings. Simplify the expression.
<STRMISMATCH>	There has been an internal error handling big strings over the network.
<SUBSCRIPT>	A subscript has an illegal value or a global reference is too long. If a global reference is too long, then either (1) the encoded length of an individual subscript exceeds 255 bytes or (2) the sum of the encoded lengths of the following items exceeds 512 bytes: the length of the global name; the lengths of all subscripts; the number of subscripts; the number of numeric subscripts; and the number of negative numeric subscripts.
<SYNTAX>	There is a syntax error (an error in the formation of a language constructs, such as a misspelled or missing keyword).
<SYSTEM>	Either there has been an attempt to do something not allowed by the operating system, or there is an error condition in Caché, in which case you should notify your support center with as much information as possible.
<TERMINATOR>	There has been an attempt to read on a terminal or device in image mode with no terminator and it was not a fixed-length read.

Caché Error	Description
<TOO MANY CLASSES>	This process has attempted to access too many active classes.
<TOO MANY OREFS>	This process has attempted to create too many simultaneously open objects.
<TOO MANY USERS OF CLASS>	Too many processes are trying to use a particular class simultaneously (more than 65561).
<TOO MANY USERS>	Too many users are attempting to use the system at the same time.
<TOOMANYFILES>	Caché is unable to open a file because the underlying operating system has run out of file descriptors.
<TRANSLATE>	Caché has read an input value for which it has no translation value. It therefore carries out the Default Action defined on the Translation tab of the Caché NLS utility.
<TRANSLOST>	A distributed transaction initiated by this job has been asynchronously rolled back by the server.
<UNDEFINED>	There has been a reference to an undefined variable.
<UNIMPLEMENTED>	There has been an attempt to use either an unimplemented function or an unimplemented argument of a legitimate command or function.
<UNKNOWN ERROR>	An unexpected error has occurred. Call your support center to resolve this serious error.
<VALUE OUT OF RANGE>	The value is outside the maximum or minimum permissible range.
<VOLUME IS NOT FORMATTED>	The volume does not have the required formatting.
<VOLUME SET ALREADY CREATED>	There has been an attempt to format a Caché database that is already formatted.
<WIDE CHAR>	Caché read a multi-byte character where a 1-byte character was expected.
<WRITE DEMON FAILED>	The write daemon is unable to continue. Call your support center to resolve this serious error.
<WRITE>	The record cannot be written.
<ZFILE>	There has been a file-related problem, such as that a file cannot be found or cannot be opened; consult the error log for more details.
<ZMISC>	An error has occurred of a different type than the others listed in this table. More detail follows the <ZMISC> code itself.
<ZOPEN>	There has been an error opening a device with TCP; the database server being shadowed is down.
<ZREAD>	There has been a failure to read from the database server that is being shadowed, indicating that the server is disconnected or not currently running.
<ZSTOP>	Shadowing has been stopped by system shutdown or the user.
<ZTRAP>	There has been an attempt to issue a ZTrap command with no argument.

1.12 DICOM error messages

1.12.1 Error message: \$ZE=<STKOV>BTXT+2^DIALOG:::6:6:

Problem: This error may occur while attempting to correct images, and while trying to input data into a field that does not exist in the FileMan Data Dictionary on the VistA side.

Solution: Have site install the latest KIDS application and restart the process images procedure.

Note: This error will only occur during the process of corrected images because the exclusive kill does not take place during this process.

1.12.2 Error message “Port in Use” or “Port Unavailable”

This message indicates that multiple processes are contending for the same socket. The resolution of this problem is usually fairly simple.

First try to close the C-Store window that displays the error message. Check that there is no other C-Store session active that is attempting to run the same program for the same instrument. If there is, that other session is probably successfully using the port that the former one could not obtain. If there is no other session attempting to run the same task, try to start it by double-clicking on its icon. This usually resolves the problem.

If the problem persists, wait until the Image Gateway is not actively processing images (all active windows show the “idle” indicator). Close the windows that display the “idle” indicator one by one. When no more telnet sessions are active, shut down and restart Caché. Once Caché is running again, restart all telnet sessions.

If, after this, the problem still persists, wait until the Image Gateway is not actively processing images (all active windows show the “idle” indicator). Close the windows that display the “idle” indicator one by one. When no more telnet sessions are active, shut down Caché. Once Caché is shut down, shut down Windows as well. Power down the PC, wait 20 seconds, and power it up again. Restart Windows. Log in using the appropriate username and password and restart Caché. Once Caché is running again, restart all telnet sessions.

If, after this, the problem still persists, log a Remedy call.

1.13 Signing on to the GUI

When a user signs on into the GUI version of the VistA Imaging software and the user has a Radiology classification as a resident and/or staff, the Kernel Broker login software will execute a Radiology routine that counts the number of “unverified reports”. In hospitals where it is not common practice that all staff members verify their radiology reports, the number of reports with a status of “waiting for action” can grow extremely large. Such a count may keep the user waiting for a successful login for several minutes. Routine ^RAUTL3 is the called routine that counts the number of radiology reports requiring action by the user who is currently logging on. Some sites have modified this routine to quit for GUI sign-on. Consult the Radiology developers

if this problem exists or have user verify all outstanding reports, which is the standard procedure.

1.14 Dumping Files that Cannot be Processed by the Gateway Application

Occasionally, input files cannot be processed properly. The distribution kit contains a public domain “octal dump” utility program (**OD.EXE**), that may be of help in diagnosing such problems. This utility program is intended to be started from the DOS-prompt, and will accept the following parameters:

```
> od -b -A {doxn} -N number -j number -t {acdfoux} filename
```

The meaning of the various parameters is described below.

-b

This parameter means that the content of the file is to be displayed in “bytes” rather than in “words” (default, when this parameter is omitted).

-A {d,o,x,n}

This parameter indicates how the counter in the left-most column is to be displayed. Only one of the four possible sub-options may be specified:

- A d means: display the counter in decimal representation
- A o means: display the counter in octal representation (default)
- A x means: display the counter in hexadecimal representation
- A n means: do not display the counter

-N number

This parameter indicates how many bytes need to be displayed. Only the first “so many” bytes of the file will be displayed, where “so many” is the number indicated by this parameter.

Note: The value of this parameter is written in decimal representation.

-j number

This parameter indicates where in the input file, the program is to start displaying characters. The program will start at the character indicated by the value of “number” in this parameter. Note: the value of this parameter is written in decimal representation.

-t {a,c,d,f,o,u,x}

When this parameter is omitted, the utility program will only display the numeric values of the various characters in the file. When additional information is desired, this parameter may be used to indicate which additional information should be displayed:

- t a means: display all ASCII characters
- t c means: display only printable ASCII characters
- t d means: display decimal representation

- t f means: display floating point representation
- t o means: display octal representation
- t u means: display unsigned decimal representation
- t x means: display hexadecimal representation

Example:

```
C:\DICOM\DICT>od -N 80 -j 50 -b -A d -t c instrument.dic <Enter>
0000000050      162 151 160 164 151 157 156 174 111 156 163 164 151 164 165 164
                r   i   p   t   i   o   n   |   I   n   s   t   i   t   u   t
0000000066      151 157 156 040 116 141 155 145 174 111 155 141 147 151 156 147
                i   o   n   N   a   m   e   |   I   m   a   g   i   n   g
0000000082      040 123 145 162 166 151 143 145 174 120 157 162 164 174 015 012
                S   e   r   v   i   c   e   |   P   o   r   t   |   \r   \n
0000000098      043 040 111 155 141 147 151 156 147 040 163 145 162 166 151 143
                #       I   m   a   g   i   n   g       s   e   r   v   i   c
0000000114      145 163 040 141 162 145 040 144 145 146 151 156 145 144 040 141
                e   s       a   r   e       d   e   f   i   n   e   d       a
0000000130
```


C:\DICOM\DICT>

1.15 Microsoft Windows Error Code

In some cases, the system reports an error number. Such error numbers are typically returned by platform-specific system calls. Usually, the software will display both the error code and the meaning. If a code, but no meaning is displayed, the program “**ErrLook**” may be of help to find out the meaning of the code.

The program “ErrLook.exe” is stored in the following directory so that is always in an accessible “execution path”: \Program Files\Vista\Imaging\DICOM.

When this program is started, it asks for an error code number. When a number is entered, click on the button labeled “**Look Up**” to display the meaning of the code.

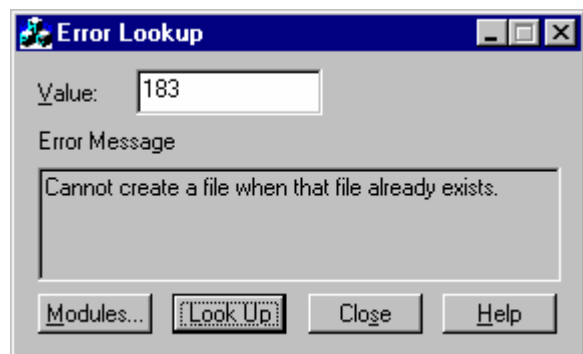

1.16 Run Time Errors Reported by the DICOM Gateway

For ease of use, all of these errors are alphabetically listed in the Index under the heading “Run Time Error”.

1.16.1 DICOM Association Errors

These errors are non-recoverable protocol violations on the part of the peer DICOM application entity. They cause both the association and the communications task to be terminated.

1.16.1.1 Unexpected “DATA” PDU

```
*****
*** DICOM GATEWAY ERROR ***
*** ***
*** PDU is out of sequence. ***
*** PDU Type is 04, data. ***
*** ***
*** Routine: ^MAGDACR1 Please Call Support Personnel ***
*****
```

This error occurs when an unexpected DICOM message is encountered during the lower-level communication “handshaking” with a peer application entity. (A PDU is a **P**rotocol **D**ata **U**nit.)

1.16.1.2 Peer application entity Requests “Abort Association”

```
*****
*** DICOM GATEWAY ERROR ***
*** ***
*** PDU Type is 07, abort association. ***
*** ***
*** Routine: ^MAGDACR1 Please Call Support Personnel ***
*****
```

The peer application entity is requesting to terminate the communication. VistA terminates the communications task.

1.16.1.3 DICOM Association – Invalid Protocol Identifier

```
*****
*** DICOM GATEWAY ERROR ***
*** ***
*** Expected protocol version 1. ***
*** Found "xxxxxxxxxx". ***
*** ***
*** Routine: ^MAGDACR1 Please Call Support Personnel ***
*****
```

The DICOM protocol version is exchanged during the initial negotiation of the association. The only legal DICOM protocol is “Version 1”.

1.16.1.4 Opcode for Unsupported Operation

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
***  No subroutine for opcode: "xxxxxxxxxx".           ***
***                                                    ***
***  Routine: ^MAGDACR1/3                               Please Call Support Personnel ***
*****

```

This error occurs when an unexpected association sub-item field is presented by a calling application entity negotiating an association with VistA. This is possibly a violation of the DICOM standard. It is a problem that a user cannot resolve and that must be reported to the VistA Imaging Project Team. The error message will include the offending opcode code.

1.16.1.5 DICOM Association – Invalid DICOM Standard UID

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
***  Expected UID for "DICOM Application Context Name", ***
***  which is "1.2.840.10008.3.1.1.1".                 ***
***  Found "xxxxxxxxxx".                               ***
***                                                    ***
***  Routine: ^MAGDACR3                               Please Call Support Personnel ***
*****

```

The DICOM Standard UID is exchanged during the initial negotiation of the association. The only legal value for DICOM is “1.2.840.10008.3.1.1.1”.

1.16.1.6 Invalid Context Identifier

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
***  Context ID should be an odd number.               ***
***  The number encountered is "xxxxxxxxxx".           ***
***                                                    ***
***  Routine: ^MAGDACR3                               Please Call Support Personnel ***
*****

```

DICOM Presentation Context IDs are odd numbers in the range 1:255. If an even number is presented, it is an error.

The error message will include the number specified.

1.16.1.7 Invalid Sub-Type

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Expected sub-type 51(hex) is missing. ***
*** Found "xxxxxxxxxx" instead. ***
***
*** Routine: ^MAGDACR3 Please Call Support Personnel ***
*****

```

The VistA Imaging DICOM Gateway found a “User Information Item Field (Item-type 50H)” in an association request. The user data field of must contain Item-type 51H, and did not. The error message will include the incorrect hexadecimal value specified.

1.16.1.8 Illegal Unique Identifier (UID)

```

*****
*** ERROR IN DICOM UNIQUE IDENTIFIER DICTIONARY ***
***
*** The UID Dictionary is not properly setup! ***
*** There is no entry in the UID dictionary for "xxxxxxxxxx" ***
***
***
*** Routine: ^MAGDUID1 Please Call Support Personnel ***
*****

```

The UID dictionary must contain the proper entries in order for the gateway applications to operate properly. This error is probably caused by an old version of the dictionary being loaded. Load the current version.

1.16.2 Storage Provider Error – Unexpected Response

```

*****
*** ERROR IN STORAGE SERVICE ***
***
*** Expected return value "STORAGE COMPLETE", ***
*** Received "xxxxxxxxxx". ***
***
*** Routine: ^MAGDCST1 Please Call Support Personnel ***
*****

```

The C-Store task running in the foreground sent an unexpected message to the MUMPS Storage Controller. It should have sent “STORAGE COMPLETE”. The erroneous text is displayed in the error message.

This error is fatal and causes the DICOM Storage Provider to terminate.

1.16.3 Errors Encountered in Reading/Writing DICOM Messages

1.16.3.1 Unsupported Value Representation (Input)

```
*****
*** DICOM GATEWAY ERROR ***
***
*** Undefined value representation: "xxxxxxxxxx". ***
***
*** Routine: ^MAGDDR3 Please Call Support Personnel ***
*****
```

This error may occur during the reading of a DICOM message. This error occurs when a data element is encountered to have a “value representation” that is not recognized by the VistA Imaging DICOM Gateway software (All legal DICOM value representations are handled by VistA Imaging DICOM Gateway software).

The error message will contain the code for the offending value representation. The erroneous message will have to be skipped.

1.16.3.2 Value Representation Mismatch

```
*****
*** DICOM TEXT GATEWAY ERROR -- WRITING DICOM FILE ***
***
*** VR Mismatch ***
*** Requested Value Representation="xxxxxxxxxx", Dictionary="xxxxxxxxxx". ***
*** Group: xxxxxxxx, Group Owner: xxxxxxxx Element: xxxxxxxx (xxxxxxxxxx) ***
***
*** Routine: ^MAGDDW3 Please Call Support Personnel ***
*****
```

This error may occur when trying to re-generate an Explicit VR image header where the VR specified for an element by a vendor is not in agreement with the DICOM Standard.

The error message will contain the offending value representation, as well as the list of value representations that is specified in the DICOM Data Dictionary. The error message will include identify the offending element.

1.16.3.3 Required Type 1 Data Item Missing

```
*****
*** DICOM TEXT GATEWAY ERROR -- WRITING DICOM FILE ***
***
*** REQUIRED TYPE 1 DATA MISSING ***
*** Group: xxxxxxxx, Group Owner: xxxxxxxx Element: xxxxxxxx (xxxxxxxxxx) ***
***
*** Routine: ^MAGDDW3 Please Call Support Personnel ***
*****
```

This error occurs when a required element’s value is not provided, while formatting a DICOM message. The error message will include identify the offending element.

This is usually a caused by a programming error that an user cannot resolve and that must be reported to the VistA Imaging Project Team.

1.16.3.4 Unsupported Value Representation (Output)

```
*****
*** DICOM TEXT GATEWAY ERROR -- WRITING DICOM FILE ***
***                                     ***
*** ERROR -- Undefined value representation: xxxxxxxxxx ***
***                                     ***
*** Routine: ^MAGDDW4                      Please Call Support Personnel ***
*****
```

This error may occur during the writing of a DICOM message. This error occurs when a data element is encountered to have a “value representation” that is not recognized by the VistA Imaging DICOM Gateway software. (All legal DICOM value representations are handled by VistA Imaging DICOM Gateway software.)

This is usually a caused by a dictionary file or programming error that a user cannot resolve and that must be reported to the VistA Imaging Project Team.

The error message will contain the offending value representation.

1.16.4 Modality Worklist – Invalid SOP Class (C-Find)

```
*****
*** ERROR IN C-FIND ***
***                                     ***
*** Expected SOP Class "Modality Worklist Information Model - FIND ***
*** Found "xxxxxxxxxx". ***
***                                     ***
*** Routine: ^MAGDFND1                      Please Call Support Personnel ***
*****
```

This error may occur while processing a C-Find request. As part of the message, the client specifies the name of the SOP class, which defines the context for the entire transaction. This error occurs when the name of this SOP class is not recognized. (Currently, only the SOP class “Modality Worklist Information Model – FIND” is recognized.) This is an error on the part of the application entity that issued the C-Find request.

The error message will include the name of the class that is specified by the client. This is a fatal error that terminates the association.

1.16.5 HL7 Error Messages

HL7 messages are generated on the main VistA hospital system in the ^MAGDHL7 global. They are read by the Text Gateway and used to update the Modality Worklist database ^MAGDWLST global and to create the DICOM messages that are sent to commercial PACS.

All of these errors are detected while reading and processing them. Bad HL7 message can be bypassed by incrementing the HL7 pointer.

1.16.5.1 Bad HL7 Message Header

```

*****
*** DICOM TEXT GATEWAY -- HL7 DATA ERROR ***
***                                     ***
*** Error: BAD HL7 MESSAGE HEADER ***
***                                     ***
*** HL7 message header should start with "MSH", ***
*** starts with "xxx". ***
***                                     ***
*** Routine: ^MAGDHR1/1A Please Call Support Personnel ***
*****

```

```

*****
*** BAD HL7 MESSAGE HEADER ***
***                                     ***
*** HL7 message header should start with "MSH", ***
*** starts with "xxx". ***
***                                     ***
*** Routine: ^MAGDHRP Please Call Support Personnel ***
*****

```

Each HL7 message should start with the prefix “MSH”. These errors will occur when the DICOM Text Gateway detects a message that does not start with this prefix. It indicates a badly garbled message that cannot be further processed.

1.16.5.2 Incomplete HL7 Message

```

*****
*** DICOM TEXT GATEWAY -- HL7 DATA ERROR ***
***                                     ***
*** Incomplete copy of HL7 message encountered. ***
*** In message number "nnn", ***
*** segment number "mmm" is missing. ***
***                                     ***
*** Routine: ^MAGDHR2 Please Call Support Personnel ***
*****

```

A portion of the HL7 message is missing, due to a problem on the main hospital system. The incomplete message will have to be bypassed manually.

1.16.5.3 Bad Fillers Order Number

```

*****
*** DICOM TEXT GATEWAY -- HL7 DATA ERROR ***
***                                     ***
*** HL7 consult processing detected a bad "Fillers Order #" ***
*** OBR-3 contained="xxxxxxxxxx". It should have been "GMRC" ***
***                                     ***
*** Routine: ^MAGDHR4A Please Call Support Personnel ***
*****

```

This might occur if the message for the consult/procedure request is incorrectly formatted. Please call the VistA Imaging Support Desk.

1.16.6 DICOM Image Processing Errors During Input

These errors in processing a DICOM image may indicate database corruption and may require manual intervention. Maybe the image processing pointer can be incremented to bypass offending images until the problem can be resolved.

When this type of error occurs, it usually indicates a serious database inconsistency. Report it immediately to the VistA Imaging Project Team.

1.16.6.1 Undefined Imaging Service

```
*****
***  DICOM IMAGE PROCESSING ERROR                               ***
***                                                                 ***
***  Undefined Imaging Service: xxxxxxxxxx                      ***
***                                                                 ***
***  Routine: MAGDIR81                                           Please Call Support Personnel ***
*****
```

This error is due to an erroneous Imaging Service field in the INSTRUMENT.DIC file. Legal values are "RAD" and "CON".

1.16.6.2 File 2005 Corruption Problem

```
*****
***  Wed 12:09 DICOM IMAGE PROCESSING ERROR - FILE 2005 CORRUPTION ***
***                                                                 ***
***  The ^MAG(2005) file has been corrupted so that new images will ***
***  overwrite old ones and general image database inconsistency ***
***  will result.                                                ***
***                                                                 ***
***  Latest internal entry number processed: 102158 Wed SEP 20, at 11:40 ***
***  Bad ^MAG(2005,0) internal entry number: 102146            ***
***  Image Gateway: "VistA DICOM Gateway Development System" ***
***                                                                 ***
***          This is a VERY SERIOUS ERROR.  Image processing ***
***          will be halted until it is resolved.              ***
***                                                                 ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                 ***
***  Routine: MAGDIR84                                           ***
*****
```

The image processing task has determined that the most recent internal entry number contained in the ^MAG(2005,0) node has been decremented. This could have been caused by restoring previously saved copy of the ^MAG global without applying the journal files. This is a VERY SERIOUS ERROR since it would allow new images to overwrite old ones, resulting in general corruption of the VistA Imaging database. Please contact IRM and the National Help Desk immediately.

1.16.6.3 Radiology Report File Corruption Problem

```

*****
***  Wed 07:02 DICOM IMAGE PROCESSING ERROR - ^RARPT FILE CORRUPTION  ***
***                                                                    ***
***  The RAD REPORT file has been corrupted so that new reports will  ***
***  overwrite old ones and general image/report database inconsistency ***
***  will result.                                                    ***
***                                                                    ***
***  Latest internal entry number processed: 77710 SEP 11 at 09:51:39 ***
***  Bad ^RARPT(0) internal entry number: 77643                     ***
***  Image Gateway: "VistA DICOM Gateway Development System"        ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing         ***
***          will be halted until it is resolved.                   ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                    ***
***  Routine: ^MAGDIR84                                              ***
*****

```

The image processing task has determined that the most recent internal entry number contained in the ^RARPT(0) node has been decremented. This could have been caused by restoring a previously saved copy of the ^RARPT global without applying the journal files. This is a VERY SERIOUS ERROR since it would allow new images to overwrite old ones, resulting in general corruption of the VistA Imaging database. Please contact IRM and the National Help Desk immediately.

1.16.6.4 Group Entry Number Problem

```

*****
***  Wed 13:07 DICOM IMAGE PROCESSING ERROR - GROUP ENTRY NUMBER PROBLEM ***
***                                                                    ***
***  The internal entry number for this group is less than that of the ***
***  last processed image.  This will cause new images to overwrite ***
***  old ones and general image database inconsistency will result. ***
***                                                                    ***
***  Latest internal entry number processed: 102158 Wed SEP 20, at 11:40 ***
***  Bad internal entry number of new group: 102100                 ***
***  Image Gateway: "VistA DICOM Gateway Development System"        ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing         ***
***          will be halted until it is resolved.                   ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                    ***
***  Routine: MAGDIR9A and MAGDIR9E                                  ***
*****

```

The image-processing task has determined that the internal entry number of the new group is less than the most recently created entry. This could have been caused by restoring a previously saved copy of the ^MAG global without applying the journal files. This is a VERY SERIOUS ERROR since it would allow new images to overwrite old ones, resulting in general corruption of the VistA Imaging database. Please contact IRM and the National Help Desk immediately.

1.16.6.5 Image Entry Number Problem

```

*****
***  Wed 12:38 DICOM IMAGE PROCESSING ERROR - IMAGE ENTRY NUMBER PROBLEM  ***
***
***  The internal entry number for this image is less than that of the      ***
***  last processed image.  This will cause new images to overwrite         ***
***  old ones and general image database inconsistency will result.         ***
***
***  Latest internal entry number processed: 102158 Wed SEP 20, at 11:40    ***
***  Bad internal entry number of new image: 102100                        ***
***  Image Gateway: "VistA DICOM Gateway Development System"               ***
***
***          This is a VERY SERIOUS ERROR.  Image processing                ***
***          will be halted until it is resolved.                          ***
***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP      ***
***
***  Routine: MAGDIR9B and MAGDIR71                                         ***
*****

```

The image-processing task has determined that the internal entry number of the new image is less than the most recently created entry. This could have been caused by restoring a previously saved copy of the ^MAG global without applying the journal files. This is a VERY SERIOUS ERROR since it would allow new images to overwrite old ones, resulting in general corruption of the VistA Imaging database. Please contact IRM and the National Help Desk immediately.

1.16.6.6 Imaging Patient Mismatch Problem

```

*****
***  Wed 12:53 DICOM IMAGE PROCESSING ERROR - PATIENT MISMATCH PROBLEM  ***
***
***  The image and the group point to different patients.                  ***
***
***  The Image points to PATIENT file internal entry number 100            ***
***  IMAGPATIENT1,ONE | 000-05-6789 | M | Dec 25,1925                     ***
***
***  The Group points to PATIENT file internal entry number 90271          ***
***  IMAGPATIENT2,TWO | 000-05-4321 | F | Oct 31,1919                     ***
***
***  Internal entry number of group: 102138                                ***
***  Image Gateway: "VistA DICOM Gateway Development System"               ***
***
***          This is a VERY SERIOUS ERROR.  Image processing                ***
***          will be halted until it is resolved.                          ***
***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP      ***
***
***  Routine: MAGDIR9A, MAGDIR9B, and MAGDIR9E                             ***
*****

```

The image-processing task has determined that the image and the group entries would point to different patients. An unknown software problem or manual modification of the data could have caused this to occur. This is a VERY SERIOUS ERROR since it might lead to a corrupted VistA Imaging database. Please contact IRM and the National Help Desk immediately.

1.16.6.7 Radiology Patient Mismatch Problem

```

*****
*** Mon 06:50 DICOM IMAGE PROCESSING ERROR - RAD PATIENT/REPORT MISMATCH ***
***
*** The image and the radiology report point to different patients. ***
***
*** The Image points to PATIENT file internal entry number 90263 ***
*** IMAGPATIENT,ONE M. | 000-01-9676 | M | Dec 25,1925 ***
***
*** The Rad Report points to PATIENT file internal entry number #17 ***
*** IMAGPATIENT,TWO N. | 000-02-7748 | F | Oct 31,1919 ***
***
*** Internal entry number of report: ^RARPT(12345) ***
*** Image Gateway: "VistA DICOM Gateway Development System" ***
***
*** This is a VERY SERIOUS ERROR. Image processing ***
*** will be halted until it is resolved. ***
***
*** Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***
*** Routine: MAGDIR9A ***
*****

```

The image-processing task has determined that the image and the radiology report entries would point to different patients. An unknown software problem or manual modification of the data could have caused this to occur. This is a **VERY SERIOUS ERROR** since it might lead to a corrupted VistA Imaging and/or Radiology database. Please contact IRM and the National Help Desk immediately.

1.16.6.8 No Radiology Case Number

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
*** No radiology case number specified for patient xxxxxxxxxx ***
*** Routine: ^MAGDIR9A Please Call Support Personnel ***
*****

```

A database inconsistency may exist in the Radiology package. There should be a radiology case number specified for this patient. Perhaps the case has been deleted. The DICOM Image Input pointer can be incremented to temporarily get around this problem.

1.16.6.9 Radiology Case Not in ^RADPT

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
*** Radiology case xxxxx is not in ^RADPT(yyyyyyy) ***
*** Routine: ^MAGDIR9A Please Call Support Personnel ***
*****

```

The entry in the Radiology Patient File (^RADPT) for this case number appears to be missing. This may be due to a database inconsistency in the Radiology package. Perhaps the case has been deleted. The DICOM Image Input pointer can be incremented to temporarily get around this problem.

1.16.6.10 Cannot Create Group for Old Radiology Images

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** IMAGE GROUP CREATION ERROR: ***
*** Radiology Report has been archived and purged. ***
*** Patient xxxxxxxxxx, Date xxxxxxxxxx, Case xxxxxxxxxx ***
***
*** Routine: ^MAGDIR9A Please Call Support Personnel ***
*****
```

This error may occur while attempting to enter information into the VistA database to establish a group of images. This error happens when the image files are for a Radiology study whose reports already have been archived and purged, as may happen when films are scanned in for 'old' studies.

1.16.6.11 Cannot Create Group for New Radiology Images

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** IMAGE GROUP CREATION ERROR: ***
*** xxxxxxxxxx ***
***
*** Routine: ^MAGDIR9A & ^MAGDIR9E Please Call Support Personnel ***
*****
```

This error may occur while attempting to create a radiology image group in the VistA database.

1.16.6.12 Cannot Find Image Group Pointer in Radiology Report

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** IMAGE GROUP LOOKUP ERROR: ***
*** Looking for 2005 cross reference in ^RARPT(xxxxxxxx) ***
***
*** Routine: ^MAGDIR9A Please Call Support Personnel ***
*****
```

This error may occur while attempting to access a radiology image group in the VistA database.

1.16.6.13 Cannot Create Image File Entry in VistA Database

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** IMAGE FILE CREATION ERROR: ***
*** xxxxxxxxxx ***
***
*** Routine: ^MAGDIR9B Please Call Support Personnel ***
*****
```

This error may occur while attempting to create the image entry in the VistA database.

1.16.6.14 Cannot Create Subdirectory to Store Image File

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** IMAGE FILE SUBDIRECTORY CREATION ERROR ***
*** Can not create the image subdirectory "xxxxxxxxxx" ***
*** xxxxxxxxxxxx ***
***
*** Routine: ^MAGDIR71 Please Call Support Personnel ***
*****

```

Images are stored in a hierarchy of subdirectories. This error may occur when the creation of a new subdirectory. The error message will contain the name of the directory that failed to be created. Check that the network file server is still “reachable”. Try to map a network drive to it and manually create the subdirectory using the command prompt.

1.16.6.15 Cannot Write Image File

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
***
*** The writing of image file "xxxxxxxxxx" failed ***
*** The error message was "xxxxxxxxxx" ***
***
*** Routine: ^MAGDIR75 Please Call Support Personnel ***
*****

```

When an image is stored on the network file server, the software checks whether the size of the created file is greater than zero. This error message is reported if the newly created file is either not created or has a zero size. Check that the network file server is still “reachable”. Try to map a network drive to it and manually copy a file using the command prompt.

Try the following steps to correct this problem:

1. From a DOS session, copy a small file to the share using the UNC path and file name for the destination. If the write fails it is usually caused by permission problems. Fix the permission problem using the Windows™ user administration tools (if necessary, find a System Administrator who has the appropriate privileges).
2. Check the available disk space on the share. If the disk space is low, run a purge on the background processor to free up space on the magnetic shares.

1.16.6.16 TIU Note Patient/Report Mismatch

```

*****
***  Mon 06:50 DICOM IMAGE PROCESSING ERROR - TIU PATIENT/REPORT MISMATCH  ***
***                                                                    ***
***  The image and the TIU note point to different patients.              ***
***                                                                    ***
***  The Image points to PATIENT file internal entry number 90263         ***
***  IMAGPATIENT,ONE M. | 000-01-9676 | M | Dec 25,1925                  ***
***                                                                    ***
***  The TIU note points to PATIENT file internal entry number #17        ***
***  IMAGPATIENT,TWO N. | 000-02-7748 | F | Oct 31,1919                  ***
***                                                                    ***
***  Internal entry number of TIU note: ^TIU(8925,1)                      ***
***  Image Gateway: "VistA DICOM Gateway Development System"             ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing              ***
***          will be halted until it is resolved.                        ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP    ***
***                                                                    ***
***  Routine: MAGDIR9E                                                    ***
*****

```

The image-processing task has determined that the image and the TIU note entries would point to different patients. An unknown software problem or manual modification of the data could have caused this to occur. This is a **VERY SERIOUS ERROR** since it might lead to a corrupted VistA Imaging and/or TIU Note database. Please contact IRM and the National Help Desk immediately.

1.16.6.17 TIU External Data Link Association Error

```

*****
***  ERROR ASSOCIATING WITH TIU EXTERNAL DATA LINK (file 8925.91)       ***
***  <error message from TIU>                                           ***
***  Image Gateway: "VistA DICOM Gateway Development System"            ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing              ***
***          will be halted until it is resolved.                        ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP    ***
***                                                                    ***
***  Routine: MAGDIR9E and MAGDHWa                                       ***
*****

```

The image-processing task has encountered an error in attempting to create an association between the image and the TIU External Data Link. An unknown software problem or manual modification of the data could have caused this to occur. This is a **VERY SERIOUS ERROR** since it might lead to a corrupted VistA Imaging and/or TIU Note database. Please contact IRM and the National Help Desk immediately.

1.16.6.18 Cannot Find Table with Additional Data Items

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
***  Data transfer dictionary <xxxxxxxxxx> is missing    ***
***                                                    ***
***  Routine: ^MAGDIR7T                                Please Call Support Personnel ***
*****

```

The error message will include the name of the missing table.

This error may occur when the designated data transfer dictionary can not be found. This is probably either a configuration or an installation problem.

The data transfer dictionaries are **DataGECT.DIC**, **Data_CR.DIC**, and **DataMisc.DIC**, all of which should be in the **f:\DICOM\Dict** subdirectory. The **Modality.DIC** file should contain the names of the data transfer dictionaries.

The missing data transfer dictionary will be created when the information from master file **F:\DICOM\Dict\Modality.DIC** is reloaded into the database of the VistA Imaging DICOM Gateway (use menu option 4, 2, 3, **Install MODALITY.DIC**).

If the file(s) in directory **F:\DICOM\Dict** are also missing, these can be re-copied from the distribution medium (either CD-ROM or network).

1.16.6.19 DICOM Element has Too Many Values

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
***  Processing "xxxxxxxxxx", tag is "xxxxxxxxxx".        ***
***  Multiplicity is "xxxxxxxxxx".                        ***
***                                                    ***
***  Routine: ^MAGDIR7T                                Please Call Support Personnel ***
*****

```

This error may occur while building the set of data elements that is to be displayed on a diagnostic workstation. Within the DICOM protocol, each data item may have 0 (zero) or more values. How many values a specific data item may have is defined in the DICOM Data Dictionary. When a data item is encountered that has more actual values than the maximum defined in the DICOM Data Dictionary, this error will occur.

The error message will identify the data stream being processed, the “tag” for the data item in question, and the actual number of values.

1.16.6.20 Cannot Create Targa™ File

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** TARGA IMAGE FILE CREATION ERROR ***
*** Unexpected error in <host system command> ***
*** ... host system error message(s) ... ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR7C Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation of the Targa™ file for an image. The actual display will include the error message that was produced by the host operating system.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started.

1.16.6.21 Cannot Create Abstract File

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** DICOM IMAGE ABSTRACT CREATION ERROR ***
*** Unexpected error in <host system command> ***
*** ... host system error message(s) ... ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR7C Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation of the “abstract” file (also known as thumbnail or icon) for an image. The actual display will include the error message that was produced by the host operating system.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started.

1.16.6.22 No Image File Created

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (1) ***
*** Unexpected error in <host system command> ***
*** Temporary TGA file <name> not in the directory ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation of an image file (Targa™ file, abstract file, etc.). The actual display will include the error message that was produced by the host operating system.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started

This particular error is most likely the result of a “disk full” or “privilege not granted” condition. Check that the disk share is still accessible through the network, check that the disk has enough free space to store new files, and check that the current user has appropriate credentials to execute the task at hand.

1.16.6.23 Error while Creating Image File

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (2) ***
*** Unexpected error in <host system command> ***
*** Error <message> for temporary <type> file <name> ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation an image file (Targa[™] file, abstract file, etc.). The actual display will include the error message that was produced by the external program.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started

This particular error is most likely the result of a “disk full” or “privilege not granted” condition. Check that the disk share is still accessible through the network, check that the disk has enough free space to store new files, and check that the current user has appropriate credentials to execute the task at hand.

1.16.6.24 Could not Copy Image File

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** <type> IMAGE COPY ERROR ***
*** Unexpected error in <host system command> ***
*** ... host system error message(s) ... ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****
```

While executing a host operating system command to copy an image file, an error occurred that prevented the proper creation of the copied file. The actual display will include the error message that was produced by the host operating system.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started.

1.16.6.25 Could not Create Image File

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
*** The writing of the image file failed. ***
*** Path: "<name>" ***
*** The error message was "<message>" ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****

```

While executing a host operating system command, an error occurred that prevented the proper creation an image file (TargaTM file, abstract file, etc.). The actual display will include the error message that was produced by the external program.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started

This particular error is most likely the result of a “disk full” or “privilege not granted” condition. Check that the disk share is still accessible through the network, check that the disk has enough free space to store new files, and check that the current user has appropriate credentials to execute the task at hand.

1.16.6.26 Image File Created, but File Size not Greater than 0

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
*** The writing of the image file failed. ***
*** Path: "<name>" ***
*** Size of written file is <number> bytes. ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****

```

While executing a host operating system command, an error occurred that prevented the proper creation an image file (TargaTM file, abstract file, etc.). The actual display will include the error message that was produced by the external program.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started

This particular error is most likely the result of a “disk full” or “privilege not granted” condition. Check that the disk share is still accessible through the network, check that the disk has enough free space to store new files, and check that the current user has appropriate credentials to execute the task at hand.

1.16.6.27 Image File Created, but Size Differs from Expected File Size

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (3) ***
*** Unexpected error in <host system command> ***
*** Temporary <type> file <name> has ZERO LENGTH ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation an image file (Targa™ file, abstract file, etc.). The actual display will include the error message that was produced by the external program.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started.

1.16.6.28 Image File Created, but File Size Differs from Expected File Size

```
*****
*** DICOM IMAGE PROCESSING ERROR ***
*** The writing of the image file failed. ***
*** Path: "<name>" ***
*** Size of written file is <wrong number> -- should be <correct number> ***
*** Imaging entry <number> successfully removed ***
*** MAGDIR75 Please Call Support Personnel ***
*****
```

While executing a host operating system command, an error occurred that prevented the proper creation an image file (Targa™ file, abstract file, etc.). The display will include the error message that was produced by the external program. When the file size on the source (the Imaging Gateway) and destination (Imaging Server) are different, the Processing on the Gateway will halt, with error messages provided through an e-mail message sent to a designated mail group and include a recipient in the Customer Support group.

Information about the image was successfully removed from the VistA database, and a new attempt to process the image can be started

This particular error is most likely the result of a problem in the site's local area network. When the size of the newly created file is reported as being approximately equal to 57kB, it is very likely that a "gigabit Ethernet card" has failed.

1.16.6.29 Could not Delete Image File

```

*****
*** DICOM IMAGE PROCESSING ERROR ***
*** ERROR DELETING <type> ***
*** Unexpected error in <host system command> ***
*** ... host system error message(s) ... ***
*** MAGDIR71 Please Call Support Personnel ***
*****

```

While executing a host operating system command to delete an image file, an error occurred that prevented the proper deletion of that file. The actual display will include the error message that was produced by the host operating system.

1.16.6.30 Cannot find Table with Additional Data Items

```

*****
*** DICOM GATEWAY ERROR ***
*** ***
*** Data transfer dictionary <<xxxxxxxxxx>> is missing ***
*** ***
*** MAGDIR7T Please Call Support Personnel ***
*****

```

The error message will include the name of the missing table.

This error may occur when the designated data transfer dictionary can not be found. This is probably either a configuration or an installation problem.

The data transfer dictionaries are **DataGECT.DIC**, **Data_CR.DIC**, and **DataMisc.DIC**, all of which should be in the **F:\DICOM\Dict** subdirectory. The **Modality.DIC** file should contain the names of the data transfer dictionaries.

The missing data transfer dictionary will be created when the information from master file **F:\DICOM\Dict\Modality.DIC** is reloaded into the database of the VistA Imaging DICOM Gateway (use menu option 4, 2, 3, **Install MODALITY.DIC**).

If the file(s) in directory **F:\DICOM\Dict** are also missing, these can be re-copied from the distribution CD-ROM.

1.16.6.31 Image Patient Mismatch (parent file 2006.5839)

```

*****
***  DICOM IMAGE PROCESSING ERROR - IMAGE GROUP MISMATCH                ***
***  The image group does not point to PARENT FILE 2006.5839.          ***
***                                                                    ***
***  The image group should point to PARENT FILE 2006.5839.            ***
***  Instead it point to PARENT FILE #nnnnnnn                          ***
***  The image group is ^MAG(2005,nnnnnn)                               ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing            ***
***          will be halted until it is resolved.                      ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                    ***
***  Routine: MAGDIR9E                                                  ***
*****

```

1.16.6.32 Image and TIU Note Associated with different Notes

```

*****
***  DICOM IMAGE PROCESSING ERROR - TIU/IMAGE GROUP MISMATCH          ***
***  The image group and TIU point to different notes.                  ***
***                                                                    ***
***  TIU points to TUI note ien #nnnnnnn                                ***
***  The image points to TIU note ien #nnnnnnn                          ***
***  TIU External Data File (8925.91) ien #nnnnnnn                      ***
***  points to image group ien #nnnnnnn                                ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing            ***
***          will be halted until it is resolved.                      ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                    ***
***  Routine: MAGDIR9E                                                  ***
*****

```

1.16.6.33 Wrong Group Object Type

```

*****
***  DICOM IMAGE PROCESSING ERROR - WRONG GROUP OBJECT TYPE            ***
***  The group entry in ^MAG(2005) does not have the proper group      ***
***  object type.                                                        ***
***                                                                    ***
***  The expected value is 11.  The value in the group entry is nnnnnn. ***
***                                                                    ***
***  Internal entry number of incorrect group: nnnnnn                  ***
***                                                                    ***
***          This is a VERY SERIOUS ERROR.  Image processing            ***
***          will be halted until it is resolved.                      ***
***                                                                    ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                                                    ***
***  Routine: MAGDIR9B                                                  ***
*****

```


1.16.6.34 Ignore Image File

```

*****
***                                     ***
*** Image <filename> will be ignored.   ***
***                                     ***
*****

```

1.16.6.35 Unknown Instrument Mnemonic

```

*****
***                                     ***
*** Image nnnnn has an unknown instrument mnemonic "xxx" skipped ***
***                                     ***
*****

```

1.16.6.36 Rename Failed

```

*****
***                                     ***
*** Renaming oldname to newname failed. ***
***                                     ***
*****

```

1.16.6.37 Image File Creation Error

```

*****
***                                     ***
*** DICOM IMAGE PROCESSING ERROR        ***
***                                     ***
*** IMAGE FILE CREATION ERROR:          ***
*** <error message>                     ***
*** Gateway: "<name>"                    ***
***                                     ***
***           This is a VERY SERIOUS ERROR. Image processing ***
***           will be halted until it is resolved.           ***
***                                     ***
*** Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                     ***
*** Problem detected by routine MAGDIR71. ***
***                                     ***
*****

```

1.16.6.38 Image File Subdirectory Creation Error

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  IMAGE FILE SUBDIRECTORY CREATION ERROR ***
***  Can not create the image subdirectory "name" ***
***  <error message>                  ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
*****

```

1.16.6.39 New Files Might Overwrite Existing Ones (Single)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  IMAGE FILE CREATION ERROR - NEW FILES MIGHT OVERWRITE EXISTING ONES ***
***  Trying to create file:            ***
***  "filename"                        ***
***                                     ***
***  However, this file already exists in that directory: ***
***  "filename"                        ***
***  Can not overwrite this with the image files.          ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed.        ***
***                                     ***
***  Gateway: "<name>"                      ***
***                                     ***
***                                     ***
***          This is a VERY SERIOUS ERROR.  Image processing ***
***          will be halted until it is resolved.          ***
***                                     ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                     ***
***  Problem detected by routine MAGDIR71.                  ***
***                                     ***
*****

```

1.16.6.40 New Files Might Overwrite Existing Ones (Multiple)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  IMAGE FILE CREATION ERROR - NEW FILES MIGHT OVERWRITE EXISTING ONES ***
***  Trying to create file:            ***
***  "filename"                       ***
***                                     ***
***  However, these nnn files already exist in that directory: ***
***  "filename"                       ***
***  . . .                           ***
***  "filename"                       ***
***  Can not overwrite these with the image files. ***
***                                     ***
***  Gateway: "<name>"                 ***
***                                     ***
***          This is a VERY SERIOUS ERROR.  Image processing ***
***          will be halted until it is resolved. ***
***                                     ***
***  Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                     ***
***  Problem detected by routine MAGDIR71. ***
***                                     ***
*****

```

1.16.6.41 Cannot Delete File

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  ERROR DELETING xxx                ***
***  Unexpected error in <command>     ***
***  <error message>                  ***
***                                     ***
***  Routine: ^MAGDIR71                Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.42 Unknown Status (Patient Safety RPC)

```

*****
***                                     ***
*** UNKNOWN STATUS IN "PATIENT SAFETY" RPC RESULT ***
*** Status = "<status>" ***
*** Gateway: "<name>" ***
***                                     ***
***           This is a VERY SERIOUS ERROR.  Image processing ***
***           will be halted until it is resolved. ***
***                                     ***
*** Call IRM and the National VistA Support Help Desk (888) 596-HELP ***
***                                     ***
*** Problem detected by routine MAGDIR74. ***
***                                     ***
*****

```

1.16.6.43 Rollback Error (number of returned elements)

```

*****
***                                     ***
*** DICOM IMAGE PROCESSING ERROR ***
***                                     ***
*** ROLLBACK M-to-M RPC BROKER ERROR ***
*** Unexpected number of array elements returned by RPC: nnn ***
*** <result text> ***
*** . . . ***
*** <result text> ***
*** Only one element ("ROLLBACK") should be returned ***
***                                     ***
*** Routine: ^MAGDIR75 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.44 Rollback error (returncode)

```

*****
***                                     ***
*** DICOM IMAGE PROCESSING ERROR ***
***                                     ***
*** ROLLBACK M-to-M RPC BROKER ERROR ***
*** Unexpected return from RPC: "<returncode>" ***
*** Expected "ROLLBACK|STATUS|IMAGEPTR" ***
***                                     ***
*** Routine: ^MAGDIR75 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.45 Problem with Temporary File (1)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (1) ***
***  Unexpected error in "<command>"    ***
***  Temporary TGA file "<filename>" not in the directory ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR75                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.46 Problem with Temporary File (2)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (2) ***
***  Unexpected error in "<command>"    ***
***  Error "<text>" for temporary <type> file "<filename>" ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR75                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.47 Problem with Temporary File (3)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  DIRECTORY PROBLEM WITH TEMPORARY <type> IMAGE FILE (3) ***
***  Unexpected error in "<command>"    ***
***  Temporary <type> file "<filename>" has ZERO LENGTH ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR75                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.48 Image Copy Error

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  <type> IMAGE COPY ERROR           ***
***  Unexpected error in "<command>"    ***
***  <error message>                  ***
***                                     ***
***  Routine: ^MAGDIR75                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.49 Image Write Failed

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  The writing of image file failed.   ***
***  Path: "<filename>"                ***
***  <error message>                  ***
***                                     ***
***  Routine: ^MAGDIR75                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.50 Targa™ Abstract Creation Error

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  TARGA IMAGE ABSTRACT CREATION ERROR ***
***  Unexpected error in "<command>"    ***
***  <error message>                  ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR7C                 Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.51 Multiframe Image File Creation Error (number of returned elements)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR                                     ***
***                                     ***
***  MULTIFRAME IMAGE FILE CREATION ERROR                             ***
***  Unexpected number of array elements returned by RPC: nnn       ***
***  Element 1) <text>                                              ***
***  . . .                                                         ***
***  Element n) <text>                                              ***
***  Two elements (PROCESSED and STORE) should be returned         ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed.                 ***
***                                     ***
***  Routine: ^MAGDIR7C                                           Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.52 Multiframe Image File Creation Error (return code not “processed”)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR                                     ***
***                                     ***
***  MULTIFRAME IMAGE FILE CREATION ERROR                             ***
***  Unexpected return from RPC: "<text>"                           ***
***  Expected "PROCESSED|location|instname|imageptr|0"             ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed.                 ***
***                                     ***
***  Routine: ^MAGDIR7C                                           Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.53 Multiframe Image File Creation Error (return code not “store”)

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR                                     ***
***                                     ***
***  MULTIFRAME IMAGE FILE CREATION ERROR                             ***
***  Unexpected return from RPC: "<text>"                           ***
***  Expected "STORE|0|imageptr^topath^file"                       ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed.                 ***
***                                     ***
***  Routine: ^MAGDIR7C                                           Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.54 Multiframe Image File Creation Error (error in command)

```
*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  DICOM IMAGE ABSTRACT CREATION ERROR ***
***  Unexpected error in "<command>"    ***
***  1) <text>                          ***
***  . . .                             ***
***  n) <text>                          ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR7D                 Please Call Support Personnel ***
***                                     ***
*****
```

1.16.6.55 Targa™ File Creation Error

```
*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***                                     ***
***  TARGA IMAGE FILE CREATION ERROR    ***
***  Unexpected error in "<command>"    ***
***  1) <text>                          ***
***  . . .                             ***
***  n) <text>                          ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed. ***
***                                     ***
***  Routine: ^MAGDIR7G                 Please Call Support Personnel ***
***                                     ***
*****
```


1.16.6.56 Could not Open Image File for Write

```

*****
***                                     ***
***  DICOM IMAGE PROCESSING ERROR      ***
***  Could not open about image text file "<filename>" for Write.      ***
***  STATUS: <xxx>                     ***
***                                     ***
***  Imaging entry #nnnnn was successfully removed.                    ***
***                                     ***
***  Routine: ^MAGDIR7T               Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.57 Data Transfer Dictionary Missing

```

*****
***                                     ***
***  DICOM GATEWAY ERROR               ***
***                                     ***
***  Data transfer dictionary <<filename>> is missing                  ***
***                                     ***
***  Routine: ^MAGDIR7T               Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.58 Invalid Multiplicity

```

*****
***                                     ***
***  DICOM GATEWAY ERROR               ***
***                                     ***
***  Processing "<node>", tag is "<tag>".                                ***
***  Multiplicity is "nnn".                                                ***
***                                     ***
***                                     ***
***  Routine: ^MAGDIR7T               Please Call Support Personnel ***
***                                     ***
*****

```

1.16.6.59 Image Entry Number Problem

```

*****
***                                     ***
*** DICOM IMAGE PROCESSING ERROR - IMAGE ENTRY NUMBER PROBLEM                 ***
***                                     ***
*** The internal entry number for this image is less than that of the          ***
*** last processed image. This will cause new images to overwrite              ***
*** old ones and general image database inconsistency will result.            ***
***                                     ***
*** Latest internal entry number processed: <nnnnn>                           ***
*** Bad internal entry number of new image: <nnnnn>                           ***
***                                     ***
*** Gateway: "<name>"                                                         ***
***                                     ***
***           This is a VERY SERIOUS ERROR. Image processing                  ***
***           will be halted until it is resolved.                           ***
***                                     ***
*** Call IRM and the National VistA Support Help Desk (888) 596-HELP          ***
***                                     ***
*** Problem detected by routine <name>.                                         ***
***                                     ***
*****

```

1.16.6.60 Routing Not Actively Used

```

*****
*** DICOM IMAGE PROCESSING ERROR                                     ***
***                                     ***
*** More than a week has elapsed since dd-Mon-yyyy                        ***
*** when the last activity occurred that is related                       ***
*** to the processing of Routed Image files.                               ***
***                                     ***
*** The site parameter for "This is a Routing Site" is                     ***
*** currently turned ON.                                                    ***
*** If this site is no longer actively routing image files                 ***
*** this site parameter must be turned OFF.                                ***
*** This parameter needs to be turned OFF on each VistA                   ***
*** DICOM Gateway that processes incoming images.                         ***
***                                     ***
*** There are currently nnn entries                                         ***
*** waiting to be processed in the evaluation queue.                       ***
***                                     ***
*** There are currently nnn entries                                         ***
*** waiting to be processed in the transmission queue.                     ***
***                                     ***
*** If this site is still a routing site, then both the                    ***
*** Routing Rule Evaluator and the Routing Transmitter                     ***
*** must be restarted.                                                      ***
*** Problem detected by routine MAGDIR82. Error Code: -405                 ***
***                                     ***
*** MAGDIR72 Please Call Support Personnel ***
*****

```

This message will appear (once per day) when a DICOM Gateway is processing images, and the setting for “Routing Active” is turned on, while no files have been routed or no rule-evaluation has taken place for over a week.

When the setting for “Routing Active” is turned on, an entry will be created in the Routing Rule Evaluation Queue for each image that is processed. If no rule-evaluator is running, this queue will grow and eventually fill any available storage space. (When a rule-evaluator is running, queue entries will be deleted after they have been processed.)

The action to be taken will depend on the situation at the site: if it is intentional that routing rule evaluation and file transmission is no longer active, then the setting for “Routing Active” should be turned off.

When Routing is intended to be active, the rule-evaluator and file transmitter(s) should be restarted.

1.16.6.61 Cannot Create Abstract

For some image files, the DICOM Gateway may be unable to create an "abstract" file. When this happens, the message stream on the Image Gateway will contain text like:

```
...
MAG_MAKEABS c:\DICOM\IMAGE_IN\Y0000138.DCM c:\DICOM\IMAGE_IN\TMP_IMAGE.ABS JPEG 100
HDR c:\DICOM\IMAGE_IN\MAG_MAKEABS_msg.txt & TYPE c:\DICOM\IMAGE_IN\MAG_MAKEABS_msg.txt
  1) 1^Failure Processing Image - Abstract Not Created
...
```

Since these messages are likely to scroll off the monitor without being noticed, an e-mail message is also sent. The contents of such a message will look like:

```
DICOM IMAGE ABSTRACT CREATION ERROR

Error reported by MAG_MAKEABS in creating the image abstract:
Image number # nnnnnnn.
1^Failure Processing Image - Abstract Not Created
Routine: MAGDIR7D
DICOM Gateway "Non-DDP Test System (Y)".
```

1.16.7 DICOM Image Processing Errors During Output

Images are internally stored in Targa[™] file format with an “about image text file,” which contains an ASCII version of the DICOM Header. When DICOM images are to be sent to another application entity, the DICOM images are “reconstituted” from the Targa[™] image files and the text files, and are then sent.

Like the VistA DICOM Storage Provider architecture, two processes are used to send a DICOM image, the MUMPS Transmission Controller and the C executable MAG_VISTA_SEND_IMAGE.exe.

The MUMPS Transmission Controller has managerial responsibilities for the effort and spawns the MAG_VISTA_SEND_IMAGE.exe program to transmit the actual image files. The two processes communicate with each other via TCP/IP, with the MUMPS routine instructing the C program with what to do.

1.16.7.1 Bad Targa™ File

```

*****
***  DICOM GATEWAY ERROR                                     ***
***                                                                 ***
***  Bad Targa File -- xxxxxxxxxxxx                         ***
***  Neither RGB color, nor gray scale.                     ***
***  Image Type is "xxxxxxxxxx"                             ***
***                                                                 ***
***  Routine: ^MAGDIW2                                     Please Call Support Personnel ***
*****

```

This error may occur while re-constituting a DICOM file from a Targa™ file and a Text file. In the Targa™ file, a code occurs that indicates the internal “image type”. Currently, only image types “gray scale” and “RGB color” are supported. When a file with a different image type is encountered, this error will occur. ”.

1.16.7.2 Cannot Establish Connection with MAG_VISTA_SEND_IMAGE

```

*****
***  ERROR IN VISTA SEND IMAGE                               ***
***                                                                 ***
***  Expected return value "Connection established to MUMPS". ***
***  Received "xxxxxxxxxx".                                  ***
***                                                                 ***
***  Routine: ^MAGDIW5                                     Please Call Support Personnel ***
*****

```

When the **MAG_VISTA_SEND_IMAGE.exe** program is spawn by the MUMPS Transmission Controller (^MAGDIW5), it tries to establish a TCP/IP connection back to the MUMPS controller task. This error message denotes that this communications failed.

1.16.7.3 Connection to Destination Storage Provider Failed

```

*****
***  ERROR IN VISTA SEND IMAGE                               ***
***                                                                 ***
***  Expected return value "CONNECTION SUCCESSFUL".         ***
***  Received "xxxxxxxxxx".                                  ***
***                                                                 ***
***  Routine: ^MAGDIW5                                     Please Call Support Personnel ***
*****

```

The **MAG_VISTA_SEND_IMAGE.exe** program to attempts to create a TCP/IP connection to the destination Storage Provider application entity. This connection could not be created.

1.16.7.4 Could Not Establish Association with Destination SCP

```

*****
*** ERROR IN VISTA SEND IMAGE ***
***                               ***
*** Expected return value "ASSOCIATION ACKNOWLEDGE". ***
*** Received "xxxxxxxxxx". ***
***                               ***
*** Routine: ^MAGDIW6           Please Call Support Personnel ***
*****

```

The MUMPS Transmission Controller created an Association Request and had the MAG_VISTA_SEND_IMAGE.exe program send it to the destination Storage SCP. The destination did not send back an Association Response.

1.16.7.5 DICOM Network Error

```

*****
*** ERROR IN VISTA SEND IMAGE ***
***                               ***
*** Job terminated due to Network Error #n ***
***                               ***
*** Routine: ^MAGDIW5           Please Call Support Personnel ***
*****

```

This error occurs while sending a DICOM message when one of the following conditions arises:

- Connection/association dropped/aborted by peer application entity
- Unknown PDU (Protocol Data Unit) type received
- Garbled transmission

These are errors caused by the DICOM system receiving images from VistA.

1.16.7.6 Image Transmission Acknowledgement Failed

```

*****
*** ERROR IN VISTA SEND IMAGE ***
***                               ***
*** Unexpected Response ***
*** Expected "IMAGE SENT", received "xxxxxxxxxx". ***
***                               ***
*** Routine: ^MAGDIW5           Please Call Support Personnel ***
*****

```

The MUMPS Transmission Controller instructed the MAG_VISTA_SEND_IMAGE.exe program to send an image to the destination Storage SCP. An error on the Storage SCP prevented successful transmission of the image.

1.16.7.7 Image Transmission Failed

```

*****
*** ERROR IN VISTA SEND IMAGE: DESTINATION STORAGE FAILURE ***
*** Status=XXXXX          Error ID: YYYYYY ***
*** Comment: ZZZZZ ***
*** Offending Element: (gggg,eeee) <name of element> ***
*** ***
*** Routine: ^MAGDIW5 ***
*****

```

The destination Storage SCP was not able to properly handle an image transmitted to it by the MAG_VISTA_SEND_IMAGE.exe program. The complete reason for the error is given in the message.

1.16.7.8 Cannot Acknowledge Establishment of Association

```

*****
*** ERROR IN VISTA SEND IMAGE ***
*** ***
*** Expected return value "ASSOCIATION ACKNOWLEDGE". ***
*** Received "xxxxxxxxxxx". ***
*** ***
*** Routine: ^MAGDIW6 Please Call Support Personnel ***
*****

```

The MUMPS Transmission Controller had the MAG_VISTA_SEND_IMAGE.exe program send an Association Acknowledge Request to the destination Storage SCP. The error is that the Storage SCP did not respond with a proper Association Acknowledge Response.

1.16.7.9 Cannot Acknowledge Release of Association

```

*****
*** ERROR IN VISTA SEND IMAGE ***
*** ***
*** Unexpected response ***
*** Expected "ASSOCIATION ENDED", Received "xxxxxxxxxxx". ***
*** ***
*** Routine: ^MAGDIW6 Please Call Support Personnel ***
*****

```

The MUMPS Transmission Controller had the MAG_VISTA_SEND_IMAGE.exe program send an Association Release Request to the destination Storage SCP. The error is that the Storage SCP did not respond with a proper Association Release Response.

1.16.8 Error Messages Report by DICOM Message Queuing Software

Between the message processing and transmission steps, all DICOM messages are stored in intermediate files (*.DCM) on the local system. For the sake of flexibility, four different mechanism are provided:

1. First-In-First-Out (FIFO) queues for PACS messages for the Text Gateway
2. Direct mode for queries (e.g., Modality Worklist) for the Text Gateway
3. FIFO queue for image acquisition (**D:\dicom\image_in\...**)

4. FIFO queue for image transmission (**D:\dicom\image_out\...**)

The error messages in this section pertain to the intermediate file processing for the Text Gateway.

1.16.8.1 Error in Directory Lookup of FIFO Queue Subdirectory

```
*****
*** DICOM TEXT GATEWAY ERROR ***
***
*** Cannot locate directory xxxxxxxxxx ***
*** Host File System Error xxxxxxxxxx ***
***
*** Routine: ^MAGDQUE0 Please Call Support Personnel ***
*****
```

This error occurs when the Text Gateway is doing a directory lookup of a FIFO queue subdirectory (typically **D:\dicom\data1 \ANNNNN**) for writing a new file. This is a catastrophic error that prevents the application from even creating the new file's subdirectory, if it did not already exist.

The error message includes the name of the subdirectory that got the error and the error code from the directory lookup operation. The most likely problem is that the path is pointing to a non-existent **D:\dicom\data1** directory.

1.16.8.2 Cannot Create FIFO Queue Subdirectory

```
*****
*** DICOM TEXT GATEWAY ERROR ***
***
*** Cannot create directory xxxxxxxxxx ***
*** Host File System Error: xxxxxxxxxx ***
***
*** Routine: ^MAGDQUE0 Please Call Support Personnel ***
*****
```

This error occurs when the Text Gateway is creating a new FIFO queue subdirectory (typically **d:\DICOM\Data1\Annnnn**).

The error message includes the name of the subdirectory that got the error and the error code from the **MKDIR** operation. The most likely problems are either that the permissions prevent the subdirectory from being created or the disk drive is out of space.

1.16.8.3 Invalid Data in Queue Pointer File

```

*****
*** DICOM TEXT GATEWAY ERROR -- WRITING DICOM FILE ***
***
*** Bad File Pointer: ***
*** "xxxxxxxxxx" = "xxxxxxxxxx". ***
*** Value of pointer should be 7 digits. ***
***
*** Routine: ^MAGDQUE0 Please Call Support Personnel ***
*****

```

This error occurs when the Text Gateway finds Queue processor invalid data in one of the FIFO queue pointers files. (The “pointer files” are named **x_READ.PTR** and **x_WRITE.PTR**, where **x** is A, B, C, D, E, F, G, H, S, T, U, V, W, X, Y, or Z.) The values in these files must be 7-digit integer numbers with leading zeroes. When any other value is encountered in one of these files, this error will occur.

The error message will include the name of the offending file and the value that is encountered in that file.

Typically, this error occurs when a queue pointer file was modified manually, not through the software in the VistA Imaging DICOM Gateway. When the displayed value indicates that a “strange” character appears leading or trailing the queue entry number, use Notepad to correct the offending file and enter the appropriate value. When the value contains an embedded control character, it may be necessary to erase the complete line, and re-enter it.

For example, suppose the **X_READ.PTR** is null (zero bytes long), and should really contain 12345. Use Notepad and type the value **0012345<Enter>**.

1.16.8.4 Cannot Create a FIFO Queue Text File

```

*****
*** DICOM TEXT GATEWAY ERROR -- WRITING DICOM FILE ***
***
*** Renaming xxxxxxxxxxx to xxxxxxxxxxx failed ***
*** Please delete the file "xxxxxxxxxx" and restart. ***
***
*** Routine: ^MAGDQUE0 Please Call Support Personnel ***
*****

```

This error occurs when the DICOM Text Gateway is not able to create a file in a FIFO Queue. When writing a DICOM file, the file is first created with a temporary name (typically **ANNNNNNN.tmp**), and then is renamed to its permanent name (in this case, **ANNNNNNN.dcm**). The error occurs whenever the temporary file cannot be renamed to its permanent name. This may occur when there is already a file present with that same name present, with permissions set so that it can't be deleted. Manually deleting the file will allow the processing to continue.

The error message will include the temporary file name and the permanent file name.

1.16.8.5 Unrecognized Command Code in DICOM Message

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Unexpected command value processing "xxxxxxxxxx". ***
*** Expected an integer numeric value, received "xxxxxxxxxx". ***
***
*** Routine: ^MAGDQUE4 Please Call Support Personnel ***
*****

```

This error occurs when the DICOM Gateway encounters an unrecognized command code in an incoming DICOM message.

The error message will include the name of the data stream being processed and the offending command code.

This is probably an error on the configuration of the peer application entity.

1.16.8.6 Unsupported Command in DICOM Message

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Unexpected command value processing "xxxxxxxxxx". ***
*** Cannot find label "xxxxxxxxxx" in routine ^MAGDQUE4. ***
***
*** Routine: ^MAGDQUE4 Please Call Support Personnel ***
*****

```

This error occurs when the DICOM Gateway encounters an unsupported command in an incoming DICOM message.

The error message will include the name of the data stream being processed and the offending command code.

This is probably an error on the configuration of the peer application entity.

1.16.9 TCP/IP Communications Errors

1.16.9.1 Cannot Connect to TCP/IP Socket (Set-Up Error)

```

*****
*** TCP not setup correctly ***
***
*** Connecting to IP Address "xxxxxxxxxx", port "xxxxxxxxxx". ***
*** Cannot open Socket ***
***
*** Routine: ^MAGDTCP1 Please Call Support Personnel ***
*****

```

This error occurs when the DICOM Gateway attempts to connect to a TCP/IP socket, and no connection can be established.

The error message will include the values of the IP Address and the port number that were used in the attempt to open the socket.

Typically, the cause of this error is that either the IP Address or the Port Number is incorrectly specified in the **F:\DICOM\Dict\SCU_List.DIC** master file. Notepad can be used to correct this file.

1.16.9.2 Cannot Find DICOM Message File to Send

```
*****
*** DICOM GATEWAY ERROR ***
***
*** Cannot find file "xxxxxxxxxx". ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****
```

This error occurs when the file containing the DICOM message file to be transmitted cannot be located.

1.16.9.3 Invalid First 8-bytes in DICOM Message

```
*****
*** DICOM GATEWAY ERROR ***
***
*** Incorrect DICOM Message -- first 8 bytes are wrong ***
*** Group=xxxxxxxxxx, Element=xxxxxxxxxx, Length=xxxxxxxxxx ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****
```

The first eight bytes of a DICOM message file should be “0 0 0 0 4 0 0 0” hexadecimal. This represents element (0000,0000) with a length of 4. This data must always be correct.

This error most likely occurs when the system is trying to send a zero-length DICOM file that was created due to running out of disk space.

1.16.9.4 Length-to-End Missing in DICOM Message

```
*****
*** DICOM GATEWAY ERROR ***
***
*** Incorrect DICOM Message -- Length to End not present ***
*** Group=xxxxxxxxxx, Element=xxxxxxxxxx, Length=xxxxxxxxxx ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****
```

The second required element in a VistA DICOM message is (0000,0001) which specifies the length to the end of the message. This is used by the TCP routine for performance to concatenate several messages into a single file (typically for Modality Worklist responses).

If this element is missing, it probably represents a programming error.

1.16.9.5 SOP Class UID Missing in DICOM Message

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Incorrect DICOM Message -- SOP Class UID not present ***
*** Group=xxxxxxxx, Element=xxxxxxxx, Length=xxxxxxxx ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****

```

The third required element in a VistA DICOM message (0000,0002) specifies the SOP Class UID for the message. This is needed to determine the Presentation Contact ID for the message.

If this element is missing, it probably represents a programming error.

1.16.9.6 Invalid SOP Class UID in DICOM Message

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Incorrect SOP Class UID format <<xxxxxxxxxx>> ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****

```

If the SOP Class UID for the message, element (0000,0002), has the wrong format (i.e., doesn't begin with a numeric), this error will occur.

This probably represents a programming error.

1.16.9.7 Unknown Presentation Context ID

```

*****
*** DICOM GATEWAY ERROR ***
***
*** Unknown Presentation Context ID for <<xxxxxxxxxx>> ***
***
*** Routine: ^MAGDTCP2 Please Call Support Personnel ***
*****

```

If the SOP Class UID for the message, element (0000,0002), has the wrong value this error will occur.

This is most likely caused by trying to send a message for a SOP Class that was not negotiated for the association. It is typically seen during "trial and error" configuration of vendor systems.

1.16.9.8 Unexpected PDU (expected 04 = P-Data)

```

*****
***  DICOM GATEWAY ERROR                               ***
***                                                    ***
*** Expected PDU type 04(hex),                          ***
*** found "xxxxxxxxxx".                                  ***
***                                                    ***
*** Routine: ^MAGDTCP3                                Please Call Support Personnel ***
*****

```

This error occurs when the DICOM Gateway attempts to receive data from a peer application entity and the code for the PDU Type is not equal to 04 (hex), which is the code for P-Data. (A PDU is a **P**rotocol **D**ata **U**nit.)

This is caused by a DICOM implementation problem with a vendor system.

1.17 Errors that may occur in Remote Procedure Calls

The error messages listed in this chapter may occur in Remote Procedure Calls. The message texts are typically displayed on the client-station from which the Remote Procedures were invoked.

1.17.1 MAG DICOM AUDIT COUNT (COUNT^MAGDRPC7)

This Remote Procedure may return the following error messages:

- -1, No Location Specified
- -2, No Message Specified

Both messages indicate that the software on the client-station omitted to specify a parameter that is needed for proper operation of the Remote Procedure.

1.17.2 MAG DICOM FIND LOCATION (FINDLOC^MAGDRPC8)

This Remote Procedure may return the following error message:

- -1, Invalid location "*name*"

This message indicates that there is no location (or division) at the current site with the name that was specified by the client-system. Either the client-system is mis-configured, or an end-user entered an invalid name.

1.17.3 MAG DICOM GET BASIC IMAGE (IMAGE^MAGDRPC2)

This Remote Procedure may return the following error messages:

- -1, Invalid IEN (*number*)
- -2, No data for "*number*".
- -3, Image #*number* has been deleted.

The first message (-1...) indicates that the value of the parameter that should identify the image for which information is to be retrieved is not a numeric value.

The second message (-2...) indicates that the value of this parameter has a numeric value, but there is no image with that (internal) number.

The third message (-3...) indicates that information is requested about an image that has been deleted. In this case, no information other than the error message is returned.

1.17.4 MAG DICOM GET NEXT QUEUE ENTRY (NEXTIMG^MAGDRPC4)

This Remote Procedure may return the following error messages:

- -1, No Origin Specified
- -2, No studies to be transmitted
- -3, No files to be transmitted

The first message indicates that the software on the client-station omitted to specify an essential parameter. Most likely, this error condition is the result of a mis-configuration of the client-station.

The second and third messages are informational only and indicate that all requested images (or studies) have been transmitted.

1.17.5 MAG DICOM GET PLACE (GETPLACE^MAGDRPC8)

This Remote Procedure may return the following error message:

- -1, Location "*name*" not found.

This message indicates that there is no location (or division) at the current site with the name that was specified by the client-system. Either the client-system is mis-configured, or an end-user entered an invalid name.

1.17.6 MAG DICOM GET SERVICE INFO (SERVICE^MAGDRPC2)

This Remote Procedure may return the following error message:

- -1, No Imaging Site Parameters defined

This message indicates that the configuration of the Imaging Site Parameters (on the VistA) system still has to be completed. Operation on the client cannot proceed until the configuration on the server is complete.

1.17.7 MAG DICOM HL7 POINTER ACTION (HL7PTR^MAGDRPC8)

This Remote Procedure may return the following error messages:

- -1, Invalid Pointer "*value*".

- -2,Invalid Request: "*action*".

The first of these messages indicates that the value of the parameter that identifies the HL7-message, as specified by the software on the client-station, cannot be used to identify a message that is currently present in the system. Either the software on the client system is working from obsolete information (server was purged since last update of the client), or an end-user entered an erroneous pointer number.

The second message indicates that the software on the client-station specified an erroneous value for the parameter that indicates the action to be taken by the Remote Procedure. This condition can only be resolved with a correction to the client-software.

1.17.8 MAG DICOM IMAGE PROCESSING (ENTRY^MAGDIR8)

This Remote Procedure may return the following error messages:

- xxxxx
Problem detected by routine xxxxx
Error Code: xxxxx

Messages of this type are always followed by one of the messages in section 1.16.6. The subsequent message provides more detailed information about the error. See section 1.16.6 for these details.

1.17.9 MAG DICOM INCORRECT IMAGE CT (CORRECT^MAGDRPC8)

This Remote Procedure may return the following error messages:

- -1,No Location Specified
- -2,No Gateway Specified

These messages indicate that the software on the client-station omitted to specify a value for a required parameter. Most likely, these error conditions are the result of a mis-configuration of the client-station.

1.17.10 MAG DICOM LOOKUP STUDY (LOOKUP^MAGDRPC4)

This Remote Procedure may return the following error messages:

- -1,No Case or Consult Number Specified
- -2,No Study Date Specified
- -3,Invalid study date "*text*".
- -4,Consult is cancelled
- -5,Consult/procedure not on file

The first two messages indicate that the software on the client-station omitted to specify a required parameter. Such an omission may happen when an end-user just presses the enter-key, rather than provide a complete answer to a question.

The other three messages indicate cases where the provided information cannot be matched successfully with information in the VistA database.

1.17.11 MAG DICOM NETWORK STATUS (ONOFFLINE^MAGDRPC5)

This Remote Procedure may return the following error message:

- -1, No Network Location Specified

This message indicates that the client system omitted to provide a value for the parameter that indicates the location (or division) at the current site from which images are to be transmitted. Most likely, the client-system is mis-configured..

1.17.12 MAG DICOM PACS CUTOFF DATE (CUTOFF^MAGDRPC1)

This Remote Procedure may return the following error messages:

- -2, No PACS Installed
- -3, PACS Retention Parameter not defined

Both messages are informational only. The first message (-2,...) indicates that the site-(location, division) is configured as if no PACS is installed. The other message indicates that as yet no value for the retention period parameter has been entered into the site configuration.

1.17.13 MAG DICOM PACS MINIMUM SPACE (MINSPACE^MAGDRPC1)

This Remote Procedure may return the following error messages:

- -2, No PACS Installed
- -3, PACS Minimum Space Percentage Parameter not defined

Both messages are informational only. The first message (-2,...) indicates that the site-(location, division) is configured as if no PACS is installed. The other message (-3,...) indicates that as yet no value for the minimum disk space parameter has been entered into the site configuration.

1.17.14 MAG DICOM QUEUE IMAGE (QUEUE^MAGDRPC3)

This Remote Procedure may return the following error messages:

- -1, No Image specified
- -2, No Destination specified
- -3, No Origin specified
- -4, Cannot Queue Image Object Type "*type*".

The first three messages indicate that the software on the client-station omitted to specify values for required parameters. Most likely, such a situation is the result of a failure to make a selection at an earlier stage in the process.

The fourth message (-4,...) indicates that the requested image is of such a nature that it is not permitted to enter that image into this transmission queue. Only images that are of the following “types” can be transmitted using the DICOM C-Store protocol:

- Type 3 = X-Ray
- Type 11 = X-Ray Group
- Type 100 = DICOM Image

1.17.15 MAG DICOM QUEUE INIT (INIT^MAGDRPC4)

This Remote Procedure may return the following error messages:

- -1, No entries at all in queue.
- -2, No entries present for "*name*".
- -3, No origin specified.

The last message (-3,...) indicates that the client system omitted to provide a value for the parameter that indicates the location (or division) at the current site from which images are to be transmitted. Most likely, the client-system is mis-configured.

The other two messages are informational only. Both indicate that currently no images are queued up to be transmitted.

1.17.16 MAG DICOM ROUTE EVAL LOG (EVALLOG^MAGDRPC6)

This Remote Procedure may return the following error messages:

- -1, No task #*number*
- -2, MAXIMUM parameter = *number* < 1
- -3, *text*

The first message (-1, ...) indicates that the client-system specified a value for the TaskMan task number that does not correspond to any active task. Most likely, the VistA system has been restarted and some clean-up or house-keeping tasks have been performed.

The second message (-2, ...) indicates that the value for the parameter that specifies the maximum number of log-messages that are to be returned to the client has an unusable value. This condition can only be remedied by correcting the software on the client-system.

The third message (-3, ...) is informational only. This message indicates that the TaskMan task, for which messages were being relayed, has stopped. The text that is included in this message contains the reason that TaskMan provided for ending the task.

1.17.17 MAG DICOM ROUTE EVAL START (START^MAGDRPC5)

This Remote Procedure may return the following error messages:

- -1, No Location Specified
- -2, No Routing Rules Specified

- -3, A Rule Evaluator is Already Running for *location*

The first two messages indicate an error or omission in the set-up of the DICOM Gateway from which the request was issued.

The third message indicates that a rule-evaluator is already running for the requested location. Only one evaluation may be running for any given location (or division) at any time.

1.17.18 MAG DICOM ROUTE NEXT FILE (XMIT^MAGDRPC5)

This Remote Procedure may return the following error messages

- -1, No Location Specified
- -2, No Valid Destinations Specified
- -1~No routable files found for image *nnnnn*

The first message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

The second message indicates that the end-user did not select any destinations to which files might be transmitted.

The third message indicates that a queue entry is being processed for an image that has no actual files associated with it, i.e. there are no files to be transmitted for that image. Either the image is currently in a “deleted” state, or there were never any files successfully stored for this image.

1.17.19 MAG DICOM ROUTE PURGE DONE (PURGDONE^MAGDRPC6)

This Remote Procedure may return the following error message:

- -1, No Location Specified

This message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

1.17.20 MAG DICOM ROUTE REMOVE OBSO (REMOBSO^MAGDRPC6)

This Remote Procedure may return the following error message:

- -1, No Location Specified

This message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

1.17.21 MAG DICOM ROUTE REQUEUE (REQUEUE^MAGDRPC6)

This Remote Procedure may return the following error message:

- -1, No Location Specified

This message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

1.17.22 MAG DICOM ROUTE TRANSACT STS (TRANSTS^MAGDRTIM)

This Remote Procedure may return the following error messages:

- -1, No Location Specified
- -2, Invalid Transaction ID

The first message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

The second message indicates an error in the software that performed the request from the client station. A status can only be provided for a “known” transaction. The client system receives transaction identifiers as new transactions are created.

1.17.23 MAG DICOM ROUTE VALID DEST (VALDEST^MAGDRPC1)

This Remote Procedure may return the following error message:

- -1, "xxxxxx" is not a valid destination

This message indicates that an attempt to re-configure a DICOM Gateway included a name for a location that is not defined in the VistA database.

1.17.24 MAG DICOM SET PACS PARAMS (SETPACS^MAGDRPC8)

This Remote Procedure may return the following error messages:

- -1, No Place Specified.
- -2, Invalid Place Specified.

The first message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.

The second message indicates that an attempt to re-configure a DICOM Gateway included a name for a location that is not defined in the VistA database.

1.17.25 MAG DICOM TEXT PROCESSING (ENTRY^MAGDHR51)

This Remote Procedure may return the following error message:

- -1 ^MAGDHL7(2006.5,nnnnnn,mmm,...) is incomplete

This message indicates that, after waiting the maximum time, an HL7 message is still not complete. Most likely, the process that generated the message did not complete successfully.

1.17.26 MAG DICOM UPDATE GATEWAY NAME (UPDTGW^MAGDRPC8)

This Remote Procedure may return the following error messages:

- -1, No Gateway Name specified
- -2, No Gateway Location specified

These messages indicate omissions in the set-up of the DICOM Gateway from which the request was issued.

1.17.27 MAG DICOM UPDATE SCU LIST (UPDTAPP^MAGDRPC8)

This Remote Procedure may return the following error message:

- -1, Missing or Inconsistent Parameters.

This indicates an omission or inconsistency in the set-up of the DICOM Gateway from which the request was issued.

1.17.28 MAG DICOM WORKSTATION VERSION (STATION^MAGDRPC1)

This Remote Procedure may return the following error message:

- -1, No Station Identifier Specified

This message indicates an omission in the set-up of the DICOM Gateway from which the request was issued.