



PROSTHETICS VISTA SUITE (GUI)

INSTALLATION GUIDE

PATCH RMPR*3*96

Version 3.0

January 2005

Department of Veterans Affairs
VistA Health System Design and Development (HSD&D)

Table of Contents

Prosthetics VISTA Suite Installation Guide	1
Overview	1
Client Installation.....	3
Extract the file.....	6
Install the GUI Application(s)	8

Prosthetics VISTA Suite (GUI) Installation Guide

Overview

Introduction

With Patch RMPR*3*96, Prosthetics offers the **View Prosthetics Billing Information** feature, which is a GUI windows-based functionality available to users. This *Prosthetics Vista Suite (GUI) Installation Guide* will help IRM and users to download the GUI portion of the patch. See the patch description for the instructions for the Vista portion of the patch. Here are the four guides that are part of the GUI Vista suite and are available on the Vista Document Library (VDL):

- Prosthetics Main Menu (GUI) User Manual
- View Prosthetic Billing Information User Manual
- Delayed Order Reports (DOR) User Manual
- NPPD Detail Display Release Notes

The new **Prosthetics Main Menu** window allows you to now access all GUI applications for Prosthetics.

Note to IRM IMPORTANT

- RMPR_3_96.EXE, the full VISTA suite, is only to be installed in Prosthetic users' computers.
- The RMPR_3_96B.EXE is only to be installed on Billing or Finance users' desktops. **This executable restricts access to users outside of Prosthetics to only the application that is relevant to their job responsibilities.**

Note: While the majority of the installation for the two files is identical, there are certain steps that are different in this Install Guide and will be noted.

Important: **Remove the VISTA Prosthetics Suite**

You must uninstall the Vista Prosthetics suite. To do this, go to the Control Panel → Add/Remove Programs option and select the Vista Prosthetics suite item.

Click the **REMOVE** button and the program will automatically uninstall.

RPC Broker

The assumption is made that the user's PC has the required RPC Broker Client Workstation set up. If this is not the case, the user should first install this software.

The installation is beyond the scope of this installation guide, but the RPC Broker files and installation instructions can be found on this VA's RPC Broker web page: <http://vista.med.va.gov/broker/index.asp>.

These files must be installed before running the GUI application.

Continued on next page

Overview, Continued

HOSTS file

A link must be established between the client and the server via the PC's Hosts file and the PC's Windows Registry. The Hosts file (no extension) is in the following:

C:\Windows	for Windows 95 or 98
C:\Windows\System32\Drivers\Etc	for Windows 2000 / XP

If the Hosts file does not exist, it must be created. (There may be a hosts.sam sample file already in place.) Note that if you create the file with Notepad, Notepad automatically adds a .txt extension.

To remove the extension, use Explorer to rename it, leaving off the .txt extension.


Client Installation

Hosts File	If the Hosts file already exists, do not delete anything.
-------------------	---

Steps	Follow these instructions for updating the Hosts file.
--------------	--

Step	Action
1	<p>If the user needs to access only one GUI server, and the ServerPort is 9200, then add a line to the Hosts file, as shown in the example (on the next page) using the IP address provided by the server administrator. (There must be at least one space between the address and the name. There can be no spaces in the name itself.)</p> <p>Continue to the next page.</p>

Notepad – Sample HOSTS file

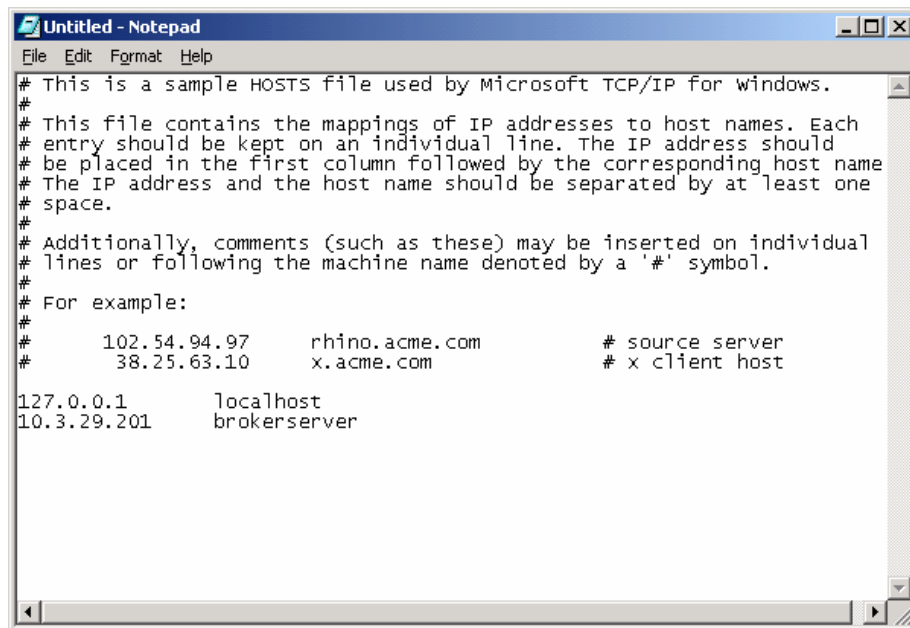

```
File Edit Format Help
# This is a sample HOSTS file used by Microsoft TCP/IP for Windows.
#
# This file contains the mappings of IP addresses to host names. Each
# entry should be kept on an individual line. The IP address should
# be placed in the first column followed by the corresponding host name
# The IP address and the host name should be separated by at least one
# space.
#
# Additionally, comments (such as these) may be inserted on individual
# lines or following the machine name denoted by a '#' symbol.
#
# For example:
#
#       102.54.94.97       rhino.acme.com       # source server
#       38.25.63.10       x.acme.com          # x client host
127.0.0.1       localhost
10.3.29.201     brokerserver
```

Continued on next page

Client Installation, Continued

Steps (continued)

Continue to follow these instructions for updating the Hosts file (if the Hosts file already exists.)

Step	Action
2	If the user needs to access only one GUI Server and the ServerPort is <u>not</u> 9200, then add one line to the Hosts file as shown in scenario 1 on the previous page, and run the ServerList utility described in scenario 3 below, or use the command line parameters shown in the next section.
3	If the user needs to access multiple GUI Servers and/or the ServerPort is NOT 9200, and the GUI Server names to the Hosts file: <ul style="list-style-type: none">• 152.111.222.333 BrokerServer1 (The names can be anything you like.)• 152.111.222.444 BrokerServer2

Continued on next page

Client Installation, Continued

Multiple servers

Then run the ServerList Utility, installed in the one of the two folders:

For Prosthetics Users: C:\Program Files\Vista\PROS folder.

For Billing Users: C:\Program Files\Vista\PROS\ViewBilling folder.

- The left windowpane of the ServerList utility shows those servers that have been set up as Host Servers in the HOSTS file.
- The right windowpane shows those servers that have been set up as BrokerServers.

Steps


Continue to follow these instructions for updating the Hosts file.

Step	Action
4	<p>To enable a BrokerServer, double-click on the name in the left window to move it to the right window, and then enter the corresponding ServerPort number.</p> <p>Note: The order of the names in the right list is the order they will be displayed in the drop-down box when starting the GUI feature. If only one server is listed on the right, then you won't see the drop down box – it will automatically select that server as the default).</p>

“Del” column

A left-click in the “**Del**” column marks and un-marks server for deletion. See the following ServerList screen print for an example.

ServerList

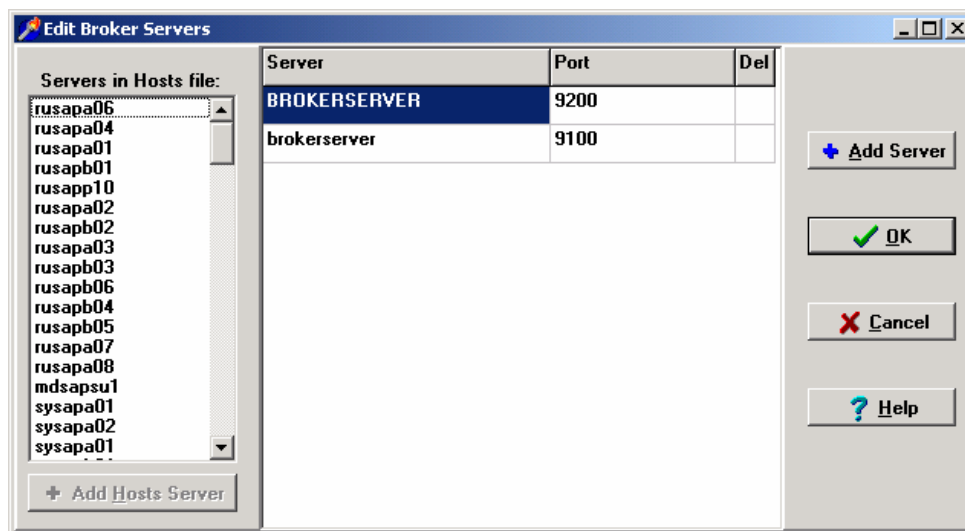

Extract the file

Anonymous directory

The file for this patch can be obtained from the ANONYMOUS.SOFTWARE directory at one of the OI Field Offices. The preferred method is to FTP the file from DOWNLOAD.VISTA.MED.VA.GOV, which will transmit the file from the first available server.

Alternatively, you may elect to retrieve the file from a specific OI Field Office. The documentation can be read via the Adobe Acrobat Reader browser program. The .EXE file is binary file and must be transferred using binary file transfer tools.

Steps

To begin the install process and unzip the files, follow these steps:

Step	Action
1	From the Anonymous Directory, download the file RMPR_3_96.EXE or RMPR_3_96B.EXE .
2	Double click the RMPR_3_96.EXE or RMPR_3_96B.EXE icon.

Warning: Wrong filename

Warning: If the filename is **RMPR_3_96.EXE;1** when you try to download it, then you **must** rename it to **RMPR_3_96.EXE** so that the filename does not have an “;1” in it.

Rename the file

To change the filename, right click on **RMPR_3_96.EXE;1**, (executable application) and select the **Rename** option. Then delete the “;1” from the filename.

Continued on next page

Extract the file, Continued

Unzipping files


The file **RMPR_3_96.EXE** or **RMPR_3_96B.EXE** is a zip file that contains the “setup” files needed for the automated installation. For setup purposes, the files will unzip to:

- C:\PROS\ folder or
- C:\PROS\ViewBilling folder

As an alternative to extracting these files to every client computer, IRM staff can unzip this file to a shared folder (directory) on a network server and then run the installation “SETUP.EXE” remotely.

Recommendation: Do not change the file folder, as future patches will need to overwrite these files.

WinZip dialog box sample (Your dialog box will contain the most current version number for the files.)

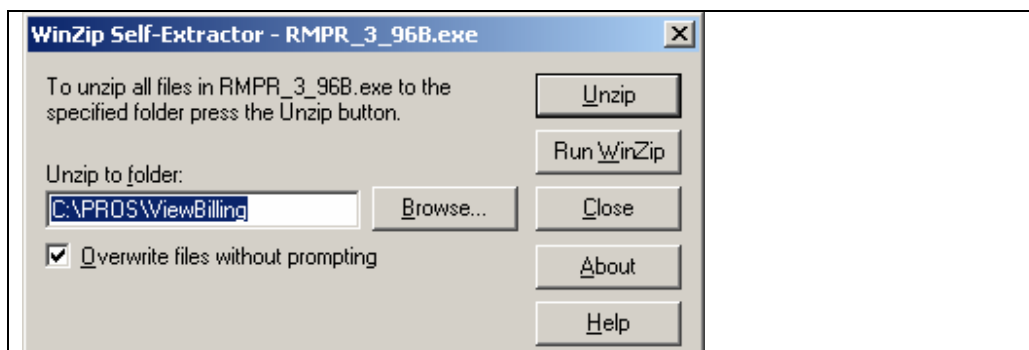

Steps

To continue to unzip the files, follow these steps:

Step	Action
3	In the WinZip Self-Extractor dialog box, click the Unzip button (as shown above).
4	A successful message box displays when the unzip process has completed. (See below.)
5	Click the OK button on the WinZip Self-Extractor confirmation box.
6	The WinZip Self-Extractor dialog box continues to display (as shown above). Click the Close button.

Successful message box



Install the GUI Application(s)

Automatic Wizard


Once you have extracted the files onto your c-drive, you can perform an installation using an automatic wizard that is provided for you. This feature has combined several steps to make installation easy.

Steps

To begin the installation, follow these steps:

Step	Action
1	Open the My Computer icon on your desktop by double clicking it.
2	Double click the C-drive in your Windows Explorer view.
3	Double click the folder that contains the files: PROS .
4	Double click the folder that contains the files: PROS or PROS\ViewBilling

Windows Explorer

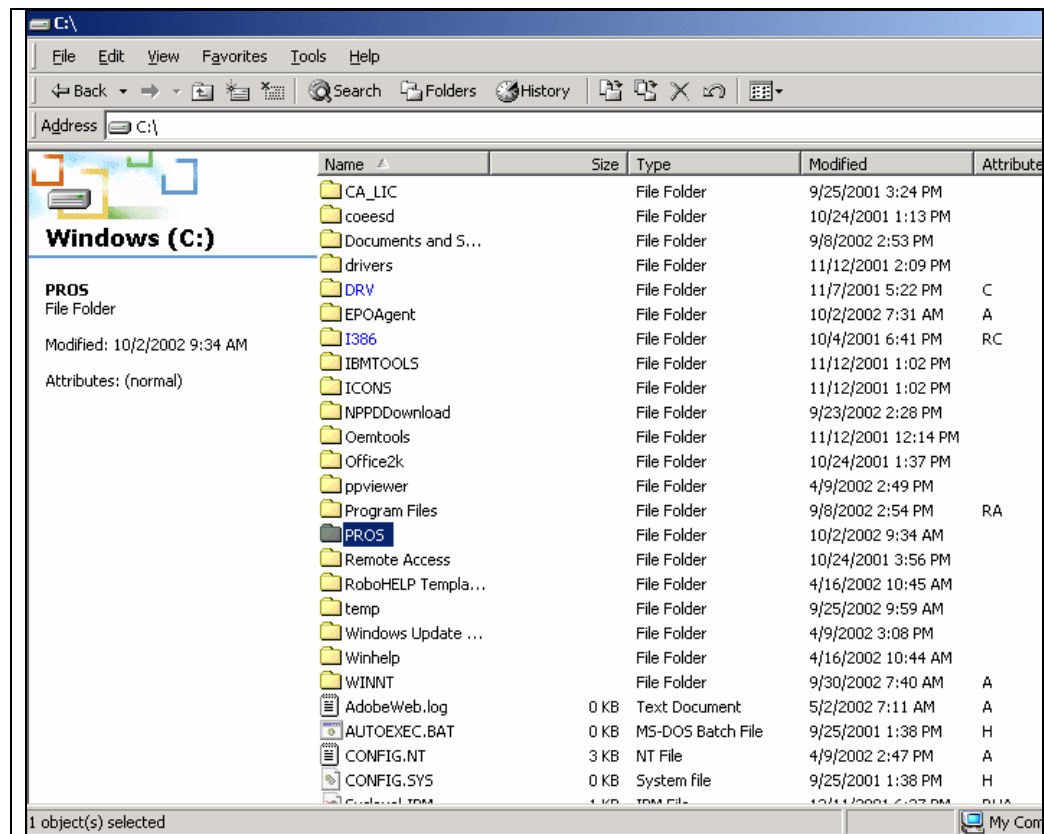

Continued on next page


Install the GUI Application(s), Continued

Steps (continued)

To run the installation program, follow these steps:

Step	Action
4	Double click the Setup.exe file in the C:\PROS or C:\PROS\ViewBilling folder shown on your Windows Explorer view.

Windows Explorer

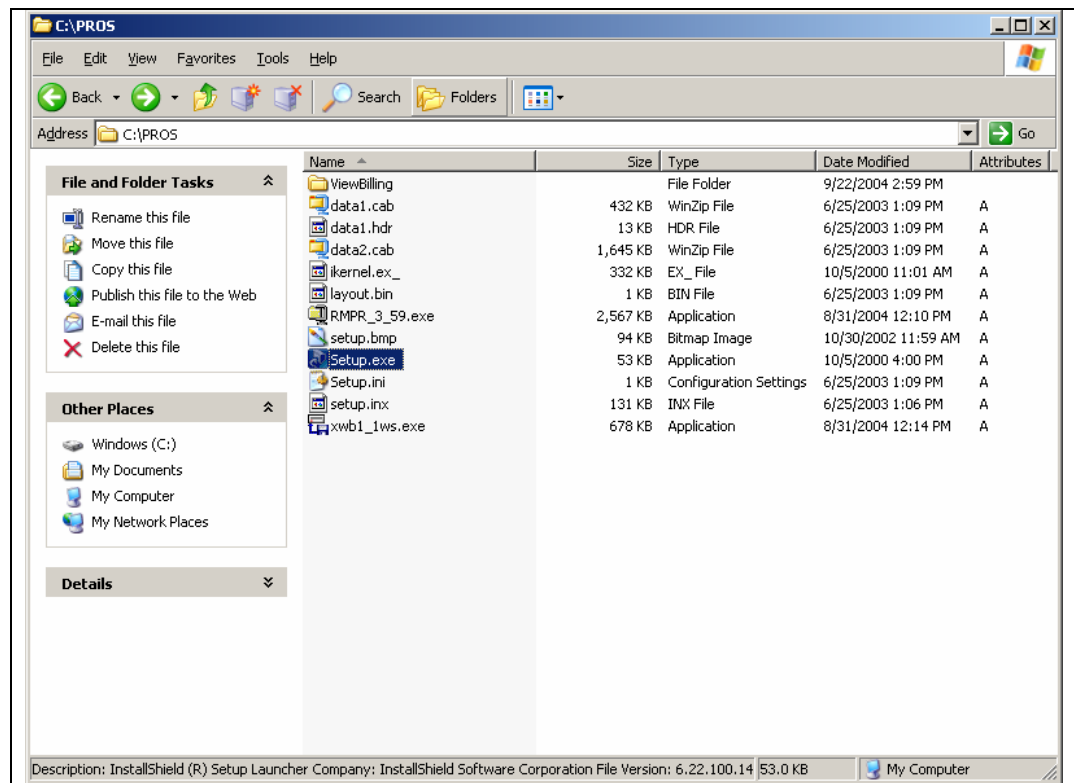

Continued on next page


Install the GUI Application(s), Continued

Steps (continued)


To run the wizard for the installation, follow these steps:

Step	Action
5	The Install Shield Wizard begins and displays the download process.
6	Click the Next button on the first window of the Install Shield Wizard for the GUI application.

Install Shield Wizard

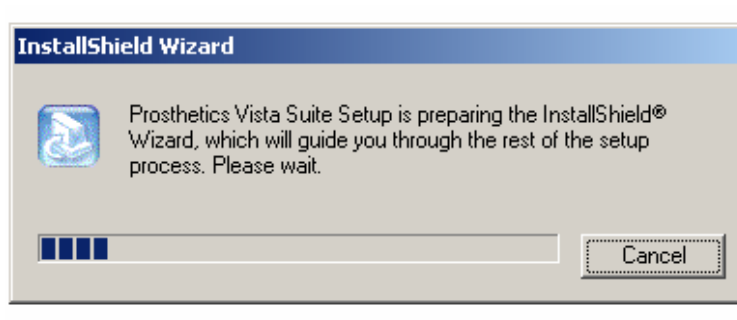

Install Shield Wizard

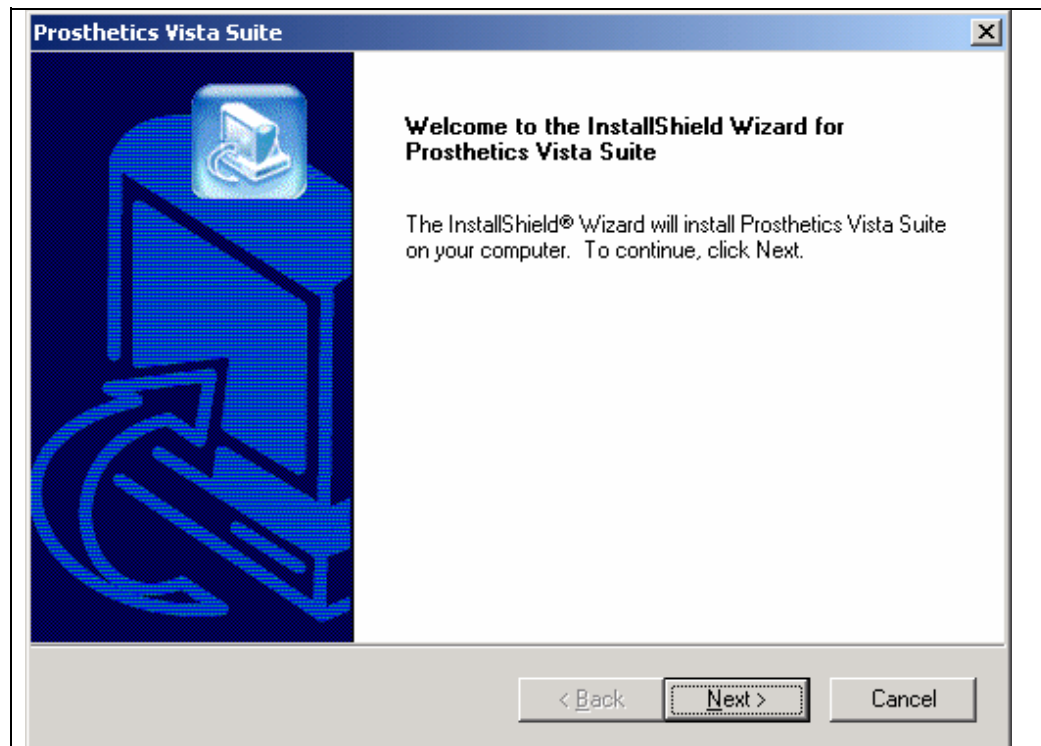

Continued on next page


Install the GUI Application(s), Continued

Steps (continued)

To continue the installation, follow these steps:

Step	Action
7	Click the Yes button on the License Agreement if you agree to the terms of the license. Note: If you choose No , the setup will automatically close.

License Agreement

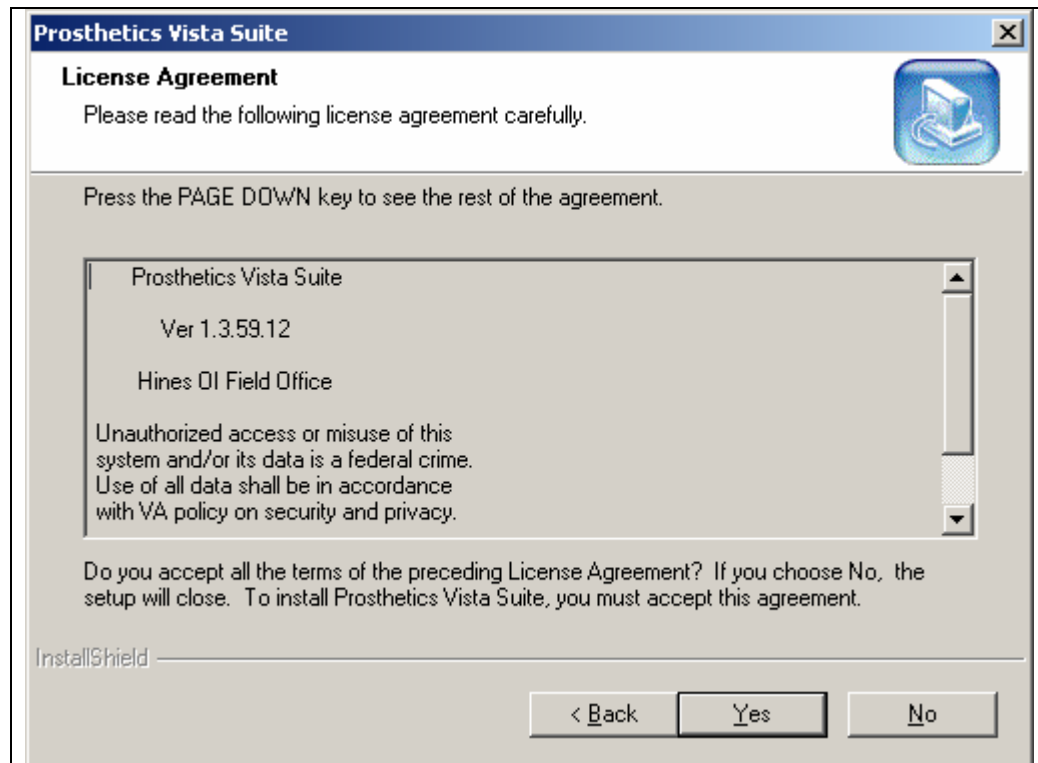

Continued on next page


Install the GUI Application(s), Continued

Recommendation The destination folder **C:\Program Files\Vista\PROS** or **C:\Program Files\Vista\PROS\ViewBilling** has been established by VA policy so as not to conflict with other VA programs. It is recommended that you do NOT change this portion of the path.


Steps To continue the installation, follow these steps:

Step	Action
8	Click the Next button on the Destination Location dialog box.
9	Click OK on the confirmation box. Notice that it states that the PROS file folder was created.

File Destination Location

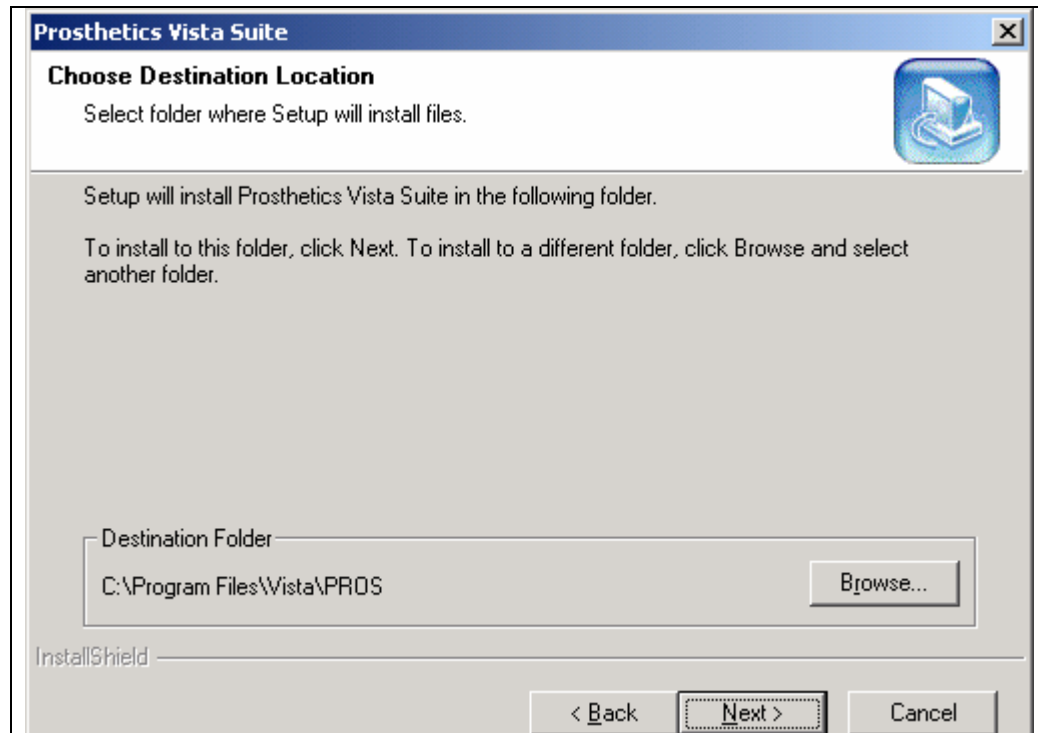

Create Directory message

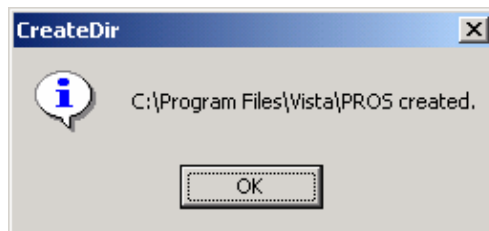

Continued on next page


Install the GUI Application(s), Continued

Steps (continued)

To continue the installation, follow these steps:

Step	Action
10	Click the Next button on the Start Copy Files dialog box.

Start Copying Files

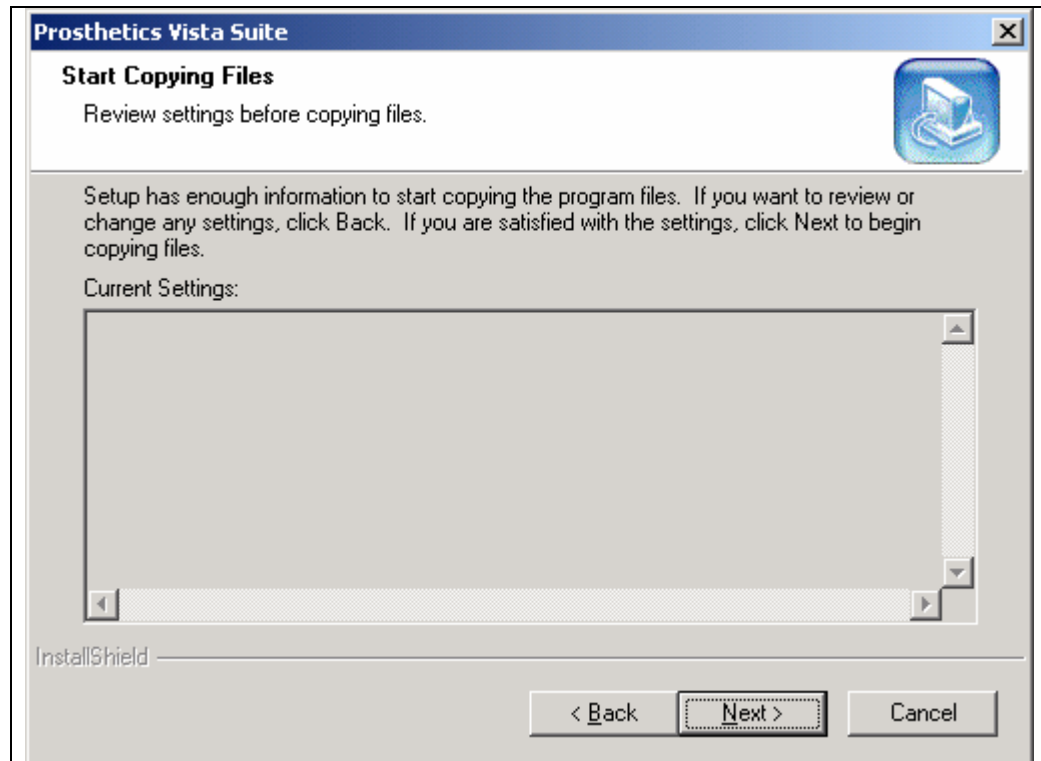

Continued on next page


Install the GUI Application(s), Continued

Steps (continued)

To continue to install the Prosthetics feature onto your desktop, follow these steps:

Step	Action
11	The Install Shield Wizard will be complete with this step.
12	Click the Finish button.

Install Shield Wizard

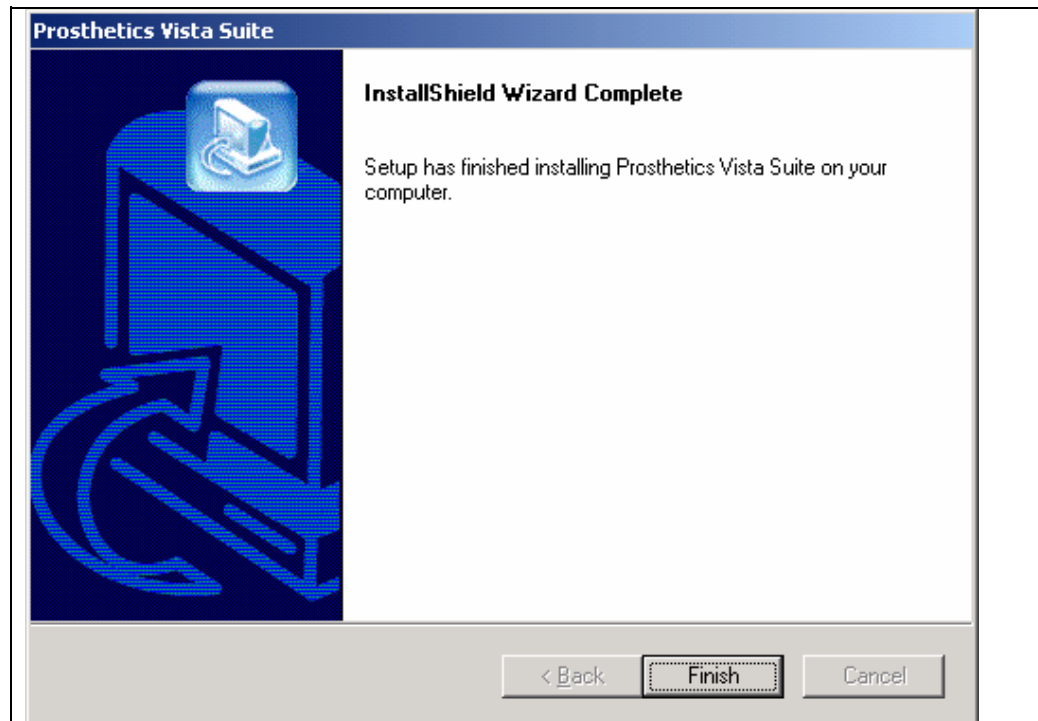

Installation Complete

The application is now installed on **C:\Program Files\Vista\PROS** or **C:\Program Files\Vista\PROS\ViewBilling**. An icon has been placed on your desktop to run the application. The application can also be run from:

- **For Prosthetics users:** Start → Programs → Prosthetics VISTA Suite
- **For Billing Users:** Start → Programs → View Prosthetics Billing Menu.

Uninstall procedures

To uninstall the application, go to the **Control Panel → Add/Remove Programs** option and select the **Prosthetics Vista Suite** item or **View Billing Menu**. Click the **REMOVE** button and the program will automatically uninstall. To complete the full removal, the setup files will need to be deleted from the file location you selected.