

Clinical Reminders

Patch PXRM*2*11 National Reminder and Dialog Updates

INSTALLATION & SETUP GUIDE

October 2008

Health Provider Systems Department of Veterans Affairs

Contents

INTRODUCTION1
Web Sites1
PRE-INSTALLATION
Required Software
Estimated Installation Time
Pre-Installation Steps2
INSTALLATION
1. Retrieve the PXRM*2.0*11 build
2. Install the build first in a training or test account
3. Load the distribution
a. Backup a Transport Global
b. Compare Transport Global to Current System
c. Verify Checksums in Transport Global
d. Print Transport Global (optional)
4. Install the build
5. Install File Print (optional)7
6. Build File Print (optional)7
7. Post-installation routine
8. Deletion of init routines
SETUP AND MAINTENANCE
Overview9
Setup Steps 12
1. Repeat the pre-installation reminder inquiry and compare the pre- and post-installation output12
2. Map Local Findings to National Terms
3. Edit Reminder Dialogs 14
Mapping a consult quick order for the VA-TBI SCREENING reminder dialog:
4. Verify that the dialogs function properly17
APPENDIX A: INSTALLATION EXAMPLE
APPENDIX B: SETTING UP A TIU/HS OBJECT FOR THE AUDIT-C TEST
APPENDIX C: COMPUTED FINDING DESCRIPTIONS

Introduction

Patch PXRM*2.0*11 includes the following modifications to the Clinical Reminders package:

- Updates eleven reminder definitions (and their dialogs if there is one attached), based on changes required by the Office of Patient Care Services and problems reported in Remedy tickets.
- Installs four new reminders.
- Installs 11 new or modified service-related computed findings.
- Updates URLs in six reminders.
- Modifies Location List exclusion functionality
- Modifies Reminder Test output
- Corrects several problems reported in Remedy tickets:
 - 229723 Nat'l reminder component "PXRM AUDC RESULT has errors
 - 233994 Discrepancy between VA-Iraq/Afghan Post Deploy reminder dialog and clinical maintenance display
 - 242238 Duplicate text in AIMS test
 - 239704 Display inconsistent for VA-IRAQ&AFGHANISTAN POST DEPLOYMENT
 - 234380. National OIF reminder is producing conflicting information in clinical maintenance vs the reminder dialog concerning section 4
 - 238209 Functionality of I/A and TBI reminders

See the Setup section of this manual for instructions about mapping terms and dialogs.

Web Sites

Site	URL	Description
National Clinical Reminders site	http://vista.med.va.gov/reminders	Contains manuals, PowerPoint presentations, and other information about Clinical Reminders
VA Mental Health Home	http://vaww.mentalhealth.va.gov/	Contains items of interest for all VA mental health clinicians
VistA Document Library	http://www.va.gov/vdl/	Contains manuals for Clinical Reminders and Mental Health YS*5.01*85

Pre-Installation

This manual describes how to install and set up Clinical Reminders patch PXRM*2.0*11.

Package/Patch	Namespace	Version	Comments
Clinical Reminders	PXRM	2.0	Fully patched
Health Summary	GMTS	2.7	Fully patched
HL7	HL	1.6	Fully patched
Kernel	XU	8.0	Fully patched
MailMan	XM	7.1	Fully patched
Mental Health	YS	5.01	
YS*5.01*85			
Pharmacy Data Management	PSS	1.0	
PSS*1.0*112			
Outpatient Pharmacy	PSO	7.0	
PSO*7.0*245			
VA FileMan	DI	22	Fully patched

Required Software

Estimated Installation Time

Setup before installation Installation Setup after installation by the Reminder Manager or CAC to address local implementation Approximately 30 minutes About 5-10 minutes Approximately one hour (depending partly on the experience level of the manager or CAC)

Pre-Installation Steps

1. Do a reminder definition inquiry on the following reminders and save for postinstallation comparisons:

VA-ALCOHOL AUDIT-C POSITIVE F/U VA-ALCOHOL USE SCREEN (AUDIT-C) VA-BL DEPRESSION SCREEN VA-BL PTSD SCREEN VA-BL OEF/OIF FEVER VA-BL OEF/OIF GI SX VA-BL OEF/OIF SKIN SX VA-DEPRESSION SCREENING VA-IRAQ & AFGHAN POST-DEPLOY SCREEN VA-MHV CERVICAL CANCER SCREEN VA-MHV COLORECTAL CANCER SCREEN VA-MHV DIABETES FOOT EXAM VA-MHV DIABETES RETINAL EXAM VA-MHV HYPERTENSION VA-MHV INFLUENZA VACCINE VA-MHV MAMMOGRAM SCREENING VA-MHV PNEUMOVAX VA-PTSD SCREENING VA-TBI SCREENING

2. Do a reminder term inquiry and save for post-installation comparisons

The following terms have changed and may need to be remapped, if previously mapped locally. **The ones in bold will almost definitely need remapping.** Please use the option, "Inquire about a reminder term" on these terms for post-installation comparison and remapping.

VA-COGNITIVE IMPAIRMENT **VA-DEPRESSION SCREEN OEF/OIF** VA-LDL VA-MHV BARIUM ENEMA VA-MHV COLONOSCOPY VA-MHV DIABETIC FOOT EXAM COMPLETE VA-MHV DIABETIC FOOT EXAM INSPECTION VA-MHV DIABETIC FOOT EXAM MONOFILAMENT VA-MHV DIABETIC FOOT EXAM PULSES VA-MHV HBA1C <8.0 VA-MHV HBA1C >7.9 VA-MHV HBA1C ALL RESULTS VA-MHV HDL – lab test and outside lab test VA-MHV INFLUENZA IMMUNIZATION VA-MHV LIPID PANEL ORDERABLE ITEMS VA-MHV OCCULT BLOOD PANEL VA-MHV OCCULT BLOOD SINGLE TEST VA-MHV OUTSIDE LDL VA-MHV OUTSIDE LDL>99 VA-MHV PNEUMOCOCCAL IMMUNIZATION VA-MHV SIGMOIDOSCOPY **VA-MHV TOTAL CHOLESTEROL** VA-MHV TRIGLYCERIDES VA-PTSD SCREEN VA-WH HYSTERECTOMY W/CERVIX REMOVED VA-WH MAMMOGRAM ORDER VA-WH MAMMOGRAM SCREEN DEFER VA-WH MAMMOGRAM SCREEN DONE VA-WH MAMMOGRAM SCREEN NOT INDICATED

VA-WH PAP SMEAR DONE VA-WH PAP SMEAR OBTAINED VA-WH PAP SMEAR ORDER HEALTH FACTOR VA-WH PAP SMEAR SCREEN NOT INDICATED

3. Create an Audit-C TIU/HS Object for the VA-Alcohol Use Screen Reminder (see instructions in <u>Appendix B.</u>

4. The reminder dialog element VA-PDIQ POLYTRAUMA CONSULT, which is used to order a TBI consult, is overwritten in the install. You will need to re-enter the finding item for this element using the menu or quick order that you use to order a TBI consult.

Installation

This build can be installed with users on the system, but it should be done during non-peak hours.

The installation needs to be done by a person with DUZ(0) set to "@."

1. Retrieve the PXRM*2.0*11 build

Use ftp to access the build – the name of the host file is PXRM_2_0_11.KID – from one of the following locations:

Albany	REDACTED	REDACTED	
Hines	REDACTED	REDACTED	
Salt Lake City	REDACTED	REDACTED	

2. Install the build first in a training or test account.

Installing in a non-production environment will give you time to get familiar with new functionality and complete the setup for reminders and dialogs prior to installing the software in production.

3. Load the distribution.

In programmer mode type, D ^XUP, select the Kernel Installation & Distribution System menu (XPD MAIN), then the Installation option, and then the option LOAD a Distribution. Enter your directory name and PXRM_2_0_11.KID at the Host File prompt.

Example

```
Select Installation Option: LOAD a Distribution
Enter a Host File: PXRM_2_0_11.KID
KIDS Distribution saved on June 14, 2008@11:53:53
Comment: PXRM*2.0*11
```

From the Installation menu, you may elect to use the following options:

a. Backup a Transport Global

This option will create a backup message of any routines exported with the patch. It will NOT back up any other changes such as DDs or templates.

b. Compare Transport Global to Current System

This option will allow you to view all changes that will be made when the patch is installed. It compares all components of the patch (routines, DDs, templates, etc.).

c. Verify Checksums in Transport Global

This option will allow you to ensure the integrity of the routines that are in the transport global. If there are any discrepancies, do not run the Install Package(s) option. Instead, run the Unload a Distribution option to remove the Transport Global from your system. Retrieve the file again from the anonymous directory (in case there was corruption in FTPing) and Load the Distribution again. If the problem still exists, log a Remedy ticket and/or call the national Help Desk (1-888-596-HELP) to report the problem.

d. Print Transport Global (optional)

This option will allow you to view the components of the KIDS build.

4. Install the build.

From the Installation menu on the Kernel Installation and Distribution System (KIDS) menu, run the option Install Package(s). Select the build PXRM*2.0*11 and proceed with the install. If you have problems with the installation, log a Remedy ticket and/or call the National Help Desk to report the problem.

Select Installation & Distribution System Option: Installation Select Installation Option: INSTALL PACKAGE(S) Select INSTALL NAME: PXRM*2.0*11

Answer "NO" to the following prompt:

Want KIDS to INHIBIT LOGONs during install? NO// NO

NOTE: DO NOT QUEUE THE INSTALLATION, because this installation may ask questions requiring responses and queuing will stop the installation. The most common are replacements for finding items or quick orders during the installation of Reminder Exchange file entries.

The reminder dialog element VA-PDIQ POLYTRAUMA CONSULT, which is used to order a TBI consult, is overwritten in the install. You will need to re-enter the finding item for this element using the menu, order dialog, or quick order that you use to order a TBI consult.

Installation Example

See Appendix A.

```
Setting up programmer environment
Terminal Type set to: C-VT220
Select OPTION NAME: XPD MAIN Kernel Installation & Distribution System menu
Select Kernel Installation & Distribution System Option: Installation
Select Installation Option: Install Package(s) [XPD INSTALL BUILD]
Select INSTALL NAME: PXRM*2.0*11
```

5. Install File Print (optional)

Use the KIDS Install File Print option to print out the results of the installation process.

```
Select Utilities Option: Install File Print
Select INSTALL NAME: PXRM*2.0*11
```

6. Build File Print (optional)

Use the KIDS Build File Print option to print out the build components.

```
Select Utilities Option: Build File Print
Select BUILD NAME: PXRM*2.0*11 CLINICAL REMINDERS
DEVICE: HOME//
```

7. Post-installation routine

The installation will place the following exchange file entries in the Reminder Exchange file #811.8:

PATCH 11 ITEMS

VA-ALCOHOL AUDIT-C POSITIVE F/U VA-ALCOHOL USE SCREEN (AUDIT-C VA-BL DEPRESSION SCREEN VA-BL PTSD SCREEN VA-BL OEF/OIF EMBEDDED FRAGMENTS VA-BL OEF/OIF FEVER VA-BL OEF/OIF GI SX VA-BL OEF/OIF SKIN SX VA-CV ELIG W/HF FOR NO SERVICE VA-CV INELG W/HF FOR SERVICE VA-DEPRESSION SCREENING VA-IRAQ & AFGHAN POST-DEPLOY SCREEN VA-MHV CERVICAL CANCER SCREEN VA-MHV COLORECTAL CANCER SCREEN VA-MHV DIABETES FOOT EXAM VA-MHV DIABETES RETINAL EXAM

VA-MHV HYPERTENSION VA-MHV INFLUENZA VACCINE VA-MHV MAMMOGRAM SCREENING VA-MHV PNEUMOVAX VA-OEF/OIF MONITOR REPORTING VA-PTSD SCREENING VA-TBI SCREENING

The post-init will install all the components of these Exchange file entries on your system. After the installation has finished, if you discover that any of these components weren't installed correctly, you can use the Reminder Exchange option on the Reminders Manager Menu to install them.

8. Deletion of init routines

After everything has been successfully installed you may delete the following init routines: PXRMP11I, PXRMP11E

Setup and Maintenance

Overview

PXRM*2.0*11 includes the following additions and modifications:

- New and modified national reminder definitions and dialogs
 - New reminders to help clerks reconcile HEC data and self-reported VA-IRAQ/AFGHAN SERVICE health factor
 - New national reminder for the monitor
- New and modified computed findings
- Modified Location List exclusion functionality
- Modified Reminder Test output
- Fixes to My HealtheVet reminders function finding pointers

New Reminder Definitions

VA-BL OEF/OIF EMBEDDED FRAGMENTS	New reminder for branching logic in the new 4 th question in the last section of the OEF/OIF reminder. Used for branching logic to indicate in the reminder dialog that an item is already completed or not needed.
VA-CV ELIG W/HF FOR NO SERVICE VA-CV INELG W/HF FOR SERVICE	These two reminders are designed to aid sites in identifying patients whose combat veteran eligibility data and data collected from the OEF/OIF Screening reminder are not in agreement. Reminder patient list functionality can be used to identify patients needing reconciliation
VA-OEF/OIF MONITOR REPORTING	New reminder for reporting this monitor. Excludes non-combat vets from the denominator; requires two visits after discharge separation.

Modified national reminder definitions

- VA-ALCOHOL AUDIT-C POSITIVE F/U EVAL
- VA-ALCOHOL USE SCREEN (AUDIT-C)
- VA-BL OEF/OIF FEVER
- VA-BL OEF/OIF GI SX
- VA-BL OEF/OIF SKIN SX

- VA-DEPRESSION SCREENING
- VA-IRAQ & AFGHANISTAN POST DEPLOYMENT SCREEN
- VA-MHV COLORECTAL CANCER SCREEN
- VA-MHV HYPERTENSION
- VA-MHV DIABETES RETINAL EXAM
- VA-PTSD SCREENING
- VA-TBI SCREENING

New and modified Service-related Computed Findings

- VA-AGENT ORANGE EXPOSURE
- VA-COMBAT SERVICE
- VA-COMBAT VET ELIGIBILITY
- VA-LAST SERVICE SEPARATION DATE (modified)
- VA-OEF SERVICE
- VA-OIF SERVICE
- VA-POW
- VA-PURPLE HEART
- VA-SERVICE BRANCH
- VA-UNKNOWN OEF/OIF SERVICE
- VA-VETERAN (modified)

Modified Location List exclusion functionality

- Exclusion location list can be defined once
- Any location lists can reference the pre-defined Exclusion location list

Modified My HealtheVet Reminders

- Updated URLs: VA-MHV CERVICAL CANCER SCREEN VA-MHV DIABETES FOOT EXAM VA-MHV DIABETES RETINAL EXAM VA-MHV HYPERTENSION VA-MHV INFLUENZA VACCINE VA-MHV MAMMOGRAM SCREENING
- VA-MHV COLORECTAL CANCER SCREEN Added reminder term for barium enema
- VA-MHV DIABETES RETINAL EXAM Updated to use only prior exam as indication of frequency of exams instead of HbA1c level and use of insulin.

- VA-MHV HYPERTENSION Converted findings to national reminder terms.
- VA-MHV INFLUENZA VACCINE Updated the date to August 2008
- VA-MHV PNEUMOVAX The reminder term VA-MHV HIGH RISK FOR PNEUMONIA was made national.

Dialog Updates

VA-DEPRESSION SCREEN	PHQ-9 not required; MH test in dialog set to '2' – required if opened; informational text for dialog edited to remove PHQ-9
VA-PTSD SCREENING	 MH test in dialog set to '2' – required if opened Added Acute Illness to dialog
VA-IRAQ & AFGHANISTAN POST DEPLOYMENT SCREEN	 Added cognitive impairment Added done elsewhere Doesn't ask the "Did you serve? " question if combat vet eligible Substituted Q4D – asks about embedded fragments

Setup Steps

1. Repeat the pre-installation reminder inquiry and compare the pre- and post-installation output.

The following reminder definitions have been changed. These may require remapping. The ones in red may require remapping as well as a check of the dialog for any orders that need to be re-linked.

Reminder	Change
VA-ALCOHOL AUDIT-C POSITIVE	Fixed branching logic
F/U EVAL	
VA-ALCOHOL USE SCREEN	Changes to dialog.
(AUDIT-C)	Removed date from term and moved to
	the health factor for "No alcohol in the
	past 1 year"
VA-BL OEF/OIF FEVER	Fixed logic
VA-BL OEF/OIF GI SX	Fixed logic
VA-BL OEF/OIF SKIN SX	Fixed logic
VA-DEPRESSION SCREENING	Changes to terms and dialog;
	verify/report any changes
VA-IRAQ & AFGHANISTAN POST	Multiple changes:
DEPLOYMENT SCREEN	Added Combat Vet to logic
	• Excluded those who did not serve
	from denominator
	• Added cognitive impairment to
	exclusions
	• Added done elsewhere to resolutions
	• Updated dates of health factors for
	depression/PTSD to 10/1/08
VA-MHV COLORECTAL CANCER	VA-MHV BARIUM ENEMA added
SCREEN	
VA-MHV HYPERTENSION	Uses reminder terms
VA-MHV DIABETES RETINAL	Frequency depends on no retinopathy
EXAM	only
VA-PTSD SCREENING	Added veteran status
	Added Acute Illness to the dialog
VA-TBI SCREENING	Two new Finding Items were added: VA-
	COMBAT VET and (FI10) VA-
	COGNITIVE IMPAIRMENT, and a
	change in the Cohort and Resolution
	logic

2. Map Local Findings to National Terms


The following terms have changed and may need to be remapped, if previously mapped locally:

The ones in bold will almost definitely need remapping.

Term	Mapping
VA-COGNITIVE IMPAIRMENT	
VA-DEPRESSION SCREEN OEF/OIF	
VA-LDL	All local lab tests for LDL
VA-MHV BARIUM ENEMA	Add outside BE
VA-MHV COLONOSCOPY	Add outside colonoscopies or other HFs
VA-MHV DIABETIC FOOT EXAM	All 3 components completed
COMPLETE	
VA-MHV DIABETIC FOOT EXAM	Local exams and/or HFs
INSPECTION	
VA-MHV DIABETIC FOOT EXAM	Local exams and/or HFs
MONOFILAMENT	
VA-MHV DIABETIC FOOT EXAM	Local exams and/or HFs
PULSES	
VA-MHV HBA1C <8.0	Local labs with condition I +V<8and
	outside HFs
VA-MHV HBA1C >7.9	
VA-MHV HBA1C ALL RESULTS	
VA-MHV HDL – lab test and outside lab	
test	
VA-MHV INFLUENZA IMMUNIZATION	Two national entries are included
VA-MHV LIPID PANEL ORDERABLE	
ITEMS	
VA-MHV OCCULT BLOOD PANEL	Lab test for panel
VA-MHV OCCULT BLOOD SINGLE	Lab test for individual tests in the occult
TEST	blood panel e.g. occult blood #1, occult
	blood #2, occult blood #3
VA-MHV OUTSIDE LDL	Local HFs for outside LDLs, national HFs
	are included
VA-MHV OUTSIDE LDL>99	Outside HFs, national HFs are included
VA-MHV PNEUMOCOCCAL	National entry is included
IMMUNIZATION	
VA-MHV SIGMOIDOSCOPY	Outside sigmoidoscopies or other HFs
VA-MHV TOTAL CHOLESTEROL	Local labs and outside HFs
VA-MHV TRIGLYCERIDES	Local labs and outside HFs
VA-PTSD SCREEN	
VA-WH HYSTERECTOMY W/CERVIX	Includes taxonomy and national HF
REMOVED	
VA-WH MAMMOGRAM ORDER	Orderable item and HF included, add local

	HFs or orders
VA-WH MAMMOGRAM SCREEN	Two national HFs included
DEFER	
VA-WH MAMMOGRAM SCREEN DONE	Taxonomies and national HF included, add
	any local HFs
VA-WH MAMMOGRAM SCREEN NOT	Includes terminal cancer and limited life
INDICATED	expectancy, add local HFs for limited life
	expectancy
VA-WH PAP SMEAR DONE	Includes taxonomy and national HF
VA-WH PAP SMEAR OBTAINED	Includes national taxonomy
VA-WH PAP SMEAR ORDER HEALTH	Already contains four national HFs
FACTOR	
VA-WH PAP SMEAR SCREEN NOT	Limited life expectancy and terminal
INDICATED	cancer taxonomy, add local HFs for limited
	life expectancy

Example: Mapping a Local Health Factor Finding to the National Reminder Term

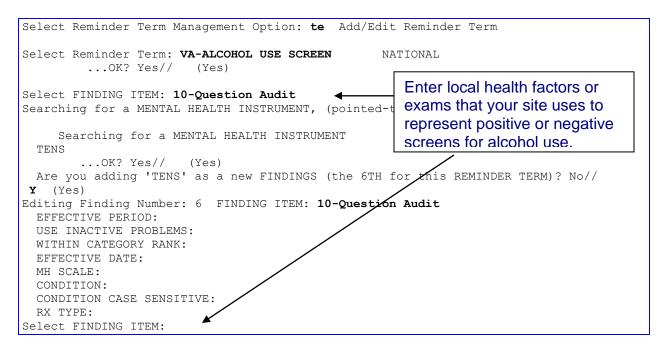

3. Edit Reminder Dialogs

Several Dialog Elements may need to be remapped after patch 11 is installed. For example, the reminder dialog element **VA-PDIQ POLYTRAUMA CONSULT**, which is used to order a TBI consult, is overwritten in the install. You will need to re-enter the finding item for this element using the menu or quick order that you use to order a TBI consult. **VA-GP ALC SCREEN POS REFERRAL** and **TBI OI CONSULT FOR KNOWN TBI** are also overwritten.

The objects in the element TXT ALC LABS & VS, in the VA-GP ALC AUDIT >7 (added after Patch 6 was installed) will also need to be remapped.

Check the dialogs for any orders that need to be re-linked for the following:

- VA-ALCOHOL AUDIT-C POSITIVE F/U EVAL
- VA-ALCOHOL USE SCREEN (AUDIT-C)
- VA-BL DEPRESSION SCREENING
- VA-IRAQ & AFGHAN POST-DEPLOY SCREEN
- VA-IRAQ/AFGHAN SERVICE
- VA-MHV HIGH RISK FOR PNEUMONIA
- VA-MHV INFLUENZA IMMUNIZATION
- VA-PTSD SCREENING

Steps to add or edit dialog elements:

a. Select Dialog management (DM) from the Reminders Manager Menu, then select Dialog (DI):

Select Reminder Managers Menu Option: DM Reminder Dialog Management DP Dialog Parameters DI Reminder Dialogs Select Reminder Dialog Management Option: DI Reminder Dialogs
Dialog ListAug 30, 2007@09:55:38Page: 1 of 8REMINDER VIEW (ALL REMINDERS BY NAME)
Item Reminder NameLinked Dialog Name & Dialog Status1AGETESTAGETEST7ALCOHOL USE SCREENALCOHOL USE SCREEN DIALOG3HIV HEPATITIS A SEROLOGIC TESTINGHIV HEPATITIS A SEROLOGIC T4HYPERTENSIONHYPERTENSION GROUPDisabled5Hypertension Screen (VHACHS)Hypertension Screen (local)
+ Enter ?? for more actions >>> AR All reminders LR Linked Reminders QU Quit CV Change View RN Name/Print Name Select Item: Next Screen//

b. Use the Change View (CV) action to see the Dialog Elements view.

Select Item: Next S	creen// CV				
Select one of	the following:				
D	Reminder Dialogs				
E	Dialog Elements				
F	Forced Values				
G	Dialog Groups				
P	Additional Prompts				
R	Reminders				
RG	Result Group (Mental Healt	h)			
RE	Result Element (Mental Hea	lth)			
TYPE OF VIEW: R// E	Dialog Elements				
Dialog List	Sep 04, 2007@10:	36	Page:	1 of	120
DIALOG VIEW (DIALOG	ELEMENTS)				
Item Dialog Name		Dialog type	!	S	tatus

1	04 AUDIT-C		Dialog	Ele	ment	
2	ABILITY FAIR (SLC)		Dialog			
	B ABILITY FAIR 1		Dialog			
4	ABILITY GOOD		Dialog	Ele	ment	
5	ABILITY POOR		Dialog	Ele	ment	
6	ADD'L INFO		Dialog	Ele	ment	
7	ADDITIONAL COMMENTS	(PXRM COMMENT)	Dialog	Ele	ment	
8	ADDITIONAL COMMENTS	(2 LINES WP)	Dialog	Ele	ment	
+	+ Next Screen	- Prev Screen	?? More Ac	tion	S	>>>
AD	Add	CV Change Vi	ew	INQ	Inquiry/Print	
CO	Copy Dialog	PT List/Prin	t All 🛛	QU	Quit	
Sele	ect Item: Next Screen/	/				

c. Use the Search List (SL) action to get to the dialog element you want the edit.

```
Select Item: Next Screen// SL SL
Search for: VA-HF DEP
...searching for 'VA-HF DEP'.....
```

e. Make edits, as needed, by entering the dialog element number.

Find Next 'VA-HF DEP'? Yes// NO	
Dialog List Sep 04, 2007@10	:25:40 Page: 77 of 120
DIALOG VIEW (DIALOG ELEMENTS)	
+Item Dialog Name	Dialog type Status
1232 VA-HF DEP 2 QUESTION NEG	Dialog Element
1233 VA-HF DEP 2 QUESTION POS	Dialog Element
1234 VA-HF DEP ACUTE ILLNESS	Dialog Element
1232 VA-HF DEP 2 QUESTION NEG 1233 VA-HF DEP 2 QUESTION POS 1234 VA-HF DEP ACUTE ILLNESS 1235 VA-HF DEP ACUTE MED CONDITION	Dialog Element
1236 VA-HF DEP CHRONIC MED CONDITION	Dialog Element
1237 VA-HF DEP REFERRAL TO MHC	Dialog Element
1238 VA-HE DEP TO BE MANAGED IN PC	Dialog Element
1239 VA-HF DEP TO BE MANAGED IN PC (2)	Dialog Element
1240 VA-HF ETOH ADVISE NEGOTIATED LEVEL	Dialog Element
1241 VA-HF ETOH ADVISE PT RESPONSE NO CHANGE	Dialog Element
1242 VA-HF ETOH ADVISE PT RESPONSE WILL PLAN	Dialog Element
1243 VA-HF ETOH DSM IV CONTINUED USE	Dialog Element
1244 VA-HF ETOH DSM IV DRINKS LARGE AMT	Dialog Element
1245 VA-HF ETOH DSM IV IMPAIRED CONTROL	Dialog Element
1246 VA-HF ETOH DSM IV NEGLECT OTHER RESPONS	Dialog Element
1247 VA-HF ETOH DSM IV TOLERANCE	Dialog Element
+ + Next Screen - Prev Screen ??	More Actions >>>
AD Add CV Change View CO Copy Dialog PT List/Print All	INQ Inquiry/Print
CO Copy Dialog PT List/Print All	QU Quit
Select Item: Next Screen// 1234	
Dialog Name: VA-HF DEP ACUTE ILLNESS	
Current dialog element/group name: VA-HF DEP ACUTE ILLNESS	
Used by: VA-GP DEP UNABLE TO SCREEN (Dialog Group)	
VA-GP PTSD UNABLE TO SCREEN (Dialog Group)	
FINDING ITEM: UNABLE TO SCREEN - ACUTE ILLNESS//	

Mapping a consult quick order for the VA-TBI SCREENING reminder dialog:

1. Type the name of the dialog element, VA-PDIQ POLYTRAUMA CONSULT (or TBI OI CONSULT FOR KNOWN TBI?)

Select Item: Next Screen// SL SL Search for: VA-TBI OI CONSULT FOR KNOWN TBI ...searching for 'VA-TBI OI CONSULT FOR KNOWN TBI'.....

- 2. At the "Find Next" prompt, type No.
- 3. Type the item number next to the element.
- 4. At the finding Item Prompt, enter Q. plus the name of the site-specific Quick Order or order menu.
- 5. Press Enter at the Additional Finding prompt.

4. Verify that the dialogs function properly

Test the reminder dialogs in CPRS, using either the exported dialogs or your locally created dialogs. Using point-and-click reminder resolution processing through CPRS GUI, verify the following:


- Correct Progress Note text is posted
- Finding Item gets sent to PCE
- Reminder is satisfied

Check the Clinical Maintenance component display in CPRS after testing dialogs to ensure that all the activities are tested are reflected in the clinical maintenance display.

Steps to test dialogs:

- 1. On the cover sheet, click on the Reminders icon.
- 2. Click on reminders in the Reminders box to see details of a reminder
- 3. Open the Notes tab and select New Note. Enter a title and visit information.
- 4. Open the Reminders drawer and review the contents.
- 5. Locate the reminder dialog and open it.
- 6. In the dialog box, check relevant actions.
- 7. Finish the reminder processing.
- 8. Review the text added to the note to assure its correctness.

- 9. Ensure that the reminder can be satisfied by the individual finding items that were mapped to the reminder terms.
- 10. Check the Clinical Maintenance component display in Health Summaries and CPRS after the reminder dialog is complete.

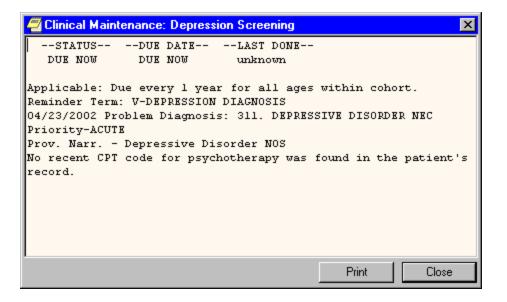

11. Repeat steps 5-10 for the Depression Assessment, Iraq & Afghan Post-Deployment Screen, and VA-PTSD SCREENING dialogs.

NOTE: Remember to "refresh" the screen after completing a dialog, if you want to see the updated status immediately. This is especially critical if you're doing the screening reminder and want to see the follow-up reminder become due because of positive results.

Appendix A: Installation Example

Select Installation Option: INSTAll Package(s) Select INSTALL NAME: PXRM*2.0*11 Loaded from Distribution 10/7/08@14:01:56 => PXRM*2.0*11 ;Created on Sept 26, 2008@11:08:46 It consisted of the following Install(s): PXRM*2.0*11 Install Started for PXRM*2.0*11 : Oct 07, 2008@05:47:49 Build Distribution Date: Sep 26, 2008 Oct 07, 2008@05:47:49 Installing Routines: Running Pre-Install Routine: PRE^PXRMP11I DISABLE options. DISABLE protocols. Removing old data dictionaries. Deleting data dictionary for file # 810.9 Installing Data Dictionaries: Oct 07, 2008@05:47:50 Installing Data: Oct 07, 2008@05:47:52 Running Post-Install Routine: POST^PXRMP11I Rebuilding Function Finding internal data structures. ENABLE options. ENABLE protocols. There are 25 Reminder Exchange entries to be installed. Installing Reminder Exchange entry PATCH 11 LOCATION LIST Installing Reminder Exchange entry VA-ALCOHOL AUDIT-C POSITIVE F/U EVAL Installing Reminder Exchange entry VA-ALCOHOL USE SCREEN (AUDIT-C) Installing Reminder Exchange entry VA-BL OEF/OIF EMBEDDED FRAGMENTS Installing Reminder Exchange entry VA-BL OEF/OIF FEVER Installing Reminder Exchange entry VA-BL OEF/OIF GI SX Installing Reminder Exchange entry VA-BL OEF/OIF SKIN SX

Installing Reminder Exchange entry VA-CV ELIG W/HF FOR NO SERVICE Installing Reminder Exchange entry VA-CV INELG W/HF FOR SERVICE Installing Reminder Exchange entry VA-DEPRESSION SCREENING Installing Reminder Exchange entry VA-IRAQ & AFGHAN POST-DEPLOY SCREEN Installing Reminder Exchange entry VA-MHV BMI > 25.0 Installing Reminder Exchange entry VA-MHV CERVICAL CANCER SCREEN Installing Reminder Exchange entry VA-MHV COLORECTAL CANCER SCREEN Installing Reminder Exchange entry VA-MHV DIABETES FOOT EXAM Installing Reminder Exchange entry VA-MHV DIABETES RETINAL EXAM Installing Reminder Exchange entry VA-MHV HYPERTENSION Installing Reminder Exchange entry VA-MHV INFLUENZA VACCINE Installing Reminder Exchange entry VA-MHV LDL CONTROL Installing Reminder Exchange entry VA-MHV MAMMOGRAM SCREENING Finding RP.MAMMOGRAM UNILAT does not exist; input a replacement or ^ to quit the install. Select RAD/NUC MED PROCEDURES NAME: MAMMOGRAPHY SCREENING (RAD Detailed) CPT:77057 Finding RP.MAMMOGRAM UNILAT does not exist; input a replacement or ^ to quit the install. Select RAD/NUC MED PROCEDURES NAME: MAMMOGRAPHY SCREENING (RAD Detailed) CPT:77057 Installing Reminder Exchange entry VA-MHV PNEUMOVAX Installing Reminder Exchange entry VA-OEF/OIF MONITOR REPORTING Installing Reminder Exchange entry VA-PTSD SCREENING Installing Reminder Exchange entry VA-TBI SCREENING Installing Reminder Exchange entry PATCH 11 ITEMS Removing PATCH 11 ITEMS transport reminder. Removing PATCH 11 DIALOG transport dialog. Updating Routine file ... Updating KIDS files... PXRM*2.0*11 Installed. Oct 07, 2008@05:51:35 Install Message sent #xxxxxxxxxxx Install Completed

Appendix B: Setting up a TIU/HS Object for the AUDIT-C Test

Before using the VA-Alcohol Use Screen reminder, a TIU/HS Object needs to be created to look for the last occurrence of an AUDIT-C test done on a patient. The TIU/HS Object must be named AUDIT-C for the object to work in the reminder dialog.

NOTE: For some reason, Health summary doesn't allow occurrence count or time period for MH tests, so you may need to enter these first, as follows:

The components MHAL and MHAS should have values of 'Y' for both fields #2 - TIME LIMITS APPLICABLE and #4 - MAXIMUM OCCURRENCES APPLICABLE.

You can edit both those fields locally by doing the following:

```
Select OPTION: ENTER OR EDIT FILE ENTRIES

INPUT TO WHAT FILE: HEALTH SUMMARY COMPONENT//

EDIT WHICH FIELD: ALL// MAXIMUM OCCURRENCES APPLICABLE

THEN EDIT FIELD:

Select HEALTH SUMMARY COMPONENT NAME: MH

1 MHA

2 MHA Administration List

3 MHA Score

4 MHV REMINDERS DETAIL DISPLAY

5 MHV REMINDERS SUMMARY DISPLAY

CHOOSE 1-5: 3 MHA Score

MAXIMUM OCCURRENCES APPLICABLE: Y yes
```

You can do this for both MHAL and MHAS components.

Create TIU/Health Summary Objects

To create the TIU/Object, use the menu option "Create TIU/Health Summary Objects" under the Manager Document Definition Menu on the TIU Maintenance Menu. An example of how to set this up is provided below. The bolded parts must be exact; the non-bolded parts are site preferences.

```
Select OPTION NAME: TIU IRM MAINTENANCE MENU
                                                TIU Maintenance Menu
                                                                         menu
   1
         TIU Parameters Menu ...
   2
        Document Definitions (Manager) ...
   3
        User Class Management ...
   4
        TIU Template Mgmt Functions ...
   5
        TIU Alert Tools
   7
         TIUHL7 Message Manager
         Title Mapping Utilities ...
   6
         Active Title Cleanup Report
```

Select TIU Maintenance Menu Option: 2 Document Definitions (Manager) --- Manager Document Definition Menu ---1 Edit Document Definitions Sort Document Definitions 2 3 Create Document Definitions Create Objects 4 5 Create TIU/Health Summary Objects Select Document Definitions (Manager) Option: 5 Create TIU/Health Summary Objects TIU Health Summary Object Jul 30, 2008@12:13:19 Page: 1 of 1 TIU Object Name Health Summary Type BRADEN SCALE 30D VA-BRADEN SCALE 30D 1 2 PRESSURE ULCER VA-PRESSURE ULCER PU INTERVENTIONS VA-PU INTERVENTIONS 3 4 TIU TPBN FUTURE APPTS TIU TPBN FUTURE APPTS Create New TIU Object Find Detailed Display/Edit TIU Object Detailed Display/Edit HS Object Quit Select Action: Quit// CRE Create New TIU Object --- Create TIU/Health Summary Object ---Enter a New TIU OBJECT NAME: AUDIT-C 🗲 this must be exact. Object Name: AUDIT-C Is this correct? YES// <Enter> Use a pre-existing Health Summary Object? NO// <Enter> Checking AUDIT-C (TIU) with Health Summary... Creating Health Summary Object 'AUDIT-C (TIU)' Select Health Summary Type: OB AUDC Are you adding 'OB AUDC' as a new HEALTH SUMMARY TYPE (the 181th)? No// YES NAME: OB AUDC// <Enter> TITLE: OB AUDC SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: <Enter> SUPPRESS SENSITIVE PRINT DATA: <Enter> LOCK: **<Enter>** OWNER: // <Enter> Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO Select COMPONENT: MHAS MHA Score MHAS SUMMARY ORDER: 5// 5 OCCURRENCE LIMIT: 1 HEADER NAME: MHA Score// <Enter> No selection items chosen. Select new items one at a time in the sequence you want them displayed. You may select up to 99 items.

Select SELECTION ITEM: AUDC Searching for a MH TESTS AND SURVEYS, (pointed-to by SELECTION ITEM) Searching for a MH TESTS AND SURVEYS AUDC 1 2 AUDCR CHOOSE 1-2: 1 AUDC Select SELECTION ITEM: <Enter> Select COMPONENT: <Enter> Do you wish to review the Summary Type structure before continuing? NO// <Enter> Please hold on while I resequence the summary order. Do you want to overwrite the TIME LIMITS in the Health Summary Type 'OB AUDC'? N//<Enter> O Print standard Health Summary Header with the Object? $\,$ N//<Enter> O Partial Header: Print Report Date? N// O Print Confidentiality Banner? N// O Print Report Header? N// O Print the standard Component Header? Y// ES Use report time/occurence limits? N// O Underline Component Header? N// O Add a Blank Line after the Component Header? N//<Enter> O Print the date a patient was deceased? N//<Enter> 0 Print a LABEL before the Health Summary Object? N//<Enter> O Suppress Components without Data? N// <Enter>O OBJECT DESCRIPTION: No existing text Edit? NO// <Enter> Create a TIU Object named: AUDIT-C Ok? YES// **<Enter>** TIU Object created successfully. Enter RETURN to continue...

Appendix C: Computed Finding Descriptions

The following are new or modified computed findings distributed by PXRM*2*11:

NAME: VA-AGENT ORANGE EXPOSURE ROUTINE: PXRMMSER ENTRY POINT: AORANGE PRINT NAME: VA-Agent Orange Exposure TYPE: MULTIPLE DESCRIPTION: This computed finding will be true if the patient has an agent orange exposure registration date in the date range specified by Beginning Date/Time and Ending Date/Time. Subscripts that can be used in a Condition are: "LOCATION" The default value is "LOCATION". CLASS: NATIONAL NAME: VA-COMBAT SERVICE ROUTINE: PXRMMSER ENTRY POINT: COMBAT PRINT NAME: VA-Combat Service TYPE: MULTIPLE DESCRIPTION: This computed finding will be true if combat service is found in the date range specified by Beginning Date/Time and Ending Date/Time. Subscripts that can be used in a Condition are: "LOCATION" The default value is "LOCATION". CLASS: NATIONAL NAME: VA-COMBAT VET ELIGIBILITY ROUTINE: PXRMMSER ENTRY POINT: CVELIG PRINT NAME: VA-Combat Vet Eligibility TYPE: MULTIPLE DESCRIPTION: This computed finding will be true if the patient qualifies as a combat vet. The possible values that can be used in a Condition are: "Eligible", "Expired", and "Not eligible". An example is: I V="Not eligible" If the patient was eligible on the evaluation date then the Value for use in a Condition will be "Eligible". If the patient is not eligible and the evaluation date is greater than the eligibility end date then the Value is "Expired", otherwise the value is "Not eligible". CLASS: NATIONAL

(modified) NAME: VA-LAST SERVICE SEPARATION DATE ROUTINE: PXRMMSER PRINT NAME: Veteran Last Separation ENTRY POINT: DISCHDT Date TYPE: SINGLE DESCRIPTION: This is a single occurrence computed finding that returns the most recent Service Separation Date for the veteran as a FileMan date. It can be used as the value in the Condition. For example: I V>3030101. The computed finding will be true if a separation date is found; otherwise it will be false. If the computed finding is true, the date of the finding will be the most recent Service Separation Date. CLASS: NATIONAL SPONSOR: Office of Quality & Performance NAME: VA-OEF SERVICE ROUTINE: PXRMMSER ENTRY POINT: OEF PRINT NAME: VA-OEF Service TYPE: MULTIPLE DESCRIPTION: This multi-occurrence computed finding will search for periods of OEF service in the date range specified by Beginning Date/Time and Ending Date/Time. Subscripts that can be used in a Condition are: "LOCATION" The default value is "LOCATION". CLASS: NATIONAL NAME: VA-OIF SERVICE ROUTINE: PXRMMSER PRINT NAME: VA-OIF Service ENTRY POINT: OIF TYPE: MULTIPLE DESCRIPTION: This multi-occurrence computed finding will search for periods of OIF service in the date range specified by Beginning Date/Time and Ending Date/Time. Subscripts that can be used in a Condition are: "LOCATION" The default value is "LOCATION". CLASS: NATIONAL NAME: VA-POW ROUTINE: PXRMMSER ENTRY POINT: POW PRINT NAME: VA-POW TYPE: MULTIPLE DESCRIPTION: This computed finding will be true if the patient was a POW in the date range specified by Beginning Date/Time and Ending Date/Time. Subscripts that can be used in a Condition are: "LOCATION" The default value is "LOCATION". CLASS: NATIONAL

NAME: VA-PURPLE HEART

ENTRY POINT: PHEART PRINT NAME: VA-Purple Heart TYPE: SINGLE DESCRIPTION: This computed finding will be true if the patient is a Purple Heart recipient.

There is no value for use in a Condition.

CLASS: NATIONAL

NAME: VA-SERVICE BRANCH

ROUTINE: PXRMMSER PRINT NAME: VA-Service Branch

ROUTINE: PXRMMSER

ENTRY POINT: SBRANCH TYPE: MULTIPLE

DESCRIPTION: This computed finding will return service branch information for a maximum of three service periods in the date range specified by Beginning Date/Time and Ending Date/Time.

Subscripts that can be used in a Condition are: "BRANCH" "DISCHARGE TYPE" "ENTRY DATE" "SEPARATION DATE" "SERVICE COMPONENT"

The default value is "DISCHARGE TYPE".

CLASS: NATIONAL

NAME: VA-UNKNOWN OEF/OIF SERVICE

ENTRY POINT: UNKOEIF PRINT NAME: VA-Unknown OEF/OIF Service TYPE: MULTIPLE DESCRIPTION: This multi-occurrence computed finding will search for periods of OEF/OIF service with an unknown location in the date range specified by Beginning Date/Time and Ending Date/Time.

Subscripts that can be used in a Condition are: "LOCATION"

The default value is "LOCATION".

CLASS: NATIONAL

(modified) NAME: VA-VETERANROUTINE: PXRMMSERENTRY POINT: VETERANPRINT NAME: Patient is a VeteranTYPE: SINGLEDESCRIPTION: This single occurrence computed finding is true if thepatient is a veteran and false otherwise.

There are no values that can be used in a Condition.

CLASS: NATIONAL

ROUTINE: PXRMMSER