

Department of Veterans Affairs
Mental Health – Suicide Prevention
Suicide Prevention Package



March 2020
Version 1.9

**Deployment, Installation, Back-Out, and Rollback Guide for
YS*5.01*150**

Submitted as CLIN 0001AX
Contract VA118-16-D-1007, Task Order VA11817F10070006

Revision History

Date	Version	Change	Project Manager	Technical Writer
March 2020	1.9	Updated version for Patch 150	REDACTED	Booz Allen SPP Project Team
June 2019	1.8	Updated version for Patch 145	REDACTED	Booz Allen SPP Project Team
May 2019	1.7	Updated version for Patch 147	REDACTED	Booz Allen SPP Project Team
February	1.6	Updated version for Patches 138, 139	REDACTED	Booz Allen SPP Project Team
September 2018	1.5	Updated version for Patch 136	REDACTED	Booz Allen SPP Project Team
August 2018	1.4	Updated version numbers for Patch 134	REDACTED	Booz Allen SPP Project Team
March 2018	1.3	Updated Title and Footers;	REDACTED	Booz Allen SPP Project Team
August 2014	1.2	Updated Title page Updated Revision History, p. i-ii Updated Table of Contents, pp. iii-iv Changed ICD9 reference to ICD	REDACTED	REDACTED
29 February 2012	1.1	Added table "MHA3 Patch YS*5.01*129	REDACTED	REDACTED

Artifact Rationale

This document describes the Deployment, Installation, Back-out, and Rollback Guide (DIBO&RG) for new products going into the Veterans Affairs (VA) Enterprise. The plan includes information about system support, issue tracking, escalation processes, and roles and responsibilities involved in all those activities. Its purpose is to provide clients, stakeholders, and support personnel with a smooth transition to the new product or software, and should be

structured appropriately, to reflect particulars of these procedures at a single or at multiple locations.

Per the Veteran-focused Integrated Process (VIP) Guide, the DIBO&RG is required to be completed prior to Critical Decision Point #2 (CD #2), with the expectation that it will be updated throughout the lifecycle of the project for each build, as needed.

Table of Contents

1	Introduction	6
1.1	Purpose	6
1.2	Dependencies	6
1.3	Constraints.....	6
2	Roles and Responsibilities.....	6
3	Deployment	7
3.1	Timeline.....	7
3.2	Site Readiness Assessment.....	7
3.2.1	Deployment Topology (Targeted Architecture).....	7
3.2.2	Site Information (Locations, Deployment Recipients).....	7
3.2.3	Site Preparation	7
3.3	Resources	8
3.3.1	Facility Specifics (optional)	8
3.3.2	Hardware	8
3.3.3	Software.....	8
3.3.4	Communications.....	8
4	Installation	8
4.1	Pre-installation and System Requirements.....	8
4.2	Platform Installation and Preparation	8
4.3	Access Requirements and Skills Needed for the Installation.....	9
4.4	Installation Procedure	9
4.5	Post Install for Instruments	10
4.5.1	Instructions using the MHA3 Utilities Option	11
4.5.2	FileMan Instructions	12
4.6	Installation Verification Procedure	13
4.7	System Configuration	15
4.8	Database Tuning.....	15
5	Back-Out Procedure	15
5.1	Back-Out Strategy	15
5.2	Back-Out Considerations.....	15
5.3	Back-Out Criteria	15
5.4	Back-Out Risks	15
5.5	Authority for Back-Out	15
5.6	Back-Out Procedure	15
5.7	Back-out Verification Procedure	16

6	Rollback Procedure	16
6.1	Rollback Considerations.....	16
6.2	Rollback Criteria	16
6.3	Rollback Risks	16
6.4	Authority for Rollback	17
6.5	Rollback Procedure	17
6.6	Rollback Verification Procedure	17
	Appendix A - Acronyms	17

1 Introduction

This document describes how to deploy and install the patch YS*5.01*150 of the Mental Health package, as well as how to back-out the product and rollback to a previous version or data set. This document is a companion to the project charter and management plan for this effort in this document.

1.1 Purpose

The purpose of this plan is to provide a single, common document that describes how, when, where, and to whom Mental Health patch YS*5.01*150 will be deployed and installed, as well as how it is to be backed out and rolled back, if necessary. The plan also identifies resources, communications plan, and rollout schedule. Specific instructions for installation, back-out, and rollback are included in this document.

1.2 Dependencies

It is assumed that this patch is being installed into a fully patched Veterans Health Information System and Technology Architecture (VistA) system. In particular patch YS*5.01*123 must be installed prior to the installation of this patch.

1.3 Constraints

There are no constraints beyond the installation into an up-to-date VistA system.

2 Roles and Responsibilities

The following describes the roles and responsibilities associated with the testing and release of YS*5.01*150. This is a VistA patch that will be deployed via the normal Mailman route.

Table 1: Deployment, Installation, Back-out, and Rollback Roles and Responsibilities

Team	Phase / Role	Tasks	Project Phase (See Schedule)
Project Manager	Deployment	Determine and document the roles and responsibilities of those involved in the deployment.	Design
Software Quality Assurance (SQA), Test Sites	Deployment	Test for operational readiness	Test
Project Manager, Release Manager	Deployment	Execute deployment	Release
Individual VistA Sites	Installation	Plan and schedule installation	Release

Team	Phase / Role	Tasks	Project Phase (See Schedule)
Release Manager	Back-out	Confirm availability of back-out instructions and back-out strategy (what are the criteria that trigger a back-out)	Release
Sustainment Team	Post Deployment	Hardware, Software and System Support	Sustain

3 Deployment

The deployment is planned as a simultaneous (National Release) rollout. Once approval has been given to nationally release, YS*5.01*150 will be available for installation and deployment at all sites.

Scheduling of test installs, testing and production deployment will be at the site's discretion. It is anticipated there will be a 30-day compliance period.

3.1 Timeline

The deployment and installation are scheduled to run from June 2019 through July 2019 as depicted in the Master Deployment Schedule in the Suicide Prevention Program (SPP) Project Management Plan.

3.2 Site Readiness Assessment

This section discusses the locations that will receive the YS*5.01*150 deployment.

3.2.1 Deployment Topology (Targeted Architecture)

MHA Update (YS*5.01*150) will be deployed to each VistA instance. This includes local sites as well as regional data centers. The executable and associated files will also be deployed to client workstations.

3.2.2 Site Information (Locations, Deployment Recipients)

The initial deployment will be to Initial Operating Capability (IOC) sites for verification of functionality. Once testing is completed and approval is given for national release, MHA Update (YS*5.01*150) will be deployed to all VistA systems.

The Production IOC testing sites are:

- REDACTED

3.2.3 Site Preparation

Other than a fully patched VistA system, there is no other preparation required.

3.3 Resources

3.3.1 Facility Specifics (optional)

N/A

3.3.2 Hardware

N/A

3.3.3 Software

N/A

3.3.4 Communications

When MHA Update (YS*5.01*150) is released, the released-patch notification will be sent from the National Patch Module to all personnel who have subscribed to notifications for the Mental Health package.

4 Installation

4.1 Pre-installation and System Requirements

Mental Health patch YS*5.01*123 and YS*5.01*147 must be installed prior to the installation of YS*5.01*150. The environment checks for YS*5.01*150 will verify this.

4.2 Platform Installation and Preparation

See previous section for required Mental Health patches.

This patch installs new instruments that will consume about 2MB of disk space in the ^XTMP global during installation. The growth of the ^YTT global will be about 1MB.

This patch can be loaded with users in the system but it is recommended that it be installed when user activity is low. Installation time will be less than 5 minutes.

It is recommended that you use the “Backup a Transport Global” option that is referenced in the installation instructions. This will be useful should it be decided to back out the installation.

To ensure the integrity of the transport global, use the “Verify Checksums in Transport Global” to compare the checksums with the list that follows:

The checksums below are new checksums and can be checked with CHECK1^XTSUMBLD.

ROUTINE NAME: YS150PST

<u>BEFORE</u> :	N/A	<u>AFTER</u> :	992
	150		<u>ROUTINE</u>

NAME: YTQAPI2D

<u>BEFORE</u> :	B7153551	<u>AFTER</u> :	B6639144	**147,150**
-----------------	----------	----------------	----------	-------------

ROUTINE NAME: YTSACE

<u>BEFORE</u> :	N/A	<u>AFTER</u> :	273
	150		<u>ROUTINE</u>

NAME: YTSBSL23

BEFORE: N/A AFTER: B3044905 **150**

ROUTINE #	YTSCMQ			
BEFORE:	N/A	AFTER:	B2552097	**150**
ROUTINE #	YTSEAT26			
BEFORE:	N/A	AFTER:	B6312697	**150**
ROUTINE #	YTSFOCI			
BEFORE:	N/A	AFTER:	B4656413	**150**
ROUTINE #	YTSMHRM			
BEFORE:	N/A	AFTER:	B2499615	**150**
ROUTINE #	YTSNUDEC			
BEFORE:	N/A	AFTER:	B13752134	**150**
ROUTINE #	YTSSIP30			
BEFORE:	N/A	AFTER:	B5052592	**150**
ROUTINE #	YTSSIPST			
BEFORE:	N/A	AFTER:	B5906306	**150**
ROUTINE #	YTSSWEM			
BEFORE:	N/A	AFTER:	B6610853	**150**

Routine list of preceding patches: 123,147

4.3 Access Requirements and Skills Needed for the Installation

Installation of MHA Update (YS*5.01*150) requires the following:

- Programmer access to VistA instance and ability to install KIDS build.
- Citrix Access Gateway (CAG) installs – access/ability to upload to the CAG.
- Network Share installs – access/ability to upload executable to the network share location.
- Individual workstation installs – access/ability to push executable and supporting files to required work stations.

4.4 Installation Procedure

1. Choose the PackMan message containing this patch and invoke the INSTALL/CHECK MESSAGE PackMan option.

2. Start up the Kernel Installation and Distribution System Menu [XPD MAIN]:

Edits and Distribution ...

Utilities ...

Installation ...

Select Kernel Installation & Distribution System Option: Installation

- | | |
|---|--|
| 1 | Load a Distribution |
| 2 | Verify Checksums in Transport Global |
| 3 | Print Transport Global |
| 4 | Compare Transport Global to Current System |
| 5 | Backup a Transport Global |
| 6 | Install Package(s) |
| | Restart Install of Package(s) |
| | Unload a Distribution |

3. From this menu, you may elect to use the following options (When prompted for the INSTALL NAME, enter YS*5.01*150):

- a. Backup a Transport Global - This option will create a backup message of any routines exported with this patch. It will not backup any other changes such as DD's or templates.
 - b. Compare Transport Global to Current System - This option will allow you to view all changes that will be made when this patch is installed. It compares all components of this patch (routines, DD's, templates, etc.).
 - c. Verify Checksums in Transport Global - This option will allow you to ensure the integrity of the routines that are in the transport global.
4. Use the Install Package(s) option and select the package YS*5.01*150
- a. When prompted "Want KIDS to INHIBIT LOGONs during the install? NO//", answer NO.
 - b. When prompted "Want to DISABLE Scheduled Options and Menu Options and Protocols? NO//", answer NO.

4.5 Post Install for Instruments

To get Progress Notes to show in CPRS Notes, use one of the following procedures - MHA3 Utilities or FileMan (see samples below).

After you have changed the PNote field, TIU Title Field or the Consult Note Title Field, make sure the: LAST EDIT DATE IN THE MH TESTS AND SURVEYS (#601.71) file is updated with the current date and time.

PROGRESS NOTE FIELD REQUIREMENTS FOR 150 INSTRUMENTS SHOULD HAVE THE WRITE FULL TEXT, PNOTE, AND R PRIVILEGE FIELDS POPULATED WITH EITHER A YES/NO, AS INDICATED IN THE TABLE BELOW:

Instrument Name	Write Full Text	Generate P-Note	R Privilege
The Florida Obsessive Compulsive Inventory (FOCI)	Yes	Yes	n/a
The Short Warwick–Edinburgh Mental Well-being Scale (SWEMWBS)	No	No	YSP
Mental Health Recovery Measure-10 (MHRM-10)	Yes	Yes	n/a
Eating Attitudes Test-26 Item (EAT-26)	Yes	Yes	n/a
Adverse Childhood Experiences (ACE)	Yes	Yes	n/a
Composite Morningness Questionnaire (CMQ)	Yes	Yes	n/a
SHORT INVENTORY OF PROBLEMS - AD (SIP-AD-30)	Yes	Yes	n/a
SHORT INVENTORY OF PROBLEMS - AD (SIP-AD-Start)	Yes	Yes	n/a
Borderline Symptom List 23 (BSL-23)	Yes	Yes	n/a
Nursing Delirium Screening Scale (NuDESC)	Yes	Yes	n/a
AD8	No	No	YSP

4.5.1 Instructions using the MHA3 Utilities Option

“Stop/Re-Start Progress Notes for an Instrument” on the MHA3 Utilities menu (YTQ MHA3 MENU) can be used to edit the PNote Field (Progress Note), TIU Title Field, and Consult Note Title Field without having to go into FileMan.

Based on what the PNote field should read, either Yes to show the Progress Note, or No to not show the Progress note, you would need to either “Start” or “Stop” the Progress Note field for the Instrument.

If a progress note is Yes and should appear, but it is NOT appearing, then you should use the START PROGRESS NOTES instructions below.

If a progress note is No and should NOT appear, but it IS appearing, then you should use the STOP PROGRESS NOTES instructions below.

START PROGRESS NOTES:

- 1 STOP/RE-START PROGRESS NOTES FOR AN INSTRUMENT
- 2 PRINT TEST FORM
- 3 DETAILED DEFINITION
- 4 DELETE PATIENT DATA
- 6 TEST USAGE
- 7 INSTRUMENT EXCHANGE
- 5 EXEMPT TEST

SELECT MHA3 UTILITIES <SPPNXT> OPTION: 1 STOP/RE-START PROGRESS FOR AN INSTRUMENT

SELECT MH TESTS AND SURVEYS NAME: [NAME HERE] GENERATE

PNOTE: NO// Y YES

TIU TITLE: MENTAL HEALTH ■ NOTE// CONSULT NOTE

TITLE: MENTAL HEALTH CONSULT NOTE//

STOP PROGRESS NOTES:

- 1 STOP/RE-START PROGRESS NOTES FOR AN INSTRUMENT
- 2 PRINT TEST FORM
- 3 DETAILED DEFINITION
- 4 DELETE PATIENT DATA
- 6 TEST USAGE
- 7 INSTRUMENT EXCHANGE
- 5 EXEMPT TEST

SELECT MHA3 UTILITIES <SPPNXT> OPTION: 1 STOP/RE-START PROGRESS FOR AN INSTRUMENT

SELECT MH TESTS AND SURVEYS NAME: [NAME HERE] GENERATE

PNOTE: YES// N NO

TIU TITLE: MENTAL HEALTH ■ NOTE// CONSULT NOTE

TITLE: MENTAL HEALTH CONSULT NOTE//

UPDATING LAST EDIT DATE IN FILEMAN IF YOU USE MHA3 UTILITIES:

SELECT OPTION: ENTER OR EDIT FILE ENTRIES

INPUT TO WHAT FILE: MH TESTS AND SURVEYS (250) EDIT
WHICH FIELD: ALL// LAST
1 LAST EDIT DATE
2 LAST EDITED CHOOSE
1-2: 1 LAST EDIT THEN EDIT
FIELD:

SELECT MH TESTS AND SURVEYS NAME: [NAME NAME HERE] LAST EDIT
DATE: FEB 5, 2019@19:31// NOW (MAR 04, 2019@06:10)

4.5.2 FileMan Instructions

1. The user needs to access FileMan and choose option 'ENTER OR EDIT FILE ENTRIES'.
2. IN INPUT TO WHAT FILE ENTER MH TESTS AND SURVEYS
3. For EDIT WHICH FIELD: – GENERATE PNOTE
4. For THEN EDIT FIELD: TIU TITLE
5. For THEN EDIT FIELD: CONSULT NOTE TITLE
6. Last Edit Date: (Enter Current Date)

Press return and at the prompt: SELECT MH TESTS AND SURVEYS NAME: ENTER INSTRUMENT
FOR FIELD GENERATE PNOTE: ENTER YES OR NO, NO THE
INSTRUCTIONS PROVIDED FOR EACH INSTRUMENT IN THIS INSTALL GUIDE.
FOR FIELD TIU TITLE: ENTER THE TIU TITLE AS DEFINED BY YOUR
SITE
FOR FIELD CONSULT NOTE TITLE: ENTER THE TITLE OF THE AS DEFINED
BY YOUR SITE
FOR LAST EDIT DATE: ENTER CURRENT DATE

EXIT FILEMAN AND YOUR NOTES SHOULD NOW SAVE TO CD

Example for Instrument: 'CES'

SELECT OPTION: ?
ANSWER WITH OPTION NAME CHOOSE
FROM:
1 ENTER OR EDIT FILE ENTRIES
2 PRINT FILE ENTRIES
3 SEARCH FILE ENTRIES
4 MODIFY FILE ATTRIBUTES
5 INQUIRE TO FILE ENTRIES
6 UTILITY FUNCTIONS
7 OTHER OPTIONS
8 DATA DICTIONARY UTILITIES
9 TRANSFER ENTRIES

SELECT OPTION: 1 ENTER OR EDIT FILE ENTRIES

```

INPUT TO WHAT FILE: MH TESTS AND  (250 EDIT
WHICH FIELD: ALL// GENERATE
THEN EDIT FIELD: TIU TITLE
THEN EDIT FIELD: CSM TITLE THEN
EDIT FIELD: LAST EDIT DATE


```

SELECT MH TESTS AND **CES** GENERATE PNOTE:
YES// YES
TIU TITLE: MENTAL HEALTH DIAGNOSTIC STUDY **NO**//
CONSULT NOTE TITLE: MENTAL HEALTH **NOTE**// LAST EDIT DATE: FEB
5, 2019@19:31// **NOW** (MAR 04,
2019@06:10)


4.6 Installation Verification Procedure

To verify that everything is installed properly, do the following:

- Launch CPRS.
- From the CPRS menu, select Tools, then Mental Health Assistant (MHA).
- As MHA starts you should see the splash screen with version 1.0.3.75 displayed in the lower right corner.

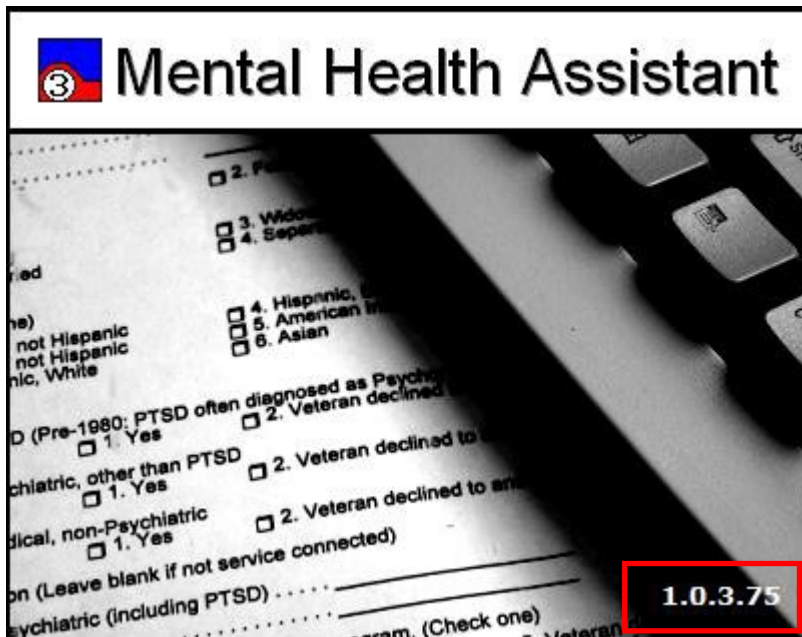

When MHA prompts for user credentials, it should use the Personal Identity Verification (PIV) Personal Identification Number (PIN) prompt, rather than access/verify codes.

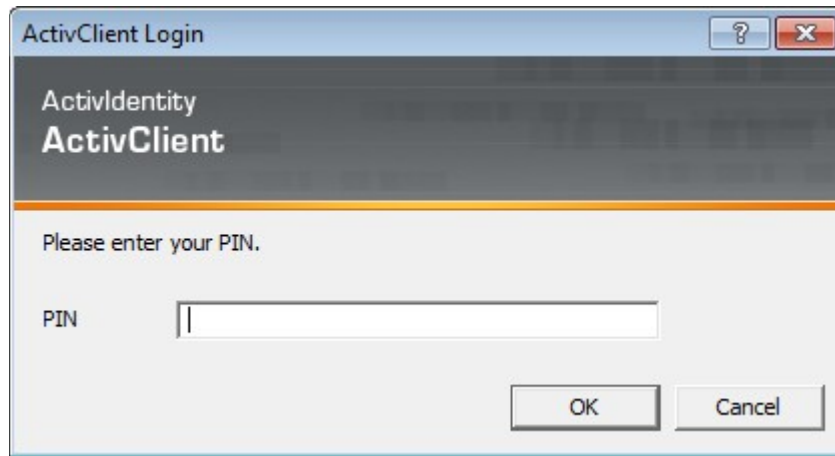

If MHA prompts for a PIV PIN as it is starting, the software has been installed properly. (If the Single Sign-On Integration (SSOi) servers are not available, you will see the usual access/verify prompt.)

After installation, notify mental health package users that the following new Suicide and Depression Screening instruments are available:

Instrument	Identifier
The Florida Obsessive Compulsive Inventory (FOCI)	FOCI
The Short Warwick-Edinburgh Mental Well-being Scale (SWEMWBS)	SWEMWBS
Mental Health Recovery Measure-10 (MHRM-10)	MHRM-10
Eating Attitudes Test-26 Item (EAT-26)	EAT-26
Adverse Childhood Experiences (ACE)	ACE
Composite Morningness Questionnaire (CMQ)	CMQ
SHORT INVENTORY OF PROBLEMS - AD (SIP-AD-30)	SIP-AD-30
SHORT INVENTORY OF PROBLEMS - AD (SIP-AD-Start)	SIP-AD-START
Borderline Symptom List 23 (BSL-23)	BSL-23
Nursing Delirium Screening Scale (NuDESC)	NuDESC
AD8	AD8

4.7 System Configuration

N/A

4.8 Database Tuning

N/A

5 Back-Out Procedure

5.1 Back-Out Strategy

It is possible to partially back-out the installation of YS*5.01*150. This would involve restoring instrument specifications to their previous state and then restoring the saved routines. The back-out of changes to the data dictionary would require a patch to a patch.

5.2 Back-Out Considerations

Please contact VistA support and the development team before attempting a back-out. The back-out procedure will still leave some changes in place. In addition, the installation of subsequent patches may be problematic if YS*5.01*150 is not installed.

5.3 Back-Out Criteria

A back-out should only be considered if there is a patient safety issue, if MHA no longer functions, or if there is some other catastrophic failure.

5.4 Back-Out Risks

The risks vary depending on what is causing the failure of the system. The main risk is that the Mental Health package would be left in an unknown configured state.

5.5 Authority for Back-Out

The VistA system manager determines if a back-out of YS*5.01*150 should be considered.

5.6 Back-Out Procedure

If you wish to restore newly installed instruments to their previous state, you must do that before any other back-out steps. See the instructions for restoring the previous instrument state in the Rollback Procedure section to do this.


To back-out routines, you must have already selected the “Backup a Transport Global” option during the installation process. To restore the previous routines:

1. Choose the PackMan message containing the backup you created during installation.
2. Invoke the INSTALL/CHECK MESSAGE PackMan option.
3. Select Kernel Installation & Distribution System Option: Installation
4. Use the Install Package(s) option to install the previously saved routines.

If you need to back-out data dictionary modifications, remove protocols, options, or templates, you will need to contact the development team for a patch.

5.7 Back-out Verification Procedure

A successful back-out may be verified by running MHA and seeing a splash screen with the highlighted version number:

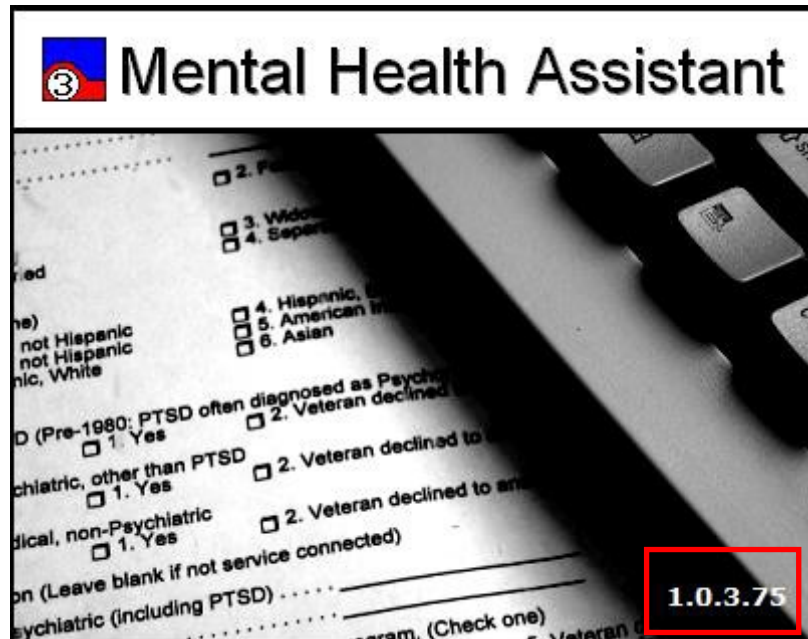

MHA should prompt for access/verify instead of PIV PIN and run successfully.

Verification of the back-out procedure would be the resolution of the problem that caused the need for the back-out.

6 Rollback Procedure

6.1 Rollback Considerations

YS*5.01*150 adds new and updates existing mental health instruments. It is possible to roll back these changes within one week of the installation.

6.2 Rollback Criteria

A rollback might be considered if the behavior of mental health instruments appears to be adversely affected after installation of YS*5.01*150. The VistA support and product development team should be contacted to determine if there is an alternative fix short of a rollback.

6.3 Rollback Risks

A rollback could adversely impact future installations of mental health instruments and cause problems with scoring existing mental health instruments.

6.4 Authority for Rollback

The VistA system manager determines if a rollback of mental health instruments distributed by YS*5.01*150 should be considered.

6.5 Rollback Procedure

These steps assume that there is a compelling reason to rollback specific instruments to their previous state.

For instruments that have been inactivated by YS*5.01*150 that need to be made active again:

- Using FileMan, edit the OPERATIONAL field (#10) and the LAST EDIT DATE field (#18) in the MH TESTS AND SURVEYS file (601.71). Select the instrument that requires re-activation.
- Change the value of the OPERATIONAL field from “Dropped” back to “Yes”
- Change the value of the LAST EDIT DATE field to ‘NOW’.

Should it be required to move instruments back to being scored in YS_MHA_AUX DLL, contact the Mental Health development team for a routine that can find the appropriate records and make the replacement.

Optionally, if you want to see how many records will be restored, choose “Trial Install” then select the number of the backup you wish to restore.

When you are ready to restore an instrument, choose “Install Exchange Entry” then select the number of the backup you want to restore.

6.6 Rollback Verification Procedure

Verify the restore by checking to see that the instrument behaves as it did prior to the install.

Appendix A - Acronyms

Acronym	Definition
CAG	Citrix Access Gateway
CD	Critical Decision
DIBO&RG	Deployment, Installation, Back-out, and Rollback Guide
IOC	Initial Operating Capability
MHA	Mental Health Assistant
PIN	Personal Identification Number
PIV	Personal Identity Verification
SPP	Suicide Prevention Package
SQA	Software Quality Assurance
SSOi	Single Sign-On Integration

Acronym	Definition
VA	Department of Veterans Affairs
VAMC	Veterans Affairs Medical Center
VIP	Veteran-focused Integration Process
VistA	Veterans Health Information System and Technology Architecture