



Update _2_0_264

Clinical Reminders

**VA-SMART DIALOG AND MAMMOGRAM REMINDER
UPDATE**

Install Guide

June 2022

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION 1
INSTALL DETAILS 9
INSTALL EXAMPLE 10
POST INSTALLATION 14

Introduction

Several revisions have been made to the VA-WH SMART BREAST IMAGING FOLLOW-UP reminder dialog and the VA-WH MAMMOGRAM SCREENING reminder definition.

1. The VA-WH NEXT BREAST PROCEDURE reminder term CONDITION has been updated. A Beginning Date/Time of 10.1.2020 has also been added to prevent the term from inappropriately evaluating procedures completed before the SMART functionality was deployed in CPRS 31b.
2. The following reminder terms have been moved from the reminder resolution and added to the reminder cohort logic as exclusions. These findings were not temporarily resolving the reminder correctly in conjunction with the Custom Due Date finding. Now, during the deferral and order period, the reminder will be N/A. It will come due again after the temporary resolution time has passed.
 - a. VA-MAMMOGRAM SCREEN DEFER
 - b. VA-WH MAMMOGRAM ORDER
 - c. VA-WH BREAST CARE ORDER HEALTH FACTOR
3. New finding FOUND TEXT has been added to display in the clinical maintenance to assist local staff understand temporary logic.
 - a. If the VA-BREAST CARE ORDER HEALTH FACTOR term is true, the following will display in clinical maintenance:

```
Cohort:
Reminder Term: VA-WH BREAST CARE ORDER HEALTH FACTOR
Health Factor: Women's Health Order Mammogram Screen HF
02/23/2022@12:49

* * * * *
*
* _
* _   A mammogram has been ordered in the past 90 days.
* _   This reminder will display as NOT APPLICABLE, even though
* _   it is applicable to this patient. If the patient
* _   does not have a mammogram within the 90 days, then the reminder
* _   will become DUE again.
* _
* _
* * * * *
```

- b. If the VA-WH MAMMOGRAM SCREEN NOT INDICATED term is true, the following will display in clinical maintenance:

```
Cohort:
Reminder Term: VA-WH MAMMOGRAM SCREEN NOT INDICATED
Health Factor: Women's Health Mammogram Screen Not Indicated
02/23/2022@12:49
Comments: TEST

* * * * *
*
* _
* _   Mammogram screening is not indicated according to documentation.
* _
* _
* * * * *
```

- c. If the VA-WH MAMMOGRAM SCREEN DEFER term is true, the following will display in clinical maintenance:

```

Cohort:
Reminder Term: VA-WH MAMMOGRAM SCREEN DEFER
Health Factor: Women's Health Mammogram Declined
02/23/2022@12:49
Comments:

* * * * *
*
* _ Mammogram Screening reminder has been deferred.
* _ During the 4 months after the deferral date, this reminder
* _ will display as NOT APPLICABLE, even though it is applicable
* _ to this patient. The reminder will become DUE again 4 months
* _ after the deferral was documented.
*
*
* * * * *

```

- d. If the VA-WH MAMMOGRAM ORDER term is true, the following will display in clinical maintenance:

```

Cohort:
Reminder Term: VA-WH MAMMOGRAM ORDER
Orderable Item: MAMMOGRAM
03/21/2022@10:41 Status: active, Start date: 03/21/2022@10:41, Stop date:
04/20/2022 Duration: 0 D

* * * * *
*
* _ A mammogram has been ordered in the past 90 days.
* _ While the patient has an active order, this reminder
* _ will display as NOT APPLICABLE, even though it is applicable
* _ to this patient. If the order is discontinued or the patient
* _ does not have a mammogram within the 90 days, then the reminder
* _ will become DUE again.
*
*
* * * * *

```

4. Two incorrect URLs have been removed from the Mammogram Screening reference information.
5. Corrected an error that caused the reminder to be not applicable on patients <40 yo who had a mammogram and a new screening date documented in the VA-WH SMART BREAST IMAGING FOLLOW-UP note.
6. Branching logic has been added to the VA-WH MAMMOGRAM SCREENING reminder dialog to hide the “Change the frequency of mammograms for this patient” option if the VA-WH NEXT BREAST PROCEDURE reminder term is true. This indicates a future screening recommendation has been documented using the SMART dialog since CPRS 31b implementation. Once the mammogram reminder due date has been set using SMART, then the original frequency change options in the mammogram screening dialog do not impact the reminder logic.
7. The option to document “Patient breast images are normal” has been revised to “Patient breast images are benign/normal”.

Before:

Please summarize result:

- Patient breast images are normal
- Patient breast images are abnormal
- Patient breast images are incomplete

After:

Please summarize result:

- Patient breast images are benign/normal
- Patient breast images are abnormal
- Patient breast images are incomplete

8. Corrected an error occurring when a user selected 'no action required' under the next steps section.
9. Corrected a reminder resolution error occurring when a user selected "Other" in the Diagnostics options section.
10. A new general finding is included that will ensure the date a procedure report was received in SMART will be stored.
11. The "Written Report Available" section has been revised to ensure required fields are easier to complete.
12. A new section has been added to document breast density per the imaging report.
13. Options for documenting outside Breast Tomosynthesis procedures have been added.

Before:

Record prior or outside mammogram:

Written Report Available

NOTE: Obtain outside written report and forward for scanning into Vista Imaging and entry of BIRADS codes by Radiology

Outside report documentation:

Outside report

Received on: * / / 2022

Location:

Select procedure to document (document next steps in next section):

Document follow-up actions:

Document how the patient was notified of the result

- Phone call
- In person
- Non-Registered Letter
- Registered Letter
- Secure Messaging
- Patient Not Notified: information forwarded to PACT/PC provider for patient notification
- Patient not contacted, chose date for tickler followup:

Additional Information:

Written Report Available

NOTE: Obtain outside written report and forward for scanning into Vista Imaging and entry of BIRADS codes by Radiology

Outside report documentation:

Outside report

Received on: * ...

Location:

Select procedure to document (document next steps in next section):

- Mammogram Screening
- Diagnostic Bilateral Mammogram:
- Diagnostic Unilateral Mammogram:
Side: *
- MRI
- Ultrasound:
Side: *
- Screening Breast Tomosynthesis
- Diagnostic Bilateral Breast Tomosynthesis
- Diagnostic Unilateral Breast Tomosynthesis:
Side: *
- Biopsy:

Document breast density per the imaging report*:

- Dense (heterogeneously dense or extremely dense)
- Not dense (predominately fatty or scattered fibroglandular tissue)
- Unknown/not documented

After:

Document how the patient was notified of the result

- Phone call
- In person
- Non-Registered Letter
- Registered Letter
- Secure Messaging
- Information forwarded to PACT/PC provider for patient notification
- Patient not contacted, chose date for tickler followup:
- Patient already notified/not needed

Additional Information:

UPDATE_2_0_264 contains 1 Reminder Exchange entry:
VA-SMART DIALOG AND MAMMOGRAM REMINDER UPDATE

The exchange file contains the following components:

TEMPLATE FIELDS

SMART WORD 40X4
VA-WH SMART OUTSIDE OTHER
VA-WH SMART ULTRASOUND TYPE
VA-WH SMART UNILAT TYPE
BLANK TEXT FIELD FOR DIALOGS
VA-WH MAMMOGRAM GUIDANCE

GENERAL FINDINGS

TICKLER REMINDER COMMENT
TICKLER REMINDER
FOLLOW UP WITH PATIENT COMMENT
FOLLOW UP WITH PATIENT
WH PROCEDURE ACTION COMMENT
RADIOLOGY REPORT RECEIVED
NOTIFICATION NOT NEEDED
EMAIL CONTACTED
SECURE MESSAGING
LETTER REGISTERED CONTACTED
LETTER (CERTIFIED)
LETTER CONTACTED
LETTER
WHO CONTACTED
IN PERSON
PHONE CALL
BREAST CARE CASCADE NAME
CLOSE CASCADE BREAST CARE
BR CLOSE AND FOLLOW COMPLETE
BR CLOSED PROCEDURE
BREAST REFER FOR BIOPSY
BIOPSY
WV NOTIFICATION PATIENT
OBTAIN PRIOR FILMS
BREAST OBTAIN PRIOR FILMS
BREAST BIOPSY ALREADY OBTAIN
BIOPSY ALREADY OBTAINED
BREAST CURRENTLY UNDER TREATMENT
CURRENTLY UNDER TREATMENT
CONSULT
BREAST CONSULT
BREAST REFER TO ONCOLOGIST
REFER TO ONCOLOGIST
BREAST REFER TO SURGEON
REFER TO SURGEON
IMAGE COMMENT

BREAST NEED MRI
RETURN FOR MRI
BREAST NEED ULTRASOUND
RETURN FOR ULTRASOUND
BREAST NEED MAMMOGRAM
RETURN FOR MAMMOGRAM
BR NO NEXT TREATMENT DATE
BR NO NEXT TREATMENT
WH PATIENT
BREAST next MAM 2Y
AGE BASED SCREENING
BREAST next MAM 1Y
BREAST next MAM AT START AGE
MAMMOGRAM SCREENING 2Y
MAMMOGRAM SCREENING 1Y
BREAST UNSATISFACTORY RESULT
BR UNSATISFACTORY RESULT
BREAST ABNORMAL RESULT
BR ABNORMAL RESULT
BREAST NORMAL RESULT
BR NEM RESULT
BI-RAD 6 RESULT
BI-RAD 6
BI-RAD 5 RESULT
BI-RAD 5
BI-RAD 4 RESULT
BI-RAD 4
BI-RAD 3 RESULT
BI-RAD 3
BI-RAD 2 RESULT
BI-RAD 2
BI-RAD 1 RESULT
BI-RAD 1
BI-RAD 0 RESULT
BI-RAD 0
BREAST UNSATISFACTORY DENSE RESULT
BREAST ABNORMAL DENSE RESULT
BREAST NORMAL DENSE RESULT
BI-RAD 6 DENSE RESULT
BI-RAD 5 DENSE RESULT
BI-RAD 4 DENSE RESULT
BI-RAD 3 DENSE RESULT
BI-RAD 2 DENSE RESULT
BI-RAD 1 DENSE RESULT
BI-RAD 0 DENSE RESULT
BREAST TOMOSYNTHESIS UNILAT REPORT
BR OUTSIDE REPORT PROCEDURE DATE
BR OUTSIDE ACCESSION
BREAST TOMOSYNTHESIS BILAT REPORT
BREAST TOMOSYNTHESIS SCREENING REPORT
OUTSIDE ULTRASOUND

OUTSIDE MRI
OUTSIDE MAMMOGRAM DX UNILAT
OUTSIDE MAMMOGRAM DX BILAT
MAMMOGRAM SCREENING REPORT
OUTSIDE REPORT
VIEW PROGRESS NOTE TEXT
ALERT ID
CLEAR ONE
WV DIAGNOSIS
WV PROCEDURE ACCESS
MAMMOGRAM TEST
VIEW LAST THREE BR TRMTS
VIEW EPISODE
SMART BI-RAD ALERT
CLEAR OTHER ALERTS ID
CLEAR ALL FOR PROVIDER
BR NEXT DUE AT AGE
BR MAMMOGRAM SCREENING NEXT TRMT

HEALTH FACTORS

VA-REMINDER UPDATES [C]
VA-UPDATE_2_0_264
WH MAMMOGRAM [C]
WH ORDER MAMMOGRAM SCREEN HF
WH OUTSIDE INCOMPLETE MAMMOGRAM
WH OUTSIDE ABNL MAMMOGRAM
WH OUTSIDE NORMAL MAMMOGRAM
VA-WH BREAST RESULT BIRAD 6
VA-WH BREAST RESULT BIRAD 5
VA-WH BREAST RESULT BIRAD 4
VA-WH BREAST RESULT BIRAD 3
VA-WH BREAST RESULT BIRAD 2
VA-WH BREAST RESULT BIRAD 1
VA-WH BREAST RESULT BIRAD 0
VA-WH BREAST OTHER
VA-WH MAMMOGRAM DX UNILATERAL
VA-WH MAMMOGRAM DX BILATERAL
VA-WH MAMMOGRAM SCREENING
WH SMART OUTSIDE REPORT INFO
WH MAMMOGRAM SCREEN NOT INDICATED
WH BR CA SCREEN N/A 5 YRS-COMORBIDITIES
WH BR CA SCREEN N/A 5 YRS-LE<5YRS
VA LIFE EXPECTANCY [C]
VA LIMITED LIFE EXPECTANCY
WH BILATERAL MASTECTOMY
WH MAMMOGRAM SCREEN FREQ - 1Y
WH MAMMOGRAM SCREEN FREQ - 6M

WH MAMMOGRAM DEFERRED
VA-WH BR CA 40-44 WANTS SCREENING
VA-WH BR CA 40-44 BEGIN AGE 45
WH MAMMOGRAM SCREEN FREQ - 2Y
WH MAMMOGRAM DECLINED
WH BR CA 40-49 WANTS SCREENING
WH MAMMOGRAM SCREEN FREQ - 4M
VA-WH BR CA DEFER 1Y
VA-WH BR CA DEFER 6M
REMINDER FACTORS [C]
INACTIVATE BREAST CANCER SCREEN
WH ORDER MAMMOGRAM UNILAT HF
WH ORDER MAMMOGRAM BILAT HF
WH MAMMOGRAM OUTSIDE

REMINDER SPONSOR

OFFICE OF PUBLIC HEALTH AND ENVIRONMENTAL HAZARDS AND THE WOMEN
VETERANS HEALTH STRATEGIC HEALTH CARE GROUP
VHA NATIONAL CENTER FOR HEALTH PROMOTION AND DISEASE PREVENTION (NCP)
AND VHA OFFICE OF WOMEN'S HEALTH
NATIONAL CENTER FOR HEALTH PROMOTION AND DISEASE PREVENTION

REMINDER TAXONOMY

VA-MASTECTOMY
VA-BREAST TUMOR
VA-TERMINAL CANCER PATIENTS
VA-MAMMOGRAM/SCREEN
VA-WH BILATERAL MASTECTOMY

REMINDER TERM

VA-REMINDER UPDATE_2_0_264
VA-WH BL ACTIVE BREAST IMAGE ALERT DATA
VA-WH OPEN PROCEDURES PER CASCADE
VA-WH BL BR ACTIVE ALERT ID
VA-PATIENT IS MALE
VA-WH MAMMOGRAM ALERT FOR MALE
VA-WH BL BR BI-RAD 6
VA-WH BL BR BI-RAD 5
VA-WH BL BR BI-RAD 4
VA-WH BL BR BI-RAD 3
VA-WH BL BR BI-RAD 0
VA-WH BL BR BI-RAD 2
VA-WH BL BR BI-RAD 1
VA-AGE >74
VA-BL AGE>74
VA-WH AGE BETWEEN 55 AND 74

VA-WH AGE BETWEEN 45 AND 54
VA-WH MAMMOGRAM AGE<55
VA-WH OPEN BR PROCEDURES WITH BIRAD CODES
VA-WH OPEN BR PROCEDURES WITHOUT BIRAD CODES AND DENSE
VA-WH OPEN BR PROCEDURES WITH BIRAD CODES DENSE
VA-WH BL GET BREAST IMAGES
VA-WH BL BREAST CASCADE EXIST
VA-WH NEXT BREAST PROCEDURE
VA-AGE >39&<45
VA-WH CHAMP VA PATIENT
VA-WH BR CA 40-44 WANTS SCREEN TERM
VA-WH BR CA 40-44 BEGIN AGE 45
VA-WH MAMMOGRAM ORDER
VA-WH HX BREAST CANCER/ABNORMAL MAM
VA-WH MAMMOGRAM SCREEN FREQ - 2Y
VA-WH MAMMOGRAM SCREEN FREQ - 1Y
VA-WH MAMMOGRAM SCREEN FREQ - 6M
VA-WH MAMMOGRAM SCREEN FREQ - 4M
VA-WH MAMMOGRAM UNSATISFACTORY IN RAD/WH PKG
VA-WH MAMMOGRAM SCREEN DEFER
VA-WH MAMMOGRAM SCREEN NOT INDICATED
VA-WH BREAST CARE ORDER HEALTH FACTOR
VA-WH MAMMOGRAM SCREEN DONE
VA-WH MAMMOGRAM SCREEN IN RAD PKG
VA-WH MAMMOGRAM SCREEN IN WH PKG
VA-TERMINAL CANCER PATIENT
VA-WH BILATERAL MASTECTOMY

REMINDER DEFINITIONS

VA-WH BL BR ALERT WITH NO OPEN PROCEDURE
VA-WH OPEN BREAST PROCEDURES PER CASCADE
VA-WH MAMMOGRAM SCREENING

REMINDER DIALOG

VA-WH SMART BREAST IMAGING FOLLOW-UP
VA-WH MAMMOGRAM SCREENING

Pre-Installation Steps

1. Record the orders mapped to the VA-WH SMART BREAST IMAGING FOLLOW-UP reminder dialog. The site will need to re-map these orders to the reminder dialog during the install.
 - a. Identify the order dialog mapped to the group VA-WH BR ADDITIONAL IMAGING. This order dialog will be used to replace the Q.BREAST TREATMENT MENU during the installation.

- b. Identify the order dialog mapped to the element VA-WH BR CONSULT. This order dialog will be used to replace the Q.GMRCOR CONSULT during the installation.
2. Record the order dialog mapped (if one is mapped) to the VA-WH GP ORDERS - MAMMOGRAM SCREEN reminder group.
3. Create a copy of the reminder inquiry output for the VA-WH MAMMOGRAM SCREENING reminder using either CPRS or Vista. Use this copy to confirm all local mapping for reminder terms are correct after post-installation.

Install Details

This update is being distributed as a web host file. The address for the host file is:

https://REDACTED/UPDATE_2_0_264.PRD

The file will be installed using Reminder Exchange, programmer access is not required.


Installation:

=====

This update can be loaded with users on the system. Installation will take less than **5** minutes.

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

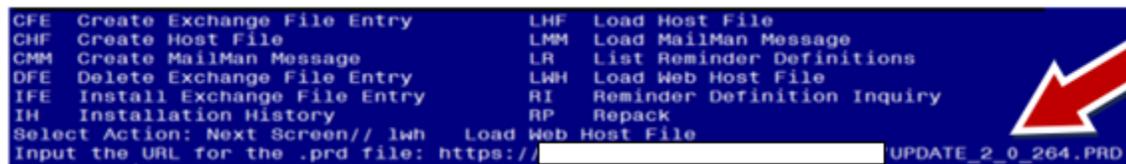

```

CFE Create Exchange File Entry      LHF Load Host File
CHF Create Host File                LMM Load MailMan Message
CMM Create MailMan Message          LR List Reminder Definitions
DFE Delete Exchange File Entry      LWH Load Web Host File
IFE Install Exchange File Entry     RI Reminder Definition Inquiry
IH Installation History              RP Repack
Select Action: Next Screen// lwh    Load Web Host File
Input the URL for the .prd file: https://[REDACTED]UPDATE_2_0_264.PRD

```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

https://vaww.va.gov/reminders/docs/UPDATE_2_0_264.PRD

You should see a message at the top of your screen that the file successfully loaded.

Search and locate an entry titled **UPDATE_2_0_264 VA-SMART DIALOG AND MAMMOGRAM REMINDER UPDATE** in reminder exchange.

```

+Item  Entry                               Source                               Date Packed
 280   UPDATE_2_0_264 VA-SMART             ←                               03/30/2022@11:54
      DIALOG AND MAMMOGRAM REMINDER
      UPDATE
 281   UPDATE_2_0_265 VA-COVID-19         03/03/2022@07:57
      TEMPLATE UPDATES
 282   UPDATE_2_0_267 GEC GERIATRIC       03/28/2022@05:59
      EVALUATION
 283   UPDATE_2_0_27  VA-CHOICE           06/14/2017@11:33
      DIALOGS
 284   UPDATE_2_0_272 VA-ATHLETES        02/22/2022@08:37
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE   Create Exchange File Entry         LHF   Load Host File
CHF   Create Host File                   LMM   Load MailMan Message
CMM   Create MailMan Message             LR    List Reminder Definitions
DFE   Delete Exchange File Entry         LWH   Load Web Host File
IFE   Install Exchange File Entry        RI    Reminder Definition Inquiry
IH    Installation History               RP    Repack
Select Action: Next Screen// IFE   Install Exchange File Entry
Enter a list or range of numbers (1-488): 280

```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry **UPDATE_2_0_264 VA-SMART DIALOG AND MAMMOGRAM REMINDER UPDATE** (in this example it is entry 280. It will vary by site). The date of the exchange file should be 03/30/2022.

```

Date Packed: 03/30/2022@11:54:28
Package Version: 2.0P71

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-WH MAMMOGRAM SCREENING
  VA-WH SMART BREAST IMAGING FOLLOW-UP

REMINDER TERM
  VA-REMINDER UPDATE_2_0_264

REMINDER DEFINITION
  VA-WH MAMMOGRAM SCREENING
Non-exchangeable order dialog(s):
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
IA    Install all Components             IS    Install Selected Component
Select Action: Next Screen// IA   Install all Components .....

```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**

- For all new entries, choose **Install**
- For REMINDER DIALOG LINK TYPE entries that exist, but the packed component is different, **OVERWRITE**
- For TERMS and TAXONOMIES that exist, but the packed component is different **SKIP**, except the following term. **Overwrite** this term if asked:
- **VA-WH NEXT BREAST PROCEDURE**
- For REMINDER DEFINITIONS that exist, but the packed component is different **OVERWRITE**.

You will be prompted to install the reminder dialog component (**VA-WH SMART BREAST IMAGING FOLLOW-UP**). You will choose to install all even though some components may have an X to indicate that they exist.

Packed reminder dialog: VA-WH SMART BREAST IMAGING FOLLOW-UP

Item	Seq.	Dialog Findings	Type	Exists		
1		VA-WH SMART BREAST IMAGING FOLLOW-UP	dialog	X		
2	10	VA-WH GP BR OPEN PROCEDURES* Finding: *NONE*	group	X		
3	10.10	VA-WH BR NO OPEN PROCEDURES* Finding: *NONE*	element	X		
4	10.15	VA-WH BR EPISODE OF CARE Finding: VIEW EPISODE (REMINDER GENERAL FINDING)	element	X		
5		PXRM GF VIEW BUTTON	prompt	X		
6	10.20	VA-WH BR LAST 3 TRMTS Finding: VIEW LAST THREE BR TRMTS (REMINDER GENERAL FINDING)	element	X		
7		PXRM GF VIEW BUTTON	prompt	X		
8	10.25	VA-WH BLANK LINE Finding: *NONE*	element	X		
9	10.30	VA-WH GP BR BI-RAD SELECTOR*	group	X		
+ + Next Screen - Prev Screen ?? More Actions						
DD		Dialog Details	DT	Dialog Text	IS	Install Selected
DF		Dialog Findings	DU	Dialog Usage	QU	Quit
DS		Dialog Summary	IA	Install All		
Select Action: Next Screen// IA Install All						

At the **Select Action** prompt, type **IA** to install the dialog: **VA-WH SMART BREAST IMAGING FOLLOW-UP**.

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

FINDING entry Q.BREAST TREATMENT MENU does not exist.

Enter a code from the list.

Select one of the following:

- D Delete
- P Replace with an existing entry
- Q Quit the install

Enter response: **P**

Select ORDER DIALOG NAME: Enter local breast treatment menu identified in Pre-Installation steps.

Item	Seq.	Dialog Findings	Type	Exists
1		VA-WH SMART BREAST IMAGING FOLLOW-UP	dialog	X
2	10	VA-WH GP BR OPEN PROCEDURES*	group	X
		Finding: *NONE*		
3	10.10	VA-WH BR NO OPEN PROCEDURES*	element	X
		Finding: *NONE*		
4	10.15	VA-WH BR EPISODE OF CARE	element	X
		Finding: VIEW EPISODE (REMINDER GENERAL FINDING)		X
5		PXRM GF VIEW BUTTON	prompt	X
6	10.20	VA-WH BR LAST 3 TRMTS	element	X
		Finding: VIEW LAST THREE BR TRMTS (REMINDER GENERAL FINDING)		X
7		PXRM GF VIEW BUTTON	prompt	X
8	10.25	VA-WH BLANK LINE	element	X
		Finding: *NONE*		
9	10.30	VA-WH GP BR BI-RAD SELECTOR*	group	X
+ + Next Screen - Prev Screen ?? More Actions				
DD	Dialog Details	DT	Dialog Text	IS Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU Quit
DS	Dialog Summary	IA	Install All	
Select Action: Next Screen// Q Quit				

After completing this dialog install, you will type **Q**.

You will be prompted to install the reminder dialog component (**VA-WH MAMMOGRAM SCREENING**) You will choose to install all even though some components may have an X to indicate that they exist.

Packed reminder dialog: VA-WH MAMMOGRAM SCREENING

Item	Seq.	Dialog Findings	Type	Exists
1		VA-WH MAMMOGRAM SCREENING	dialog	X
2	5	VA-WH MAMMOGRAM SCREENING HEADER*	group	X
		Finding: *NONE*		
3	5.5	VA-WH GP MAM SCREEN HEADER*	group	X
		Finding: *NONE*		
4	5.5.5	VA-WH MAM HEADER1 - SCREEN	element	X
		Finding: *NONE*		
5	5.5.10	VA-WH GP MAM SCREEN HEADER-EXPANDED	group	X
		Finding: *NONE*		
6	5.5.10.5	VA-WH MAM HEADER2 - SCREEN	element	X
		Finding: *NONE*		
7	5.5.15	VA-WH DOTTED LINE	element	X
		Finding: *NONE*		
8	5.5.20	VA-WH GP MAMMOGRAM SCREENING	group	X
+ + Next Screen - Prev Screen ?? More Actions				
DD	Dialog Details	DT	Dialog Text	IS Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU Quit
DS	Dialog Summary	IA	Install All	
Select Action: Next Screen// IA				

At the **Select Action** prompt, type **IA** to install the dialog: **VA-WH MAMMOGRAM SCREENING**).

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

FINDING ITEM: [Add order dialog here.](#)

- 2. Using the reminder inquiry output from the pre-installation step, compare all reminder terms in the VA-WH MAMMOGRAM SCREENING reminder definition to confirm all local mapping is correct.**
- 3. Confirm all local orders are mapped correctly in the VA-WH SMART BREAST IMAGING FOLLOW-UP reminder dialog.**