



Update_2_0_408

Clinical Reminders

**VA-PGX (PHASER) CLOPIDOGREL/PREVIOUS
RESULTS CROC UPDATE**

Install Guide

February 2024

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION 1

PRE-INSTALL..... 2

INSTALL DETAILS 2

INSTALL EXAMPLE 2

POST INSTALLATION 6

Introduction

Reminder Update 408 distributes 2 reminder order checks to evaluate pharmacogenomic influence on prescribing for clopidogrel plus 4 reminder order checks to alert providers to already available pharmacogenomic test results should they attempt to order another pharmacogenomic test panel. Below are the Reminder Order check groups that are distributed. These groups also contain the corresponding Reminder Order check rules. These order checks complement the existing portfolio of PHASER pharmacogenomic order checks previously distributed.

PHASER VA Provider information and Site Deployment Materials: [PHASER-SharePoint](#) (only accessible within the VA firewall). PHASER public website: www.cancer.va.gov/phaser.asp

UPDATE_2_0_408 contains 1 Reminder Exchange entry:

UPDATE_2_0_408 VA-PGX (PHASER) CLOPIDOGREL/PREVIOUS RESULTS CROC

The exchange file contains the following components:

REMINDER SPONSOR

- 1 VHA PHARMACY BENEFITS MANAGEMENT
- 2 VHA ONCOLOGY

HEALTH FACTORS

- 3 VA-PGX GENE TESTING [C]
- 4 VA-PGX 16 GENE PANEL
- 5 VA-PGX 11 GENE PANEL
- 6 VA-PGX 8 GENE PANEL
- 7 VA-PGX 18 GENE PANEL
- 8 VA-PGX 15 GENE PANEL
- 9 VA-PGX CYP3A5 POOR METABOLIZER
- 10 VA-PGX CYP3A5 NORMAL METABOLIZER
- 11 VA-PGX CYP3A5 INTERMED METABOLIZER
- 12 VA-PGX CYP2C19 UNINFORMATIVE
- 13 VA-PGX CYP2C19 POOR METABOLIZER
- 14 VA-PGX CYP2C19 NORMAL METABOLIZER
- 15 VA-PGX CYP2C19 RAPID METABOLIZER
- 16 VA-PGX CYP2C19 ULTRARAPID METABOLIZER
- 17 VA-PGX CYP2C19 INTERMED METABOLIZER
- 18 VA-REMINDER UPDATES [C]
- 19 VA-UPDATE_2_0_408
- 20 VA-PGX TESTING ACCEPTED

REMINDER TERM

- 21 VA-PGX CYP2C19 NOT INTERMEDIATE
- 22 VA-PGX CYP2C19 INTERMED METABOLIZER
- 23 VA-PGX CYP2C19 NOT POOR
- 24 VA-PGX CYP2C19 POOR METABOLIZER
- 25 VA-REMINDER UPDATE_2_0_408

REMINDER DEFINITION

- 26 VA-PGX PREVIOUS 16 GENE PANEL
- 27 VA-PGX PREVIOUS 11 GENE PANEL
- 28 VA-PGX PREVIOUS 8 GENE PANEL

- 29 VA-PGX PREVIOUS 18 GENE PANEL
- 30 VA-PGX CYP2C19IM MOST RECENT
- 31 VA-PGX CYP2C19PM MOST RECENT

REMINDER ORDER CHECK RULES

- 32 VA-PGX PREVIOUS 16 GENE PANEL (R)
- 33 VA-PGX PREVIOUS 11 GENE PANEL (R)
- 34 VA-PGX PREVIOUS 8 GENE PANEL (R)
- 35 VA-PGX PREVIOUS 18 GENE PANEL (R)
- 36 VA-PGX CLOPIDOGREL CYP2C19IM (R)
- 37 VA-PGX CLOPIDOGREL CYP2C19PM (R)

REMINDER ORDER CHECK ITEMS GROUP

- 38 VA-PGX PREVIOUS 16 GENE PANEL (G)
- 39 VA-PGX PREVIOUS 11 GENE PANEL (G)
- 40 VA-PGX PREVIOUS 8 GENE PANEL (G)
- 41 VA-PGX PREVIOUS 18 GENE PANEL (G)
- 42 VA-PGX CLOPIDOGREL CYP2C19IM (G)
- 43 VA-PGX CLOPIDOGREL CYP2C19PM (G)

REMINDER DIALOG

- 44 VA-EL PHASER CONSENT ONLY PT AGREES

Pre-Install

None

Install Details

This update is being distributed as a web host file. The address for the host file is:
https://REDACTED /UPDATE_2_0_408.PRD

The file will be installed using Reminder Exchange, programmer access is not required.

Installation:

=====

This update can be loaded with users on the system. Installation will take less than 5 minutes.

Install Example

If your site is unable to use LWH functionality, use VistA mail to install the exchange file.

To Load the Web Host File. Navigate to Reminder exchange in Vista

```

+      + Next Screen  - Prev Screen  ?? More Actions
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                 LMM  Load MailMan Message
CMM  Create MailMan Message           LR   List Reminder Definitions
DFE  Delete Exchange File Entry       LWH  Load Web Host File
IFE  Install Exchange File Entry      RI   Reminder Definition Inquiry
IH   Installation History             RP   Repack
Select Action: Next Screen// LWH  Load Web Host File
Input the URL for the .prd file:
https://[REDACTED]/UPDATE_2_0_408.PRD

```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File
 At the **Input the url for the .prd file:** prompt, type the following web address:
https://REDACTED /UPDATE_2_0_408.PRD

```

https://[REDACTED]/UPDATE_2_0_408.PRD successfully loaded.

```

You should see a message at the top of your screen that the file successfully loaded

Search for and locate a reminder exchange entry named **UPDATE_2_0_408 VA-PGX (PHASER) CLOPIDOGREL/PREVIOUS RESULTS CROC**

+Item	Entry	Source	Date Packed
444	UPDATE_2_0_408 VA-PGX (PHASER) CLOPIDOGREL/PREVIOUS RESULTS CROC		
445	UPDATE_2_0_41 FUNCTIONAL AND PSYCHOSOCIAL TEMPLATE		
446	UPDATE_2_0_43 VA-SUICIDE RISK ASSESSMENT		
447	UPDATE_2_0_44 VA-MHRRTP TEMPLATES		
448	UPDATE_2_0_45 VA-MATERNITY		

```

+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                 LMM  Load MailMan Message
CMM  Create MailMan Message           LR   List Reminder Definitions
DFE  Delete Exchange File Entry       LWH  Load Web Host File
IFE  Install Exchange File Entry      RI   Reminder Definition Inquiry
IH   Installation History             RP   Repack
Select Action: Next Screen// IFE  Install Exchange File Entry
Enter a list or range of numbers (1-660): 444

```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry titled **UPDATE_2_0_408 VA-PGX (PHASER) CLOPIDOGREL/PREVIOUS RESULTS CROC** (in this example it is entry 444, it will vary by site). The Date Packed should read 01/25/2024.

Description:

The following Clinical Reminder items were selected for packing:

REMINDER ORDER CHECK ITEMS GROUP

VA-PGX CLOPIDOGREL CYP2C19PM (G)

VA-PGX CLOPIDOGREL CYP2C19IM (G)

VA-PGX PREVIOUS 18 GENE PANEL (G)

VA-PGX PREVIOUS 8 GENE PANEL (G)

VA-PGX PREVIOUS 11 GENE PANEL (G)

VA-PGX PREVIOUS 16 GENE PANEL (G)

REMINDER TERM

VA-REMINDER UPDATE 2 0 408

+ Enter ?? for more actions

IA Install all Components

IS Install Selected Component

Select Action: Next Screen// IA

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**

During the installation, you will see several prompts:

You will **Install (I)** all new components

Overwrite (O) any other existing components

Once you have installed the Order Check Groups and Rules, you will be prompted to install the dialog element VA-EL PHASER CONSENT ONLY PT AGREES. You do not need to install this piece and may type Q to quit this part of the install.

```
Packed reminder dialog: VA-EL PHASER CONSENT ONLY PT AGREES
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-EL PHASER CONSENT ONLY PT AGREES	element	X
		Finding: VA-PGX TESTING ACCEPTED (HEALTH FACTOR)		X
2		PXRM VISIT DATE FORCED TODAY	forced	X

```

+ Next Screen  - Prev Screen  ?? More Actions
DD Dialog Details  DT Dialog Text    IS Install Selected
DF Dialog Findings DU Dialog Usage   QU Quit
DS Dialog Summary  IA Install All
Select Action: Quit// Q

```

Type Q to quit the dialog install

```
Description:
The following Clinical Reminder items were selected for packing:
REMINDER ORDER CHECK ITEMS GROUP
  VA-PGX CLOPIDOGREL CYP2C19PM (G)
  VA-PGX CLOPIDOGREL CYP2C19IM (G)
  VA-PGX PREVIOUS 18 GENE PANEL (G)
  VA-PGX PREVIOUS 8 GENE PANEL (G)
  VA-PGX PREVIOUS 11 GENE PANEL (G)
  VA-PGX PREVIOUS 16 GENE PANEL (G)

REMINDER DIALOG
  VA-EL PHASER CONSENT ONLY PT AGREES
+ + Next Screen  - Prev Screen  ?? More Actions  >>>
IA Install all Components      IS Install Selected Component
Select Action: Next Screen// Q
```

Install is complete and you may type "Q" to Quit


Note: If your site is NOT a participating PHASER site, you do not need to proceed to the Post Install steps. No further action is needed. Go to VA PHASER SITE LISTING REDACTED to check your site's PHASER participation status.

Post Installation

1. For the 4 new Reminder Order Check Groups, the PHASER Panel orderable item(s) will need to be added to the below Order Check Item List. When entering the Orderable Items, if your standard is to have your station number as a prefix, then use that to search. If you do not know the name of your Orderable Item, please contact your LIM to obtain the name.

VA-PGX PREVIOUS 18 GENE PANEL (G)
VA-PGX PREVIOUS 8 GENE PANEL (G)
VA-PGX PREVIOUS 11 GENE PANEL (G)
VA-PGX PREVIOUS 16 GENE PANEL (G)

Example of how to add the Orderable Item(s)

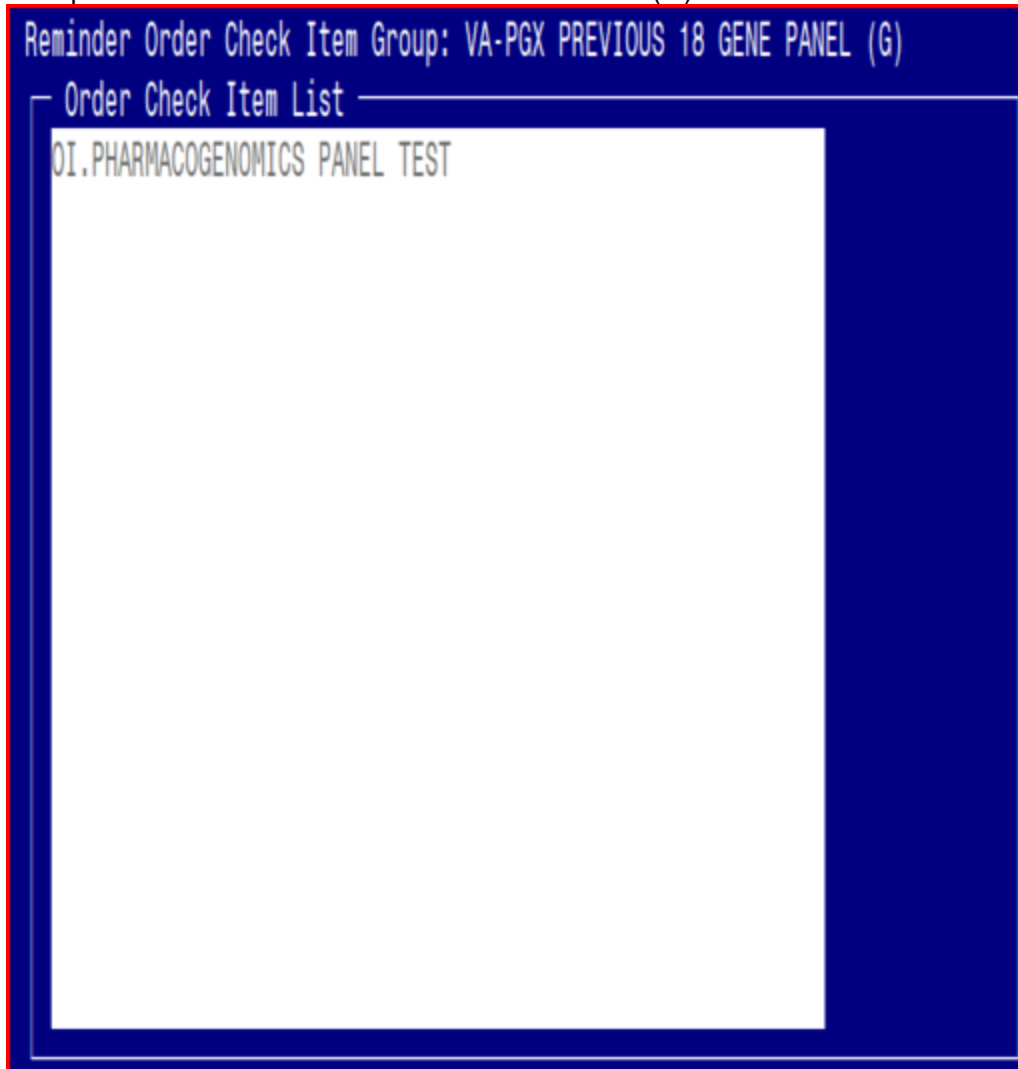

The screenshot shows a terminal window with a dark blue background and white text. At the top, it displays "Reminder Order Check Item Group: VA-PGX PREVIOUS 18 GENE PANEL (G)". Below this, there is a section titled "Order Check Item List" which is enclosed in a white rectangular box. Inside this box, the text "OI.PHARMACOGENOMICS PANEL TEST" is visible. The terminal window is framed by a red border.

Are you adding 'PHARMACOGENOMICS PANEL TEST' as a new ORDER CHECK ITEMS? No// **Y**

2. Once all Orderable Items have been added to the new Reminder Order Check groups, you will need to change the STATUS from INAC to PROD for all 6 Reminder Order Check Rules

VA-PGX CLOPIDOGREL CYP2C19PM (R)
VA-PGX CLOPIDOGREL CYP2C19IM (R)
VA-PGX PREVIOUS 18 GENE PANEL (R)
VA-PGX PREVIOUS 8 GENE PANEL (R)
VA-PGX PREVIOUS 11 GENE PANEL (R)
VA-PGX PREVIOUS 16 GENE PANEL (R)

RULE NAME: VA-PGX PREVIOUS 18 GENE PANEL (R)
DISPLAY NAME: Pharmacogenomics Warning - PGx 18-Gene Panel Results Available
STATUS: INAC