



Airborne Hazard and Open Burn Pit Initial Evaluation
Reminder Dialog

PXRM*2.0*39

INSTALLATION and SETUP GUIDE

November 2014

Product Development
Department of Veterans Affairs

Contents

| | |
|---|---------------|
| INTRODUCTION..... | 3 |
| PRE-INSTALLATION | 6 |
| Required Software for patch PXR..... | 6 |
| Related Documentation | 6 |
| Web Sites | 6 |
| INSTALLATION..... | 6 |
| 1. Retrieve the host file containing the patch..... | 7 |
| 2. Install the build first in a training or test account..... | 7 |
| 3. Load the distribution. | 7 |
| a. Verify Checksums in Transport Global | 7 |
| 4. Install the build. | 8 |
| 5. Install File Print..... | 8 |
| 6. Build File Print..... | 9 |
| 7. Post-installation routines | 9 |
| POST-INSTALLATION SETUP | 6 |
| APPENDIX A INSTALLATION EXAMPLE | 30 |
| APPENDIX B HEALTH FACTOR LIST..... | 32 |
| ACRONYMS..... | 35 |

Introduction

This patch releases the Airborne Hazard and Open Burn Pit Registry Initial Evaluation reminder dialog without any changes to routines, data dictionaries, or other package functions – “content” only. The reminder dialog is VA-AH-BPR INITIAL EVALUATION NOTE

Airborne Hazard and Open Burn Pit Registry Initial Evaluation

The Airborne Hazard and Burn Pit Registry Clinical Reminder Dialog Template was developed by the War-Related Illness and Injury Study Center (WRIISC) and Office of Public Health - Post-Deployment Health. The purpose of this template is to support the provider in a focused evaluation of a Veteran who requests an examination after completing the Airborne Hazards and Open Burn Pit Registry Self-Assessment. The template may also be useful for encounters where deployment health concerns are expressed regardless of whether the Veteran is participating in the registry. The template may be completed by a primary care provider, environmental health clinician, post-deployment champion, or a pulmonologist. This template is critical to improving the process and documentation of clinical care, and facilitating high quality care and programmatic quality improvement.

1. Recommendations

It is recommended that the reminder dialog first be installed into TEST account. After installation into TEST account the CAC/HIS will need to provide the appropriate stakeholder(s) access to TEST account so the dialog can be reviewed. The stakeholders should be familiar with the functionality and content of the dialog before implementing into Production.

This dialog should be linked to the progress note title, AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE. The CAC/HIS will need to create the note title and link it to the VHA Enterprise Standard Title, PERSIAN GULF REGISTRY E & M NOTE. [Instructions are included.](#)

2. Component Inventory

The VA-AH-BPR INITIAL EVALUATION NOTE Reminder Dialog contains:

- 1 dialog
- 3 branching terms
- 21 TIU objects
- 7 Health Factor categories
- 68 Health Factors
- 2 TIU template fields
- 1 Sponsor

2.1 Reminder Dialog Name:

VA-AH-BPR INITIAL EVALUATION NOTE

2.2 Terms:

The following three branching reminder terms are exported with the reminder. **These terms require local modification and mapping.** Instructions are included.

VA-AH-BPR AUDIT C DUE (BRANCHING LOGIC)
 VA-AH-BPR PTSD SCREEN DUE (BRANCHING LOGIC)
 VA-AH-BPR DEPRESSION SCREEN DUE (BRANCHING LOGIC)

2.3 Local Order Dialogs

The dialog has eight elements used for ordering appropriate interventions through the note. The “Additional findings” field in the elements will need to be mapped to the correct local order dialog (quick order, generic order, order set, or order menu). Instructions are included below.

2.4 Local Groups/Elements

The dialog has several “Local” groups and elements, so each site can configure the section as necessary. The local groups/elements begin with the namespace AH-BPR and are not prefixed with VA-.

A. The dialog has a group named AH-BPR MEDS LOCAL RECONCILIATION designated as local”, so each site can add local groups/elements to meet medication reconciliation documentation requirements. The group is active with installation. If it is not needed it must be inactivated, so it doesn’t display in the dialog.

B. All groups and elements that include a TIU object are designated “local” in case a site can’t map to their documentation methods. An example is the element AH-BPR PFT RESULTS. If a site captures PFT results in a method that proves impossible to map a TIU Object, then this element can be inactivated at the site, so that option does not display in CPRS. Another option is for the site to include information in the element directing the provider to the PFT results.

2.5 Patient data objects:

Twenty-one TIU/health summary objects/types will install with this patch; however, only twelve are new. **The new health summary types will need to be mapped to local items.** Instructions are included. The nine existing TIU objects should not require local modification, but are included below for reference.

NOTICE: The GMTSMGR security key is necessary for the staff responsible for mapping the health summary types.

| TIU Object Name | Health Summary Object Name | Health Summary Type Name |
|----------------------------------|----------------------------------|---------------------------|
| AH-BPR PULM CONSULT RESULTS;1;2Y | AH-BPR PULM CONSULT RESULTS;1;2Y | AH-BPR PULM CONSULT;1;2Y |
| AH-BPR ENT CONSULT RESULTS;1;2Y | AH-BPR ENT CONSULT RESULTS;1;2Y | AH-BPR ENT CONSULT;1;2Y |
| AH-BPR ABG;1;2Y | AH-BPR ABG;1;2Y | AH-BPR ABG;1;2Y |
| AH-BPR ECHO;1;2Y | AH-BPR ECHO;1;2Y | AH-BPR ECHO;1;2Y |
| AH-BPR CHEST CT;1;2Y | AH-BPR CHEST CT;1;2Y | AH-BPR CHEST CT;1;2Y |
| AH-BPR PFT RESULTS;1;2Y | AH-BPR PFT RESULTS;1;2Y | AH-BPR PFT RESULTS;1;2Y |
| AH-BPR CHEST 2 VIEWS;1;2Y | AH-BPR CHEST 2 VIEWS;1;2Y | AH-BPR CHEST 2 VIEWS;1;2Y |
| AH-BPR CBC W/DIFF;1;2Y | AH-BPR CBC W/DIFF;1;2Y | AH-BPR CBC W/DIFF;1;2Y |
| AH-BPR PULSE OX | AH-BPR PULSE OX | AH-BPR PULSE OX |

| | | |
|--------------------------|--------------------------|--------------------------|
| AH-BPR AUDIT-C | AH-BPR AUDIT-C | AH-BPR AUDIT-C |
| AH-BPR PC-PTSD PCL | AH-BPR PC-PTSD PCL | AH-BPR PC-PTSD PCL |
| AH-BPR DEPRESSION SCREEN | AH-BPR DEPRESSION SCREEN | AH-BPR DEPRESSION SCREEN |

New TIU objects:

Existing TIU objects: These should require no local configuration.

PATIENT WEIGHT
PATIENT HEIGHT
PAIN
BLOOD PRESSURE
RESPIRATION
PULSE
TEMPERATURE
ACTIVE MEDICATIONS
VA-WRIISC ACTIVE PROBLEMS

2.6 Health Factors:

There are 7 new health factors categories containing 68 new health factors.

Three health factors should already exist; CURRENT SMOKER, PREVIOUS SMOKER, and LIFETIME NON-SMOKER.

Pre-Installation

*Required Software for PXR*2*39*

| Package/Patch | Namespace | Version | Comments |
|--------------------------|-----------|---------|---------------|
| Clinical Reminders | PXR | 2.0 | Fully patched |
| Health Summary | GMTS | 2.7 | Fully patched |
| Kernel | XU | 8.0 | Fully patched |
| NATIONAL DRUG FILE | PSN | 4.0 | Fully patched |
| Pharmacy Data Management | PSS | 1.0 | Fully patched |
| Outpatient Pharmacy | PSO | 7.0 | Fully patched |
| VA FileMan | DI | 22 | Fully patched |

Related Documentation

| Documentation | Documentation File name |
|------------------------------|-------------------------|
| Installation and Setup Guide | PXR_2_0_39_IG.PDF |
| | |

Web Sites

| Site | URL | Description |
|---------------------------------------|---|--|
| National Clinical Reminders site | http://vista.med.va.gov/reminders | Contains manuals, PowerPoint presentations, and other information about Clinical Reminders |
| National Clinical Reminders Committee | http://vaww.portal.va.gov/sites/ncrc/public/default.aspx | This committee directs the development of new and revised national reminders |
| VistA Document Library | http://www.va.gov/vdl/ | Contains manuals for Clinical Reminders and |

Installation

This patch can be installed with users on the system, but it should be done during non-peak hours. Estimated Installation Time: 10-15 minutes

1. Retrieve the host file from one of the following locations (with the ASCII file type):

| | | |
|----------------|----------|----------|
| Albany | REDACTED | REDACTED |
| Hines | REDACTED | REDACTED |
| Salt Lake City | REDACTED | REDACTED |

2. Install the patch first in a training or test account.

Installing in a non-production environment will give you time to get familiar with new functionality and complete the setup for reminders and dialogs prior to installing the software in production.

3. Load the distribution.

In programmer mode type, D ^XUP, select the Kernel Installation & Distribution System menu (XPD MAIN), then the Installation option, and then the option LOAD a Distribution. Enter your directory name. KID at the Host File prompt.

Example

```
Select Installation <TEST ACCOUNT> Option: 1  Load a Distribution
Enter a Host File: <your directory>PXRM_2_0_39.KID

KIDS Distribution saved on Jan 21, 2014@10:07:41
Comment: AIRBORNE HAZARD

This Distribution contains Transport Globals for the following Package(s):
  PXRM*2.0*39
Distribution OK!

Want to Continue with Load? YES//
Loading Distribution...

  PXRM*2.0*39
Use INSTALL NAME: PXRM*2.0*39 to install this Distribution.
```

From the Installation menu, you may elect to use the following options:

4. Backup a Transport Global

This option will create a backup message of any routines exported with the patch. It will NOT back up any other changes such as DDs or templates.

a. Compare Transport Global to Current System

This option will allow you to view all changes that will be made when the patch is installed. It compares all components of the patch (routines, DDs, templates, etc.).

5. Install the build.

From the Installation menu on the Kernel Installation and Distribution System (KIDS) menu, run the option Install Package(s). Select the build PXRM*2.0*36 and proceed with

the install. If you have problems with the installation, log a Remedy ticket and/or call the National Help Desk to report the problem.

```
Select Installation & Distribution System Option: Installation
Select Installation Option: INSTALL PACKAGE(S)
Select INSTALL NAME: PXRM*2.0*39
```

Answer "NO" to the following prompt:

Want KIDS to INHIBIT LOGONs during install? NO// **NO**

Want to DISABLE Scheduled Options, Menu Options, and Protocols? NO// **NO**

Installation Example

See [Appendix A](#).

6. Install File Print

Use the KIDS Install File Print option to print out the results of the installation process. You can select the multi-package build or any of the individual builds included in the multi-package build.

```
Select Utilities Option: Install File Print
Select INSTALL NAME: PXRM*2.0*39
```


7. Build File Print

Use the KIDS Build File Print option to print out the build components.

```
Select Utilities Option:  Build File Print
Select BUILD NAME: PXRMP39I
DEVICE: HOME//
```

8. Post-installation routines

After successful installation, the following init routines may be deleted:

```
PXRMP39E
PXRMP39I
```

Post-Installation Setup

1. Post Install Details

The CAC/HIS has several steps to complete after the patch is installed before the dialog should be deployed. Detailed instructions for each step are below. PXRMP39I_CAC_CKLIST.PDF is available as part of the documentation for this patch. This is a checklist that will help ensure you have done everything correctly.

2. Configure three reminder terms for branching logic

Three reminder terms need to be modified to local practice before deployment.

Before beginning, identify the local reminder definitions in use for Depression screening (PHQ-2), PTSD screening, and Alcohol Use screening (AUDIT-C). All have been deployed nationally; however, sites may have created local versions.

```
CF      Reminder Computed Finding... RM
        Reminder Definition Management ...
SM      Reminder Sponsor TXM
        Reminder Taxonomy ... TRM
        Reminder Term Management ...
LM      Reminder Location List RX
        Reminder Exchange
RT      Reminder Test
OS      Other Supporting Menus ...
INFO    Reminder Information... DM
        Reminder Dialog Management ...
CP      CPRS Reminder RP
        Reminder Reports ...
        Select Reminder Managers Menu Option: TRM
```

TL List Reminder Terms
TI Inquire about ■TE Add/Edit
Reminder Term
TC Copy Reminder Term

Select Reminder Term Management <TEST ACCOUNT> Option: E

Select Reminder Term: VA-AH-BPR AUDIT C DUE (BRANCHING)

Reminder Term has no findings!

Select Finding: CF.VA-REMINDER DEFINITION

Searching for a REMINDER COMPUTED FINDING, (pointed-to ■ITEM)

Searching for a REMINDER ■ VA-REMINDER
DEFINITION NATIONAL
...OK? Yes// Yes

Are you adding 'VA-REMINDER DEFINITION' as a ■(the 1ST for this REMINDER TERM)?
No// Yes

Computed Finding Description: *Description removed from this document*

Editing Finding Number: 1

FINDING ITEM: VA-REMINDER BEGINNING

DATE/TIME:

ENDING DATE/TIME:

OCCURRENCE COUNT:

COMPUTED FINDING PARAMETER: Enter exact name of your site's definition use
for alcohol use screening using the AUDIT-C

CONDITION: I ■ CONDITION CASE ■
USE STATUS/COND IN SEARCH:

Choose from:

CF VA-REMINDER DEFINITION

Finding # 1

Select Finding:

Input your edit ■dit?

NO//

Select Reminder Term: VA-AH-BPR DEPRESSION SCREEN DE ■ LOGIC)

Reminder Term has no findings!

Select Finding: CF.VA-REMINDER DEFINITION

Searching for a REMINDER COMPUTED FINDING, (pointed-to ■ITEM)

Searching for a REMINDER ■ VA-REMINDER
DEFINITION NATIONAL
...OK? Yes// Yes

Are you adding 'VA-REMINDER DEFINITION' as a ■ (the 1ST for this REMINDER TERM)?
No// Yes

Computed Finding Description: *Description removed from this document*

Editing Finding Number: 1

FINDING ITEM: VA-REMINDER BEGINNING DATE/TIME:

ENDING DATE/TIME:

OCCURRENCE COUNT:

COMPUTED FINDING PARAMETER: *Enter exact name of your site's definition use
for depression screening using the PHQ-2 or PHQ9
depending on which one your site uses for
depression screening*

CONDITION: I ■ CONDITION CASE ■

USE STATUS/COND IN SEARCH:

Choose from:

CF VA-REMINDER DEFINITION

Finding # 1

Select Finding:

Input your edit ■ edit?

NO//

Select Reminder Term: VA-AH-BPR PTSD SCREEN DUE ■ ■

Reminder Term has no findings!

Select Finding: CF.VA-REMINDER DEFINITION

Searching for a REMINDER COMPUTED FINDING, (pointed-to ■ ITEM)

Searching for a REMINDER ■ VA-REMINDER
DEFINITION NATIONAL
...OK? Yes// Yes

Are you adding 'VA-REMINDER DEFINITION' as a ■ ■ (the 1ST for this
REMINDER TERM)? No// Yes

Computed Finding Description: *Description removed from this document*

Editing Finding Number: 1

FINDING ITEM: VA-REMINDER BEGINNING DATE/TIME:

ENDING DATE/TIME:

OCCURRENCE COUNT:

COMPUTED FINDING PARAMETER: *Enter exact name of your site's definition use for PTSD screening.*

CONDITION: ☐ CONDITION CASE ☐

USE STATUS/COND IN SEARCH:

Choose from:

CF VA-REMINDER DEFINITION

Finding # 1

Select Finding:

Input your edit ☐ edit?

NO//

3. Configure the 12 new TIU objects/Health Summary objects/types

NOTICE: The GMTSMGR security key is necessary for the staff responsible for mapping the health summary types.

If the following message is displayed when accessing the health summary types you will need to have the GMTSMGR security key added to your profile before beginning:

“Alteration of this summary report is restricted to its owner.
See the Clinical Coordinator if you need additional help.”

All of the TIU objects are already embedded in the dialog, but will not display information until the health summary type for each one is mapped to local items. **Before beginning, identify the following local information for each object/type.** Detailed instructions for editing the health summary types follow.

1. HEALTH SUMMARY TYPE: AH-BPR PULM CONSULT;1;2Y
 - a. What note title is used to result Pulmonary ☐
2. HEALTH SUMMARY TYPE: AH-BPR ENT CONSULT;1;2Y
 - a. What note title is used to result ENT ☐
3. HEALTH SUMMARY TYPE: AH-BPR ABG;1;2Y
 - a. What are the local lab test name(s) for an ☐
4. HEALTH SUMMARY TYPE: AH-BPR ECHO;1;2Y
 - a. How are echo results stored locally?
 - b. If a note title, then what is the title?
5. HEALTH SUMMARY TYPE: AH-BPR CHEST CT;1;2Y
 - a. What are the local imaging test names ☐T THORAX with and/or without contrast?

6. HEALTH SUMMARY TYPE: AH-BPR PFT RESULTS;1;2Y
 - a. How are PFT results stored locally?
 - b. If a note title, then what is the title?
7. HEALTH SUMMARY TYPE: AH-BPR CHEST 2 VIEWS;1;2Y
 - a. What are the local imaging test names for LAT x-rays?
8. HEALTH SUMMARY TYPE: AH-BPR CBC W/DIFF;1;2Y
 - a. What are the local lab test name(s) for a CBC?
9. HEALTH SUMMARY TYPE: AH-BPR PULSE OX
 - a. Nothing needed
10. HEALTH SUMMARY TYPE: AH-BPR AUDIT-C
 - a. Nothing needed
11. HEALTH SUMMARY TYPE: AH-BPR PC-PTSD PCL
 - a. Nothing needed
12. HEALTH SUMMARY TYPE: AH-BPR DEPRESSION SCREEN
 - a. What is the local name of the reminder definition?

Health Summary Maintenance Menu

- 1 Disable/Enable Health Summary Component
- 2 Create/Modify Health Summary Components
- 3 Edit Ad Hoc Health Summary Type
- 4 Rebuild Ad Hoc Health Summary Type
- 5 Resequence a Health Summary Type
- 6 Create/Modify Health Summary Type**
- 7 Edit Health Summary Site Parameters
- 8 Health Summary Objects Menu ...
- 9 CPRS Reports Tab 'Health Summary Types List' M
- 10 CPRS Health Summary Display/Edit Site Defaults.

Select Health Summary Maintenance Menu Option: **6**

Select Health Summary Type: **AH-BPR PULM CONSULT;1;2Y**

AH-Bpr Pulm Consult;1;2y (Pulmonary Consult Results, B Owner **REDACTED**)

OK? **YES**

WARNING: You are about to edit a **Health Summary Type** that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue?NO// **YES**

NAME: AH-BPR PULM CONSULT;1;2YReplace
 TITLE: Pulmonary Consult ~~Re~~ Replace
 SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
 SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ☒ Summary Type? YES//
NO Select COMPONENT: **SPN** PROGRESS NOTES SELECTED SPN
 SUMMARY ~~ORDER~~ 5
 OCCURRENCE LIMIT: **1**
 TIME LIMIT: **2Y**
 HEADER NAME: Selected Prog ~~Notes~~ **Consult** No selection items

chosen.

Select new items one at a time in the ~~summary~~ you want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of progress note title. Can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR PULM CONSULT;1;2Y
Title: Pulmonary ☒ Owner:
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord  Component  Name  Max  Hos  ICD  Pro  CPT
Occ  Time b Text  Nar  Mod  Selection
-----
SPN   5      Pulmonary  consult  1   2Y
                                           PULMONARY ☒ RESULT NOTE
  
```

Select Health Summary Type: **AH-BPR ENT CONSULT;1;2Y**

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// **YES**
 NAME: AH-BPR ENT CONSULT;1;2YReplace
 TITLE: ENT Consult Results//
 SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
 SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ☒ Summary Type? YES// **NO**
 Select COMPONENT: **SPN** PROGRESS NOTES SELECTED SPN
 SUMMARY ORDER: 5// 5

OCCURRENCE ☐ TIME

LIMIT: 2Y

HEADER NAME: Selected Prog ☐ Consult No selection

items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of progress note title. Can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR ENT CONSULT;1;2Y
Title: ENT Consult Owner:
SUPPRESS PRINT OF COMPONENTS ☐ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord  Component  Name  Max  Hos  ICD  Pro  CPT
      Occ  Time b Text  Nar  Mod  Selection
-----
SPN   5      ENT consult info  1   2Y
                                           ENT CONSULT ☐ NOTE
```

Select Health Summary Type: AH-BPR ABG;1;2Y

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// ☒

NAME: AH-BPR ABG;1;2Y//

TITLE: ABG//

SUPPRESS PRINT OF COMPONENTS ☐ yes// SUPPRESS

SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ☐ Summary Type? YES//

☒ NO Select COMPONENT: SLT LAB TESTS SELECTED SLT

SUMMARY ~~01P~~5

OCCURRENCE LIMIT: 1

TIME LIMIT: 2Y

HEADER NAME: Lab Tests ☐ ABG No selection

items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select up to 99 items.

Select SELECTION ITEM: *Enter exact name of local ABG lab test. Can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

Type Name: AH-BPR ABG;1;2Y
 Title: ~~A~~
 Owner:
 SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
 SENSITIVE PRINT DATA: NO SSN

| Abb | Ord | Component | Name | Max Occ | Hos Time | ICD Text | Pro Nar | CPT Mod | Selection |
|-----|-----|-----------|------|------------|-------------|-------------|------------|------------|-----------|
| SLT | 5 | Last | ABG | 1 | 2Y | | | | ABG |

Select Health Summary Type: AH-BPR ECHO;1;2Y

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// ☒
 NAME: AH-BPR ECHO;1;2Y//
 TITLE: Echo//
 SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
 SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ☒ Summary Type? YES// NO
 Select COMPONENT: SPN PROGRESS NOTES SELECTED SPN
If local results are not stored with a progress note title, then a different component should be chosen.

SUMMARY ~~ORD~~ 5
 OCCURRENCE LIMIT: 1
 TIME LIMIT: 2Y
 HEADER NAME: Selected Prog ~~SPN~~ Echo No selection

items chosen.

Select new items one at a time in the ~~spn~~ you want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of progress note title. Can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

Type Name: AH-BPR ECHO;1;2Y
 Title: ~~E~~
 Owner:
 SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
 SENSITIVE PRINT DATA: NO SSN

| Abb | Ord | Component | Name | Max Occ | Hos Time | ICD Loc | Pro Text | CPT Nar | Selection |
|-----|-----|-----------|------|------------|-------------|------------|-------------|------------|-----------|
| SPN | 5 | Last | Echo | 1 | 2Y | | | | CP ECHO |

Select Health Summary Type: AH-BPR CHEST CT;1;2Y

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue?NO// YES NAME:

AH-BPR CHEST ~~7~~ Replace

TITLE: Chest CT//

SUPPRESS PRINT OF COMPONENTS ~~1~~ yes// SUPPRESS

SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ~~1~~ Summary Type? YES//

NO Select COMPONENT: SII IMAGING IMPRESSION SELECTED SII

SUMMARY ~~01P35~~

OCCURRENCE LIMIT: 1

TIME LIMIT: 2Y

HEADER NAME: Sel Image Impression Replace ... With ~~1~~ CT No

selection items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of Thorax Chest ~~01P35~~ and/or without contrast imaging tests. You can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR CHEST CT;1;2Y
Title: 1
Owner:
SUPPRESS PRINT OF COMPONENTS 1 yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord      Component Name      Max      Hos  ICD  Pro  CPT
      Occ      Time b Text  Nar  Mod  Selection
-----
SII   5      Last Chest CT      5      2Y
                                           CT THORAX 1 CT THORAX 1
                                           CT THORAX W/O CONT
```

Select Health Summary Type: AH-BPR PFT RESULTS;1;2Y

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue?NO// YES

NAME: AH-BPR PFT RESULTS;1;2Y Replace

TITLE: PFT Results//

SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: yes//

SUPPRESS SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing Summary Type? YES//

NO Select COMPONENT: **SPN** PROGRESS NOTES SELECTED SPN

If local results are not stored with a progress note title, then a different component should be chosen.

SUMMARY ~~0005~~

OCCURRENCE LIMIT: **1**

TIME LIMIT: **2Y**

HEADER NAME: Selected Prog **Last PFT** No

selection items chosen.

Select new items one at a time in the ~~summary~~ you want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of progress note title. Can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR PFT RESULTS;1;2Y
Title: PFT Owner:
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord  Component  Name  Max  Hos  ICD  Pro  CPT
Occ  Time b Text  Nar  Mod  Selection
-----
SPN   5    Last PFT      1   2Y
                                     CP PFT
```

Select Health Summary Type: **AH-BPR CHEST 2 VIEWS;1;2Y**

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// **YES**

NAME: AH-BPR CHEST 2 VIEWS;1;2Y Replace

TITLE: Chest x-ray 2 views//

SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS

SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing Summary Type? YES//

NO Select COMPONENT: **SII** IMAGING IMPRESSION SELECTED SII

SUMMARY ~~0005~~

OCCURRENCE LIMIT: **1**

TIME LIMIT: 2Y

HEADER NAME: Sel Image Impression Replace ... With Last x-ray No

selection items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select any number of items.

Select SELECTION ITEM: *Enter exact name of Chest x-ray PA & AT imaging tests. You can enter more than one at the next prompt if needed.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR CHEST 2 VIEWS;1;2Y
Title: Chest 2 views Owner:
SUPPRESS PRINT OF COMPONENTS ■ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord  Component Name      Max      Hos  ICD  Pro  CPT
      Occ      Time b Text   Nar  Mod  Selection
-----
SII   5      Last Chest x-ray      1      2Y
                                           CHEST 2 VIEWS CHEST 4
                                           VIEWS
```

Select Health Summary Type: AH-BPR CBC W/DIFF;1;2Y

WARNING: You are about to edit a ~~lab~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// YES

NAME: AH-BPR CBC W/DIFF;1;2Y Replace

TITLE: CBC w/diff//

SUPPRESS PRINT OF COMPONENTS ■ yes// SUPPRESS

SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing Summary Type? YES//

NO Select COMPONENT: : SLT LAB TESTS SELECTED SLT

SUMMARY 0155

OCCURRENCE LIMIT: 1

TIME LIMIT: 2Y

HEADER NAME: Lab Tests ~~Selected~~ CBC W/DIFF No

selection items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select up to 99 items.

Select SELECTION ITEM: *Enter exact name of local CBC w/diff tests. You may need to enter the individual lab tests separately.*

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR CBC W/DIFF;1;2Y
Title: CBC Owner:
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

```

| Abb | Ord | Component | Max Occ | Time | Hos | ICD Text | Pro Nar | CPT Mod | Selection |
|-----|-----|-----------|------------|------|-----|-------------|------------|------------|---|
| SLT | 5 | Last CBC | 1 | 2Y | | | | | WBC RBC HGB HCT MCV MCH MCHC PLT RDW SEGS BANDS LYMPHS MONOS EOSINO BASO PLT |

Select Health Summary Type: **AH-BPR PULSE OX**

WARNING: You are about to edit a ~~HEALTH~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do you want to continue? NO// ☒
NAME: AH-BPR PULSE OX//
TITLE: Pulse Oximetry//
SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ~~HEALTH~~ Summary Type? YES//
NO Select COMPONENT: **VITAL SIGNS SELECTED**

SUMMARY ~~OBJECT~~ 5
OCCURRENCE LIMIT: **1**
TIME LIMIT:
HEADER NAME: Vital Signs Selected ~~HEALTH~~. With **P02**: No

selection items chosen.

Select new items one at a time in the ~~summary~~ want them displayed. You may select any number of items.

Select SELECTION ITEM: **PULSE OXIMETRY**

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR PULSE OX
Title: Pulse Oximeter
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

```

| Abb | Ord | Component | Name | Max Occ | Hos Time | ICD b Text | Pro Nar | CPT Mod | Selection |
|-----|-----|-----------|------|------------|-------------|---------------|------------|------------|----------------|
| SVS | 5 | P02: | | 1 | | | | | PULSE OXIMETRY |

Select Health Summary Type: **AH-BPR AUDIT-C**

WARNING: You are about to edit a **Health Summary Type** that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do you want to continue? NO// ☒
NAME: AH-BPR AUDIT-C//
TITLE: AUDIT-C//
SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
SENSITIVE PRINT DATA: NO SSN//

Select COMPONENT: **MHAS** ~~MHA~~ MHAS
SUMMARY ORDER: 5// 5
HEADER NAME: MHA Score// **AUDIT-C:**

Select SELECTION ITEM: **AUDC**

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR AUDIT-C
Title: Owner
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

```

| Abb | Ord | Component | Name | Max Occ | Hos Time | ICD b Text | Pro Nar | CPT Mod | Selection |
|------|-----|-----------|------|------------|-------------|---------------|------------|------------|-----------|
| MHAS | 5 | AUDIT-C: | | 1 | | | | | AUDC |

Select Health Summary Type: **AH-BPR DEPRESSION SCREEN**

WARNING: You are about to edit a **Health Summary Type** that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do you want to continue? NO// **YES**
NAME: AH-BPR DEPRESSION SCREEN Replace
TITLE: Depression Screen//

SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
SENSITIVE PRINT DATA: NO SSN//

Do you wish to copy COMPONENTS from an existing ☒ Summary Type? YES//
NO Select COMPONENT: **CF** CLINICAL REMINDERS FINDINGS CF
SUMMARY ORDER: 5// 5
HEADER NAME: Reminders Findings// **PT Screen** No selection items

chosen.

Select new items one at a time in the ~~summary~~ you want them displayed. You may select any number of items.

Select SELECTION ITEM: **VA-DEPRESSION SCREENING**

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR DEPRESSION SCREEN
Title: Depression Owner:
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

Abb  Ord  Component Name      Max      Hos  ICD  Pro  CPT
      Occ      Time b Text   Nar  Mod Sh
-----
CF    5      Reminders Findings
                                           VA-DEPRESSION SCREENING
```

Select Health Summary Type: ☒ **PC-PTSD PCL**

WARNING: You are about to edit a ~~Health~~ Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do you want to continue? NO// ☒
NAME: AH-BPR PC-PTSD PCL//
TITLE: PTSD Screen:
SUPPRESS PRINT OF COMPONENTS ☒ yes// SUPPRESS
SENSITIVE PRINT DATA: NO SSN//

Select COMPONENT: **MHAS** ~~MHA~~
MHAS SUMMARY ORDER:

5// 5
HEADER NAME: MHA Score// **PTSD Screen:**

Select SELECTION ITEM: **PC PTSD**

The Health Summary Type structure should look similar to this:

```

Type Name: AH-BPR PC-PTSD PCL
Title:
Owner: FAHNER,JEFF
SUPPRESS PRINT OF COMPONENTS ☒ yes SUPPRESS
SENSITIVE PRINT DATA: NO SSN

```

| Abb | Ord | Component | Max Occ | Hos Time | ICD Text | Pro Nar | CPT Mod |
|------|-----|-----------|------------|-------------|-------------|------------|------------|
| MHAS | 5 | MHA Score | 2 | | | | |

PC ☒ CLC PCLM

4. Add appropriate local order dialogs to order section

Eight local orders need to be added to eight separate elements in the dialog. . **Before beginning, identify the following local information for each element.** Detailed instructions for editing the health summary types follow.

1. Identify local order dialog name for an outpatient Chest X-ray PA & LAT.
2. Identify local order dialog name for an outpatient Pulmonary Function Test.
3. Identify local order dialog name for an outpatient Chest CT w & w/o contrast.
4. Identify local order dialog name for an outpatient echo.
5. Identify local order dialog name for an outpatient ABG.
6. Identify local order dialog name for an outpatient CBC w/diff.
7. Identify local order dialog name for an outpatient ENT Consult.
8. Identify local order dialog name for an outpatient Pulmonary Consult

From the Reminder Manager's Menu

```

CF      Reminder Computed Findings... RM
        Reminder Definition Management ...
SM      Reminder Sponsor ☒ TXM
        Reminder Taxonomy ☒ ... TRM
        Reminder Term Management ...
LM      Reminder Location List ☒ RX
        Reminder Exchange
RT      Reminder Test
OS      Other Supporting Menus ...
INFO    Reminder Information ☒ ... DM
        Reminder Dialog Management ...
CP      CPRS Reminder ☒ RP
        Reminder Reports ...
MST     Reminders MST Synchronization ☒ ... PL      Reminder
Patient List Menu ...
PAR     Reminder Parameters ...
ROC     Reminder Order ☒ ... XM
        Reminder Extract Menu ...
GEC     GEC Referral Report

```

Select Reminder Managers Menu <TEST ACCOUNT> Option: **D**

DP Dialog Reminders
 DR Dialog Reports ...
 IA Inactive Codes Mail Message

Select Reminder Dialog Management Option: **DI** REMINDER VIEW

(ALL REMINDERS BY NAME)

| Item | Reminder Name | Linked Dialogs | Dialog Status |
|------|----------------------|----------------|---------------|
| 1 | ACUTE MYOCARDIAL | | |
| 2 | AGP TEST | AGP TEST ALL | |
| 3 | CF ALLERGY TEST | | |
| 4 | CHEST XRAY | | |
| 5 | CHOLESTEROL SCREEN | | |
| 6 | CRITICAL DRUG | | |
| 7 | DEPRESSION SCREEN | | |
| 8 | FAHNER TEST REMINDER | | |
| 9 | LONG TERM STEROID | | |
| 10 | MAMMOGRAM | | |
| 11 | MAMMOGRAMS | | |
| 12 | MAMMOGRAPHY | | |
| 13 | PAP SMEAR | | |
| 14 | PATCH 11 LOCATION | | |
| 15 | PPD | | |
| 16 | PPD1 | | |

+ Enter? for more >>>

AR All reminders LR Linked
 Quit CV Change View QU
 Name/Print Name RN

Select Item: Next Select one of the

following:

D Reminder Dialogs
 E Dialog Elements
 F Forced Values
 G Dialog Groups
 P Additional Prompts
 R Reminders
 RG Result Group (Mental Health)
 RE Result Element (Mental Health)

TYPE OF VIEW: R// **D**

DIALOG VIEW (REMINDER DIALOGS - SOURCE REMINDER NAME)

| Item | Reminder Dialog Name | Source | Status |
|-----------|------------------------|---------------------------------------|--------|
| 1/21/2014 | Clinical Reminders PXR | M*2.0*39 Installation and Setup Guide | 24 |

Current dialog element/group name: VA-AH-BPR XRAY- PA AND LATERAL
Used by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: AH-BPR ORDER CHEST XRAY //

Select ADDITIONAL FINDINGS: *Q.Enter exact name of local order dialog*

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION

+Item Seq. Dialog Summary

| | | |
|-----|-----------|--|
| 184 | 35.25.5 | Element: VA-AH-BPR ASSESSMENT AND PLAN ADDL INFO |
| 185 | 35.25.10 | Element: VA-AH-BPR WORD PROCESSING 5 LINES |
| 186 | 35.30 | Element: VA-AH-BPR SPACE (NOTE ONLY) |
| 187 | 35.35 | Group: VA-AH-BPR ORDERS |
| 188 | 35.35.2 | Group: VA-AH-BPR DIAGNOSTIC EVAL DISPLAY ONLY |
| 189 | 35.35.2.5 | Group: VA-AH-BPR DIAGNOSTIC EVAL DISPLAY ONLY |
| 190 | 35.35.5 | Element: VA-AH-BPR ORDER CHEST XRAY- PA AND LATERAL |
| 191 | 35.35.10 | Element: VA-AH-BPR ORDER FULL PULMONARY FUNCTION TESTS |

+ + Next Screen - Prev Screen ?? More Actions

>>>

| | | | | | |
|-----|--------------------|----|--------------------|-----|---------------|
| ADD | Add Element/Group | DS | Dialog | INQ | Inquiry/Print |
| CO | Copy Dialog | DO | Dialog Overview | QU | Quit |
| DD | Detailed Display | DT | Dialog Text | | |
| DP | Progress Note Text | ED | Edit/Delete Dialog | | |

Select Item: Next Screen// *Find VA-AH-BPR ORDER FULL PULMONARY FUNCTION TESTS*

and select the appropriate item number to open the element

Current dialog element/group name: VA-AH-BPR PULMONARY FUNCTION TESTS

Used by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: AH-BPR ORDER PFT //

Select ADDITIONAL FINDINGS: *Q.Enter exact name of local order dialog*

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION

+Item Seq. Dialog Summary

| | | |
|-----|----------|---|
| 192 | 35.35.15 | Element: VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST |
| 193 | 35.35.20 | Element: VA-AH-BPR ORDER ECHOCARDIOGRAM |
| 194 | 35.35.25 | Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS |
| 195 | 35.35.30 | Element: VA-AH-BPR ORDER CBC W/DIFF |

196 35.35.35 Element: VA-AH-BPR ORDER ENT CONSULT

197 35.35.40 Element: VA-AH-BPR ORDER PULMONARY CONSULT

198 40 Element: VA-AH-BPR SPACE (DIALOG & NOTE)

```

+          + Next Screen - Prev Screen ?? More Actions
>>>
ADD Add Element/Group DS Dialog INQ
    Inquiry/Print      CO      Copy Dialog
    DO      Dialog Overview    QU      Quit
DD Detailed Display DT Dialog Text
DP Progress Note Text ED Edit/Delete Dialog
Select Item: Next Screen// Find VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST
and select the appropriate item number to open the element

```

Current dialog element/group name: VA-AH-BPR RESOLUTION CT CHEST
 Used by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER CHEST CT//
 Select ADDITIONAL FINDINGS: Q. Enter exact name of local order dialog

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION
 +Item Seq. Dialog Summary

192 35.35.15 Element: VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST

193 35.35.20 Element: VA-AH-BPR ORDER ECHOCARDIOGRAM

194 35.35.25 Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS

195 35.35.30 Element: VA-AH-BPR ORDER CBC W/DIFF

196 35.35.35 Element: VA-AH-BPR ORDER ENT CONSULT

197 35.35.40 Element: VA-AH-BPR ORDER PULMONARY CONSULT

198 40 Element: VA-AH-BPR SPACE (DIALOG & NOTE)

```

+          + Next Screen - Prev Screen ?? More Actions
>>>
ADD Add Element/Group DS Dialog INQ
    Inquiry/Print      CO      Copy Dialog
    DO      Dialog Overview    QU      Quit
DD Detailed Display DT Dialog Text
DP Progress Note Text ED Edit/Delete Dialog
Select Item: Next Screen// Find VA-AH-BPR ORDER ECHOCARDIOGRAM
and select the appropriate item number to open the element

```

Current dialog element/group name: VA-AH-BPR ECHOCARDIOGRAM
 Used by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER ECHO //

Select ADDITIONAL FINDINGS: *Q.Enter exact name of local order dialog*

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION ■

+Item Seq. Dialog Summary

| | | |
|-----|----------|--|
| 192 | 35.35.15 | Element: VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST |
| 193 | 35.35.20 | Element: VA-AH-BPR ORDER ECHOCARDIOGRAM |
| 194 | 35.35.25 | Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS |
| 195 | 35.35.30 | Element: VA-AH-BPR ORDER CBC W/DIFF |
| 196 | 35.35.35 | Element: VA-AH-BPR ORDER ENT CONSULT |
| 197 | 35.35.40 | Element: VA-AH-BPR ORDER PULMONARY CONSULT |
| 198 | 40 | Element: VA-AH-BPR SPACE (DIALOG & NOTE) |

| + | + Next ■ | - Prev | Screen ?? ■ | Actions |
|-----|---------------|--------|----------------|-------------|
| >>> | | | | |
| ADD | Add Element | DS | Dialog Summary | INQ Inquiry |
| CO | Copy Dialog | DO | Dialog Only | QU Quit |
| DD | Detailed ■ | DT | Dialog Text | |
| DP | Progress Note | ED | Edit/Delete ■ | |

Select Item: Next Screen// *Find VA-AH-BPR ORDER ARTERIAL BLOOD GAS and select the appropriate item number to open the element*

Current dialog element/group name: ~~VA-AH-BPR~~ ARTERIAL BLOOD GAS Used by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER ABG //

Select ADDITIONAL FINDINGS: *Q.Enter exact name of local order dialog*

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION ■

+Item Seq. Dialog Summary

| | | |
|-----|----------|--|
| 192 | 35.35.15 | Element: VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST |
| 193 | 35.35.20 | Element: VA-AH-BPR ORDER ECHOCARDIOGRAM |
| 194 | 35.35.25 | Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS |
| 195 | 35.35.30 | Element: VA-AH-BPR ORDER CBC W/DIFF |
| 196 | 35.35.35 | Element: VA-AH-BPR ORDER ENT CONSULT |
| 197 | 35.35.40 | Element: VA-AH-BPR ORDER PULMONARY CONSULT |

198 40 Element: VA-AH-BPR SPACE (DIALOG & NOTE)

```
+      + Next Screen - Prev Screen ?? More Actions
>>>
ADD  Add Element/Group  DS  Dialog  INQ
    Inquiry/Print      CO      Copy Dialog
    DO      Dialog Overview  QU  Quit
DD   Detailed Display  DT   Dialog Text
DP   Progress Note Text ED   Edit/Delete Dialog
Select Item: Next Screen// Find VA-AH-BPR ORDER CBC
and select the appropriate item number to open the element
```

Current dialog element/group name: CBC W/DIFF Used by:
VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER CBC //

Select ADDITIONAL FINDINGS: Q. Enter exact name of local order dialog

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION
+Item Seq. Dialog Summary

```
192 35.35.15 Element: VA-AH-BPR ORDER HIGH RESOLUTION CT
193 35.35.20 Element: VA-AH-BPR ORDER ECHOCARDIOGRAM
194 35.35.25 Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS
195 35.35.30 Element: VA-AH-BPR ORDER CBC W/DIFF
196 35.35.35 Element: VA-AH-BPR ORDER ENT CONSULT
197 35.35.40 Element: VA-AH-BPR ORDER PULMONARY CONSULT
198 40 Element: VA-AH-BPR SPACE (DIALOG & NOTE)
```

```
+      + Next - Prev Screen ?? Actions
>>>
ADD  Add Element  DS  Dialog  INQ  Inquiry/Print
CO   Copy Dialog  DO  Dialog Overview  QU  Quit
DD   Detailed Display  DT  Dialog Text
DP   Progress Note Text ED  Edit/Delete
Select Item: Next Screen// Find VA-AH-BPR ORDER ENT
and select the appropriate item number to open the element
```

Current dialog element/group name: ENT CONSULT Used
by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER ENT CONSULT //

Select ADDITIONAL FINDINGS: Q. Enter exact name of local order dialog

REMINDER DIALOG NAME: VA-AH-BPR INITIAL EVALUATION

+Item Seq. Dialog Summary

| | | |
|-----|----------|---|
| 192 | 35.35.15 | Element: VA-AH-BPR ORDER HIGH RESOLUTION CT CHEST |
| 193 | 35.35.20 | Element: VA-AH-BPR ORDER ECHOCARDIOGRAM |
| 194 | 35.35.25 | Element: VA-AH-BPR ORDER ARTERIAL BLOOD GAS |
| 195 | 35.35.30 | Element: VA-AH-BPR ORDER CBC W/DIFF |
| 196 | 35.35.35 | Element: VA-AH-BPR ORDER ENT CONSULT |
| 197 | 35.35.40 | Element: VA-AH-BPR ORDER PULMONARY CONSULT |
| 198 | 40 | Element: VA-AH-BPR SPACE (DIALOG & NOTE) |

```

+          + Next  - Prev Screen ??  Actions
>>>
ADD Add Element    DS Dialog Summary    INQ Inquiry
CO Copy Dialog     DO Dialog Order      QU Quit
DD Detailed Note   DT Dialog Text
DP Progress Note   ED Edit/Delete
Select Item: Next Screen// Find VA-AH-BPR ORDER PULMONARY
and select the appropriate item number to open the element

```

Current dialog element/group name: VA-AH-BPR PULMONARY CONSULT Used
by: VA-AH-BPR ORDERS (Current Reminder Dialog)

FINDING ITEM: HF AH-BPR ORDER PULM CONSULT //

Select ADDITIONAL FINDINGS: Q. Enter exact name of local order dialog

5. Enabling the dialog:

Before linking the dialog to a note title you must first add into the TIU Template Reminder Dialog Parameter list.

Adding dialog to the TIU Reminder Dialog Parameter (starting from the Reminder Manager Menu)

Select Reminder Managers Menu <TEST ACCOUNT> Option: CP CPRS Reminder Configuration

| | |
|-----|--|
| CA | Add/Edit Reminder Categories |
| CL | CPRS Lookup Categories |
| CS | CPRS Cover Sheet Reminder List |
| MH | Mental Health Dialogs Active |
| PN | Progress Note Headers |
| RA | Reminder GUI Resolution Active |
| TIU | TIU Template Reminder Dialog Parameter |

DL Default Outside Location
 PT Position Reminder Text at Cursor
 NP New Reminder Parameters
 GEC GEC Status Check Active
 WH WH Print Now Active

Select CPRS Reminder Configuration <TEST ACCOUNT> Option: **TIU**

Reminder Dialogs allows as Templates may be set for the following:

1 User USR [choose from NEW PERSON]
 3 Service SRV [choose from SERVICE/SECTION]
 4 Division DIV [NAME OF YOUR DIVISION]
 5 System SYS [YOURSERVER.YOURSITE.MED.VA.GOV]

Enter selection: **5** System [YOURSERVER.YOURSITE.MED.VA.GOV]

Select Display Sequence: **?** **← enter ? to determine last number in sequence**

Display Sequence Value

 1 TOBACCO CESSATION NOTE
 2 PATIENT EDUCATION ASSESSMENT
 3 TEST CLINIC DIALOG


This list is
condensed for
space reasons.

....
 412 NON VA CARE COORDINATION NOTE (D) **← last sequence number**
on my list yours will be different (con't next page)

Select Display Sequence: **Select appropriate Display Sequence number**

Are you adding 413 as a new Display Sequence? Yes// **YES**



Display Sequence: 413// 413

Clinical Reminder Dialog: **VA-AH-BPR INITIAL EVALUATION**  reminder dialog
 NATIONAL

...OK? Yes// **Yes**

a. Remove “disabled” from Reminder Dialog

Open reminder dialog manager, change view to ‘DIALOG’ locate ~~VA-BPR~~
 INITIAL EVALUATION NOTE dialog on the list.

- Edit the VA-AH-BPR INITIAL EVALUATION  and remove the DISABLE field.
- Remember VA-AH-BPR INITIAL EVALUATION  will not have an associated source reminder.

Select Item: Next Screen// **ED** Edit/Delete Dialog

NAME: VA-AH-BPR INITIAL EVALUATION █
Replace DISABLE: DISABLED IN EXCHANGE Replace @
SURE YOU WANT TO █ CLASS:
NATIONAL//
SPONSOR: Office of Public Health-Post-Deployment █
// REVIEW █
SOURCE REMINDER://

=====

While using the TIU editors if at any time the settings do not look correct NEVER click the apply button, but instead click CANCEL, then start over

=====

b. Create note title, AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE

From the Document Definitions (Manager) Menu: Create document definitions
Type Next Level at prompt
Enter corresponding number for Progress Notes at prompt
Hit enter until you find Medicine
Type Next Level and enter corresponding number for Medicine at prompt
Type Title at prompt
Type **AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE** at prompt
Map title to **PERSIAN GULF REGISTRY E & M NOTE**
Enter A at status prompt

- 1 TIU Parameters Menu ...
- 2 Document Definitions (Manager) ...
- 3 User Class Management ...
- 4 TIU Template Mgmt Functions ...
- 5 TIU Alert Tools
- 7 TIUHL7 Message Title Mapping █
- 6 Active Title Cleanup Report

Select TIU Maintenance Menu <TEST ACCOUNT> █ 2 Document
Definitions (Manager)

--- Manager Document Definition Menu ---

- 1 Edit Document Definitions
- 2 Sort Document Definitions
- 3 Create Document Definitions
- 4 Create Objects
- 5 Create TIU/Health Summary Objects

Select Document Definitions (Manager) <TEST>Option: 3
 Create Document Definitions.....

| | Name | Type |
|---|-----------------------|------|
| 1 | CLINICAL DOCUMENTS | CL |
| 2 | PROGRESS NOTES | CL |
| 3 | ADDENDUM | DC |
| 4 | DISCHARGE SUMMARY | CL |
| 5 | CLINICAL PROCEDURES | CL |
| 6 | LR LABORATORY REPORTS | CL |
| 7 | SURGICAL REPORTS | CL |

| + | ?Help | >ScrollRight | PS/PL PrintScrn/List | +/- | >>> |
|---------------------|-------|------------------|----------------------|-----------------------|-----|
| Class/DocumentClass | | Next Level | | Detailed Display/Edit | |
| (Title) | | Restart | | Status... | |
| (Component) | | Boilerplate Text | | Delete | |

Select CLINICAL DOCUMENTS Item (Line 2-7): 2

| | Name | Type |
|---|--|------|
| 2 | PROGRESS NOTES | CL |
| 3 | MEDICINE | DC |
| 4 | MEDICAL OBSERVATION NOTE | TL |
| 5 | CANCER CARE | TL |
| 6 | CARDIAC CATH LAB PROCEDURE REPORT (T) | TL |
| 7 | CARDIAC ELECTROPHYSIOLOGY PROCEDURE REPORT (T) | TL |
| 8 | CARDIOLOGY ACS TEAM FOLLOW-UP NOTE | TL |
| 9 | CARDIOLOGY AICD IMPLANTATION REPORT | TL |

| + | ?Help | >ScrollRight | PS/PL PrintScrn/List | +/- | >>> |
|-----------------------|-------|------------------|----------------------|-----------------------|-----|
| (Class/DocumentClass) | | Next Level | | Detailed Display/Edit | |
| Title | | Restart | | Status... | |
| (Component) | | Boilerplate Text | | Delete | |

Select Action: Next Screen// title Title
 Enter the Name of a new MEDICINE: AIRBORNE HAZARD/BURN PIT INITIAL EVALUATION
 NOTE
 CLASS OWNER: CLINICAL COORDINATOR// CLINICAL COORDINATOR

EVERY Local Title must be mapped to a VHA Enterprise Title.

Direct Mapping to Enterprise Standard Title...
 Your LOCAL Title is: AIRBORNE HAZARD/BURN PIT INITIAL EVALUATION NOTE

NOTE: Only ACTIVE Titles may be selected...

Select VHA ENTERPRISE STANDARD TITLE: PERSIAN GULF E & M NOTE


I found a match of: PERSIAN GULF REGISTRY E & M NOTE
 ... OK? Yes// YES

Ready to map LOCAL Title: AIRBORNE HAZARD/BURN PIT REGISTRY N:
NOTE to
VHA Enterprise Standard Title: **PERSIAN GULF REGISTRY E & M NOTE.**
... OK? Yes//**YES**


STATUS: (A/I/T): INACTIVE// **A** ACTIVE Entry Activated.

4.1 Linking Progress Note Title to Dialog

From Template Editor
Edit Shared Template
Click + Document Titles

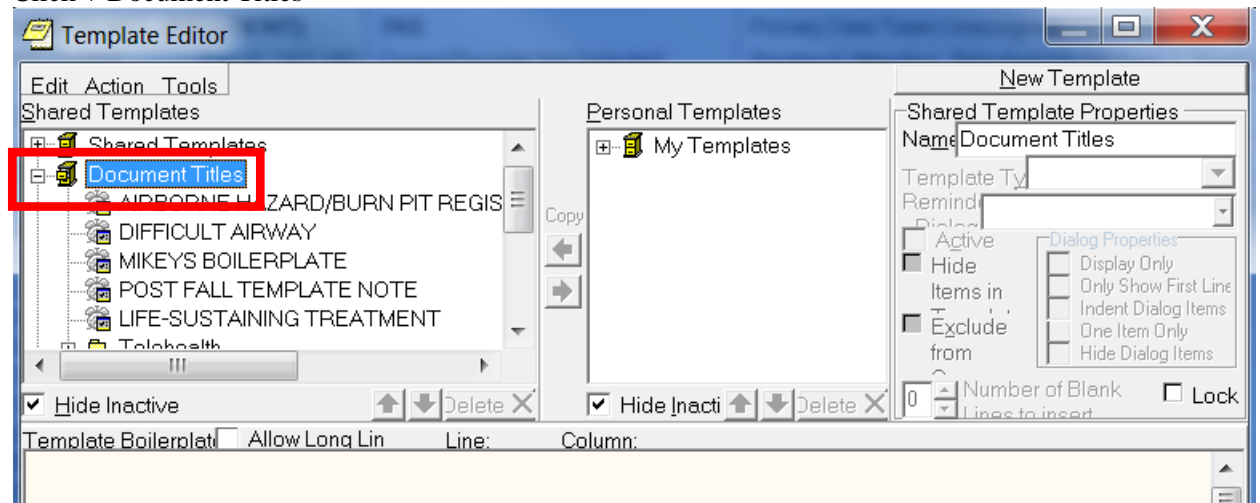


Click New Template button

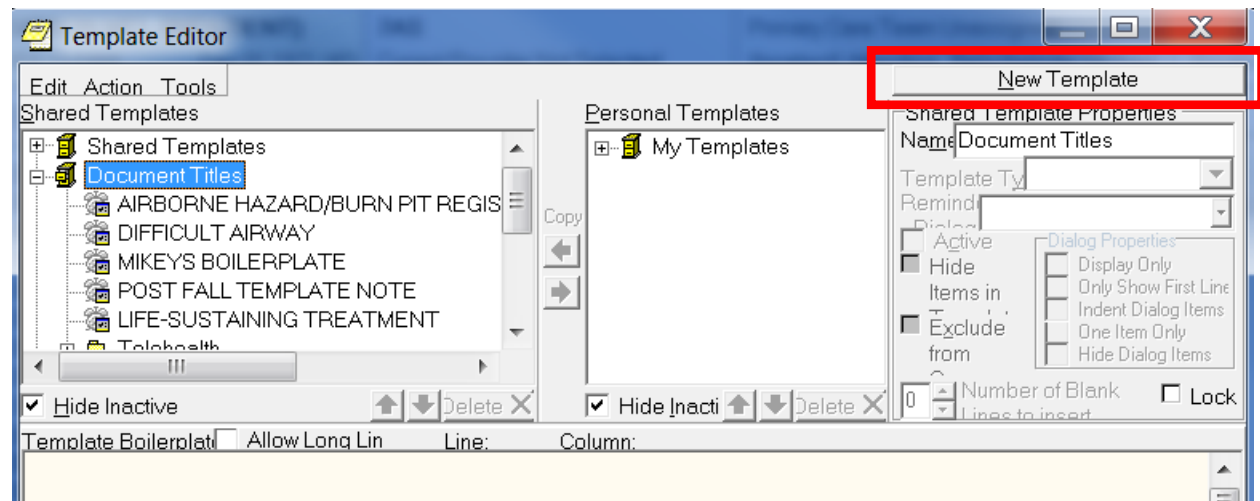


Type AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION in Name box


Click down arrow and Click on Reminder Dialog in Template Type Box

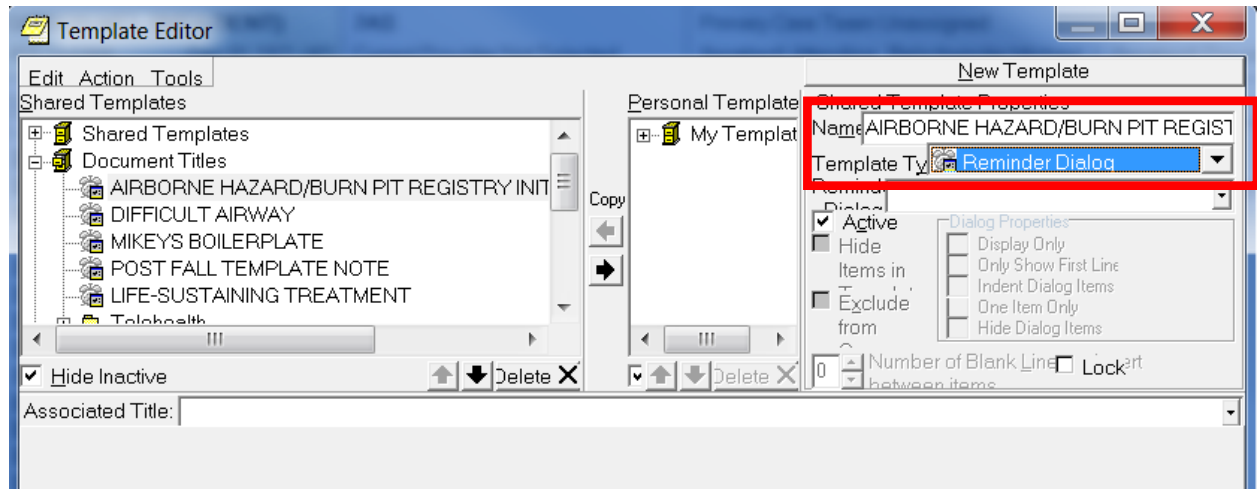


Type VA-AH-BPR INITIAL EVALUATION NOTE in Reminder Dialog box

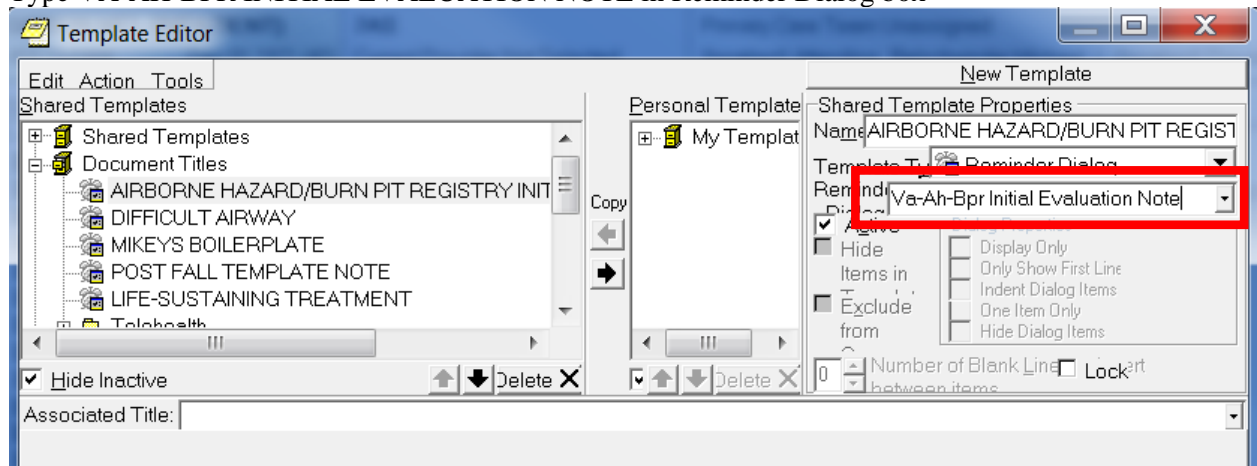

Type AIRBORNE HAZAED/BURN PIT REGISTRY INITIAL EVALUATION NOTE in Associated Title Box

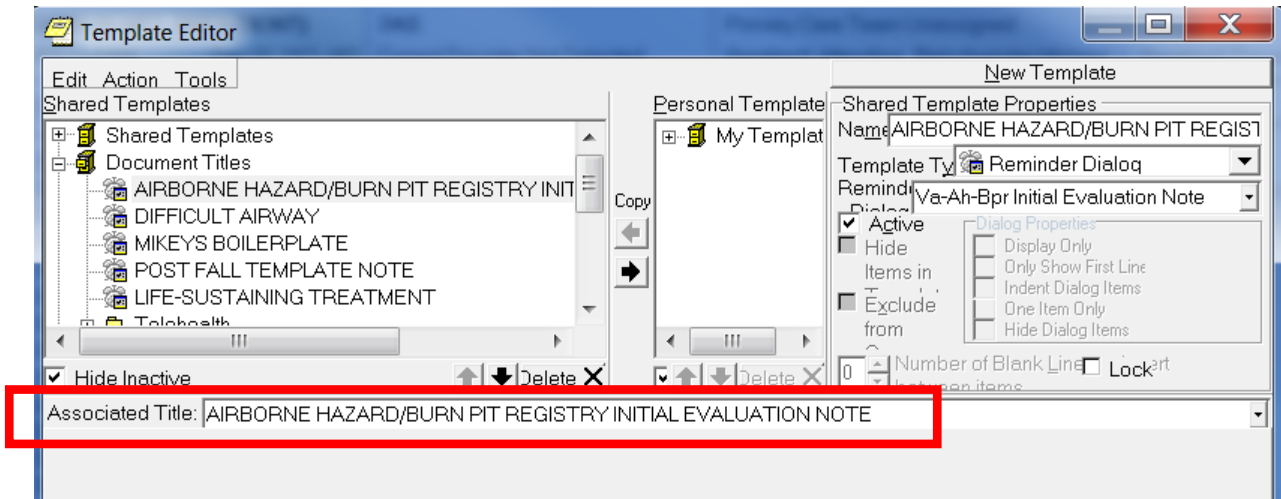

Click Apply in lower right corner

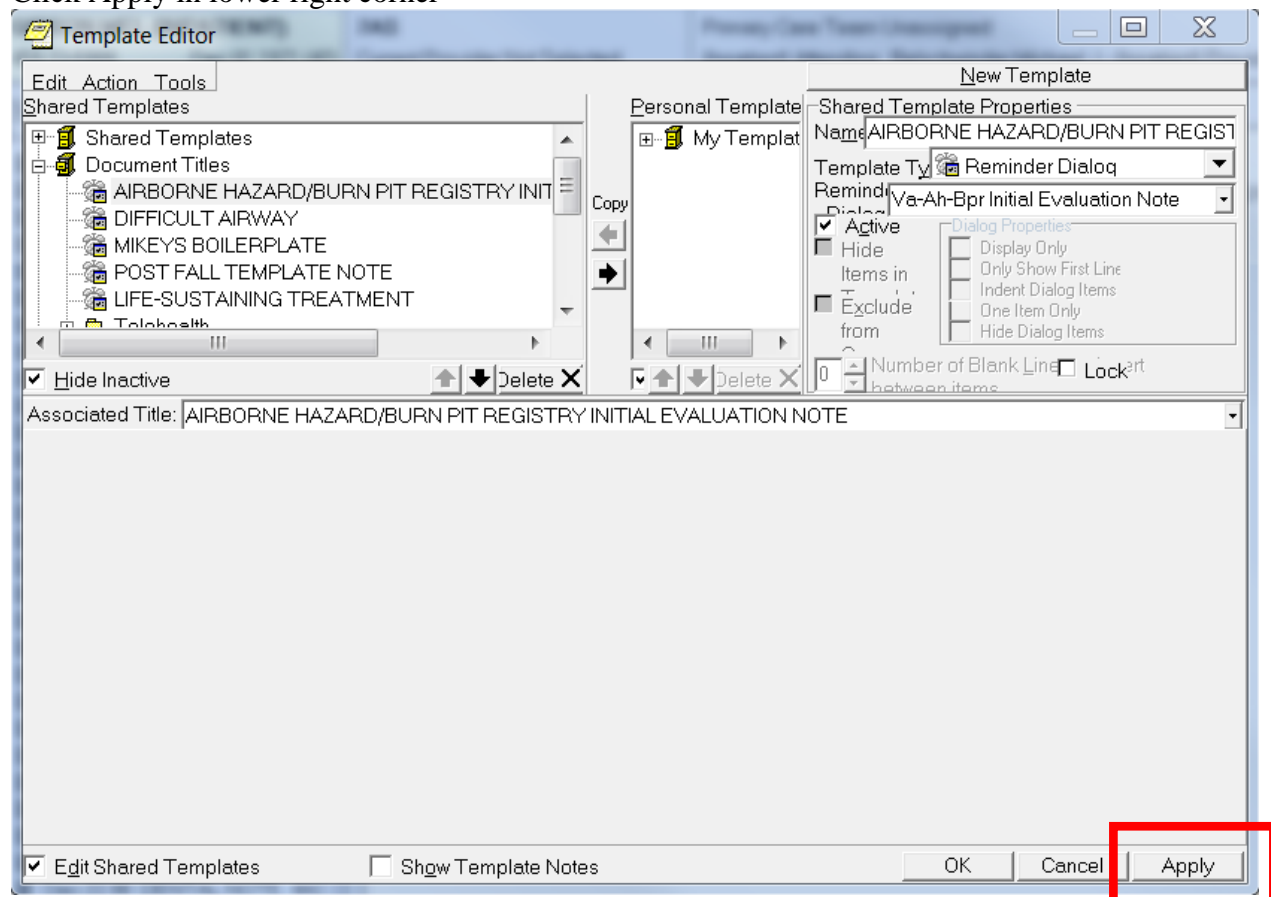

Confirm title is active in CPRS and dialog is attached

Progress Note Properties

Progress Note Title: AIRBORNE <AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE>

AIRBORNE <AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE>
AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE
AIRBORNE HAZARD/BURN PIT REGISTRY INITIAL EVALUATION NOTE <AIRBORNE HA
AIRWAY <DIFFICULT AIRWAY>
ALERT <INFECTION CONTROL ALERT NOTE>
ALLERGY <ADVERSE REACTION/ALLERGY>
ALLERGY <ALLERGY/IMMUNIZATION>

Date/Time of Note: Jan 21, 2014@07:28 ...

Author: Fahner, Jeff

OK Cancel

Appendix A: Installation Example

XPD Main Menu

- 1 Load a Distribution
- 2 Verify Checksums in Transport Global
- 3 Print Transport Global
- 4 Compare Transport Global to Current System
- 5 Backup a Transport Global
- 6 Install Package(s)
Restart Install of Package(s)
Unload a Distribution

Select Installation <TEST ACCOUNT> Option: 1 Load a Distribution
Enter a Host File: <your directory>PXRM_2_0_39.KID

KIDS Distribution saved on Sep 11, 2014@09:18:02
Comment: Airborne hazard

This Distribution contains Transport Globals for the following Package(s):
Build PXRM*2.0*39 has been loaded before, here is when:

PXRM*2.0*39 Install Completed
was loaded on Sep 11, 2014@10:42:09

OK to continue with Load? NO// YES

Distribution OK!

Want to Continue with Load? YES//
Loading Distribution...

PXRM*2.0*39
Use INSTALL NAME: PXRM*2.0*39 to install this Distribution.

- 1 Load a Distribution
- 2 Verify Checksums in Transport Global
- 3 Print Transport Global
- 4 Compare Transport Global to Current System
- 5 Backup a Transport Global
- 6 Install Package(s)
Restart Install of Package(s)
Unload a Distribution

Select Installation <TEST ACCOUNT> Option: INstall Package(s)
Select INSTALL NAME: PXRM*2.0*39 10/9/14@12:21:25
=> Airborne hazard ;Created on Sep 11, 2014@09:18:02

This Distribution was loaded on Oct 09, 2014@12:21:25 with header of
Airborne hazard ;Created on Sep 11, 2014@09:18:02
It consisted of the following Install(s):
PXRM*2.0*39
Checking Install for Package PXRM*2.0*39

Install Questions for PXRM*2.0*39

Incoming Files:

811.8 REMINDER EXCHANGE (including data)
Note: You already have the 'REMINDER EXCHANGE' File.
I will OVERWRITE your data with mine.

Want KIDS to INHIBIT LOGONS during the install? NO//
Want to DISABLE Scheduled Options, Menu Options, and Protocols? NO//

Enter the Device you want to print the Install messages.
You can queue the install by enter a 'Q' at the device prompt.
Enter a '^' to abort the install.

DEVICE: HOME// TELNET PORT

Install Started for PXRM*2.0*39 :
Oct 09, 2014@12:21:36

Build Distribution Date: Sep 11, 2014

Installing Routines:
Oct 09, 2014@12:21:36

Running Pre-Install Routine: PRE^PXRMP39I

DISABLE options.

DISABLE protocols.

Installing Data Dictionaries:
Oct 09, 2014@12:21:36

Installing Data:
Oct 09, 2014@12:21:38

Running Post-Install Routine: POST^PXRMP39I

ENABLE options.

ENABLE protocols.

There are 2 Reminder Exchange entries to be installed.
1. Installing Reminder Exchange entry VA-AH/BPR PXRM*2*39
2. Installing Reminder Exchange entry VA-PATCH 39 POST COMPONENTS

Updating Routine file...

Updating KIDS files...

PXRM*2.0*39 Installed.
Oct 09, 2014@12:22:13

Not a production UCI

Install Completed PXRM*2.0*39 I

Appendix B: Health Factor List

CATEGORY: AH-BPR CHIEF COMPLAINT

AH-BPR CC CANCER
AH-BPR CC CHEST PAIN
AH-BPR CC CHILDREN
AH-BPR CC CHRONIC SINUS INFECTION
AH-BPR CC CONCERNED
AH-BPR CC COUGH > 3 WKS
AH-BPR CC DEC EXERCISE ABILITY
AH-BPR CC EYE PROBLEM
AH-BPR CC GI PROBLEM
AH-BPR CC HAY FEVER/RESP ALLERGY
AH-BPR CC HEART PROBLEM
AH-BPR CC IMMUNE PROBLEM
AH-BPR CC NEURO PROBLEM
AH-BPR CC OTHER
AH-BPR CC RUNNY NOSE
AH-BPR CC SKIN PROBLEM
AH-BPR CC SOB
AH-BPR CC SORE THROAT
AH-BPR CC SPUTUM >3 WKS
AH-BPR CC WHEEZING/WHISTLING

CATEGORY: AH-BPR EXPOSURE CONCERNS

AH-BPR CONCERN HOBBIES
AH-BPR CONCERN MILITARY JOB
AH-BPR CONCERN OFF BASE AIR
AH-BPR CONCERN ON BASE AIR
AH-BPR CONCERN OTHER
AH-BPR CONCERN SMOKING
AH-BPR CONCERN UNKNOWN

CATEGORY: AH-BPR ASTHMA

AH-BPR ASTHMA NO
AH-BPR ASTHMA UNSURE
AH-BPR ASTHMA YES

CATEGORY: AH-BPR SMOKING DEPLOYMENT

AH-BPR SMOKING DEPLOYMENT NO
AH-BPR SMOKING DEPLOYMENT UNKNOWN
AH-BPR SMOKING DEPLOYMENT YES

CATEGORY: TOBACCO

CURRENT SMOKER
PREVIOUS SMOKER
LIFETIME NON-SMOKER

CATEGORY: AH-BPR BIRTH DEFECTS

AH-BPR BIRTH DEFECTS N/A
AH-BPR BIRTH DEFECTS NO
AH-BPR BIRTH DEFECTS UNCLEAR
AH-BPR BIRTH DEFECTS YES

CATEGORY: AH-BPR ABNORMALITIES

AH-BPR ABN CARDIO
AH-BPR ABN COGNITIVE/MEMORY
AH-BPR ABN DERM
AH-BPR ABN ENT
AH-BPR ABN FATIGUE
AH-BPR ABN FEVER
AH-BPR ABN GI
AH-BPR ABN GU
AH-BPR ABN HEME/ONC
AH-BPR ABN LYMPHADENOPATHY
AH-BPR ABN MUSCULOSKELETAL
AH-BPR ABN NEURO
AH-BPR ABN NONE
AH-BPR ABN OTHER
AH-BPR ABN PSYCH
AH-BPR ABN RESPIRATORY
AH-BPR ABN TEMP SENS
AH-BPR ABN WEIGHT GAIN
AH-BPR ABN WEIGHT LOSS

CATEGORY: AH-BPR ORDERS

AH-BPR ORDER CHEST XRAY

AH-BPR ORDER ABG

AH-BPR ORDER CBC

AH-BPR ORDER CHEST CT

AH-BPR ORDER ECHO

AH-BPR ORDER ENT CONSULT

AH-BPR ORDER PULM CONSULT

AH-BPR ORDER PFT

Acronyms

The OIT Master Glossary is available at

http://vaww.oed.wss.va.gov/process/Library/master_glossary/masterglossary.htm

| Term | Definition |
|------------|---|
| AH-BPR | Airborne Hazard-Open Burn Pit |
| ASU | Authorization/Subscription Utility |
| Clin4 | National Customer Support team that supports Clinical Reminders |
| CPRS | Computerized Patient Record System |
| DBA | Database Administration |
| DG | Registration and Enrollment Package namespace |
| ESM | Enterprise Systems Management (ESM) |
| FIM | Functional Independence Measure |
| GMTS | Health Summary namespace (also HSUM) |
| GUI | Graphic User Interface |
| HRMH/HRMHP | High Risk Mental Health Patient |
| IAB | Initial Assessment & Briefing |
| ICD-10 | International Classification of Diseases, 10th Edition |
| ICR | Internal Control Number |
| IOC | Initial Operating Capabilities |
| LSSD | Last Service Separation Date |
| MH | Mental Health |
| MHTC | Mental Health Treatment Coordinator |
| OHI | Office of Health Information |
| OI | Office of Information |
| OIF/OEF | Operation Iraqi Freedom/Operation Enduring Freedom |
| OIT/OI&T | Office of Information Technology |
| OMHS | Office of Mental Health Services |
| ORR | Operational Readiness Review |
| PCS | Patient Care Services |
| PD | Product Development |
| PIMS | Patient Information Management System |

| Term | Definition |
|-------------|--|
| PMAS | Program Management Accountability System |
| PTM | Patch Tracker Message |
| PXRM | Clinical Reminder Package namespace |
| RSD | Requirements Specification Document |
| SD | Scheduling Package Namespace |
| SQA | Software Quality Assurance |
| USR | ASU package namespace |
| VA | Department of Veteran Affairs |
| VHA | Veterans Health Administration |
| VISN | Veterans Integrated Service Network |
| VistA | Veterans Health Information System and Technology Architecture |