

Enterprise Health Management Platform (eHMP)

USER GUIDE

eHMP Version 1.2

Viewer Edition

August 2016

VA



U.S. Department
of Veterans Affairs

Contents

Chapter 1: Introduction to eHMP	1
About eHMP	1
System Requirements.....	1
eHMP User Interface (UI)	1
eHMP Help.....	2
Online Help Feature	2
National Service Desk (NSD).....	2
eHMP Help Resources.....	2
The Organization of this User Guide	2
Chapter 2: Getting Started	3
Logging into eHMP via IE 11	3
Adding eHMP to Favorites in IE11	4
Accessing eHMP using IE11 Favorites	5
Accessing eHMP through CPRS	6
Logging Out of eHMP	7
Manual Sign Out	7
Auto Sign Out	7
Chapter 3: Patient Search and Selection.....	9
About Patient Search and Selection	9
Patient Search	9
General Patient Search	9
My Site Search	10
Nationwide Search.....	11
Selecting a Patient.....	12
Accessing an Unrestricted Patient Record	12
Accessing a Restricted Patient Record	13
Accessing a Flagged Patient Record	14
Main Application Window	15
Chapter 4: The eHMP Header	16

About the eHMP Header.....	16
Patient Selection Button.....	16
eHMP Banner Information	16
User Information Dropdown	17
Patient Photo.....	18
Patient Demographics	18
Patient Postings.....	19
Visit Information	20
Care Team Information.....	20
Chapter 5: Global Timeline Date Filter	22
About the Global Timeline Date Filter	22
Using the Preset Date Ranges.....	23
Using the Custom Date Range	24
Using the All Events Timeline	24
Using the Timeline Summary.....	25
Chapter 6: Applet Views	27
Active & Recent Medications	27
Active & Recent Medications: Trend View	27
Active & Recent Medications: Summary View	29
Active & Recent Medications: Detail View	30
Allergies.....	31
Allergies: Trend View	31
Allergies: Summary View	31
Allergies: Expanded View.....	32
Allergies: Detail View	32
Appointments & Visits	33
Appointments & Visits: Summary View	33
Appointments & Visits: Expanded View	33
Appointments & Visits: Detail View	34
Clinical Reminders	35
Clinical Reminders: Summary View	35

Clinical Reminders: Expanded View	35
Clinical Reminders: Detail View.....	36
Community Health Summaries.....	36
Community Health Summaries: Summary View	36
Community Health Summaries: Expanded View	37
Community Health Summaries: Detail View	37
Documents	38
Documents: Summary View	38
Documents: Expanded View	38
Documents: Detail View	39
Encounters.....	40
Encounters: Trend View	40
Encounters: Detail View	41
Immunizations	42
Immunizations: Trend View	42
Immunizations: Summary View	42
Immunizations: Expanded View.....	43
Immunizations: Detail View	43
Medications Review	44
Medications Review: Expanded View.....	44
Medications Review: Detail View	48
Military History	49
Military History: Summary View	49
Military History: Expanded View	49
Military History: Detail View	50
Military History: Edit Form.....	50
Narrative Lab Results.....	51
Narrative Lab Results: Summary View	51
Narrative Lab Results: Expanded View.....	52
Narrative Lab Results: Detail View	52
Numeric Lab Results	53

Numeric Lab Results: Trend View	53
Numeric Lab Results: Summary View	54
Numeric Lab Results: Expanded View	55
Numeric Lab Results: Detail View	55
Orders.....	56
Orders: Summary View	56
Orders: Expanded View.....	57
Orders: Detail View	57
Problems	58
Problems: Trend View	58
Problems: Summary View	59
Problems: Expanded View	60
Problems: Detail View	61
Reports	62
Reports: Summary View	62
Reports: Detail View	62
Stacked Graphs.....	63
Stacked Graphs: Expanded View	63
Add a Graph to Stacked Graphs Applet.....	64
Edit the Stacked Graphs Applet	65
Timeline	66
Timeline: Summary View	66
Timeline: Expanded View.....	66
Timeline: Detail View	67
VistA Health Summaries	67
VistA Health Summaries: Summary View	68
VistA Health Summaries: Detail View	68
Vitals.....	69
Vitals: Trend View	69
Vitals: Summary View	70
Vitals: Expanded View.....	71

Vitals: Detail View	72
Chapter 7: Applet Features	73
About Applets.....	73
Refresh Button.....	73
Online Help Button	73
Filter Button.....	74
User-Defined Applet Filters	74
Maximize Applet Button	75
Minimize Applet Button.....	76
Sortable Column Headers	76
Manual Tile Sort.....	76
Infobutton.....	77
Details Form Button.....	78
Quicklook Button.....	79
Submenu Button.....	79
Comment Indicators	80
Chapter 8: Workspaces	81
About Workspaces.....	81
Predefined Workspaces	81
Workspace Manager	81
Add a New Workspace	82
Workspace Filter	82
Default View	83
Workspace Title	83
Associated Conditions	83
Workspace Description	84
Duplicate Workspaces.....	84
Rearrange Workspace Listings.....	85
Delete Workspaces	85
Preview Workspace.....	86
Customize a New Workspace.....	86

Launch Workspaces 88

Chapter 9: Other eHMP Features..... 90

Search Record Field 90

Patient Record Refresh and Status Bar 92

Table of Figures

Figure 2-1 eHMP Login Screen.....	3
Figure 2-2 Log In Error Message	4
Figure 2-3 Star Button	4
Figure 2-4 Add a Favorite.....	5
Figure 2-5 Star Button	5
Figure 2-6 Selecting the eHMP Link.....	6
Figure 2-7 eHMP Link on the CPRS Tools Menu.....	6
Figure 2-8 Sign Out Button on Header Dropdown	7
Figure 2-9 Auto Logout Warning Message	7
Figure 2-10 Inactivity Logout Message	8
Figure 3-1 Patient Selection Screen	9
Figure 3-2 General Patient Search	10
Figure 3-3 My Site/Clinics.....	10
Figure 3-4 My Site/Wards.....	11
Figure 3-5 Nationwide Search.....	12
Figure 3-6 Too Many Results Have Returned	12
Figure 3-7 Patient Detail Dialog	13
Figure 3-8 Restricted Record Acknowledgement.....	13
Figure 3-9 Confirm Selection of Restricted Patient Record	14
Figure 3-10 Confirm Selection of Flagged Patient Record.....	14
Figure 3-11 Main Application Window.....	15
Figure 4-1 The eHMP Header.....	16
Figure 4-2 Detailed eHMP Banner Information.....	17
Figure 4-3 Sign Out Option	17
Figure 4-4 Patient Photo (Example).....	18
Figure 4-5 Patient Demographics	18
Figure 4-6 Patient Demographics Bar with Broken Chain Link	18
Figure 4-7 Expanded Patient Information	19
Figure 4-8 Postings Icons.....	19

Figure 4-9 Expanded Postings Information.....	20
Figure 4-10 Visit Information	20
Figure 4-11 Care Team Information	21
Figure 4-12 Detailed Care Team Information.....	21
Figure 5-1 Global Timeline Date Filter	22
Figure 5-2 Expanded Global Timeline Date Filter.....	23
Figure 5-3 Preset Date Ranges	23
Figure 5-4 Calendar Dropdown	24
Figure 5-5 All Events Timeline	24
Figure 5-6 Filter Button.....	25
Figure 5-7 Patient Activity Detail Dialog Box	26
Figure 6-1 Active & Recent Medications Applet Trend View	27
Figure 6-2 Active & Recent Medications Applet Toolbar Menu – Trend View	28
Figure 6-3 Active & Recent Medications Quick View.....	29
Figure 6-4 Active & Recent Medications Applet Summary View.....	29
Figure 6-5 Active & Recent Medications Applet Toolbar Menu – Summary View	29
Figure 6-6 Active & Recent Medications Details Form Button	30
Figure 6-7 Active & Recent Medications Detail Dialog Box	30
Figure 6-8 Allergies Applet Trend View	31
Figure 6-9 Allergies Applet Summary View	31
Figure 6-10 Allergies Applet Expanded View	32
Figure 6-11 Allergy Detail Dialog Box.....	32
Figure 6-12 Appointments & Visits Applet Summary View	33
Figure 6-13 Appointments & Visits Applet Expanded View	34
Figure 6-14 Appointments & Visits Detail Dialog Box	34
Figure 6-15 Clinical Reminders Applet Summary View	35
Figure 6-16 Clinical Reminders Applet Expanded View	35
Figure 6-17 Clinical Reminders Detail Dialog Box.....	36
Figure 6-18 Community Health Summaries Applet Summary View.....	36
Figure 6-19 Community Health Summaries Applet Expanded View	37
Figure 6-20 Community Health Summaries Detail Dialog Box	37

Figure 6-21 Documents Applet Summary View	38
Figure 6-22 Documents Applet Expanded View.....	39
Figure 6-23 Documents Detail Dialog Box	39
Figure 6-24 Encounters Applet Trend View	40
Figure 6-25 Encounters List in Trend View (Expanded).....	40
Figure 6-26 Encounters Applet Toolbar Menu	41
Figure 6-27 Encounters Detail Dialog Box.....	41
Figure 6-28 Immunizations Applet Trend View.....	42
Figure 6-29 Detail of an Immunization in Trend View	42
Figure 6-30 Immunizations Applet Summary View	42
Figure 6-31 Immunizations Applet Expanded View.....	43
Figure 6-32 Immunizations Detail Dialog Box.....	43
Figure 6-33 Medications Review Applet Expanded View.....	44
Figure 6-34 Medications Review Outpatient Meds	45
Figure 6-35 Outpatient Medications Graph Indicators	46
Figure 6-36 Medications Review Inpatient Meds	47
Figure 6-37 Medications Review Clinic Order Meds	48
Figure 6-38 Medications Review Applet Detail View	48
Figure 6-39 Military History Applet Summary View	49
Figure 6-40 Military History Applet Expanded View	50
Figure 6-41 Military History Detail Dialog Box	50
Figure 6-42 Edit Military History	51
Figure 6-43 Narrative Lab Results Applet Summary View.....	51
Figure 6-44 Narrative Lab Results Applet Expanded View	52
Figure 6-45 Narrative Lab Results Detail Dialog Box.....	53
Figure 6-46 Numeric Lab Results Applet Trend View.....	54
Figure 6-47 Numeric Lab Results Quick View.....	54
Figure 6-48 Numeric Lab Results Applet Summary View.....	54
Figure 6-49 Numeric Lab Results Applet Expanded View	55
Figure 6-50 Numeric Lab Results Detail Dialog Box.....	56
Figure 6-51 Orders Applet Summary View.....	56

Figure 6-52 Orders Applet Expanded View	57
Figure 6-53 Orders Detail Dialog Box	58
Figure 6-54 Problems Applet Trend View	58
Figure 6-55 Problems Applet Trend View Toolbar Menu	59
Figure 6-56 Problems Applet Quick View	59
Figure 6-57 Problems Applet Summary View.....	60
Figure 6-58 Problems Applet Summary View Toolbar Menu	60
Figure 6-59 Problems Applet Expanded View	61
Figure 6-60 Problems Detail Dialog Box	61
Figure 6-61 Reports Applet Summary View	62
Figure 6-62 Reports Detail Dialog Box	63
Figure 6-63 Stacked Graphs Applet Expanded View	64
Figure 6-64 Add Graph to Stacked Graphs Applet	64
Figure 6-65 Graph Added to Stacked Graphs Applet	64
Figure 6-66 Stacked Graphs Applet Toolbar Menu	65
Figure 6-67 Delete Stacked Graphs Warning	65
Figure 6-68 Timeline Summary View.....	66
Figure 6-69 Timeline Expanded View	66
Figure 6-70 Timeline Detail Dialog Box.....	67
Figure 6-71 VistA Health Summaries Summary View	68
Figure 6-72 VistA Health Summaries Detail Dialog Box.....	68
Figure 6-73 Vitals Applet Trend View	69
Figure 6-74 Vitals Applet Trend View Toolbar Menu.....	70
Figure 6-75 Vitals Applet Quick View	70
Figure 6-76 Vitals Applet Summary View.....	70
Figure 6-77 Vitals Detail Dialog Box	71
Figure 6-78 Vitals Applet Expanded View	72
Figure 6-79 Vitals Detail Dialog Box.....	72
Figure 7-1 Applet Feature Buttons	73
Figure 7-2 PDF Version of the eHMP User Guide	74
Figure 7-3 Filtered Applet	75

Figure 7-4 Filter Name	75
Figure 7-5 Sortable Column Header	76
Figure 7-6 Tile Sorting.....	77
Figure 7-7 Manual Tile Sorting	77
Figure 7-8 Infobutton in the Problems Applet.....	78
Figure 7-9 Detail Dialog Box for Selected Item	78
Figure 7-10 Quicklook Button	79
Figure 7-11 Selecting Submenu Button to List Associated Workspaces	79
Figure 7-12 Comment Indicator	80
Figure 7-13 Detailed Comment Information.....	80
Figure 8-1 Workspace Manager Window	82
Figure 8-2 Filter Workspaces	82
Figure 8-3 Workspace Title Field.....	83
Figure 8-4 Search Problems Field	83
Figure 8-5 Problems Search Results.....	84
Figure 8-6 Number of Associated Problems	84
Figure 8-7 Workspace Add Description Field	84
Figure 8-8 Duplicated Workspace	85
Figure 8-9 Rearrange Workspaces	85
Figure 8-10 Delete User-Defined Workspace Confirmation Message	86
Figure 8-11 Workspace Preview	86
Figure 8-12 Workspace Editor Screen	87
Figure 8-13 Select a View	87
Figure 8-14 Using Double Arrows to Modify Applet Size	88
Figure 8-15 Launch Workspaces	88
Figure 9-1 Search Record Field Dropdown.....	90
Figure 9-2 Search Record Results	90
Figure 9-3 Selected Keyword	91
Figure 9-4 Selected Keyword Detail Dialog Box	91
Figure 9-5 eHMP Status Bar.....	92
Figure 9-6 eHMP Data Sources.....	92

Chapter 1: Introduction to eHMP

About eHMP

The Enterprise Health Management Platform (eHMP) project is a multi-year effort to evolve a modern, service-oriented platform which provides a web-based user interface (UI), clinical data services (CDS), and assembles patient clinical data from federated Veterans Health Information Systems and Technology Architecture (VistA) repositories, Department of Defense (DoD), and private partner data sources, reflective of each location providing care to the patient. This federated data is aggregated into an enterprise patient record. eHMP service components will span all application layers, including presentation, business and core services, and data access.

Release 1.2 introduces critical viewer edition enhancements to provide new capabilities to the Department of Veterans Affairs (VA) beyond what is available today via Computerized Patient Record System (CPRS), Joint Legacy Viewer (JLV), and VistAWeb. The system will provide enhanced presentations of clinical data that will range from trend views that provide a quick snapshot of easily understandable data, to detailed views that provide the user with a full range of options for examining longitudinal patient medical records. Users will be able to configure these views into a limitless number of custom workspaces in order to support a variety of clinical workflows. There will also be multiple pre-configured workspaces available to the user, which are filtered for specific conditions. The workspaces will provide the appropriate clinical information for a selected condition (e.g., COPD, Diabetes). Further enhancements will include the Military History applet, improved text search, and online application help screens.

A significant step toward reaching interoperability goals is the incorporation of discrete, clinical information received from community health partners via C32 (Continuity of Care Document (CCD), which focuses on patient summary information) and Consolidated-Clinical Document Architecture (C-CDA) documents. The narrative portions of these documents are available in eHMP.


System Requirements

eHMP User Interface (UI)

eHMP is a web-based application, and is designed to run in a web browser. eHMP release v1.2 is optimized for use with Internet Explorer 11 (IE 11).

eHMP Help

Online Help Feature

There are online help buttons throughout the eHMP application. The online help feature provides detailed information about specific system topics. Select the **online help** button  to open context-specific help in a separate window.

National Service Desk (NSD)

The National Service Desk (NSD) is available to report any issues while using eHMP. To contact the NSD for eHMP support, call 1-855-NSD-Help (1-855-673-4357).

Please note the following:

- The NSD sets priorities for tickets based on an accepted NSD Priority Matrix.
- Requesting support via an email is not as expedient as placing a phone call, as emails are deemed to be non-urgent. Emails are processed for support within 24 hours of receipt.
- Password resets cannot be achieved via email; you must contact NSD by phone.

eHMP Help Resources

The VA eHMP Deployment SharePoint site¹ and the eHMP Deployment page on VA Pulse² provides resources for more reference information on the eHMP application. Some of the available resources include frequently asked questions (FAQs), quick reference guides, and system demonstrations.

In addition, the VA Software Document Library³ (VDL) houses a variety of eHMP documents.

The Organization of this User Guide

This User Guide is organized to mimic the use of the application, beginning with how to log in to eHMP, and progresses by giving instructions for each of the features available within the application.

¹ http://vawww.oed.portal.va.gov/pm/iehr/vista_evolution/ehmpdeployment/SitePages/Home.aspx

² <https://www.vapulse.net/community/vista-evolution/enterprise-availability/overview>

³ <http://www.va.gov/vdl/>

Chapter 2: Getting Started

eHMP can be accessed in two ways: through IE11, or through CPRS. The steps for accessing eHMP are described below.


Note: Not all sites have eHMP set up on the Tools Menu in CPRS.

Logging into eHMP via IE 11

1. Open **IE11**. Enter the **eHMP v1.2 web address**⁴ into the address field, and the login screen displays.

The login screen displayed in Figure 2-1 is comprised of disclosure and warning language, with application information on the left side of the window, and credential entry fields on the right side of the window.

Figure 2-1 eHMP Login Screen

2. Select a **facility** from the *Select a facility...* dropdown menu.
 3. Enter your **VistA access code** into the *Access Code* field, and then press the **Tab** key to advance to the following field.
 4. Enter your **VistA verify code** into the *Verify Code* field.
- Note:** Select the **Show** icon  to display the entered text in the *Access Code* and *Verify Code* fields.
5. Review the acknowledgement and then select the **checkbox**.


⁴ <https://ehmp.va.gov>

Note: eHMP does not allow you to sign in unless the checkbox is selected.

6. Select the **Sign In** button to display the Patient Selection screen.

Note: If you receive the message below (Figure 2-2) when trying to log in, you need to change the expired Verify Code using CPRS. The current version of eHMP does not allow you to change an expired Verify Code.

Figure 2-2 Log In Error Message

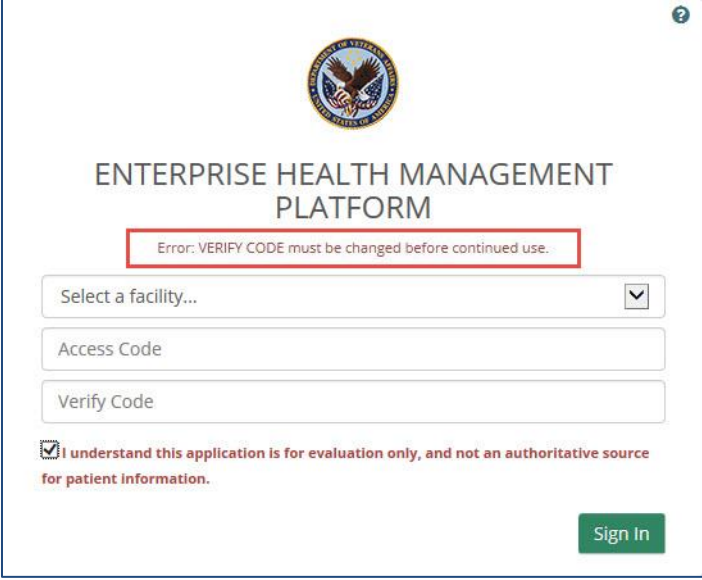


The screenshot shows the login interface for the Enterprise Health Management Platform. At the top center is the Department of Veterans Affairs seal. Below it, the text "ENTERPRISE HEALTH MANAGEMENT PLATFORM" is displayed. A red-bordered error message box states: "Error: VERIFY CODE must be changed before continued use." Below the error message are three input fields: "Select a facility..." (a dropdown menu), "Access Code", and "Verify Code". At the bottom left, there is a checkbox that is checked, with the text "I understand this application is for evaluation only, and not an authoritative source for patient information." At the bottom right is a green "Sign In" button.

Adding eHMP to Favorites in IE11


1. Access **eHMP**.
2. Select the **star** button (Figure 2-3) in the upper right corner of the application window.

Figure 2-3 Star Button

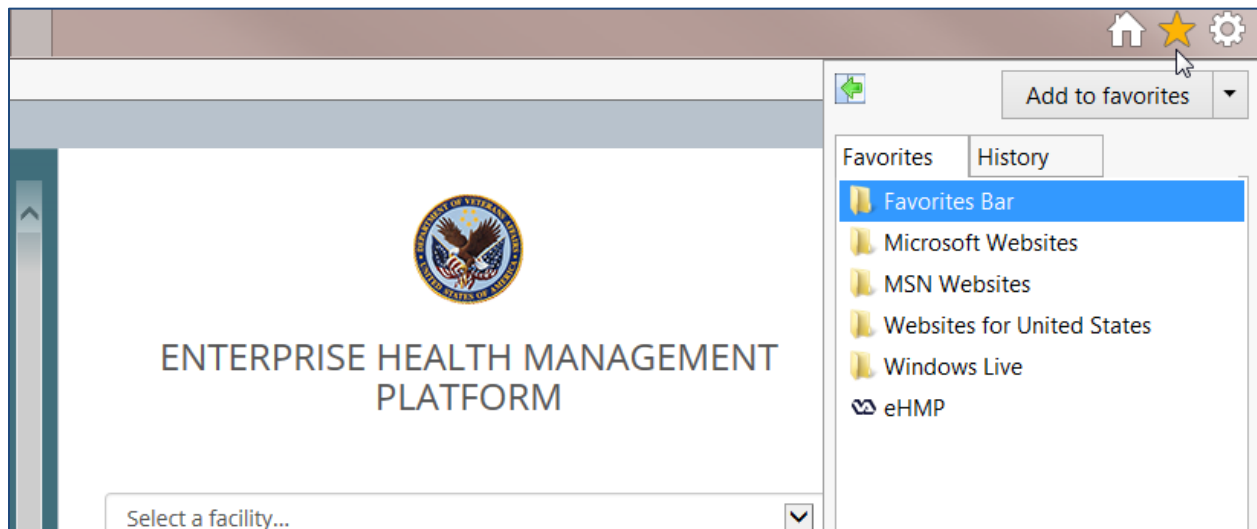

3. Select **Add to Favorites**. The Add a Favorite dialog (Figure 2-4) pops up.

Figure 2-4 Add a Favorite

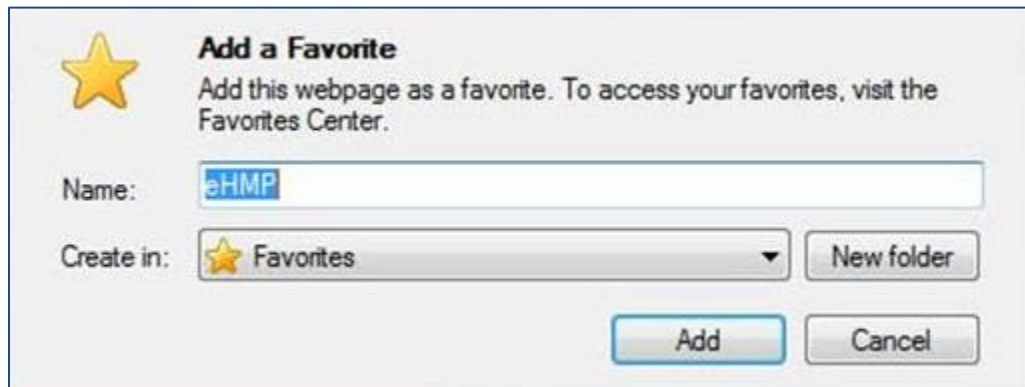


4. Select **Add**.

Accessing eHMP using IE11 Favorites

1. Open **IE11**.
2. Select the **star** button (Figure 2-5) in the upper right corner of the application window.

Figure 2-5 Star Button

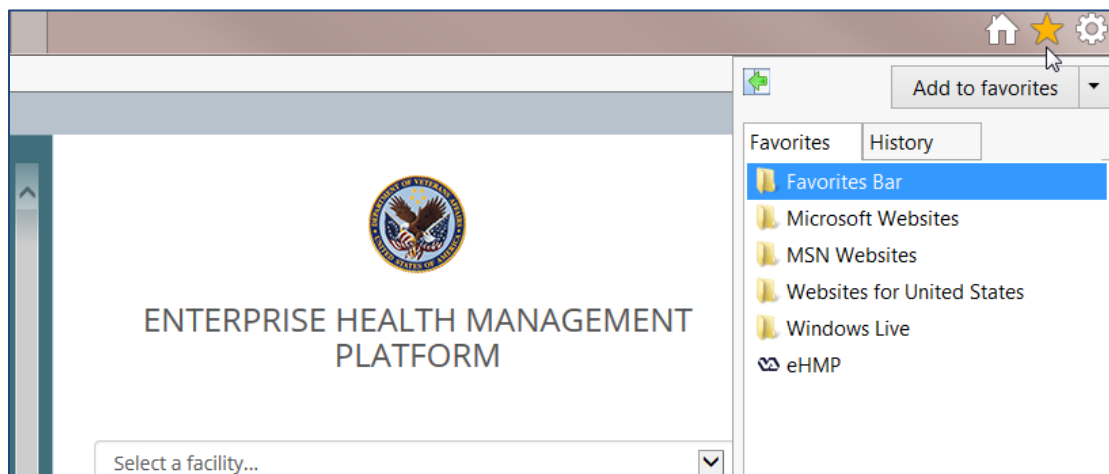

3. Select the **eHMP link** (Figure 2-6).

Figure 2-6 Selecting the eHMP Link




4. The eHMP application displays in the browser window.

Accessing eHMP through CPRS

1. Access **CPRS** as usual.
2. Select **eHMP** from the CPRS Tools dropdown menu (Figure 2-7).

Note: Not all sites have eHMP set up on the Tools Menu in CPRS.

Figure 2-7 eHMP Link on the CPRS Tools Menu

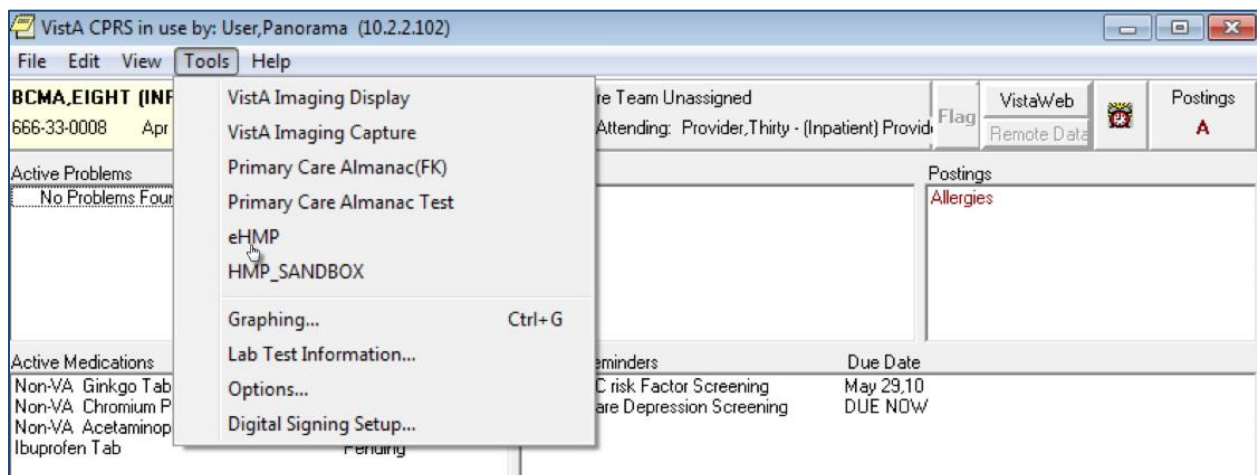

3. The eHMP application opens.

Note: If you accessed a patient's record in CPRS prior to logging into eHMP, that same patient's record is the default patient view displayed when you are automatically logged into eHMP. Likewise, if you access a particular patient in eHMP, that patient's record is the default view presented when you switch back to CPRS. This is available only if Sentillion client or CCOW has been installed on your workstation.


Logging Out of eHMP

There are two ways to log out of eHMP: signing out manually, or being automatically logged out due to inactivity.

Manual Sign Out

1. Select your **user name** on the header in the upper right corner, and a dropdown displays.
2. Select **Sign Out** (Figure 2-8) to log out of eHMP.


Figure 2-8 Sign Out Button on Header Dropdown

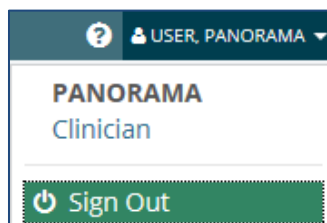

Auto Sign Out


When you have been inactive for 12 minutes, a warning message informing you your session ends in 3 minutes (Figure 2-9) displays. You are given the option to stay logged in by selecting the **Continue** button, or to log out of the application by selecting the **Logout** button.

Figure 2-9 Auto Logout Warning Message

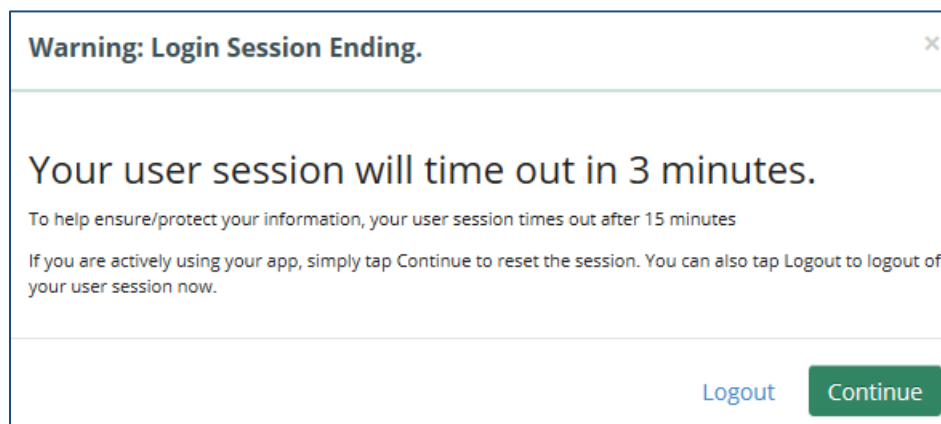

When you have been inactive in the eHMP application for 15 minutes or more, you are automatically logged out of the application, and a warning message (Figure 2-10) displays.

Figure 2-10 Inactivity Logout Message

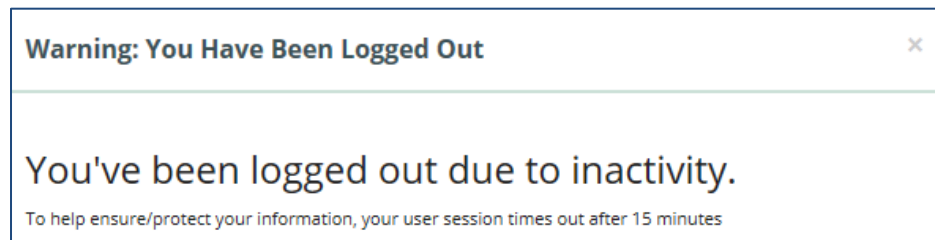


Chapter 3: Patient Search and Selection

About Patient Search and Selection

If you accessed eHMP through IE11, the Patient Selection screen displays (Figure 3-1), where you have the option to search for a patient.

Figure 3-1 Patient Selection Screen

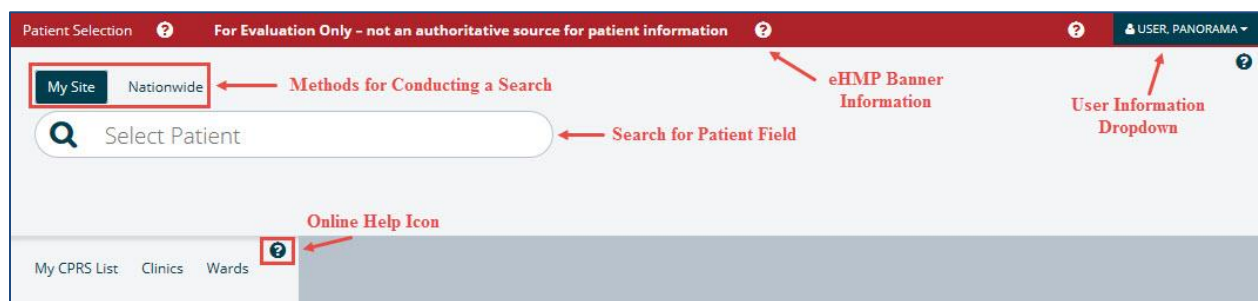

Patient Search

In order to access patient information and the detailed views of the eHMP application, you must conduct a patient search.

There are a couple of methods for conducting a patient search in eHMP:

- General Patient search
- My Site search
- Nationwide search

General Patient Search

To run a general patient search:

1. Enter the **patient's first initial of their last name followed by the last four digits of their social security number (SSN)**, or enter a **generic name** in the *Select Patient* field.
2. Select the **magnifying glass** or press **enter** to run the search (Figure 3-2).

Figure 3-2 General Patient Search

Patient Name	SSN	Date of Birth	Gender
EIGHT, PATIENT	***-**-0008	04/07/1935 (80y)	Male

My Site Search

The My Site search allows you to search for a patient from several perspectives, including a comprehensive list of all patients, from your defined default patient selection list in CPRS, patients from a specific clinic, or patients in a specific ward at your local VA facility.

My CPRS List

If you have defined a default patient selection list in CPRS, that same list of patients is presented in eHMP.

1. Select the **My CPRS List** tab. A list of patients displays.
2. Select a **patient** from the list to display the Patient Detail Dialog.

Clinics

To search for a patient using the **My Site/Clinics** method:

1. Select the **My Site** tab.
2. Select the **Clinics** tab. A list of clinics at the local VA facility displays.
3. Filter the list of results by entering the **name of the clinic** in the *Filter clinics* field. The results populate as you type.
4. Select the **clinic name** from the list. The list of patients from that clinic are displayed (Figure 3-3).

Figure 3-3 My Site/Clinics

Patient Name	SSN	Date/Time	DOB	Gender
No results found.				

5. Use the **preset date range** buttons to filter through the list of patients at the selected clinic.
6. Select a **patient's name** to display the Patient Detail Dialog.

Wards

To search for a patient (within a specific ward) using the **My Site/Wards** method:

1. Select the **My Site** tab.
2. Select the **Wards** tab. The list of wards at the local VA facility are displayed.
3. Filter the list of results by entering the **ward name** in the *Filter wards* field.
4. Select the **ward name** from the list. The list of patients in that ward are displayed (Figure 3-4).

Figure 3-4 My Site/Wards



Patient Name	SSN	Date of Birth	Gender	Room-Bed
EIGHT,PATIENT	***,**-0008	04/07/1935 (80y)	Male	722-B
EIGHTEEN,INPATIENT	***,**-0818	03/09/1945 (70y)	Male	730-B
ELEVEN,INPATIENT	***,**-0811	03/09/1945 (70y)	Male	727-A
FIFTEEN,INPATIENT	***,**-0815	03/09/1945 (70y)	Male	729-A
FIFTY,OUTPATIENT	***,**-0650	03/09/1945 (70y)	Male	735-A
FIVE,INPATIENT	***,**-0805	03/09/1945 (70y)	Male	724-C
FOUR,INPATIENT	***,**-0804	03/09/1945 (70y)	Male	724-B
FOURTEEN,INPATIENT	***,**-0814	03/09/1945 (70y)	Male	728-B
NINE,INPATIENT	***,**-0809	03/09/1945 (70y)	Male	726-C
NINETEEN,INPATIENT	***,**-0819	06/01/1995 (20y)	Female	731-A
ONE,INAPTIENT ONE	***,**-0801	03/09/1945 (70y)	Male	722-A
SEVEN,INPATIENT	***,**-0807	03/09/1945 (70y)	Male	726-A

5. Select a **patient's name** to display the Patient Detail Dialog.

Nationwide Search

The Nationwide search feature is used to run a patient search across the enterprise, including all VistA systems, as well as external systems, such as DoD and the Nationwide Health Information Network (NwHIN). This search allows you to enter the patient's last name (required) and one or more fields to enable search (first name, date of birth (DOB), and SSN).

To run a **Nationwide Patient** search:

1. Select the **Nationwide** tab (Figure 3-5).

Figure 3-5 Nationwide Search

My Site **Nationwide**

Nationwide Search
Enter the patient's first name, last name, and SSN to enable search.

2 Last Name (req) 3 First Name (req) DOB: MM/DD/YYYY 4 SSN: ###-##-#### (req) 5 Search

2. Enter the **patient's last name** in the *Last Name (required)* field.
3. Enter the **patient's first name** in the *First Name (req)* field.
4. Enter the **patient's SSN** in the *SSN: ###-##-#### (req)* field.
5. Select the **Search** button. The search results display below the search bar.

Note: Depending on the search results, the request may result in no patient found, 1 to 10 matching results found, or more than 10 records found, which returns a notification message in place of the records (Figure 3-6).

Figure 3-6 Too Many Results Have Returned

Patient Name	SSN	Date of Birth	Gender
Too many results have returned. Please be more specific in your search criteria.			

Selecting a Patient

Once the correct patient has been found, select that patient's name from the search results list. If a patient record is not restricted, their record can be displayed after confirming the selection (as demonstrated in the Accessing an Unrestricted Patient Record section) of that patient.

If a patient record is restricted, a warning appears prior to accessing the record, requiring acknowledgement of the restricted nature of the record before accessing it (as demonstrated in the Accessing a Restricted Patient Record section).

Accessing an Unrestricted Patient Record

To access an **unrestricted patient** record:

1. Select the **name of the patient** from the search results list. The Patient Detail Dialog appears (Figure 3-7).

Figure 3-7 Patient Detail Dialog

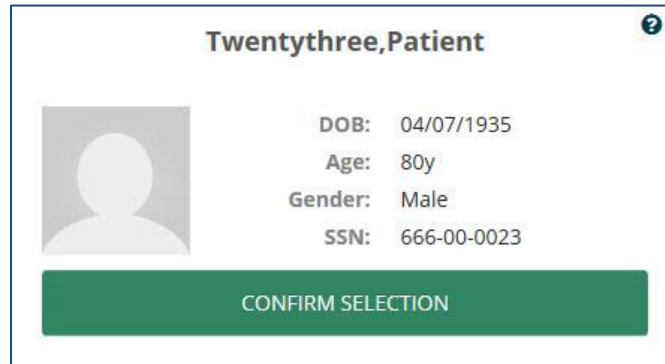


Figure 3-7 Patient Detail Dialog

The dialog displays the patient's name, "Twentythree, Patient", at the top. Below the name is a placeholder for a patient photo. To the right of the photo, the following information is displayed:

- DOB: 04/07/1935
- Age: 80y
- Gender: Male
- SSN: 666-00-0023

At the bottom of the dialog is a green button labeled "CONFIRM SELECTION".

2. Select the **CONFIRM SELECTION** button.
3. The patient's record displays in the Main Application Window.

Accessing a Restricted Patient Record

To access a **restricted patient** record:

1. Select the **name of the patient** from the search results list.

Note: The word *SENSITIVE* displays instead of the patient's SSN and DOB.

2. The Restricted Record Dialog appears (Figure 3-8).

Figure 3-8 Restricted Record Acknowledgement



Figure 3-8 Restricted Record Acknowledgement

The dialog displays the patient's name, "Zzzretfiveeightyfive, Patient", at the top. Below the name is a dark blue header bar labeled "RESTRICTED RECORD". Below the header bar, the following text is displayed:

This record is protected by the Privacy Act of 1974 and the Health Insurance Portability and Accountability Act of 1996. If you elect to proceed, you will be required to prove you have a need to know. Accessing this patient is tracked, and your station Security Officer will contact you for your justification.

At the bottom of the dialog is a green button labeled "ACKNOWLEDGE RESTRICTED RECORD".

The notification that appears in the dialog makes you aware that the patient's record is restricted, and advises that if you proceed to "break the glass" and open the record, your activity is tracked. In addition, your station Security Officer will contact you for your justification in accessing the restricted patient record.


3. Review the notification dialog, and select the **ACKNOWLEDGE RESTRICTED RECORD** button. The Restricted Record Patient Detail Dialog appears (Figure 3-9), displaying your acknowledgement of the restricted record.

Note: Select **ACKNOWLEDGED** ☒ to display the warning again.

Figure 3-9 Confirm Selection of Restricted Patient Record

Zzzretfiveeightyfive, Patient

RESTRICTED RECORD **ACKNOWLEDGED** ☒


DOB: 04/07/1935
Age: 80y
Gender: Male
SSN: ***-**-6886

CONFIRM SELECTION

4. Select the **CONFIRM SELECTION** button.
5. The patient's record displays in the Main Application Window.

Accessing a Flagged Patient Record

There are patients in the system with flagged records. The purpose of a patient record flag is to alert providers about certain patients, such as patients with disruptive, threatening, and/or violent behaviors, or research patients. The patient record flags can be national (Category I), which are shared among the VA facilities, or local (Category II), which are only shared at the local site. You must review and acknowledge the patient record flag during patient selection before accessing the patient record.

To access a **flagged patient** record:

1. Select the **name of the patient** from the search results list. The Patient Detail Dialog appears.
2. Select the **CONFIRM SELECTION** button. The Patient Flag Dialog appears (Figure 3-10).

Note: Some patients may have more than one patient record flag. Scroll through the patient flag dialog to review all of the notifications associated with a patient.

Figure 3-10 Confirm Selection of Flagged Patient Record

Eight, Patient

WANDERER **CATEGORY II (LOCAL)**

This Patient Likes To Wander Around The Hospital. Please Notify The Security Office

Initial Assigned Date: 12/23/2014 16:00
Approved by: PROGRAMMER,EIGHT
Assignment Status: Active
Owner Site: CAMP MASTER
Originating Site: CAMP MASTER

CONFIRM

3. Review the notification dialog, and Select the **CONFIRM** button.
4. The Main Application Window appears, displaying the patient record.

Main Application Window

Once the patient has been selected and confirmed, the patient record opens and you are navigated to the Main Application Window (Figure 3-11). The Main Application Window is comprised of the eHMP Header, Global Timeline Date Filter, Workspace Manager and listing, search record field, and the user workspace. Each of these areas is described in subsequent chapters.

Figure 3-11 Main Application Window

The screenshot displays the eHMP Main Application Window for a patient named Twentythree.Patient. The interface is organized into several sections:

- Header:** Includes patient selection, evaluation status, patient details (DOB, SSN, Gender), and provider information.
- Global Timeline Date Filter:** Set to 04/27/2014 - 04/27/2016.
- Workspace Manager:** Contains tabs for Overview, Search Record, and a search icon.
- Clinical Reminders:** A table listing reminders for Hepatitis C risk factor screening, Primary Care Depression Screening, Hypertension, Iraq/Afghan Post-Deployment Screen, and TBI Screening, all due on 10/27/2015.
- Encounters:** A table showing visits, appointments, admissions, and procedures, all with a last occurrence of None and 0 occurrences.
- Reports:** A table for reports with no records found.
- Conditions:** A table listing conditions such as Essential Hypertension (Disorder), Hyperlipidemia, and Acute Myocardial Infarction, with their acuity, last occurrence, and occurrence count.
- Active & Recent Medications:** A table showing medications like Methocarbamol 500 MG Oral Tablet and Aspirin 81 MG Enteric Coated Tablet, with their refills, status, and frequency.
- Immunizations:** A table showing immunizations like PNEUMOCOCCAL 15y.
- Allergies:** A table showing allergies like CHOCOLATE, CHOCOLATE, CHOCOLATE, PENICILLIN, PENICILLIN, and MILK.
- Vitals:** A table showing vital signs like BPS, BPD, Pulse, RR, and Temp, with no records found.
- Numeric Lab Results:** A table for lab results with no records found.

The footer shows the eHMP version 1.2.2.49391 and a Refresh All Data button.

Chapter 4: The eHMP Header

About the eHMP Header

The eHMP Header (Figure 4-1) displays abbreviated information. It lists the logged in user, allows for another patient search, and presents demographics about the displayed patient, any postings associated with the patient, the visit information for a patient with an inpatient status, and the patient's care team information.


The following elements make up the eHMP Header:

- Patient Selection button
- eHMP Banner Information
- User information dropdown
- Patient photo
- Patient demographics, inpatient/outpatient status, and clinical context object workgroup (CCOW) status
- Patient postings
- Visit Information (inpatient status only)
- Care team information

Figure 4-1 The eHMP Header



Patient Selection Button

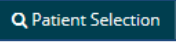The Patient Selection button  is located in the top left corner of the eHMP Header.

1. Select the **button** to navigate back to the Patient Selection screen to search for another patient.
2. Select the **X** in the upper right corner of the Patient Selection screen to cancel a search, and return to the current patient record.

eHMP Banner Information

The eHMP banner informs you that the application is for evaluation only, and should not be used as an authoritative source for patient information.

To display detailed information:


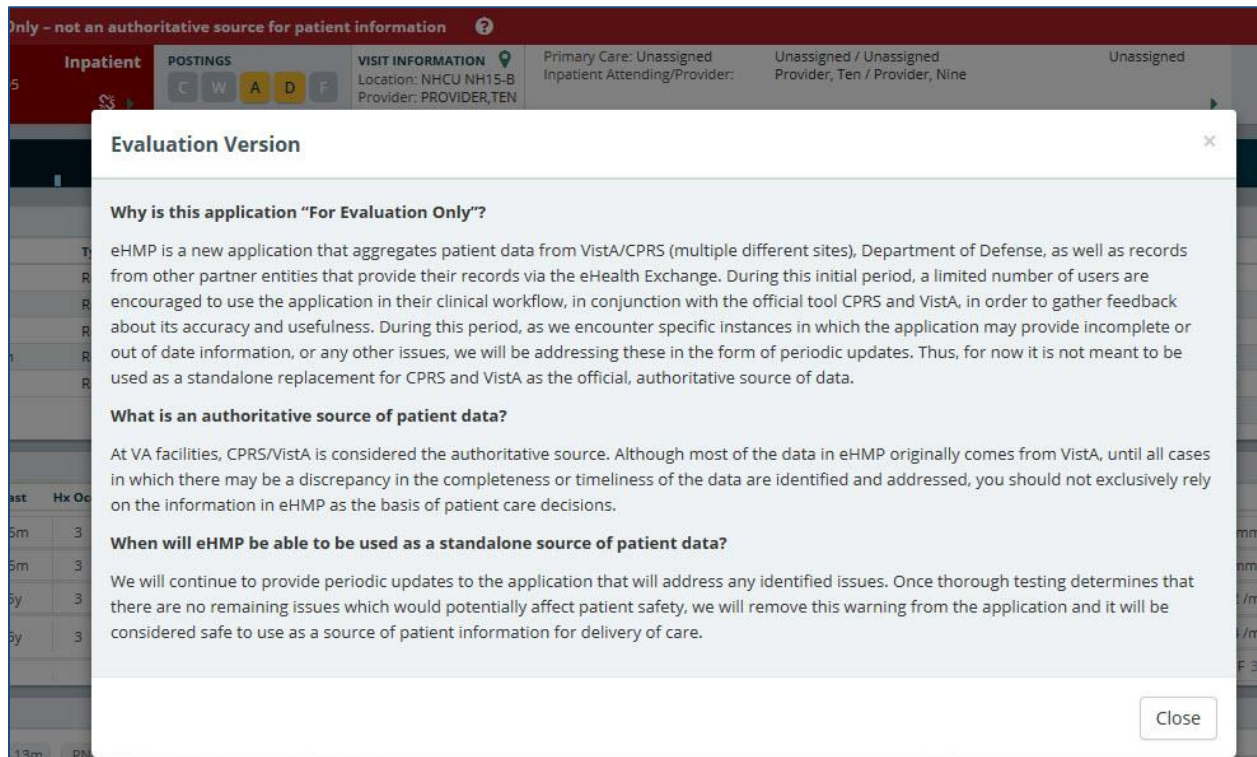
1. Select the online help button  next to the banner text. A Detail Dialog box displays with information regarding the eHMP banner (Figure 4-2).

Figure 4-2 Detailed eHMP Banner Information



2. Select the **Close** button to close the Detail Dialog box and return to the patient record.

User Information Dropdown

Your user information displays in the upper right corner of the eHMP application. As depicted in Figure 4-3, selecting **your name** results in the display of a dropdown menu with the *Sign Out* option.

Figure 4-3 Sign Out Option



Patient Photo

eHMP displays the patient image as found in VistA Imaging. The image is retrieved from the Veteran Health Identification Card (VHIC) system and is displayed in the header. If an image is unavailable for that patient, then a gender-neutral image displays (Figure 4-4).

Figure 4-4 Patient Photo (Example)




Patient Demographics

The selected patient is listed on the header in the upper left corner. By default, the patient's name, DOB, age, SSN, and Gender display (Figure 4-5). In addition, a status of inpatient or outpatient, and the status of CCOW is indicated. When eHMP is in patient context, a green chain link icon displays as displayed in Figure 4-5. A red demographics bar with a broken chain link indicates that eHMP is not in patient context as displayed in Figure 4-6.

Figure 4-5 Patient Demographics



Figure 4-6 Patient Demographics Bar with Broken Chain Link

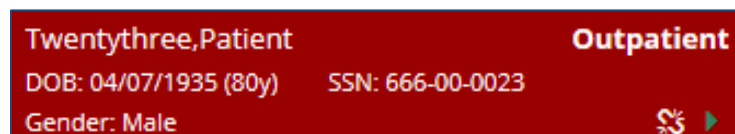

To display expanded patient information:

1. Select the **dropdown arrow** to open the Patient Demographics Dialog.
The dialog contains sections with additional patient information.
2. When the icon is displayed, non-local demographic data is available and can be displayed by selecting the **section** (Figure 4-7).

Figure 4-7 Expanded Patient Information

Twentythree, Patient **Outpatient** **POSTINGS** Primary Care: Green Mental Health: Mh Team Provider, Fifteen Vehu, One

DOB: 04/07/1935 (81y) SSN: 666-00-0023 Gender: Male

Phone

Home (843) 555-1234

Cell (843) 555-5678

Work (843) 555-2345

Addresses

Home Home Address Charleston, Sc, 29492

Temporary No Record Found

Email

23@EXAMPLE.COM

Emergency Contact

Sister Veteran, Sister

Home (843) 555-0987

Work (843) 555-9876

Address Sist Address Mount Pleasant, Sc, 29464

Health Benefits And Insurance

Service Connected Yes

Service Connected Conditions

DYSLEXIA (100%)

DEMENTIA (50%)

MULTIPLE SCLEROSIS (50%)

INSOMNIA (30%)

AUTISM (60%)

NARCOLEPSY (70%)

Insurance

Insurance PRIVATE

Name INSURANCE CO

INC

Service And Social History

Veteran Status Yes

Marital Status Legally Separated

Religion Baptist

Next Of Kin

Vista Site	Name	Relationship	Home	Work	Address
POLAROID	Veteran, Brother	Brother	(843) 555-8765	(843) 555-7654	Brot Address West Columbia, Sc, 29169
KODAK	Veteran, Brother	Relationship Unknown	No Record Found	No Record Found	No Record Found

Address Brot Address West Columbia, Sc, 29169

3. Select the **dropdown arrow** again or anywhere outside of the dialog box to exit the dialog and return to the previous view.


Patient Postings

The Postings section of the eHMP Header presents shortcuts to the following pieces of information:

- Crisis notes
- Warnings
- Allergies and Adverse Drug Reactions
- Directives
- Patient Record Flags

The Postings area of the eHMP header is depicted in Figure 4-8.

Figure 4-8 Postings Icons

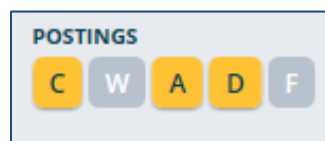

The highlighted letters in the Postings area indicate that detailed information can be accessed.

1. Selecting a **highlighted letter** opens a dialog with detailed information (Figure 4-9).

Figure 4-9 Expanded Postings Information

POSTING: Allergies

VISIT INFORMATION

Location: NHCU NH15-B
Provider: PROVIDER, TEN

Primary Care: Unassigned
Inpatient Attending/Provider: Unassigned
Provider: Unassigned

Allergies

MILK NAUSEA, VOMITING

Drug Classes

Nature of reaction	ALLERGY	Severity
Entered By:	PROGRAMMER TWENTY	Originated 01/19/1994 08:22
Verified	01/19/1994	Observed/Historical
Facility	New Jersey HCS	Observed Date
Comment	TESTING DIETETICS	

2. Selecting the **highlighted letter** again, or anywhere outside of the box, closes the dialog and returns to the previous view.

Visit Information

If the patient has an inpatient status, then the visit information displays with the ward location and the Attending Inpatient Provider (Figure 4-10). If the Attending Inpatient Provider does not exist, then it displays the Inpatient Provider.

Figure 4-10 Visit Information

VISIT INFORMATION

Location: NHCU NH15-B
Provider: PROVIDER, TEN

Care Team Information

Abbreviated Care Team Information for the selected patient is displayed in the eHMP Header (Figure 4-11). The information displayed here includes:

- The primary care team and provider
- The associate name (if available)

- The team phone number
- The mental health treatment team and treatment coordinator (if available). If the patient has an inpatient status, the inpatient attending and provider are listed here instead of the mental health treatment team and treatment coordinator.

Figure 4-11 Care Team Information

Primary Care: Green	Provider, Fifteen / Pcomm-resident, One	(555) 555-5858
Mental Health: Mh Team	Vehu, One	

To access detailed Care Team information:

1. Select the **dropdown arrow**. Detailed information displays.

Note: If the patient has visited other VistA sites, you are able to review the listing for each care team by the site for which the patient has a VistA record (Figure 4-12).

2. Select the **dropdown arrow** again to close the dialog.

Figure 4-12 Detailed Care Team Information

Primary Care: Unassigned		Unassigned / Unassigned		Unassigned	
Inpatient Attending/Provider:		Provider, Ten / Provider, Nine			
Provider Title		Name	Analog Pager	Digital Pager	Office Phone
Primary Care Provider		Unassigned	not specified	not specified	not specified
Primary Care Assoc Provider		Unassigned	not specified	not specified	not specified
Inpatient Attending Provider		Provider, Ten	not specified	not specified	not specified
Inpatient Provider		Provider, Nine	not specified	(555) 555-9998	not specified
MH Treatment T		Unassigned	not specified	not specified	not specified
MH Treatment C		Other Site Care Providers			
Facility		Name	Analog Pager	Digital Pager	Office Phone
15m	POLAROID	Provider, Nine	not specified	(555) 555-9998	not specified
15m	KODAK	Provider, Nine	not specified	(555) 555-9998	not specified
	HDR	Provider, Nine	not specified	(555) 555-9998	not specified
	VLER	Provider, Nine	not specified	(555) 555-9998	not specified


Chapter 5: Global Timeline Date Filter

About the Global Timeline Date Filter

The Global Timeline Date Filter controls the amount of historical data displayed throughout the application. It displays two years of patient information by default, and allows you to modify the date range to allow for easier investigation of records.

The summary view of the Global Timeline Date Filter displays the date range selected with a graphical representation. The bars (blue denotes outpatient and green denotes inpatient) depict the number of patient activities recorded for a specific date (Figure 5-1). Hovering your cursor over the date displays a tool tip with detailed information for the specified date. The red line on the graph represents the current date.


Figure 5-1 Global Timeline Date Filter

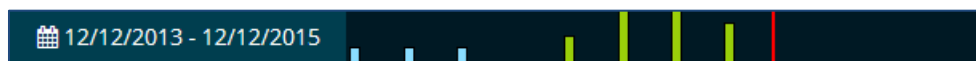

The list below describes the elements of the expanded Global Timeline Date Filter as displayed in Figure 5-2:

1. Preset date ranges
2. Custom date range
3. All Events timeline graph
4. Timeline Summary

Figure 5-2 Expanded Global Timeline Date Filter

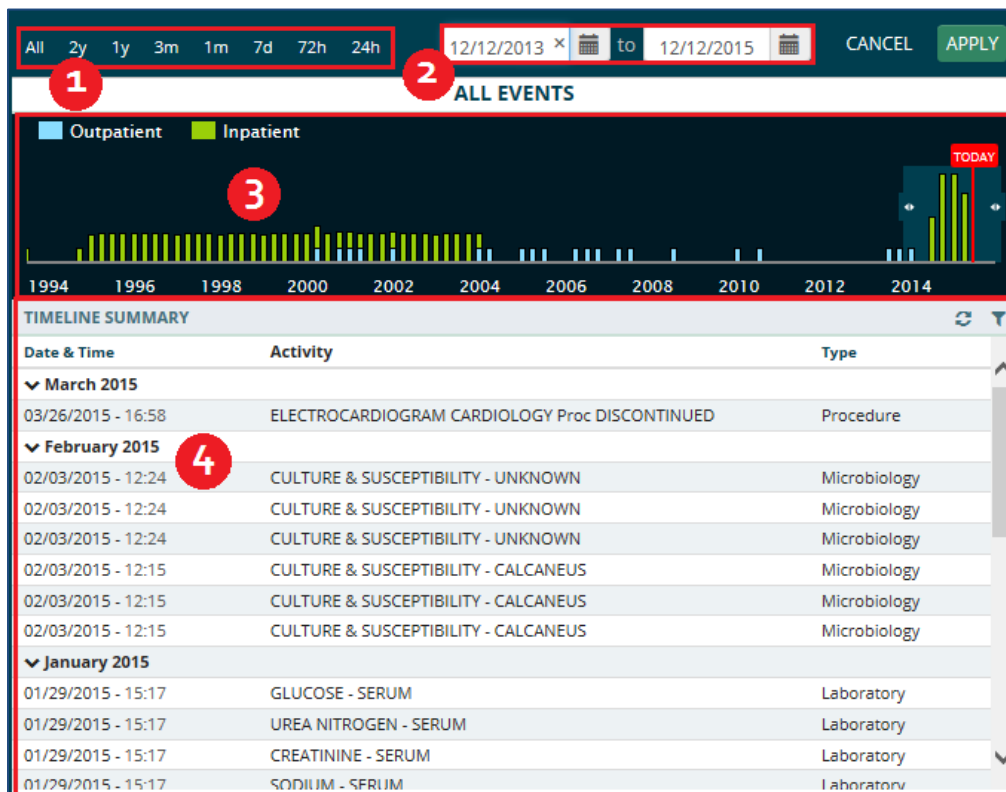


Using the Preset Date Ranges

To use preset date ranges:

1. Select the **date range** to open the Global Timeline Date Filter.
2. Select the desired **preset date range** (Figure 5-3). The dropdown calendar, All Events timeline graph, and Timeline Summary update with the selected date range.

Figure 5-3 Preset Date Ranges

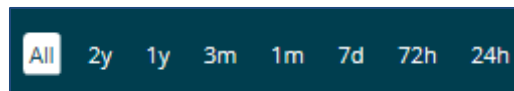

3. Select **APPLY** to set the date range or select **CANCEL** to exit the Global Timeline Date Filter without changing the dates. The patient information for the selected date range displays.

Note: If you select **outside** of the Global Timeline Date Filter box prior to selecting **APPLY**, the box closes and the selected dates are not applied.


Using the Custom Date Range

To use custom data ranges:

1. Select the **date range** to open the Global Timeline Date Filter.
2. Use the **calendar dropdown** (Figure 5-4) or enter the **dates** using the MM/DD/YYYY format to choose a custom date range.

Note: You can only select a **start date** that is in the past and an **end date** that is the current date or a date in the future.

Figure 5-4 Calendar Dropdown

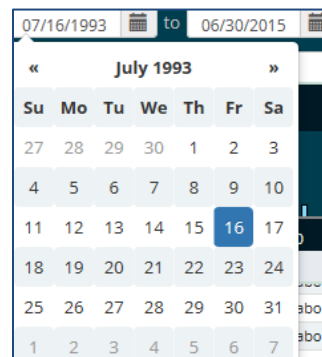


3. Select **APPLY** to set the date range, or select **CANCEL** to exit the Global Timeline Date Filter without changing the dates. The patient information for the selected date range displays.

Note: If you select **outside** of the Global Timeline Date Filter box prior to selecting **APPLY**, the box closes and the selected dates are not applied.

Using the All Events Timeline

The All Events Timeline (Figure 5-5) displays a graphical representation of events from the full patient historical record that includes both inpatient and outpatient information. It allows you to limit the data set more accurately to perform a quick analysis.

Figure 5-5 All Events Timeline

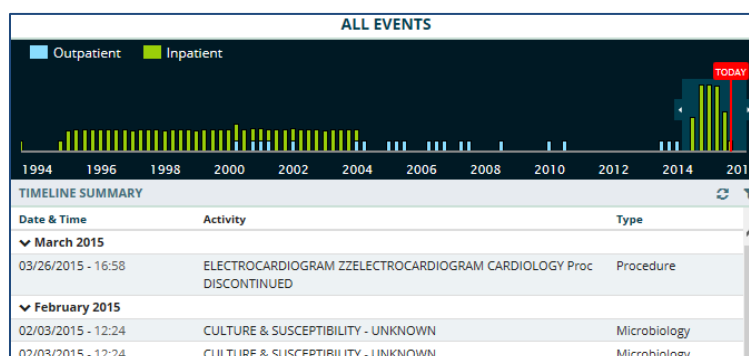

To use the all events timeline:

1. Select the **date range** to open the Global Timeline Date Filter.
2. Select and drag the **double arrows** on the left of the shaded timeline area to select the *start* date.
3. Select and drag the **double arrows** on the right of the shaded timeline area to select the *end* date. The dates in the calendar boxes above, as well as the Timeline Summary, populate as you select the date range.
4. Select **APPLY** to set the date range or select **CANCEL** to exit the Global Timeline Date Filter without changing the dates. The patient information for the selected date range displays.

Note: If there is no recent data, only the red **Today** reference line displays.

Note: If you select **outside** of the Global Timeline Date Filter box prior to selecting **APPLY**, the box closes and the selected dates are not applied.

Using the Timeline Summary

As the date range for the patient record changes, the Timeline Summary refreshes to display patient activities within the selected date range. The Timeline Summary is a list view grouped by month. The groups are collapsible and can be sorted by date and time, or type.

To use the timeline summary:

1. Select the **Filter** button. The *Enter your text filter* field displays.
2. Enter the **filter text**. The results populate as you enter the text (**Error! Reference source not found.**).

Figure 5-6 Filter Button

TIMELINE SUMMARY		
<input type="text" value="cholesterol"/>		
Date & Time	Activity	Type
▼ March 2010		
03/05/2010 - 10:00	LDL CHOLESTEROL - SERUM	Laboratory
03/05/2010 - 10:00	CHOLESTEROL - SERUM	Laboratory
03/05/2010 - 10:00	LDL CHOLESTEROL - SERUM	Laboratory
03/05/2010 - 10:00	CHOLESTEROL - SERUM	Laboratory
03/05/2010 - 10:00	LDL CHOLESTEROL - SERUM	Laboratory
03/05/2010 - 10:00	CHOLESTEROL - SERUM	Laboratory
▼ December 2009		
12/01/2009 - 12:00	CHOLESTEROL - SERUM	Laboratory
12/01/2009 - 12:00	CHOLESTEROL - SERUM	Laboratory
12/01/2009 - 12:00	CHOLESTEROL - SERUM	Laboratory
12/01/2009 - 10:00	LDL CHOLESTEROL - SERUM	Laboratory

3. Select **an activity** in the Timeline Summary to display the Detail Dialog box (**Error! Reference source not found.**).

Figure 5-7 Patient Activity Detail Dialog Box

GENERAL INTERNAL MEDICINE [X]

Date	10/10/2013 - 13:00
Type	Compensation & Pension
Category	Outpatient Visit
Patient Class	Ambulatory
Appointment Status	NO ACTION TAKEN
Location	General Medicine
Stop Code	GENERAL INTERNAL MEDICINE
Facility	CAMP MASTER

Providers

Additional Provider: Provider, Eight

Close

4. Select **Close** to close the Detail Dialog box and return to the previous view.

Chapter 6: Applet Views

Applets are small applications that run within eHMP and sorts patient data into segments throughout the patient record. Five applet views are available, which allows you to customize the display of patient data. The five applet views are described below:

- Trend View – a concise presentation of information pertinent to helping a clinician decide what to do next, often incorporating graphics or symbols to communicate at high level information.
- Summary View – a table format of the applet content.
- Expanded View – displays more columns and details than the trend and summary views.
- Detailed View – displays a popup window with detailed information for the selected item.
- Maximized View – displays in full screen with all information for the applet.

This chapter describes the different applets within eHMP and the available views for each applet. The applets are listed in alphabetical order.

Active & Recent Medications

The Active & Recent Medications applet provides you with a quick view of a patient's active and recent medications. If the patient does not have any medications for the defined date range, then a message indicating that no records have been found displays.

Note: If a patient has an inpatient status, then only active inpatient medications and IVs display. If a patient has an outpatient status, then only active outpatient and non-VA medications display.

Active & Recent Medications: Trend View

Figure 6-1 displays the trend view that lists the active and recent medication and dosage, any refills, the status of the medication, and when the medication was last filled. It also includes all medications that have expired within the last six months.

Figure 6-1 Active & Recent Medications Applet Trend View

ACTIVE & RECENT MEDICATIONS			
Medication	Refills	Status/Fillable	
gabapentin 300 MG Oral Capsule Give: 300MG PO Q6H	NA	Active	2m
gabapentin 600 MG Oral Tablet Give: 1200MG PO Q6H	NA	Active	2m
Methadone Hydrochloride 10 MG Oral Tablet Give: 15MG PO Q12H	NA	Active	2m
gabapentin 300 MG Oral Capsule Give: 300MG PO Q6H	NA	Active	2m

Refills Column

The refills column indicates the refill status for an active or recent medication.

- When there is no refill count information available, 'NA' displays (usually on DoD medications).
- When a medication has zero refills remaining, a '0', surrounded by a red highlighted box, displays.
- If a medication has one refill left, a '1', surrounded by an orange highlighted box, displays.
- Medications that have two or more refills remaining displays the number of refills left.

Status/Fillable Column

The status/fillable column indicates the status of an active or recent medication.

- 'Active' displays when the medication is active.
- 'Expired' displays in red for a medication that is about to expire, or has expired.
- 'New' displays if the medication is relatively recent and pending.

Note: Hover your cursor over the status button to display a tooltip with more information on the status of the medication.

To display more information for a medication in the trend view:

1. Select the **left side of the medication tile** to display the applet toolbar menu (Figure 6-2).

Figure 6-2 Active & Recent Medications Applet Toolbar Menu – Trend View






2. Select the **Infobutton**  to open a new browser window that displays specific medication information.
3. Select the **Details form** button  to open the Detail Dialog box for the selected medication.
4. Select the **Quicklook** button  to display expanded information for the selected medication.
5. Select the **right side of the medication tile** to display a quick view containing up to the last five fills for the selected medication consisting of the last update date, medication name with dosage, sig, and time since last change (Figure 6-3).

Figure 6-3 Active & Recent Medications Quick View

ACTIVE & RECENT MEDICATIONS				VITALS									
Medication	Refills	Status/Fillable		Type	Result								
gabapentin 300 MG Oral Capsule Give: 300MG PO Q6H	NA	Active	2m	BPS	158 mm[Hg]								
gabapentin 600 MG Oral Tablet Give: 1200MG PO Q6H	NA	Active	2m	BPD	91 mm[Hg]								
Methadone Hydrochloride 10 MG Oral Tablet Give: 15MG PO Q12H	NA	Active	2m	Pulse	82 /min								
gabapentin 300 MG Oral Capsule Give: 300MG PO Q6H	NA			BP	144 /min								
				<table><tr><th>Last Update</th><th>Medication</th><th>Sig</th><th>Since</th></tr><tr><td>03/15/2016</td><td>Methadone Hydrochloride 10 MG Oral Tablet</td><td>Give: 15MG PO Q12H</td><td>2m</td></tr></table>		Last Update	Medication	Sig	Since	03/15/2016	Methadone Hydrochloride 10 MG Oral Tablet	Give: 15MG PO Q12H	2m
Last Update	Medication	Sig	Since										
03/15/2016	Methadone Hydrochloride 10 MG Oral Tablet	Give: 15MG PO Q12H	2m										
ALLERGIES													

Active & Recent Medications: Summary View

Figure 6-4 depicts the summary view listing the medication and dosage, and the facility where the medication was prescribed.

Figure 6-4 Active & Recent Medications Applet Summary View



ACTIVE & RECENT MEDICATIONS	
Medication	Facility
DOCUSATE NA CAP,ORAL (ACTIVE) Give: 100MG PO QDAY	BAY
DOCUSATE NA CAP,ORAL (ACTIVE) Give: 100MG PO QDAY	BAY
DOCUSATE NA CAP,ORAL (EXPIRED) Give: 100MG PO QDAY	BAY
GABAPENTIN TAB (DISCONTINUED) Give: 1200MG PO Q6H	BAY
GABAPENTIN TAB (DISCONTINUED) Give: 1200MG PO Q6H	BAY

To display more information for a medication in the summary view:

1. Select a **medication** to display the applet toolbar menu (Figure 6-5).

Figure 6-5 Active & Recent Medications Applet Toolbar Menu – Summary View

ACTIVE & RECENT MEDICATIONS	
<div> <div>i</div> <div>Details</div> </div>	Facility
GABAPENTIN CAP,ORAL (ACTIVE) Give: 300MG PO Q6H	BAY
GABAPENTIN TAB (ACTIVE) Give: 1200MG PO Q6H	BAY
METHADONE TAB (ACTIVE) Give: 15MG PO Q12H	BAY
GABAPENTIN CAP,ORAL (ACTIVE) Give: 300MG PO Q6H	BAY
GABAPENTIN TAB (ACTIVE) Give: 1200MG PO Q6H	BAY
METHADONE TAB (ACTIVE) Give: 15MG PO Q12H	BAY


2. Select the **Info** button  to open a new browser window that displays specific medication information.
3. Select the **Details form** button  to open the Detail Dialog box for the selected medication.

Active & Recent Medications: Detail View

To display a detailed view of an active or recent medication from the trend or summary view:

1. Select an **active or recent medication** from the list to display the Details form button (Figure 6-6).

Figure 6-6 Active & Recent Medications Details Form Button

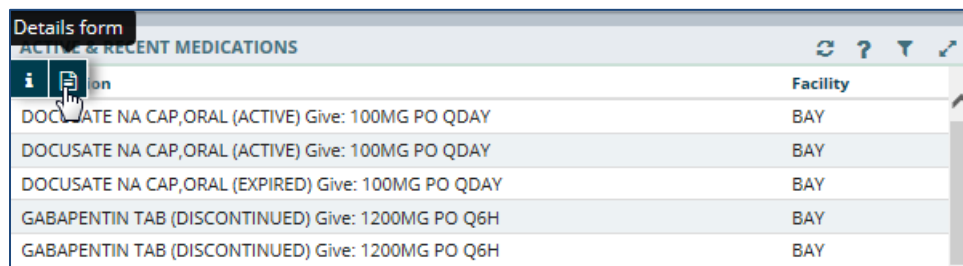

The screenshot shows a table titled 'ACTIVE & RECENT MEDICATIONS'. A tooltip labeled 'Details form' is pointing to an icon in the first column. The table has two columns: 'Medication' and 'Facility'.

Medication	Facility
DOCUSATE NA CAP,ORAL (ACTIVE) Give: 100MG PO QDAY	BAY
DOCUSATE NA CAP,ORAL (ACTIVE) Give: 100MG PO QDAY	BAY
DOCUSATE NA CAP,ORAL (EXPIRED) Give: 100MG PO QDAY	BAY
GABAPENTIN TAB (DISCONTINUED) Give: 1200MG PO Q6H	BAY
GABAPENTIN TAB (DISCONTINUED) Give: 1200MG PO Q6H	BAY


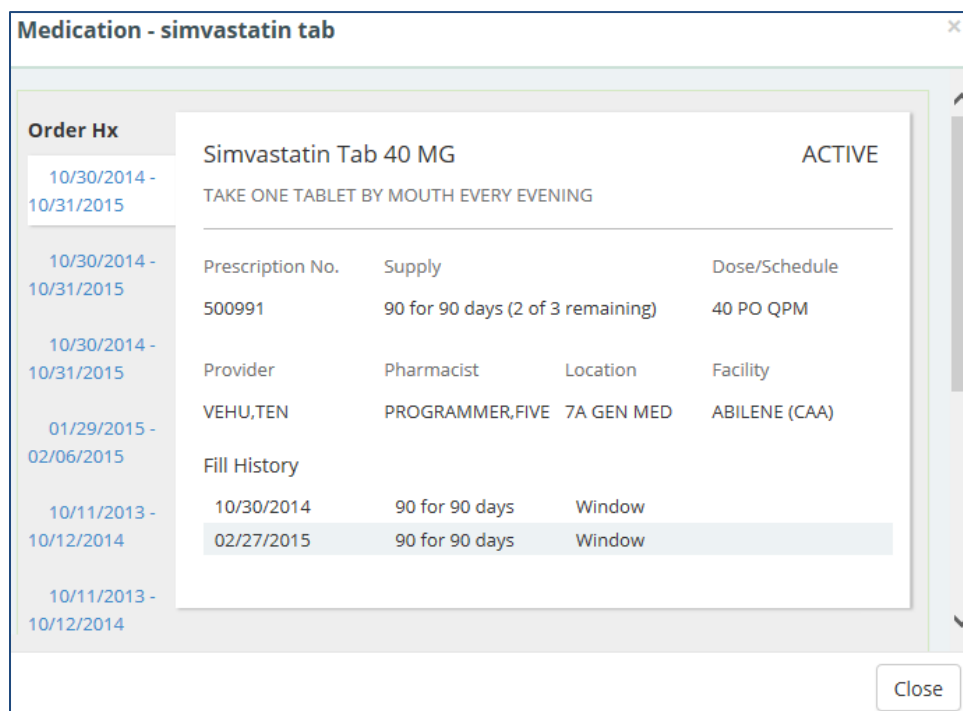
2. Select the **Details form** button  to open the Active & Recent Medication Detail Dialog box (Figure 6-7).

Figure 6-7 Active & Recent Medications Detail Dialog Box



The screenshot shows a dialog box titled 'Medication - simvastatin tab'. It displays detailed information for 'Simvastatin Tab 40 MG' which is 'ACTIVE'. The dialog includes sections for 'Order Hx', 'Prescription No.', 'Supply', 'Dose/Schedule', 'Provider', 'Pharmacist', 'Location', 'Facility', and 'Fill History'.

Order Hx		Simvastatin Tab 40 MG		ACTIVE
10/30/2014 - 10/31/2015		TAKE ONE TABLET BY MOUTH EVERY EVENING		
10/30/2014 - 10/31/2015	Prescription No.	Supply	Dose/Schedule	
	500991	90 for 90 days (2 of 3 remaining)	40 PO QPM	
10/30/2014 - 10/31/2015	Provider	Pharmacist	Location	Facility
	VEHU,TEN	PROGRAMMER,FIVE	7A GEN MED	ABILENE (CAA)
01/29/2015 - 02/06/2015	Fill History			
10/11/2013 - 10/12/2014	10/30/2014	90 for 90 days	Window	
	02/27/2015	90 for 90 days	Window	
10/11/2013 - 10/12/2014				

Close

All available orders for the medication display on tabs to the left under Order Hx.

1. Select a **date tab** below the *Order Hx* section to display detailed medication information for the selected date range.
2. Select a **resource link** below the *Links* and *Patient Education* sections, to open a new window with more information on the selected medication.

3. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Active & Recent Medication Dialog box, to return to the applet.

Allergies

The Allergies applet provides a trend, summary, and expanded view of identified patient allergies and adverse drug reactions, with standardized coding references from all sources. If the patient does not have any allergies, then a message indicating that no records have been found displays.

Note: *It is common for a patient to have more than one instance of the same allergy listed.*

Allergies: Trend View

Top allergy information for a given patient is displayed in pill format (Figure 6-8). A red pill indicates a severe allergy, an orange pill indicates a moderate allergy, a white pill indicates a mild allergy, and a gray pill indicates that the severity is not noted.

Figure 6-8 Allergies Applet Trend View



Allergies: Summary View

The summary view (Figure 6-9) lists the allergen name, reaction, and level of severity for each allergen. The default summary view lists the allergies sorted first by severity, and then by the date the allergy was entered into the VA system.

Figure 6-9 Allergies Applet Summary View

ALLERGIES			🔄 ? ⚙️ ↗️
Allergen Name	Reaction	Severity	
NUTS	NAUSEA,VOMITING	Severe	
CHOCOLATE	ANXIETY; ITCHING,WATERING EYES; DROWSINESS; NAUSEA,VOMITING; DIARRHEA; RASH	Moderate	
DUST	DROWSINESS; ITCHING,WATERING EYES	Moderate	
GRASS POLLEN	ITCHING,WATERING EYES	Mild	

Allergies: Expanded View

The default expanded view lists the allergies sorted first by severity, and then by the date the allergy was entered into the VA system.

The expanded view of Allergies (Figure 6-10) displays the following information in sortable columns:

- Allergen Name
- Standardized Allergen
- Reaction
- Severity
- Drug class
- Entered By
- Facility

Figure 6-10 Allergies Applet Expanded View

ALLERGIES						
Allergen Name	Standardized Allergen	Reaction	Severity	Drug Class	Entered By	Facility
NUTS	Nuts	NAUSEA,VOMITING	Severe		USER,PANORAMA	CAMP MASTER
CHOCOLATE	Chocolate	ANXIETY; ITCHING,WATERING EYES; DROWSINESS; NAUSEA,VOMITING; DIARRHEA; RASH	Moderate		USER,PANORAMA	CAMP MASTER
DUST	House dust (Greer Labs) extract	DROWSINESS; ITCHING,WATERING EYES	Moderate		USER,PANORAMA	CAMP MASTER
GRASS POLLEN	Grass pollen	ITCHING,WATERING EYES	Mild		USER,PANORAMA	CAMP MASTER

Allergies: Detail View

To display a detailed view of an allergy from the trend, summary, or expanded views:


1. Select an **allergy** to display the applet toolbar menu.
2. Select the **Details form** button  to open the Allergy Detail Dialog box (Figure 6-11).

Figure 6-11 Allergy Detail Dialog Box

Allergen - GRASS POLLEN

Previous
Next

GRASS POLLEN

Symptoms: ITCHING,WATERING EYES

Mild

Drug Classes:
Nature of Reaction: Adverse Reaction
Entered By: USER,PANORAMA
Originated: 08/24/2015 - 10:18
Verified: 08/24/2015 - 10:19
Observed/Historical: Observed
Observed Date: 08/24/2015 - 00:00
Site: CAMP MASTER

Comments:

Close

3. Select the **Previous** and **Next** buttons to navigate through the allergies.
4. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Allergy Detail Dialog box, to return to the applet.

Appointments & Visits

The Appointments & Visits applet lists any future outpatient or specialty care appointments, and past clinic visits for a selected patient at VA and DoD facilities.

Appointments & Visits: Summary View

Figure 6-12 displays the summary view listing the appointment dates, encounter descriptions, the locations where the patient was seen, the status of the appointment, and the facility they visited.

Figure 6-12 Appointments & Visits Applet Summary View

APPOINTMENTS & VISITS				
Date	Description	Location	Status	Facility
03/14/2015 - 16:58	Visit	NHCU		TST1
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE		TST1
08/31/2014 - 14:00	Visit	PRIMARY CARE		TST2
08/31/2014 - 14:00	Visit	PRIMARY CARE		TST1
08/06/2014 - 11:47	Visit	DERMATOLOGY		TST1
08/06/2014 - 11:47	Visit	DERMATOLOGY		TST2
08/06/2014 - 11:47	Visit	DERMATOLOGY		TST1

Appointments & Visits: Expanded View

Figure 6-13 depicts the expanded view of the Appointments & Visits applet. The following information is displayed in sortable columns:

- Date
- Description
- Location
- Status
- Type
- Provider
- Reason
- Facility

Figure 6-13 Appointments & Visits Applet Expanded View

APPOINTMENTS & VISITS

All2yr1yr3mo1mo7d72hr24hr02/26/2014to02/26/2016Apply

Q Enter your text filterAdd

Source: All VA + DOD

Date	Description	Location	Status	Type	Provider	Reason	Facility
03/14/2015 - 16:58	Visit	NHCU		Daily Hospitalization Data	Nurse, Eighteen		TST1
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	Regular	Provider, Fifteen		BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	Regular	Provider, Fifteen		BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE	SCHEDULED/KEPT	Regular	Provider, Fifteen		BAY
08/31/2014 - 14:00	Visit	PRIMARY CARE		Office/outpatient Visit Est	Provider, Fifteen	UNSPECIFIED INFECTION OF LOWER LEG BONE	TST1

The default view contains both VA and DoD data, but you can change the Source to filter the results. The Source filtering labels include Local VA, All VA, and All VA + DoD.

To filter the data:

1. Select the **Source** dropdown menu to review a list of data filtering labels.
2. Select the **desired source**.
3. Enter **filter text** in the *Enter your text filter* field.
4. Use the **preset date ranges** or **calendar dropdown** to select a specific date range.

Appointments & Visits: Detail View

To display the details of an appointment or visit (Figure 6-14) from the summary or expanded applet:

1. Select an **appointment or visit** from the list to display the Appointments & Visits Detail Dialog box.
2. Select the **Next** and **Previous** buttons to navigate between appointments and visits.
3. Select the **Close** button, the **X** button in the upper right-hand corner, or anywhere outside of the Appointments & Visits Detail Dialog box, to return to the applet.

Figure 6-14 Appointments & Visits Detail Dialog Box

Visit

↑ Previous

↓ Next

X

Date: 08/31/2014 - 14:00

Type: Regular

Description: Visit

Patient Class: Ambulatory

Location: PRIMARY CARE

Status: SCHEDULED/KEPT

Stop Code: PRIMARY CARE/MEDICINE

Provider: Provider, Fifteen

Facility: BAY

Reason:

Close

Clinical Reminders

The Clinical Reminders applet provides a summary and expanded view of the clinical reminders in effect for a patient from all sites. These reminders advise you when upcoming events are due for the patient.

Clinical Reminders: Summary View

Figure 6-15 displays the summary view that lists the priority (if any), title, type, and due date for any clinical reminders.

Figure 6-15 Clinical Reminders Applet Summary View

CLINICAL REMINDERS			
Priority	Title	Type	Due Date
None	Hepatitis C risk Factor Screening	Reminder	6/1/2015
None	Primary Care Depression Screening	Reminder	6/1/2015
None	Hypertension	Reminder	6/1/2015
None	Hypertension and BP>140/90	Reminder	6/1/2015
None	Iraq&Afghan Post-Deployment Screen	Reminder	6/1/2015

Clinical Reminders: Expanded View

Figure 6-16 depicts the expanded view of the Clinical Reminders applet. The following information is displayed in sortable columns:

- Priority
- Title
- Type
- Due Date
- Done Date

Figure 6-16 Clinical Reminders Applet Expanded View

CLINICAL REMINDERS				
Priority	Title	Type	Due Date	Done Date
None	Hepatitis C risk Factor Screening	Reminder	8/25/2015	
None	Primary Care Depression Screening	Reminder	8/25/2015	
None	Hypertension	Reminder	8/25/2015	
None	Hypertension and BP>140/90	Reminder	8/25/2015	
None	Iraq&Afghan Post-Deployment Screen	Reminder	8/25/2015	
None	TBI Screening	Reminder	8/25/2015	

Clinical Reminders: Detail View

To display the details of a clinical reminder from the summary or expanded applet:

1. Select a **clinical reminder** from the list to display the Clinical Reminders Detail Dialog box (Figure 6-17).

Figure 6-17 Clinical Reminders Detail Dialog Box

2. Select the **Close** button, the **X** button in the upper right-hand corner, or anywhere outside of the Clinical Reminders Detail Dialog box, to return to the applet.

Community Health Summaries

The Community Health Summaries applet displays the selected patient's clinical health care summaries received from VA's external Health Information Exchange (HIE) partners who participate in the NwHIN.

Community Health Summaries: Summary View

The summary view (Figure 6-18) lists the date of service and authoring institution(s) for each care summary.

Figure 6-18 Community Health Summaries Applet Summary View

COMMUNITY HEALTH SUMMARIES	
Date	Authoring Institution
06/17/2014	Kaiser Permanente Mid-Atlantic STSTMA2
06/17/2014	Inland Northwest Health Services
06/17/2014	Regenstrief Institute Clinic
06/17/2014	HAWAII PACIFIC HEALTH SA
06/17/2014	Kaiser Permanente Southern California - RESC
06/17/2014	Conemaugh Health System
03/11/2014	Allscripts CCDA Example
11/16/2014	Epic CCDA Example 1
05/17/2014	Cerner CCDA Example

Community Health Summaries: Expanded View

2. Scroll through the document or use the **Table of Contents** hyperlinks to go directly to a section.
3. Select the **Next** and **Previous** buttons to navigate between Community Health Summaries.
4. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Community Health Summaries Detail Dialog box, to return to the expanded applet.

Documents

The Documents applet lists multiple categories of documentation from various sources. The document categories include clinical notes, discharge summaries, advanced directives, crisis notes, warnings, lab results, and imaging and radiology reports.

Documents: Summary View

The summary view (Figure 6-21) lists the document date, type of document, and who entered the document in the system.

Figure 6-21 Documents Applet Summary View

DOCUMENTS		
Date	Type	Entered By
▼ August 2014		
08/04/2014	Consult	Provider,Seven
08/02/2014	Consult	Provider,Twenty
08/09/2014	Procedure	Provider,Seven
08/08/2014	Consult	Provider,Twenty
08/04/2014	Consult	Provider,Seven
08/02/2014	Consult	Provider,Twenty

Documents: Expanded View

The expanded view of the Documents applet (Figure 6-22) displays the following information in sortable columns:

- Date
- Description
- Type
- Entered By
- Facility

Figure 6-22 Documents Applet Expanded View

DOCUMENTS				
Date	Description	Type	Entered By	Facility
▼ August 2014				
08/04/2014	DERMATOLOGY Cons	Consult	Provider,Seven	CAMP MASTER
08/02/2014	ORTHOPEDIC SURGERY Cons	Consult	Provider,Twenty	CAMP MASTER
08/09/2014	EEG CP ELECTROENCEPHALOGRAPHY Cons	Procedure	Provider,Seven	CAMP MASTER
08/08/2014	OCCUPATIONAL THERAPY Cons	Consult	Provider,Twenty	CAMP MASTER
08/04/2014	DERMATOLOGY Cons	Consult	Provider,Seven	CAMP BEE
08/02/2014	ORTHOPEDIC SURGERY Cons	Consult	Provider,Twenty	CAMP BEE
08/09/2014	EEG CP ELECTROENCEPHALOGRAPHY Cons	Procedure	Provider,Seven	CAMP BEE
08/08/2014	OCCUPATIONAL THERAPY Cons	Consult	Provider,Twenty	CAMP BEE
08/04/2014	DERMATOLOGY Cons	Consult	Provider,Seven	CAMP MASTER
08/02/2014	ORTHOPEDIC SURGERY Cons	Consult	Provider,Twenty	CAMP MASTER
08/09/2014	EEG CP ELECTROENCEPHALOGRAPHY Cons	Procedure	Provider,Seven	CAMP MASTER

Documents: Detail View

To display the details of an item in the summary and expanded list views:

1. Select a **document** in the list. The Document Detail Dialog box (Figure 6-23) opens.

Figure 6-23 Documents Detail Dialog Box

Ir microbiology report Details

Facility CAMP MASTER

Author None

Status Completed

Date/Time 02/03/2015 - 12:24

Accession [UID]: MI 15 4 [1315000004]

Received: Feb 03, 2015@12:24

Collection sample: UNKNOWN

Collection date: Feb 03, 2015 12:24

Provider: PROVIDER,TWENTY

Test(s) ordered: CULTURE & SUSCEPTIBILITY..... completed: Feb 03, 2015

* BACTERIOLOGY FINAL REPORT => Feb 03, 2015 14:20 TECH CODE: 119

GRAM STAIN:

CULTURE RESULTS: ACINETOBACTER ANITRATUS - Quantity: >25,000 - <50,000 CFU/ML

ANTIBIOTIC SUSCEPTIBILITY TEST RESULTS:

ACINETOBACTER ANITRATUS

:

AMIKACAN..... R

CEFAZOLIN..... S

CEFOXITIN..... S

CEFOTAXIME..... S

GENTAMICIN..... R

TOBRMCM..... R

TRMSULF..... R

VANCOMYCIN..... R

RESTRICTED

Close

2. Select the **Close** button, the **X** in the upper right-hand corner, to close the Documents Detail Dialog box, and return to the document list.

Encounters

The Encounters applet provides a high-level view of the patients' encounters for outpatient visits, appointments, admissions, and procedures to quickly understand the types of care the patient has received.

Encounters: Trend View

Figure 6-24 displays the trend view of the Encounters applet, with data grouped by the following encounter type: visits, appointments, admissions, and procedures. In addition, it displays when the patient was last seen for the encounter, and the history for each encounter type, with a graphical representation.

Figure 6-24 Encounters Applet Trend View

ENCOUNTERS			
Encounter	Last	Hx Occurrence	
▸ Visits	17m	4	
Appointments	None	0	
▸ Admissions	10m	3	
▸ Procedures	2m	1	

Note: Select the **dropdown arrow** to display a list of encounters for the selected group (Figure 6-25).

Figure 6-25 Encounters List in Trend View (Expanded)


ENCOUNTERS											
Encounter	Last	Hx Occurrence									
▸ Visits	17m	16									
▾ Appointments	18m	3									
<table> <tr> <th>Appointment Type</th><th>Last</th><th>Hx Occurrence</th><th></th></tr> <tr> <td>GENERAL INTERNAL MEDICINE</td><td>18m</td><td>3</td><td></td></tr> </table>				Appointment Type	Last	Hx Occurrence		GENERAL INTERNAL MEDICINE	18m	3	
Appointment Type	Last	Hx Occurrence									
GENERAL INTERNAL MEDICINE	18m	3									
▸ Admissions	10m	3									
▸ Procedures	2m	1									

Encounters: Detail View

To display the details of an item in the trend list view:

1. Select the **Encounter group dropdown** arrow.
2. A list of encounter types for that group display.
3. Select the desired **encounter type** and the applet toolbar menu displays (Figure 6-26).

Figure 6-26 Encounters Applet Toolbar Menu

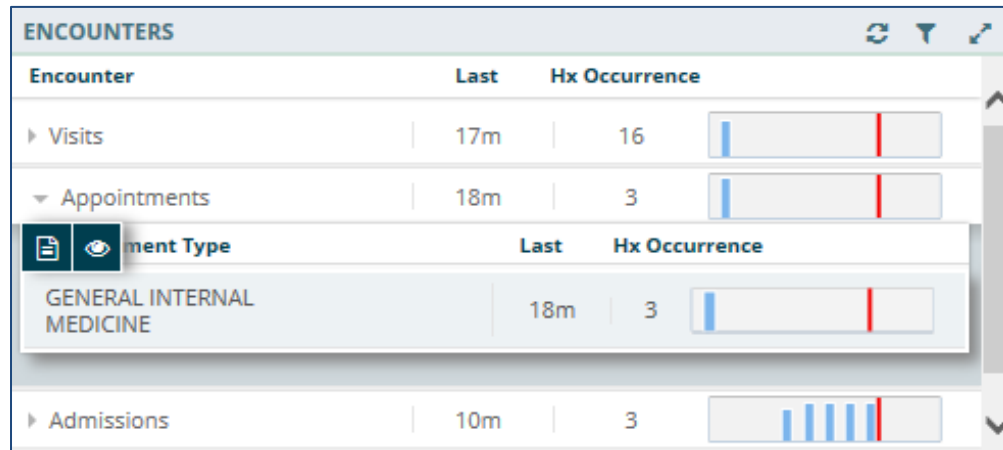

4. Select the **Details form** button to open the Encounters Detail Dialog box (Figure 6-27), or the **Quicklook** button to open the Encounters Detail list of the last five occurrences for the selected item.

Figure 6-27 Encounters Detail Dialog Box

Hospitalization

Date: 08/14/2014 - 13:07
 Type: Hospitalization
 Category: Admission
 Patient Class: Inpatient
 Location: 7A Gen Med
 Stop Code:
 Facility: ABILENE (CAA)

Providers

Additional Provider: Provider,Thirty
 Primary: Provider,Twenty

Reason

R/O MI

Close

5. Select the **Close** button, the **X** in the upper right-hand corner, to close the Encounters Detail Dialog box and return to the previous view.


Immunizations

The Immunizations applet presents a list of vaccines from all sources for a given patient.

Immunizations: Trend View


Figure 6-28 displays the trend view of immunizations a patient has received. All immunizations are represented in pill format. Each pill displays the immunization name, the series number (if available), and the date the immunization was last administered.

Figure 6-28 Immunizations Applet Trend View

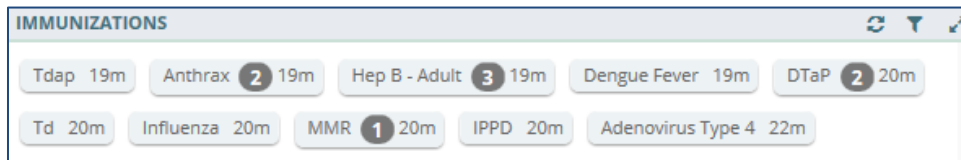

Detailed information for an immunization can be displayed by hovering your cursor over a pill (Figure 6-29).

Figure 6-29 Detail of an Immunization in Trend View

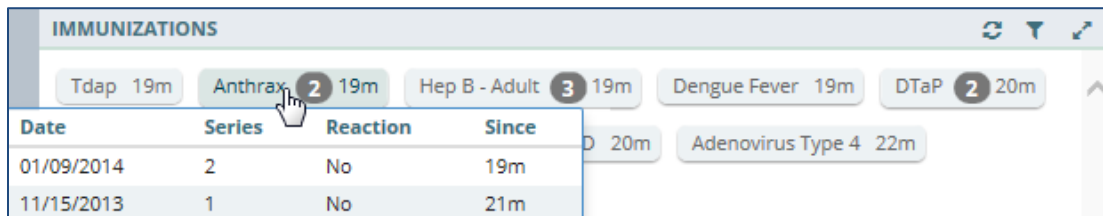

Immunizations: Summary View

Figure 6-30 displays the summary view that lists the vaccine name, reaction (if any), date administered, and facility.

Figure 6-30 Immunizations Applet Summary View

IMMUNIZATIONS			
Vaccine Name	Reaction	Date	Facility
Tdap		01/13/2014	DOD
Anthrax		01/09/2014	DOD
Hep B - Adult		01/09/2014	DOD
Dengue Fever		01/09/2014	DOD
DTaP		12/05/2013	DOD
Td		12/03/2013	DOD
Influenza		12/03/2013	DOD
Hep B - Adult		12/03/2013	DOD

Immunizations: Expanded View

The expanded view of Immunizations (Figure 6-31) displays the following information in sortable columns:

- Vaccine Name
- Standardized Name
- Reaction
- Series
- Repeat Contraindicated
- Date
- Facility

Figure 6-31 Immunizations Applet Expanded View

IMMUNIZATIONS							↺ ? ↻ ×
Vaccine Name	Standardized Name	Reaction	Series	Repeat Contraindicated	Date	Facility	
Tdap	tetanus toxoid, reduced diphtheria toxoid, and acellular pertussis vaccine		0	No	01/13/2014	DOD	
Anthrax	anthrax vaccine		2	No	01/09/2014	DOD	
Hep B - Adult	hepatitis B vaccine, adult dosage		3	No	01/09/2014	DOD	
Dengue Fever			0	No	01/09/2014	DOD	
DTaP	diphtheria, tetanus toxoids and acellular pertussis vaccine		2	No	12/05/2013	DOD	
Td			0	No	12/03/2013	DOD	

Immunizations: Detail View

To display the details of an Immunization from the trend, summary, or expanded views:


1. Select an **Immunization** from the list to display the applet toolbar menu.
2. Select the **Details form** button  to display the Vaccine Detail Dialog box (Figure 6-32).

Figure 6-32 Immunizations Detail Dialog Box

Vaccine - Anthrax

↑ Previous

↓ Next

×

Name	Reaction	Series	Repeat Contraindicated	Date	Facility	Site
Anthrax		2	No	01/09/2014	DOD	DOD

Viewing 12/01/2013 to 12/01/2015

All

2yr

1yr

3mo

1mo

7d

72hr

24hr

12/01/2013

📅

to

12/01/2015

📅

Apply

ANTHRAX

Immunizations: 1

Date	Summary	Reaction	Series	Repeat Contraindicated	Facility
01/09/2014	Anthrax		2	No	DOD

Close

3. To change the date range of the displayed vaccine, select either one of the **preset date ranges**, **OR**.
4. Use the **calendar dropdowns** and select **Apply** to choose a custom date range.
5. Select the **Next** and **Previous** buttons to navigate between immunizations.
6. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Immunization Detail Dialog box, to return to the applet.

Medications Review

The Medications Review applet provides a single view of the medication list with additional graphical data to support analysis of the existing patient medical data in order to facilitate better clinical decisions.

Legend:

- ' = Minutes
- h = Hours
- d = Days
- m = Months
- y = Years
- N/A = Data is Not Available

Medications Review: Expanded View

The expanded view of the Medications Review applet (Figure 6-34) displays the selected patient's medication history in both list and graph formats grouped by Inpatient, Clinic Order, and Outpatient medications. When a category is selected, the group is expanded to display the list and graph of items.

Note: If a patient has an inpatient status, then the Inpatient group displays and the Outpatient group is collapsed. If a patient has an outpatient status, then the Outpatient group displays and the Inpatient group is collapsed.

Figure 6-33 Medications Review Applet Expanded View

MEDICATION REVIEW			
▼ INPATIENT MEDS			
Name	Sig	Status/Next	
Metformin	500MG PO Q12H	Active	Graphs for Inpatient Medications are not available at this time
	500MG PO Q12H	Expired 10m	Graphs for Inpatient Medications are not available at this time
Metoprolol	50MG PO BID	Active	Graphs for Inpatient Medications are not available at this time
	50MG PO BID	Expired 10m	Graphs for Inpatient Medications are not available at this time
Simvastatin	40MG PO QPM	Active	Graphs for Inpatient Medications are not available at this time
	40MG PO QPM	Expired 10m	Graphs for Inpatient Medications are not available at this time
Warfarin	5MG PO QDAY-WARF	Active	Graphs for Inpatient Medications are not available at this time
	5MG PO QDAY-WARF	Expired 10m	Graphs for Inpatient Medications are not available at this time
► CLINIC ORDER MEDS			
► OUTPATIENT MEDS			

Outpatient Medications

The outpatient medication types include Outpatient, supply, and non-VA (external). The data listed for each outpatient medication order consists of the following, as displayed in Figure 6-34:

- Medication name: represents the name of the medication ingredient
- Non-local indicator: an indication that the medication was not ordered at your home facility
- Sig column
 - First preference: Dose + Units + Route + Frequency
 - If any of the above are missing, then concatenation: Strength (i.e. size of tablet) + Product form + Sig
- Status/Fillable column
 - Active order – displays 'Fillable' and the amount of time the medication can be filled
 - Zero refills – displays '0 Refills'
 - Pending order – displays 'Pending'
 - Expired order – displays the 'Expired' and the amount of time from the expiration date
 - Discontinued order – displays 'Discontinued' and the amount of time from the discontinued date
 - Time Annotations
 - >60 days to 24 months = 'm' (ex. 14m)
 - >47 hours to 60 days = 'd' (ex. 24d)
 - >24 months = 'y' (ex. 5y)
 - 60 minutes to 47 hours = 'h' (ex. 37h)
 - 0-59 minutes = The actual minutes expressed with quotations (ex. "48")


Figure 6-34 Medications Review Outpatient Meds



OUTPATIENT MEDICATIONS REVIEW INDICATORS

The Outpatient Medications graph displays a medication's start date, stop date, and the dispensing dates. The graph indicators for outpatient medications are displayed in Figure 6-35.

Figure 6-35 Outpatient Medications Graph Indicators

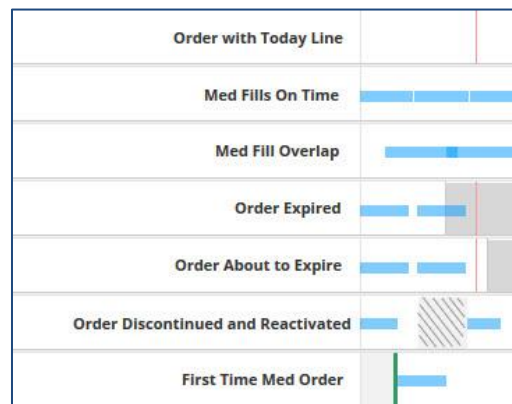

Some tips for the Outpatient Medications graph:

- A white background depicts an active order
- A gray background depicts when a medication has expired
- A hashed gray background depicts when a medication has been discontinued
- If two orders overlap, an active order overrides a discontinued or expired order
- Active/Hold and Active/Suspend orders count as being active
- A green bar depicts when a medication has been filled for the first time

Inpatient Medications

The data listed for each inpatient medication order (Figure 6-36) consists of the following:


- Name Column: represents the name of the medication ingredient
- Non-local indicator : an indication that the medication was not ordered at your home facility
- Sig Column
 - First preference Dose + Units + Route + Frequency
 - If any of the above are missing, then concatenation: Strength (i.e. size of tablet) + Product form + Sig
- Status/Next Column: the medication status or the time of the next administration
 - Pending order – displays 'pending'
 - Expired order – displays the 'Expired' and the amount of time from the expiration date
 - Discontinued order – displays 'Discontinued' and the amount of time from the discontinued date
 - Time Annotations
 - >60 days to 24 months = 'm' (ex. 14m)
 - >47 hours to 60 days = 'd' (ex. 24d)
 - >24 months = 'y' (ex. 5y)
 - 60 minutes to 47 hours = 'h' (ex. 37h)
 - 0-59 minutes = The actual minutes expressed with quotations (ex. "48")

Figure 6-36 Medications Review Inpatient Meds

MEDICATION REVIEW		
▼ INPATIENT MEDS		
Name	Sig	Status/Next
Bupropion	100MG PO TID	Active
Docusate	100MG PO QDAY	Active
	100MG PO QDAY	Active
	100MG PO QDAY	Expired 3m
Gabapentin	300MG PO Q6H	Active
	1200MG PO Q6H	Active
	300MG PO Q6H	Active
	1200MG PO Q6H	Active
	1200MG PO Q6H	Active
	300MG PO Q6H	Discontinued 23h
	1200MG PO Q6H	Discontinued 6m
	1200MG PO Q6H	Discontinued 6m

Clinic Order Meds

The data listed for each clinic order medication (Figure 6-37) consists of the following:

- Name Column: represents the name of the medication ingredient
- Non-local indicator: an indication that the medication was not ordered at your home facility
- Sig Column
 - First preference Dose + Units + Route + Frequency
 - If any of the above are missing, then concatenation: Strength (i.e. size of tablet) + Product form + Sig
- Status/Next Column: the medication status or the time of the next administration
 - Active order – displays 'Fillable' and the amount of time the medication can be filled
 - Zero refills – displays '0 Refills'
 - Pending order – displays 'pending'
 - Expired order – displays the 'Expired' and the amount of time from the expiration date
 - Discontinued order – displays 'Discontinued' and the amount of time from the discontinued date
 - Time Annotations
 - >60 days to 24 months = 'm' (ex. 14m)
 - >47 hours to 60 days = 'd' (ex. 24d)
 - >24 months = 'y' (ex. 5y)
 - 60 minutes to 47 hours = 'h' (ex. 37h)
 - 0-59 minutes = The actual minutes expressed with quotations (ex. "48")

Figure 6-37 Medications Review Clinic Order Meds

MEDICATION REVIEW		
▶ INPATIENT MEDS		
▼ CLINIC ORDER MEDS		
Name	Sig	Status/Fillable
Gabapentin	800MG PO Q6H	Discontinued 6y
	800MG PO Q6H	Discontinued 6y
	800MG PO Q6H	Discontinued 6y
▶ OUTPATIENT MEDS		

Medications Review: Detail View

What displays for the detail view of a medication depends on the status of the patient and if a VA or a Non-VA provider prescribed the medication.

To display a detailed view of a medication:


1. Select the **Inpatient Meds**, **Clinic Order Meds**, or the **Outpatient Meds** category. The category expands to display a list of medications.
2. Select a **medication**. The Infobutton and Details form button display.
3. Select the **Details form** button  and the detail view for the selected medication displays (Figure 6-38).

Figure 6-38 Medications Review Applet Detail View

MEDICATION REVIEW																				
<div> <div>Details form</div> <div> </div> </div>																				
	Sig	Status/Next																		
Metformin	500MG PO Q12H	Expired 7m																		
	500MG PO Q12H	Active																		
<div> <div> Order Hx <div> 01/29/2015 - 02/06/2015 02/01/2010 - 03/31/2010 02/01/2010 - 03/31/2010 </div> </div> <div> Links <div> Clinical Pharmacology MDConsult UpToDate VisualDx </div> </div> <div> Patient Education <div> Krames StayWell UpToDate VisualDx </div> </div> </div> <div> <div> Metformin Tab, Oral 500 MG <div>EXPIRED</div> </div> <div> Give: 500MG PO Q12H </div> <div> <table border="1"> <thead> <tr> <th>Prescription No.</th> <th>Supply</th> <th colspan="2">Dose/Schedule</th> </tr> </thead> <tbody> <tr> <td></td> <td></td> <td colspan="2">500 PO Q12H</td> </tr> <tr> <td>Provider</td> <td>Pharmacist</td> <td>Location</td> <td>Facility</td> </tr> <tr> <td>PROVIDER, TEN</td> <td>PROGRAMMER, FIVE</td> <td>7A GEN MED</td> <td>ABILENE (CAA)</td> </tr> </tbody> </table> </div> <div> Fill History </div> <div> No Fill History </div> </div>					Prescription No.	Supply	Dose/Schedule				500 PO Q12H		Provider	Pharmacist	Location	Facility	PROVIDER, TEN	PROGRAMMER, FIVE	7A GEN MED	ABILENE (CAA)
Prescription No.	Supply	Dose/Schedule																		
		500 PO Q12H																		
Provider	Pharmacist	Location	Facility																	
PROVIDER, TEN	PROGRAMMER, FIVE	7A GEN MED	ABILENE (CAA)																	

In addition to the detailed information for the selected medication, the *Links* and *Patient Education* sections provide access to external information and resources. This information is found on the bottom left-hand side of the Medications Detail Dialog box.

To use a resource link in the Links and Patient Education sections of the Medications Detail Dialog box:

1. Select one of the **resources** under the *Links* or *Patient Education* sections to launch the external resource. The external resource opens in a new browser tab.
2. Select the **X** on the external resource's browser tab to close it, or select the **VA eHMP** browser tab to return to the eHMP application without closing the external resource tab.

Military History

The Military History applet displays the following military history items: branches of service, years of service, areas served, and military occupational specialties. It allows you to edit the text field information that is written back to the selected patient's record and stored in eHMP.

Military History: Summary View

Figure 6-39 displays the summary view that lists the name and description for each section.

Figure 6-39 Military History Applet Summary View

MILITARY HISTORY	
Name	Description
Branch(s) of Service	Navy
Service Date(s)	Date Entered: 2008 NOV 27 - Date Separated: 2009 MAY 12 /// Date Entered: 2010 AUG 20 - Date Separated: 2011 FEB 09 ///
Areas of Service	AEROMED EVAC SQ (AMC)
Occupational Specialties	4A071, Health Service Mgt Craftsman

Military History: Expanded View

The expanded view of the Military History applet (Figure 6-40) displays the following information in sortable columns:

- Name
- Description
- Last Modified
- Location (VA site that modified military history information)
- Modified By

Figure 6-40 Military History Applet Expanded View

MILITARY HISTORY ↺ ? ✕				
Name	Description	Last Modified	Location	Modified By
Branch(s) of Service	Navy	12/09/2015	PANORAMA	VIHAAN KHAN
Service Date(s)	Date Entered: 2008 NOV 27 - Date Separated: 2009 MAY 12 /// Date Entered: 2010 AUG 20 - Date Separated: 2011 FEB 09 ///	12/02/2015	KODAK	PANORAMA USER
Areas of Service	AEROMED EVAC SQ (AMC)	12/02/2015	KODAK	PANORAMA USER
Occupational Specialties	4A071, Health Service Mgt Craftsman	12/02/2015	KODAK	PANORAMA USER

Military History: Detail View

To display a detailed view of an available military history item from the summary or expanded applet:


1. Select a **list item** and the applet toolbar menu displays.
2. Select the **Details form** button  to open the Military History Detail Dialog box (Figure 6-41).

Figure 6-41 Military History Detail Dialog Box



The dialog box has a title bar with a close button (X) in the top right corner. The main content area is titled "Occupational Specialties" and contains the text "30CX Commander, Mission Support /// 33SX Communications & Information ///". A "Close" button is located in the bottom right corner.

3. Select the **Close** button or the **X** in the upper right-hand corner to return to the applet.

Military History: Edit Form

To edit a military history item from the summary or expanded applet:


1. Select a **list item** and a set of buttons display.
2. Select the **Edit form** button  to open the Edit Military History box (Figure 6-42).

Figure 6-42 Edit Military History

3. Edit the **information** in the text box.
4. Select either **Save** to save the information or **Cancel** to return to the applet without saving changes.

Narrative Lab Results

The Narrative Lab Results applet lists the reports of the patient's recorded laboratory results.

Narrative Lab Results: Summary View

Figure 6-43 displays the summary view that lists the date and time taken, description, type, and facility.

Figure 6-43 Narrative Lab Results Applet Summary View

NARRATIVE LAB RESULTS				
Date	Description	Type	Facility	
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	TST1	
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	TST1	
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	TST2	
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	TST1	
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	TST1	
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	TST2	

Narrative Lab Results: Expanded View

The expanded view of the Lab Results applet (Figure 6-44) displays the following information in sortable columns:

- Date
- Description
- Type
- Author or Verifier
- Facility

Figure 6-44 Narrative Lab Results Applet Expanded View

NARRATIVE LAB RESULTS					?	Y	⚙	↗
Date	Description	Type	Author or Verifier	Facility				
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST2				
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST2				
03/15/2010 - 09:40	LR SURGICAL PATHOLOGY REPORT	Laboratory Report	None	TST1				
03/15/2010 - 09:40	LR SURGICAL PATHOLOGY REPORT	Laboratory Report	None	TST1				
03/15/2010 - 09:40	LR SURGICAL PATHOLOGY REPORT	Laboratory Report	None	TST2				
03/04/2010 - 12:34	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				
03/04/2010 - 12:34	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1				

Narrative Lab Results: Detail View

To display a detailed view of a lab result from the summary or expanded applet:


1. Select a **list item** to display the applet toolbar menu.
2. Select the **Details form** button  to open the Narrative Lab Results Detail Dialog box (Figure 6-45).

Figure 6-45 Narrative Lab Results Detail Dialog Box

Ir microbiology report Details Previous Next ×

Facility CAMP MASTER
 Author None
 Status Completed
 Date/Time 02/03/2015 - 12:24

Accession [UID]: MI 15 4 [1315000004] Received: Feb 03, 2015@12:24
 Collection sample: UNKNOWN Collection date: Feb 03, 2015 12:24
 Provider: PROVIDER, TWENTY

Test(s) ordered: CULTURE & SUSCEPTIBILITY..... completed: Feb 03, 2015

* BACTERIOLOGY FINAL REPORT => Feb 03, 2015 14:20 TECH CODE: 119
 GRAM STAIN:

CULTURE RESULTS: ACINETOBACTER ANITRATUS - Quantity: >25,000 - <50,000 CFU/ML

ANTIBIOTIC SUSCEPTIBILITY TEST RESULTS:

	ACINETOBACTER ANITRATUS	
AMIKACAN.....	R	
CEFAZOLIN.....	S	RESTRICTED
CEFOXITIN.....	S	
CEFOTAXIME.....	S	
GENTAMICIN.....	R	
TOBRAMICIN.....	R	
TRMSULF.....	R	
VANCOMYCIN.....	R	

Close






3. Select the **Next** and **Previous** buttons to navigate between lab results.
4. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Narrative Lab Results Detail Dialog box, to return to the applet.

Numeric Lab Results

The Numeric Lab Results applet lists a patient's recorded laboratory results.

Numeric Lab Results: Trend View

Figure 6-46 displays the trend view that lists the following:

- Lab Test
- Results (numeric data value)
- When the lab result was last documented
- Data Range Graph that includes the following rules:
 - Blue diamond  - current normal
 - Yellow diamond  - current abnormal low/high
 - Red star  - current critical low/high
 - Black dot with line attached to current flag  - previous value and indicates low or high value in comparison to the current value
 - Black dot inside current flag  - no change between current and previous values

- White background depicts the reference range
- Gray background depicts outside of the reference range

Figure 6-46 Numeric Lab Results Applet Trend View



Select the **data range graph** to display a list of the five previous results of a given lab, including the test value, reference range, observed date and time, and facility (Figure 6-47).

Figure 6-47 Numeric Lab Results Quick View

NUMERIC LAB RESULTS				
Lab Test	Result	Last		
GLUCOSE	221 mg/dL	16m		
UREA NITROGEN	11 mg/dL	16m		
CREATININE				
SODIUM				
TROPONIN				

	Value	Ref. Range	Observed	Facility
CREATININE	11 mg/dL	5-24 mg/dL	01/29/2015 15:1700	TST1
SODIUM	11 mg/dL	5-24 mg/dL	01/29/2015 15:1700	TST2
TROPONIN	11 mg/dL	5-24 mg/dL	01/29/2015 15:17	TST1

Numeric Lab Results: Summary View

Figure 6-48 displays the summary view that lists the date and time taken, lab test (includes a Panel button that opens a test drawer where applicable), flag (yellow indicates abnormal results and red indicates critical), and result of the patient's lab results.

Figure 6-48 Numeric Lab Results Applet Summary View

NUMERIC LAB RESULTS			
Date	Lab Test	Flag	Result
01/29/2015 - 15:17	Panel CBC BLOOD SP LB #18415	H	
01/29/2015 - 15:17	Panel CBC BLOOD SP LB #18415	H	
01/29/2015 - 15:17	Panel CBC BLOOD SP LB #18415	H	
01/29/2015 - 15:17	TROPONIN - SERUM	H	1 ug/mL




Numeric Lab Results: Expanded View

The expanded view of the Numeric Lab Results applet (Figure 6-49) displays the following information in sortable columns:

- Date
- Lab Test
- Flag
- Result
- Unit
- Ref Range
- Facility


Note: The H, H+, L, and L+ icons indicate abnormal (yellow) and critical (red) highs and lows.

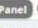Figure 6-49 Numeric Lab Results Applet Expanded View

NUMERIC LAB RESULTS							
Date	Lab Test	Flag	Result	Unit	Ref Range	Facility	
01/29/2015 - 15:17	 CBC BLOOD SP LB #18415	H				TST1	
01/29/2015 - 15:17	 CBC BLOOD SP LB #18415	H				TST2	
01/29/2015 - 15:17	 CBC BLOOD SP LB #18415	H				TST1	
01/29/2015 - 15:17	TROPONIN - SERUM	H	1	ug/mL	0-0.5	TST1	
01/29/2015 - 15:17	 CHEM 7 BLOOD SERUM SP LB #18415	H				TST1	
01/29/2015 - 15:17	TROPONIN - SERUM	H	1	ug/mL	0-0.5	TST2	
01/29/2015 - 15:17	 CHEM 7 BLOOD SERUM SP LB #18415	H				TST2	
01/29/2015 - 15:17	TROPONIN - SERUM	H	1	ug/mL	0-0.5	TST1	
01/29/2015 - 15:17	 CHEM 7 BLOOD SERUM SP LB #18415	H				TST1	
05/07/2013 - 10:43	Sodium, Blood Quantitative - PLASMA		139	mmol/L	134-146	DOD	
05/05/2013 - 14:10	Potassium, Serum or Plasma Quantitative - PLASMA	H	5.4	mmol/L	3.5-4.7	DOD	


Numeric Lab Results: Detail View

To display a detailed view of a numeric lab result from the trend, summary or expanded applet:

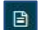1. Select a **list item** and the applet toolbar menu displays.
2. Select the **Details form** button  to open the Numeric Lab Results Detail Dialog box (Figure 6-50).

Note: If the Panel icon is displayed next to a list item, the Detail Dialog box does not open. Instead, select the **item** and the applet toolbar menu displays. Select the **Details form** button, and the specific tests run within the panel display, along with the results. Select a **lab test** under the panel and then the **Details form** button. The Lab Results Detail Dialog box displays a summary list and a historic graphical representation for the selected test that can be filtered by using the date filter.

Figure 6-50 Numeric Lab Results Detail Dialog Box

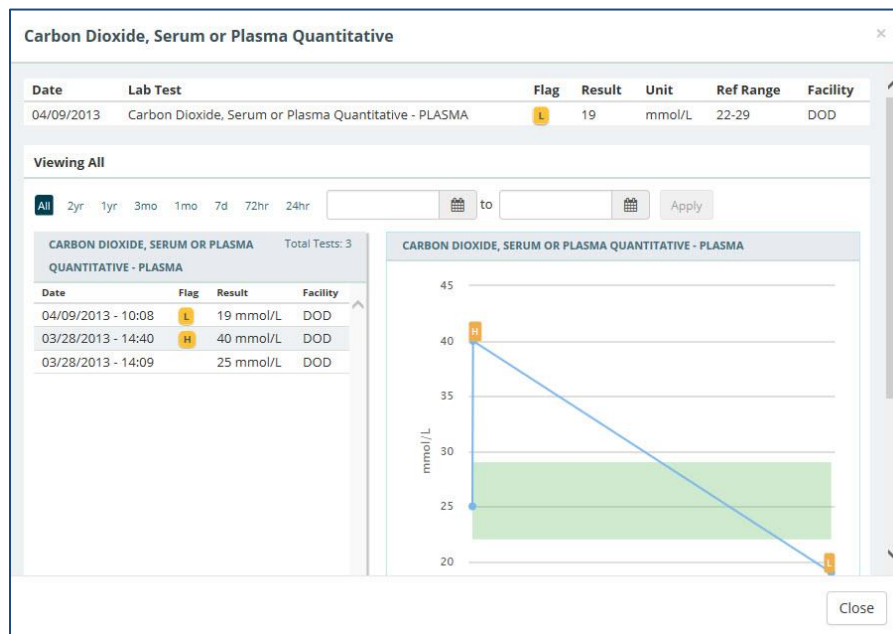

3. Use the **preset date ranges** or the **calendar dropdowns** and select **Apply** to change the date range for the displayed lab result.
4. Hover your cursor over **different areas of the graph** to display a tooltip with the Lab Result and Reference Range information.
5. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Numeric Lab Results Detail Dialog box, to return to the applet.

Orders

The Orders applet displays all orders for the selected patient submitted from all sites.

Orders: Summary View

Figure 6-51 displays the Orders applet summary view that lists the order date, status, order, and facility. Orders can be filtered by type using the *Order Type* dropdown box within the applet.

Figure 6-51 Orders Applet Summary View

ORDERS			
Order Type: All			
Order Date	Status	Order	Facility
03/26/2015	DISCONTINUED	ELECTROCARDIOGRAM CARDIOLOGY Proc Bedside <DISCONTINUED BY SERVICE>	BAY
08/14/2014	ACTIVE	REGULAR Diet	BAY
08/14/2014	ACTIVE	REGULAR Diet	BAY
08/14/2014	ACTIVE	REGULAR Diet	BAY
02/03/2015	COMPLETE	CULTURE & SUSCEPTIBILITY UNKNOWN WC LB #18424	BAY

Orders: Expanded View

Similar to the Orders summary view, orders can be filtered by type using the *Order Type* dropdown in expanded view. The expanded view of Orders (Figure 6-52) displays the following information in sortable columns:

- Order Date
- Status
- Order
- Type
- Provider Name
- Start Date
- Stop Date
- Facility

Figure 6-52 Orders Applet Expanded View

Order Date	Status	Order	Type	Provider Name	Start Date	Stop Date	Facility
03/26/2015	DISCONTINUED	ELECTROCARDIOGRAM CARDIOLOGY Proc Bedside <DISCONTINUED BY SERVICE>	Consult	Provider, Five	03/26/2015	03/26/2015	BAY
08/14/2014	ACTIVE	REGULAR Diet	Dietetics Order	Programmer, One	08/14/2014		BAY
08/14/2014	ACTIVE	REGULAR Diet	Dietetics Order	Programmer, One	08/14/2014		BAY
08/14/2014	ACTIVE	REGULAR Diet	Dietetics Order	Programmer, One	08/14/2014		BAY
02/03/2015	COMPLETE	CULTURE & SUSCEPTIBILITY UNKNOWN WC LB #18424	Laboratory	Provider, Twenty	02/03/2015	02/03/2015	BAY
02/03/2015	COMPLETE	CULTURE & SUSCEPTIBILITY UNKNOWN WC LB #18424	Laboratory	Provider, Twenty	02/03/2015	02/03/2015	BAY

Orders: Detail View

To display a detailed view of an order from the summary or expanded applet:

1. Select an order from the list within the applet. The Orders Detail Dialog box opens (Figure 6-53).

Figure 6-53 Orders Detail Dialog Box

2. Select the **Next** and **Previous** buttons to navigate between orders.
3. Select the **Close** button, the **X** in the upper right corner, or anywhere outside of the Orders Detail Dialog box, to return to the applet.

Problems

The Problems applet displays a list of problems, with standardized coding references, that has been compiled by the patient's providers.

Problems: Trend View

Figure 6-54 displays the trend view that lists the problem, acuity of problem, status of the problem, and the facility that documented the problem. If applicable, the problems are grouped by their standardization codes Systemized Nomenclature of Medicine Clinical Terms (SNOMED CT).

Figure 6-54 Problems Applet Trend View

PROBLEMS				
Problem	Acuity	Status	Facility	
Hand Joint Pain (Finding)	Chronic	Active	DOD	
Shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD	
Bone Pain (Finding)	Chronic	Active	DOD	
Swelling Of Limb (Finding)	Chronic	Active	DOD	

To display more information for a problem in the trend view:

1. Select the **left side of the problems tile** to display the applet toolbar menu (Figure 6-55).

Figure 6-55 Problems Applet Trend View Toolbar Menu



PROBLEMS			
	Acuity	Status	Facility
Hand Joint Pain (Finding)	Chronic	Active	DOD
Shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD
Bone Pain (Finding)	Chronic	Active	DOD
Swelling Of Limb (Finding)	Chronic	Active	DOD





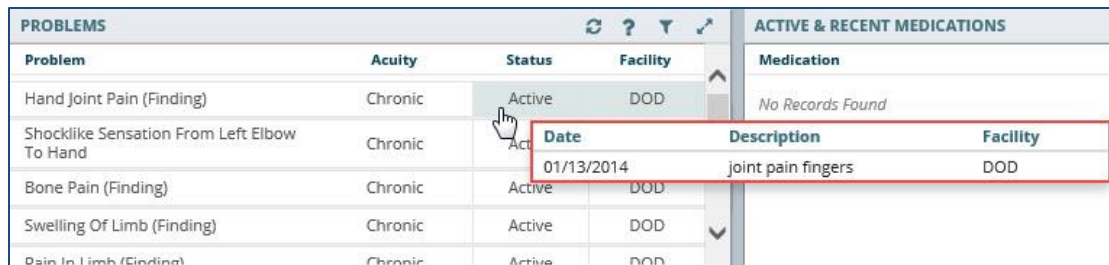
2. Select the **Info** button  to open a new browser window that displays information for the selected problem.
3. Select the **Details form** button  to open the Problems Detail Dialog box.
4. Select the **Quicklook** button  to display expanded information for the selected problem.
5. Select the **Submenu** button  (if applicable) to open a new workspace that is associated with the selected problem.
6. Select the **right side of the problem tile** to display a quick view containing up to the last five instances for the selected problem (Figure 6-56).

Figure 6-56 Problems Applet Quick View



PROBLEMS				ACTIVE & RECENT MEDICATIONS
Problem	Acuity	Status	Facility	Medication
Hand Joint Pain (Finding)	Chronic	Active	DOD	No Records Found
Shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD	
Bone Pain (Finding)	Chronic	Active	DOD	
Swelling Of Limb (Finding)	Chronic	Active	DOD	
Pain In Limb (Finding)	Chronic	Active	DOD	

Date	Description	Facility
01/13/2014	joint pain fingers	DOD

Problems: Summary View

Figure 6-57 displays the summary view that lists the description, the acuity of the patient's problem, and the status of the problem.

Figure 6-57 Problems Applet Summary View

PROBLEMS			↺ ? ↻ ⚙ ↗
Description	Acuity	Status	
Diabetes Mellitus Type II or unspecified	Chronic	Active	⬆
Chronic Systolic Heart failure	Chronic	Active	
Acute myocardial infarction, unspecified site, episode of care unspecified	Unknown	Active	
Hypertension	Chronic	Active	
Hyperlipidemia	Chronic	Active	⬇

To display more information for a problem in the summary view:

1. Select a **problem** to display the applet toolbar menu (Figure 6-58).

Figure 6-58 Problems Applet Summary View Toolbar Menu

PROBLEMS			↺ ? ↻ ⚙ ↗
Description	Acuity	Status	
i D S itus Type II or unspecified	Chronic	Active	⬆
Chronic Systolic Heart failure	Chronic	Active	
Acute myocardial infarction, unspecified site, episode of care unspecified	Unknown	Active	
Hypertension	Chronic	Active	
Hyperlipidemia	Chronic	Active	
Diabetes Mellitus Type II or unspecified	Chronic	Active	⬇

2. Select the **Info** button **i** to open a new browser window that displays information for the selected problem.
3. Select the **Details form** button **D** to open the Problems Detail Dialog box.
4. Select the **Submenu** button **S** (if applicable) to open a new workspace that is associated with the selected problem.

Problems: Expanded View

The expanded view of the Problems applet (Figure 6-59) displays the following information in sortable columns:

- Description
- Standardized Description
- Acuity
- Status
- Onset Date
- Last Updated
- Provider
- Facility

Figure 6-59 Problems Applet Expanded View

PROBLEMS							
Description	Standardized Description	Acuity	Status	Onset Date	Last Updated	Provider	Facility
Diabetes Mellitus Type II or unspecified		Chronic	Active	05/02/1998	03/30/2004	Vehu, Eight	TST1
Chronic Systolic Heart failure	Chronic systolic heart failure (disorder)	Chronic	Active	03/09/2004	03/09/2004	Labtech, Special	TST1
Acute myocardial infarction, unspecified site, episode of care unspecified		Unknown	Active	03/17/2005	03/17/2005	Vehu, Eight	TST1
Hypertension	Essential hypertension (disorder)	Chronic	Active	04/07/2005	04/10/2007	Vehu, Onehundred	TST1
Hyperlipidemia		Chronic	Active	04/07/2005	04/10/2007	Vehu, Onehundred	TST1
Occasional, uncontrolled chest pain	Impending infarction (disorder)	Acute	Active	03/15/1996	05/14/1996	Programmer, Twenty	NJS
Diabetes Mellitus Type II or unspecified		Chronic	Active	05/02/1998	03/30/2004	Vehu, Eight	TST2
Chronic Systolic Heart failure	Chronic systolic heart failure (disorder)	Chronic	Active	03/09/2004	03/09/2004	Labtech, Special	TST2
Acute myocardial infarction, unspecified site, episode of care unspecified		Unknown	Active	03/17/2005	03/17/2005	Vehu, Eight	TST2

Problems: Detail View

To display the details of a patient's problem from a trend, summary, or expanded applet:


1. Select a **problem** from the list and applet toolbar menu displays.
2. Select the **Details form** button  to display the Problems Detail Dialog box (Figure 6-60).

Figure 6-60 Problems Detail Dialog Box

Acute myocardial infarction, unspecified site, episode of care unspecified (ICD-9-CM 410.90)

↑ Previous

↓ Next

×

Primary ICD-9-CM: 410.90

SNOMED CT:

Onset: 03/17/2005

Acuity: Unknown

Provider: Vehu, Eight

Facility: CAMP MASTER

Location: General Medicine

Status: Active

Entered: 03/17/2005

Updated: 03/17/2005

Comments:

Close


3. Select the **Next** and **Previous** buttons to navigate between problems.
4. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Problems Detail Dialog box, to return to the applet.

Reports

Reports: Summary View

Figure 6-61 displays the summary view of the Reports applet that lists available reports by date, type, and entered by. Reports are listed in reverse chronological order so that providers are able to find the most recent report.

Figure 6-61 Reports Applet Summary View

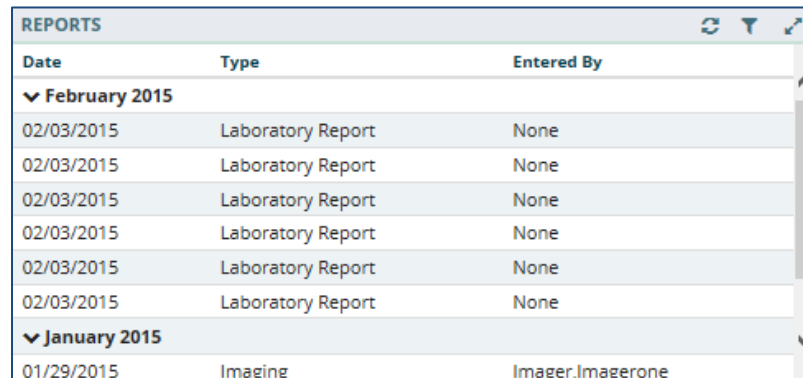

REPORTS		
Date	Type	Entered By
▼ February 2015		
02/03/2015	Laboratory Report	None
02/03/2015	Laboratory Report	None
02/03/2015	Laboratory Report	None
02/03/2015	Laboratory Report	None
02/03/2015	Laboratory Report	None
02/03/2015	Laboratory Report	None
▼ January 2015		
01/29/2015	Imaging	Imager.Imagerone

Reports: Detail View

To display a detailed view of a report from the summary applet:

1. Select a **report** from the list within the applet. The Reports Detail Dialog box opens (Figure 6-62).

Figure 6-62 Reports Detail Dialog Box

computed tomography, head or brain; without contrast material, followed by contrast material(s) and further sections Details ✕


Facility CAMP MASTER
Type Imaging
Status COMPLETE
Date/Time 01/29/2015 - 15:37
Providers Imager,Imagerone

Order Information

Requesting Provider IMAGER,IMAGERONE
 Orderable Item CT HEAD W&WO CONT
 Reason TBI

Results

[CT HEAD W&WO CONT](#)

Local Title: CT HEAD W&WO CONT
 Date Of Note: 01/29/2015 - 15:37
 Author: Imager,Imagerone 
 Status: Verified

EIGHT,PATIENT 666-00-0008 DOB-APR 07, 1935 M
 Exm Date: JAN 29, 2015@15:37
 Req Phys: VEHU,TEN

Pat Loc: 7A GEN MED/06-03-2015@16:30
 Img Loc: CT SCAN
 Service: MEDICINE

(Case 5 COMPLETE) CT HEAD W&WO CONT (CT Detailed) CPT:70470
 Reason for Study: TBI

Clinical History:
 CAR ACCIDENT

Close


2. Select the **hyperlink** (if available) or scroll down for more information.
3. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Reports Detail Dialog box, to return to the applet.

Stacked Graphs

Stacked Graphs: Expanded View

The Stacked Graphs applet (Figure 6-63) is user defined, and provides graphing functionality. The applet enables you to graph different types of data to a standardized, x-axis timeline, similar to functionality currently available in CPRS.

Figure 6-63 Stacked Graphs Applet Expanded View

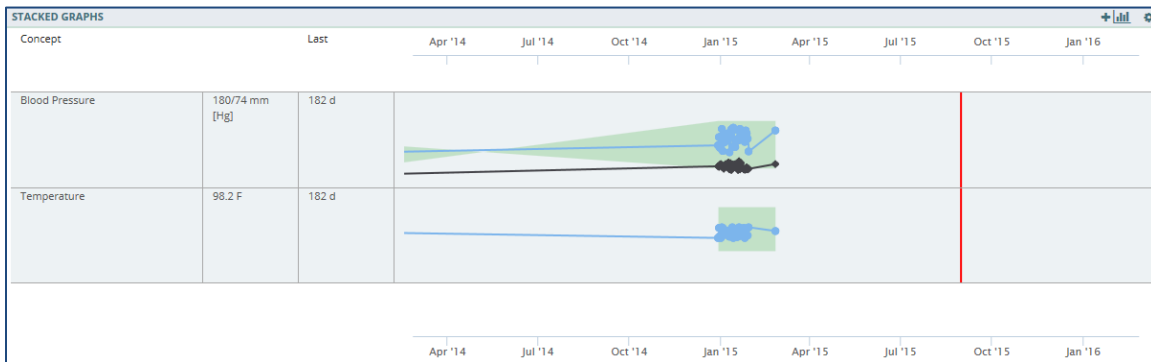

Add a Graph to Stacked Graphs Applet

To add graphs to the stacked graphs applet:


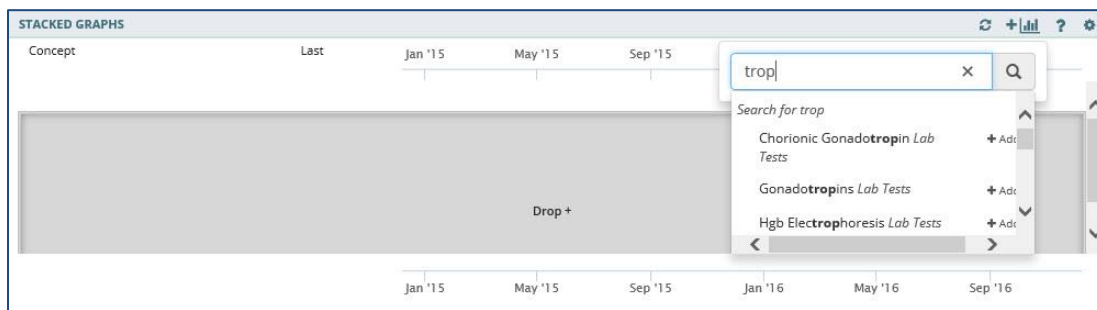
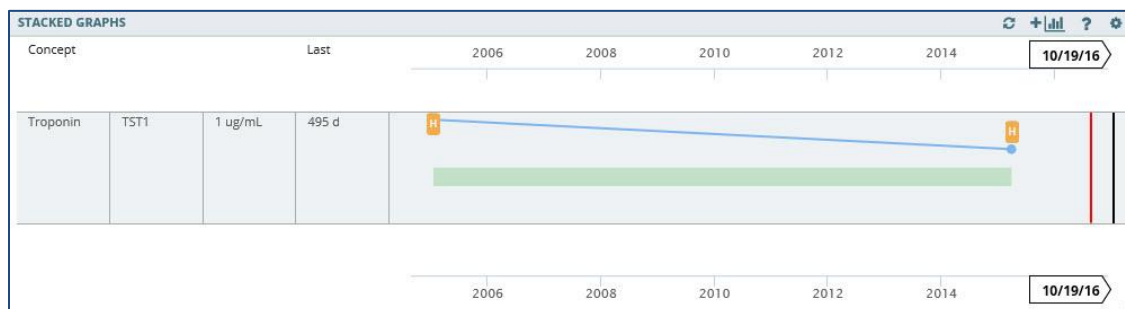
1. Create a **new workspace**. Refer to *Workspaces*.
2. Add the **stacked graphs** applet to the workspace.
3. Open the **workspace** with the stacked graph applet.
4. Select the **Add a graph** button . A Search field displays.
5. Enter the **name** of the desired graph type (i.e., temperature, blood pressure, etc.). The results populate as you enter text (Figure 6-64).

Figure 6-64 Add Graph to Stacked Graphs Applet



6. Select the desired **graph**. The graph displays in the Stacked Graphs applet (Figure 6-65).

Figure 6-65 Graph Added to Stacked Graphs Applet



7. Repeat **steps 4-6** to continue adding graphs to the applet.


Edit the Stacked Graphs Applet

Once a graph is added to the Stacked Graphs applet, you have the option to rearrange the graphs in the applet, review additional information for the graph, or delete the graph from the applet.

To edit the graphs:

1. Select the **graph** you want to edit and the applet toolbar menu displays (Figure 6-66).

Figure 6-66 Stacked Graphs Applet Toolbar Menu

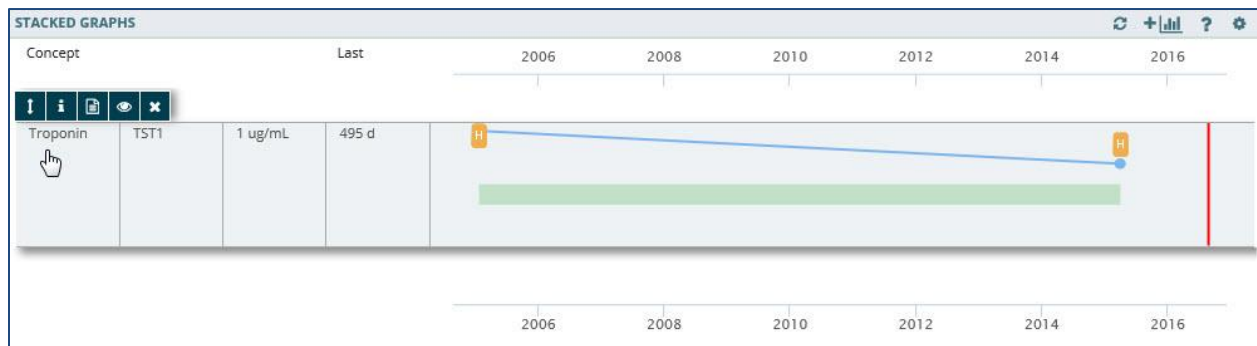






2. Select the **Tile sort** button  to sort the order of the graphs in the applet.
3. Select the **Infobutton**  to open a new browser window that displays information for the selected graph.
4. Select the **Details form** button  to open a Detail Dialog box.
5. Select the **Quicklook** button  to display expanded information for the selected graph.
6. Select the **Delete stacked graph** button  to remove the graph from the Stacked Graphs applet. A warning displays asking if you are sure you want to remove the graph (Figure 6-67).

Figure 6-67 Delete Stacked Graphs Warning



7. Select **Delete** to remove the graph or **Cancel** to keep the graph in the Stacked Graphs applet.

Timeline

Timeline view provides detailed information of a patient's visit history.

Timeline: Summary View

Figure 6-68 displays the Timeline applet in summary view. It lists the date and time, activity name, and type for each event.

Figure 6-68 Timeline Summary View

TIMELINE			↺ ? ⌵ ✕
Date & Time	Activity	Type	
▼ March 2015			
03/26/2015 - 16:58	ELECTROCARDIOGRAM ZZELECTROCARDIOGRAM CARDIOLOGY Proc DISCONTINUED	Procedure	
▼ February 2015			
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology	
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology	
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology	

Timeline: Expanded View

Figure 6-69 displays the expanded view of Timeline that lists the date and time, activity, type, entered by, and facility for the event.

To display the expanded view for the Timeline applet:

1. Select **Timeline** from the Workspace dropdown menu. The patient's historical visit data displays, grouped by month. The default view displays the date in reverse chronological order.

Figure 6-69 Timeline Expanded View

TIMELINE					↺ ? ⌵ ✕
Date & Time	Activity	Type	Entered By	Facility	
▼ March 2015					
03/26/2015 - 16:58	ELECTROCARDIOGRAM ZZELECTROCARDIOGRAM CARDIOLOGY Proc DISCONTINUED	Procedure		CAMP MASTER	
▼ February 2015					
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology		CAMP MASTER	
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology		CAMP BEE	
02/03/2015 - 12:24	CULTURE & SUSCEPTIBILITY - UNKNOWN	Microbiology		CAMP MASTER	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology		CAMP MASTER	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology		CAMP BEE	
02/03/2015 - 12:15	CULTURE & SUSCEPTIBILITY - CALCANEUS	Microbiology		CAMP MASTER	

Timeline: Detail View

To display a detailed view of a specific activity:

1. Select an **item** from the list, and the Timeline Detail Dialog box displays (Figure 6-70).

Figure 6-70 Timeline Detail Dialog Box

The screenshot shows a dialog box titled "electrocardiogram zzelectrocardiogram cardiology proc Details". The dialog contains the following information:

- Facility:** CAMP MASTER
- Type:** Procedure
- Status:** DISCONTINUED
- Date/Time:** 03/26/2015 - 16:58
- Order Information:**
 - To Service: CARDIOLOGY
 - From Service: 7A GEN MED
 - Requesting Provider: PROVIDER,FIVE
 - Place: Bedside
 - Urgency: Routine
 - Orderable Item: ELECTROCARDIOGRAM
 - Procedure: ZZELECTROCARDIOGRAM
 - Reason: !!! TEST
- Activity:**
 - CPRS RELEASED ORDER
 - Date/Time: 03/26/2015 - 16:58
 - Responsible: PROVIDER,FIVE
 - DISCONTINUED
 - Date/Time: 03/26/2015 - 16:59
 - Responsible: PROVIDER,FIVE

A "Close" button is located at the bottom right of the dialog box.

2. Select the **Close** button, the **X** located in the upper right-hand corner, or anywhere outside of the Timeline Detail Dialog box, to close the detailed visit information and return to the default Timeline view.

VistA Health Summaries

The VistA Health Summaries applet provides Health Summary Report functionality in eHMP. This is a re-creation of the reports available in the CPRS Reports tab under 'Health Summaries', a feature that is also provided in VistA Web, listed as 'Health Summaries', and representing each facility where the patient has a record.

A Health Summary is a clinically oriented, structured report that extracts many kinds of data from VistA web and displays it in a standard format. The individual patient is the focus of health summaries. The data displayed covers a wide range of health-related information such as demographic data, allergies, current active medical problems, and laboratory results.

VistA Health Summaries: Summary View

Figure 6-71 displays the summary view of the VistA Health Summaries applet that groups each facility where the patient has a record. The reports are listed under the collapsible facility name, with the number of reports available indicated.

Figure 6-71 VistA Health Summaries Summary View

VISTA HEALTH SUMMARIES	
Facility	Report
^ TST1 28	
^ TST2 28	
^ 500 28	

VistA Health Summaries: Detail View

To display a detailed view of a VistA Health Summary:

1. Select the **dropdown arrow** next to the facility name. A list of reports for that facility display.
2. Select a **report name** and the Report Detail Dialog box opens (Figure 6-72).

Figure 6-72 VistA Health Summaries Detail Dialog Box

TST2 - CARDIOLOGY REPORTS

↑ Previous

↓ Next

×

EIGHT,PATIENT, 04/07/1935, 80y, 666-00-0008

07/10/2015 11:06

***** CONFIDENTIAL CARDIOLOGY REPORTS SUMMARY pg. 1 *****

EIGHT,PATIENT

666-00-0008

7A GEN MED 722-B

D08: 04/07/1935

----- SPN - Selected Prog Notes -----

No data available for CARD HEART ABN CNS

----- MEDB - Med Brief Report -----

CONSULT	COMPLETED PROCEDURES	DATE/TIME	PROCEDURE CODE
NUMBER		PERFORMED	
629	EEG	MAR 14, 2010@15:31	

*** END ***** CONFIDENTIAL CARDIOLOGY REPORTS SUMMARY pg. 1 *****

Close

3. Select the **Next** and **Previous** buttons to navigate between reports.
4. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Report Detail Dialog box, to return to the applet.

Vitals

The Vitals applet displays the patient's most recently recorded vitals. The information can be viewed both numerically and in graph form.

Vitals: Trend View

Figure 6-73 displays the trend view for Vitals with the following data:

- Type – type of vital collected, i.e. blood pressure, pain, weight
- Result – data captured in relation to the vital measured
- Last – timeframe vital was last collected
- Data Range Graph, which includes:
 - A diamond represents the last (current) value
 - Blue diamond indicates a current normal value
 - Orange diamond indicates a high/low value
 - A dot represents previously recorded value
 - A dot within a diamond indicates no change between previous entries. The current value is the same as previous.
 - If a diamond is to the left of the dot, then the value of the diamond is lower than the dot
 - If a diamond is to the right of the dot, then the value of the diamond is higher than the dot




Figure 6-73 Vitals Applet Trend View

VITALS			
Type	Result	Last	
BPS	180 mm[Hg]	3m	
BPD	74 mm[Hg]	3m	
Pulse	80 /min	3m	

To display more information for a vital in the trend view:

1. Select the **left side of the vitals tile** to display the applet toolbar menu (Figure 6-74).

Figure 6-74 Vitals Applet Trend View Toolbar Menu

VITALS				
				
	Result	Last		
BPS	180 mm[Hg]	15m		
BPD	74 mm[Hg]	15m		
Pulse	80 /min	15m		
RR	15 /min	15m		




2. Select the **Info** button  to open a new browser window that displays information for the selected vital.
3. Select the **Details form** button  to open the Vitals Detail Dialog box.
4. Select the **Quicklook** button  to display expanded information for the selected vital.
5. Select the **right side of the vitals tile** to display a quick view containing up to the last five instances for the selected vital (Figure 6-75).

Figure 6-75 Vitals Applet Quick View

VITALS				
Type	Result	Last		
BPS	180 mm[Hg]	15m		
BPD				
Pulse				
RR				
Temp				
	Date	Result	Ref. Range	Facility
	02/24/2015 - 12:17	180 mm[Hg]	100mm[Hg] - 210mm[Hg]	BAY
	02/24/2015 - 12:17	180 mm[Hg]	100mm[Hg] - 210mm[Hg]	BAY
	02/24/2015 - 12:17	180 mm[Hg]	100mm[Hg] - 210mm[Hg]	BAY
	01/29/2015 - 17:05	114 mm[Hg]	100mm[Hg] - 210mm[Hg]	TST2
	01/29/2015 - 17:05	114 mm[Hg]	100mm[Hg] - 210mm[Hg]	TST1

Vitals: Summary View

Figure 6-76 displays the summary view that lists the patient's blood pressure (BP), pulse (P), respiration (R), temperature (T), pulse oximetry (PO₂), pain (PN), weight (WT), and body mass index (BMI).


Figure 6-76 Vitals Applet Summary View

VITALS					
BP	180/74 mm[Hg]	02/24/2015	PO ₂	99 %	02/24/2015
P	80 /min	02/24/2015	PN	1	02/24/2015
R	15 /min	02/24/2015	WT	205 lb / 93.18 kg	02/24/2015
T	98.2 F / 36.8 C	02/24/2015	BMI	28.6	02/24/2015

To display more information for a vital in the summary view:

1. Select a **vital** from the applet to open the Vitals Detail Dialog box (Figure 6-77).

Figure 6-77 Vitals Detail Dialog Box

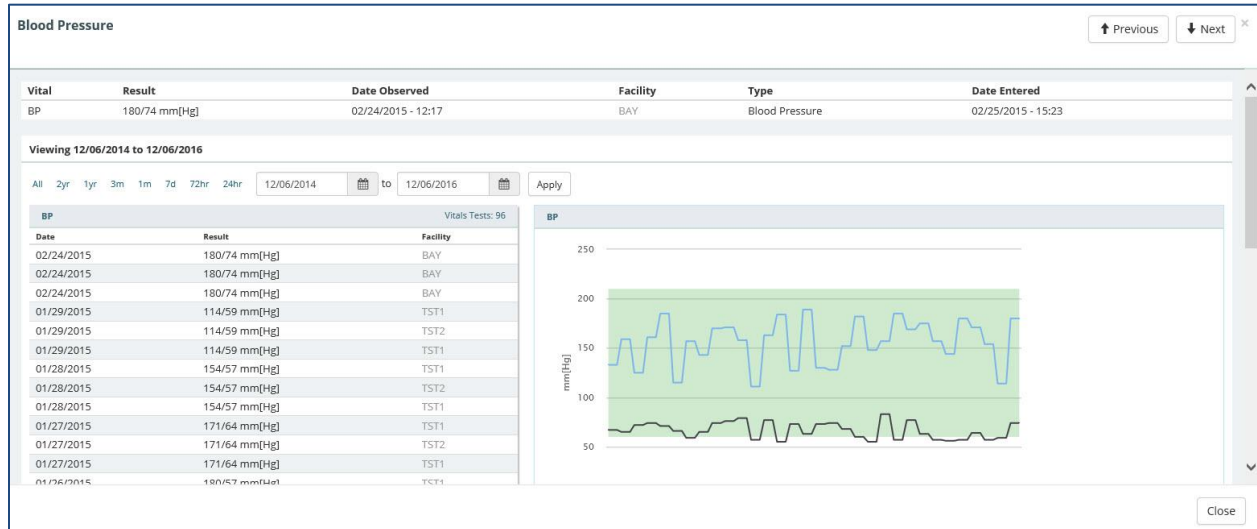

2. Select a **preset date** or use the **calendar dropdowns** and select **Apply** to filter the information.
3. Hover your cursor over **different areas of the graph** to display a tooltip with the Vital Result and Reference Range information.
4. Select the **Next** and **Previous** buttons to navigate between vitals.
5. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Vitals Detail Dialog box, to return to the applet.

Vitals: Expanded View

The expanded view of the Vitals applet (Figure 6-78) displays the following information in sortable columns:

- Date Observed
- Type
- Result
- Date Entered
- Qualifiers
- Facility

Figure 6-78 Vitals Applet Expanded View

VITALS

All2yr1yr3mo1mo7d72hr24hr

02/25/2014to02/25/2016

Apply

Q Enter your text filter

Add

Date Observed	Type	Result	Date Entered	Qualifiers	Facility
02/24/2015 - 12:17	BMI	28.6			BAY
02/24/2015 - 12:17	Pain	1	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Pulse Oximetry	99 %	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Weight	205 lb (93.18 kg)	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Height	71 in (180.34 cm)	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Pulse	80 /min	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Respiration	15 /min	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Temperature	98.2 F (36.8 C)	02/25/2015 - 15:23		BAY
02/24/2015 - 12:17	Blood Pressure	180/74 mm[Hg]	02/25/2015 - 15:23		BAY

Use the **preset dates** or the **calendar dropdowns** to filter the vitals by date. Enter **text filter** in the *Enter your text filter* field to filter the information using a keyword.

Vitals: Detail View

To display a detailed view of vitals from the trend, summary or expanded applet:


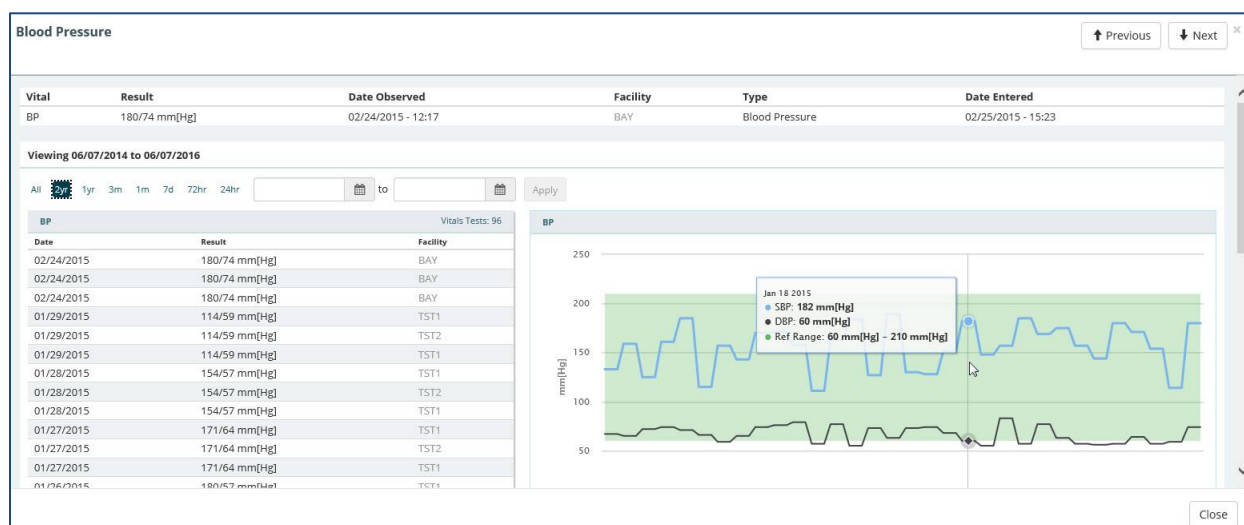
1. Select a **vital** (e.g., Blood Pressure) and the applet toolbar displays.
2. Select the **Details form** button  and the Vitals Detail Dialog box opens.
3. A historical, detailed view of that vital is presented in list view and graph view. The shaded area on the graph depicts the reference range, if available, for that vital.
4. Hovering over a **time on the graph** displays a tool tip with more information (Figure 6-79).

Figure 6-79 Vitals Detail Dialog Box



5. Select the **Close** button, the **X** in the upper right-hand corner, or anywhere outside of the Vitals Detail Dialog box, to return to the applet.

Chapter 7: Applet Features

About Applets

The applets in the eHMP application are widgets that sort patient data into segments throughout the patient record.

Most of the applets can be refreshed, filtered, maximized, and minimized (Figure 7-1) to display various levels of detail, provide different views of an applet's information, and open search capability.


Figure 7-1 Applet Feature Buttons



Refresh Button

The **Refresh** button found on each applet in eHMP updates patient data, and should be used after entering new patient information in CPRS.

To refresh patient data:

1. After entering any new patient information in CPRS, access eHMP again.
2. Select the **Refresh** button  for an applet. The new information displays.

Online Help Button

The **Online Help** button allows you to easily access context-driven, application-oriented help information.

To access Online Help:


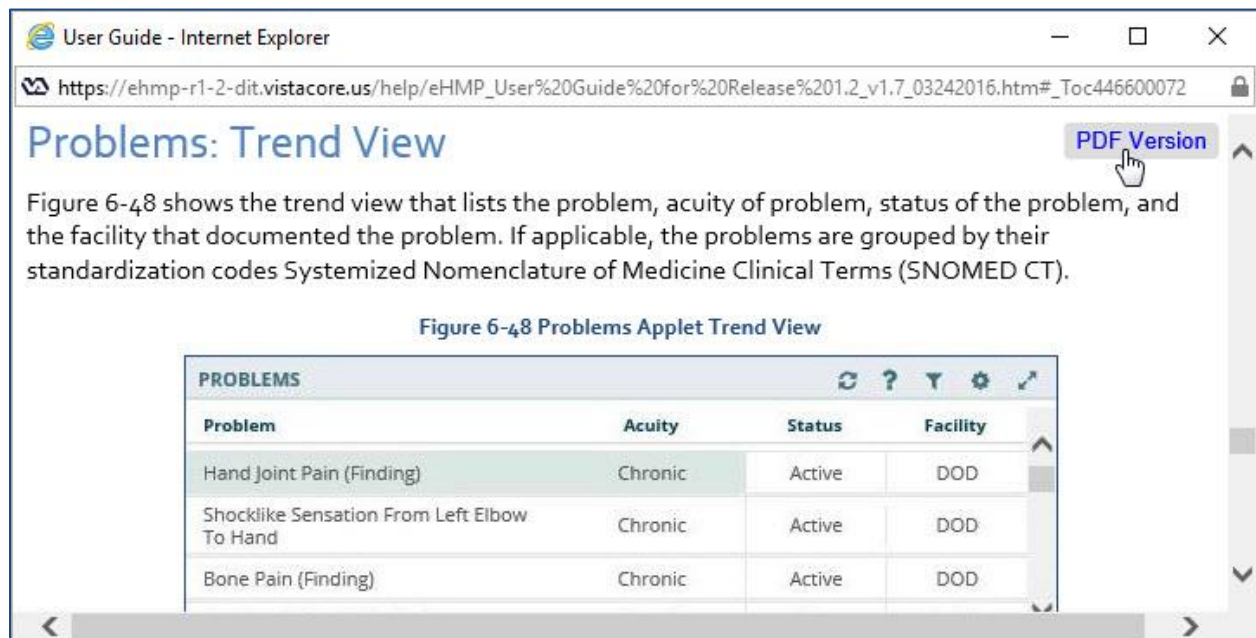
1. Select the **Online Help** button  for more information on a specific topic, and a new window displays with information on that topic.
2. Select the **PDF Version** link to download the eHMP User Guide (Figure 7-2).

Figure 7-2 PDF Version of the eHMP User Guide




3. Close out of the **Help window** to return to eHMP.

Filter Button

The **Filter** button allows you to filter data, limiting the items displayed within an applet.

To filter data for a particular applet:

1. Select the **Filter** button  to filter the data by entering key terms in the **Enter your text filter** field.
***Note:** If more than one word is used, the search results contain all of the filter words.*
2. Select the **X** next to the *Enter your text filter* field to return to the default list of items for that applet.
3. Select the **Filter** button again to remove the filter field.


User-Defined Applet Filters

You can add more than one filter text to user-defined workspaces that persist from patient-to-patient and session-to-session.

To add filters to a user-defined workspace:


1. Open the **user-defined workspace**.
2. Select the **Filter** button and the *Enter Filter* field displays.
3. Enter the **filter text**. The results populate as you enter text.
4. Select **Add** or press the **Enter** key on the keyboard to add the filter text. The filter text displays with an x next to the text and the applet header becomes darker with the word Filtered added (Figure 7-3).

Figure 7-3 Filtered Applet

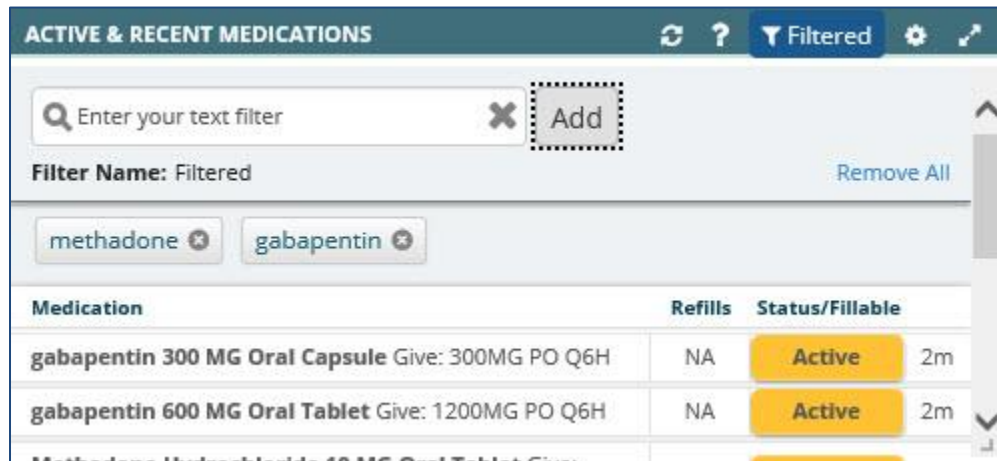


5. Repeat **steps 3 and 4** to continue adding filter text to the applet.
6. Select **Filtered** to the right of *Filter Name* to change the filter name.
7. Enter the **new filter name** and press **Enter** on the keyboard. The new filter name displays to the right of Filter Name and in the applet header (Figure 7-4).

Figure 7-4 Filter Name

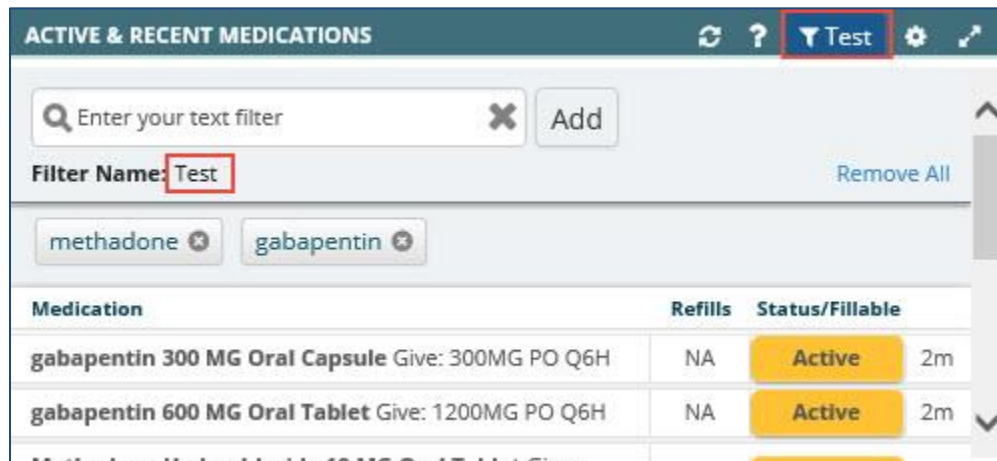


8. Select either the **x** to the right of the filtered text to remove it from the filtered data or **Remove All** to remove all of the filtered text.

Maximize Applet Button

The **Maximize Applet** button allows you to expand the applet. The maximized view provides more detailed information for the applet. Select the **Maximize Applet** button  located in the upper right corner of the applet to maximize it.

Minimize Applet Button

When viewing an applet in the maximized view, select the **Minimize** button  to return to the previous view.

Sortable Column Headers

Most applet column headers can be selected in the trend, summary, maximized, or expanded views to sort the information in ascending or descending order.

To sort a column in an applet:

1. Select the **column name** (e.g., Description or Facility). The list sorts in ascending order (Figure 7-5).

Figure 7-5 Sortable Column Header



Problem ▲	Acuity	Status	Facility
Acute Myocardial Infarction, Unspecified Site, Episode Of Care Unspecified	Unknown	Active	TST1
Aneurysm Of Gastroduodenal Artery (Disorder)	Chronic	Active	DOD
Ankle Pain (Finding)	Chronic	Active	DOD
Bone Pain (Finding)	Chronic	Active	DOD

2. Select the **column name** again and the list sorts in descending order.

Manual Tile Sort

Tiles are list items displayed when working in the trend view of an applet. Some applets provide the ability to manually sort tiles by dragging and dropping them to a new location in the list. Tiles can be sorted so that the order of display is more clinically relevant to the workspace. The sorted data within the workspace persists from session to session.

Note: The tile-sorting feature is available only for user-defined workspaces.

Tile sorting is available in the trend view for the following applets:

- Active & Recent Medications
- Numeric Lab Results
- Problems
- Vitals

To sort a tile in the Problems applet (example):

1. Select the **tile** to be moved (Figure 7-6). The selected tile is highlighted and the applet toolbar menu displays.


Figure 7-6 Tile Sorting

PROBLEMS			
Problem	Acuity	Status	Facility
Hand Joint Pain (Finding)	Chronic	Active	DOD
Shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD
Swelling Of Limb (Finding)	Chronic	Active	DOD

2. Drag and drop the **tile** to the desired location in the list. Once the tile has been moved, the word *Manual* appears in the header (Figure 7-7).

Figure 7-7 Manual Tile Sorting

PROBLEMS			
Problem	Acuity	Status	Facility
Manual	Chronic	Active	DOD
Swelling Of Limb (Finding)	Chronic	Active	DOD
Shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD
Bone Pain (Finding)	Chronic	Active	DOD

3. Repeat **steps 1 and 2** to reorder the remaining tiles in the applet.
4. To delete the manual sort and revert to the default view, select the **Clear your manual sort** button  next to *Manual*.

Infobutton

Some applets provide context-sensitive medical information specific to certain concepts, such as medications, immunizations, and labs. The **Infobutton** feature allows you to better assess and treat patients.

To access the Infobutton within the Problems applet (example):


1. Select **an item** from the Problems applet. The applet toolbar menu displays.
2. Select the **Infobutton**  (Figure 7-8).

Figure 7-8 Infobutton in the Problems Applet



Problem	Acuity	Status	Facility
shocklike Sensation From Left Elbow To Hand	Chronic	Active	DOD
Bone Pain (Finding)	Chronic	Active	DOD
Swelling Of Limb (Finding)	Chronic	Active	DOD

3. A new browser window opens that displays specific immunization information.
4. Close the **browser window** to return to the previous view.

Details Form Button

The **Details form** button is an alternate shortcut to a detail dialog box for a selected item.

To access the detailed dialog box using the Details form button:


1. Select the **row of the item** and the applet toolbar menu displays.
2. Select the **Details form** button . The detail dialog box for the selected item opens (Figure 7-9).

Figure 7-9 Detail Dialog Box for Selected Item



shocklike sensation from left elbow to hand
↑ Previous
↓ Next
×

Primary ICD-9-CM:
 SNOMED CT:

 Onset:
 Acuity: Chronic
 Provider: Bhie, Userone
 Facility: DOD
 Location: NH Great Lakes II
 Status: Active

 Entered: 01/10/2014
 Updated: 01/10/2014

 Comments:

Close

3. Select **Close** to close the dialog and return to the previous view.

Comment Indicators

A Comment Indicator displays when there is additional information about an item. (Figure 7-12).

Figure 7-12 Comment Indicator



To display a comment:


1. Select an **item with the Comment Indicator**.
2. Select the **Details form** button  and the detail dialog box displays. The comment displays at the bottom of the dialog (Figure 7-13).

Figure 7-13 Detailed Comment Information

Occasional, uncontrolled chest pain (ICD-9-CM 411.1)

↑ Previous ↓ Next ×

Primary ICD-9-CM:	411.1
SNOMED CT:	25106000
Onset:	03/15/1996
Acuity:	Acute
Provider:	Programmer, Twenty
Facility:	New Jersey HCS
Location:	
Status:	Active
Entered:	05/14/1996
Updated:	05/14/1996

Comments:

05/14/1996 PROGRAMMER, TWENTY
SHERIDAN PROBLEM

Close

Chapter 8: Workspaces

About Workspaces

Workspaces represent an area of activity with specific applets that support that activity. There are two types of workspaces for eHMP v1.2: predefined and user-defined workspaces.

Predefined Workspaces

Predefined workspaces have been created to provide you with different views of patient information. These workspaces cannot be edited or deleted. The Coversheet, Timeline, Overview, Meds Review, and Documents are all predefined workspaces.

The Depression, Diabetes Mellitus, Hypertension, and Pre-Procedure workspaces are also predefined and cannot be edited or deleted. These were developed for common conditions and expose a combination of concepts, graphing, and items to support clinical management of these conditions in an efficient and safe manner.

Workspace Manager

The Workspace Manager feature allows you to add, edit, arrange and delete user-defined workspaces within the patient record.

Select the **Workspace Manager** button  to open the Workspace Manager window.

The following are the elements of the Workspace Manager (Figure 8-1):

- Add New Workspace
- Workspace Filter
- Default View
- Title
- Associated Conditions
- Description
- Author
- Duplicate
- Rearrange
- Locked/Delete
- Preview
- Customize/Launch

Figure 8-1 Workspace Manager Window

Workspace Manager									
	Title	Assoc. Conditions	Description	Author					
☆	🔒 Coversheet	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Overview	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Timeline	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Documents	👁 0				I	🔒	Preview	Launch
☆	🔒 Depression CBW	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Diabetes Mellitus CBW	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Hypertension CBW	👁 0			📄	I	🔒	Preview	Launch
☆	🔒 Pre-Procedure CBW	👁 0			📄	I	🔒	Preview	Launch
☆	Timeline Copy	👁 0	Add description	PANORAMA USER	📄	I	🗑	Preview	Customize

Note: A shaded row with a lock icon beside the title indicates that workspace cannot be edited or deleted.

Add a New Workspace

The Workspace Manager allows you to create a new user-defined workspace that persists from patient to patient, and session to session.

To create a user-defined workspace:

1. Select the **Add New Workspace** button  to add a user-defined workspace to the workspace manager listing.

Note: Once a workspace has been added to the Workspace Manager, you can customize the new workspace.


Workspace Filter

The Workspace Filter allows you to filter the workspaces by title or description in the Workspace Manager.

To filter through the list of workspaces:


1. Select the **Filter** button . The *Filter by title or description* field displays (Figure 8-2).

Figure 8-2 Filter Workspaces

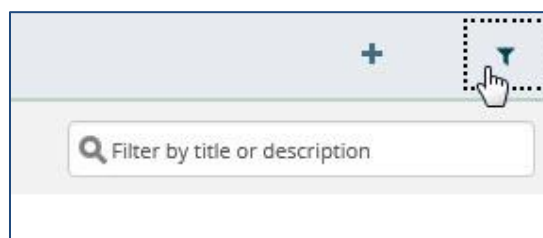

2. Enter the desired **workspace title or description**. The results populate as you enter text.
3. Select the **x** to remove the text and start a new filter.
4. Select the **filter** button again, to remove the *Filter by title or description* field.

Default View






The Default View button allows you to select the workspace you want to set as your default view when you open eHMP. Select the **Default View** button  to set the desired default view.

Note: *The Overview workspace is the default view until another workspace has been selected.*

Workspace Title

The Workspace Title field allows you to create or edit the title for a user-defined workspace. Enter or edit a **workspace title** in the *Title* field (Figure 8-3).

Figure 8-3 Workspace Title Field

Workspace Manager			
Title	Assoc. Conditions	Description	Author
☆  Coversheet	 0		
☆ <input type="text" value="User Defined Workspace 25"/>	 0	Add description	PANORAMA USER
★  Overview	 0		

Note: *The shaded rows cannot be edited.*

Associated Conditions

The Associated Conditions feature allows you to associate a problem with a user-defined workspace. When a problem is associated to a workspace, you are able to select the problem in the Problems applet and access associated workspaces for treatment of that problem.

Note: *You can only associate a problem to user-defined workspaces.*

To add an associated condition to a workspace:


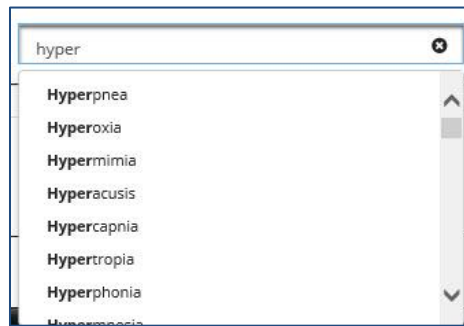
1. Select the **Assoc. Conditions** button . The *Search Problems* field displays (Figure 8-4).

Figure 8-4 Search Problems Field

No Associations Added

2. Enter a **problem** in the field. The results populate as you enter text (Figure 8-5).

Figure 8-5 Problems Search Results



3. Select the **problem** from the list.
4. You can add another problem by entering it in the *Search Problems* field, or exit by selecting the **CLOSE** button.
5. The problem is now associated with that workspace and a number displays next to Assoc. Conditions button (Figure 8-6).

Figure 8-6 Number of Associated Problems



Workspace Description

The Description field allows you to describe a user-defined workspace. This feature is associated with the workspace filter function. Select in the **Add description** field (Figure 8-7) to add a workspace description.

Note: *The workspace description is optional and not required.*

Figure 8-7 Workspace Add Description Field


Workspace Manager		
Assoc. Conditions	Description	Author
0		
0	Add description	PANORAMA USER

Note: *A description cannot be added to the shaded rows.*

Duplicate Workspaces

You are able to duplicate any workspace so that you can modify the cloned workspace.

To create a copy of an existing workspace:

1. Select the **Duplicate** button  to create a duplicated workspace. A copy of that workspace displays in the Workspace Manager list, under the duplicated workspace (Figure 8-8).

Note: Refer to *Customize a New Workspace* for guidance on customizing a duplicated workspace.

Figure 8-8 Duplicated Workspace

	Title	Assoc. Conditions	Description	Author					
☆	Coversheet	0				1		Preview	Launch
☆	Coversheet Copy	0	Add description	PANORAMA USER		1		Preview	Customize

Rearrange Workspace Listings

You are able to rearrange the workspaces so that you can order which workspaces you would like listed first in the Workspace Manager and the workspace dropdown listing.

To rearrange the Workspace Manager listing:











1. Open the **Workspace Manager**.
2. Select the **Rearrange** button . The selected row is highlighted (Figure 8-9).

Figure 8-9 Rearrange Workspaces

Workspace Manager									
<div>  Filter by title or description </div>									
	Title	Assoc. Conditions	Description	Author					
☆	Coversheet	0				1		Preview	Launch
☆	Timeline	0				1		Preview	Launch
★	Overview	0				1		Preview	Launch
☆	Meds Review	0				1		Preview	Launch

3. Select and drag, or use the up and down arrow keys on the keyboard, to move the **workspace row** to the desired spot.
4. Repeat **steps 2 and 3** to continue rearranging the workspace listing.


Delete Workspaces

You are able to delete any user-defined or duplicated workspace. Once deleted, these workspaces are no longer available.

To delete a user-defined or duplicated workspace:

1. Select the **Delete** button . A confirmation message displays (Figure 8-10).

Figure 8-10 Delete User-Defined Workspace Confirmation Message

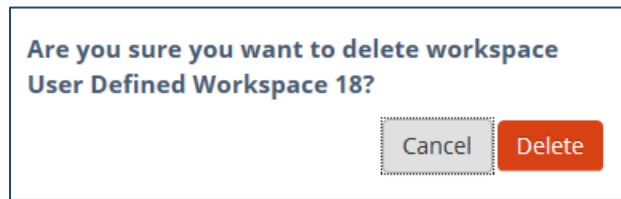

2. Select **Delete** to delete the workspace, or select **Cancel** to return to the Workspace Manager without deleting the workspace.

Note: You are not able to delete workspaces in shaded rows.

Preview Workspace

You can preview any workspace in the Workspace Manager. The preview link displays the workspace so you can review the layout prior to launching it in the Main Application Window.


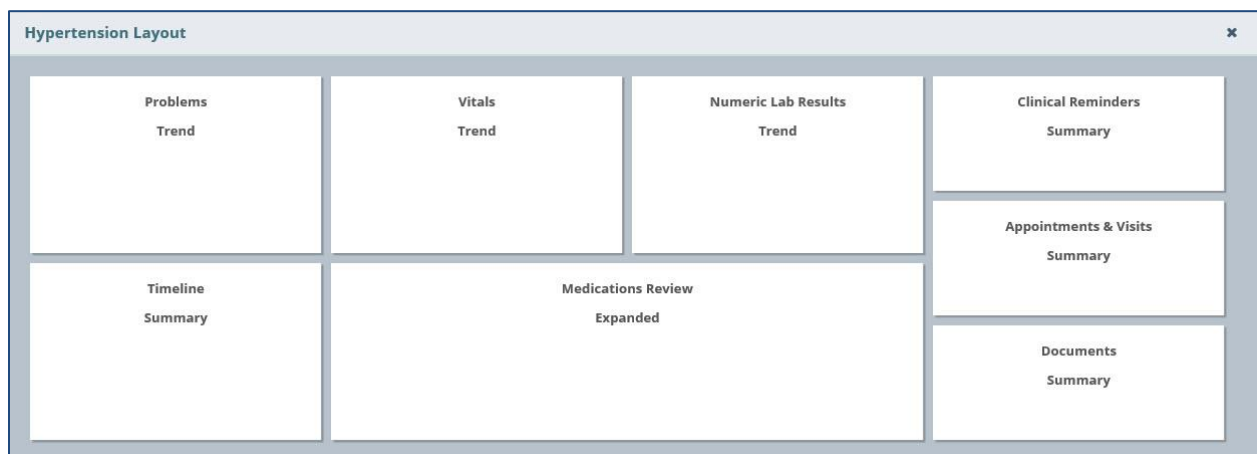
1. Select the **Preview** link  in the Workspace Manager to open a preview of the selected workspace in a new window (Figure 8-11).

Figure 8-11 Workspace Preview



Note: The Preview functionality is not available for user-defined workspaces that have not been customized. The Preview is shaded out to indicate the link is unavailable.

2. Select the **X** to return to the Workspace Manager.

Customize a New Workspace

You can customize workspaces you create to enable unique experiences within the application.

To implement the customization of a new workspace:


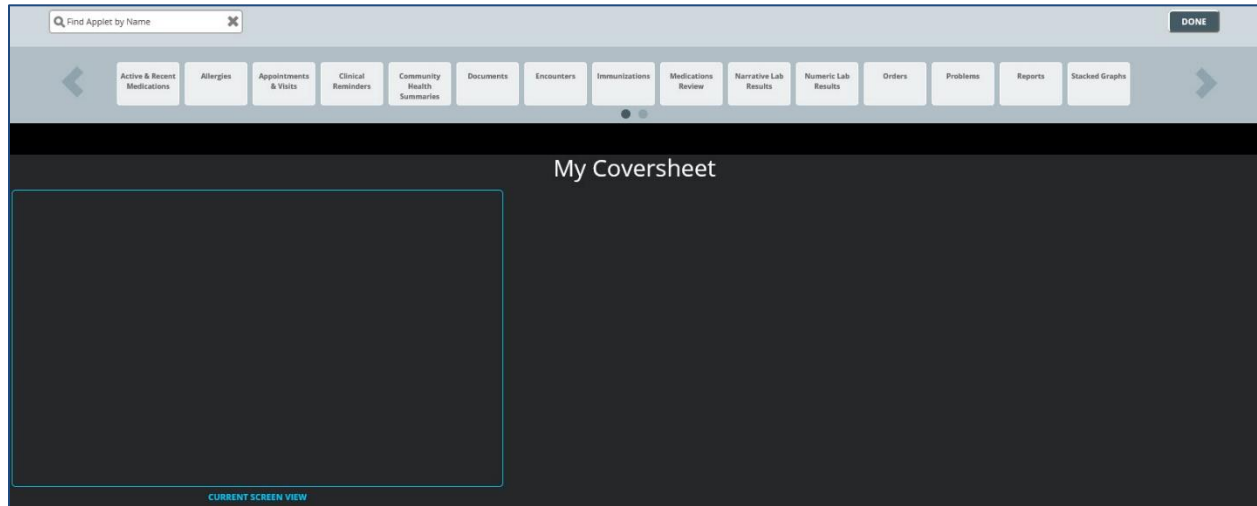
1. Select the **Customize** link  in the Workspace Manager to open the Workspace Editor screen. The Workspace Editor screen displays with a listing of the applets in alphabetical order (Figure 8-12).

Figure 8-12 Workspace Editor Screen



2. Use the arrows to scroll through the **applets**, or enter the **applet name** in the *Find Applet by Name* field above the applets.
3. Drag and drop the desired **applet** in the **CURRENT SCREEN VIEW** box.
Note: Any applets outside the **CURRENT SCREEN VIEW** box are included in the workspace. Use the scrollbar at the bottom of the screen to view applets outside of the box when viewing the workspace in the Main Application Window.
4. Select a View option displays (Figure 8-13). Select the **desired view** for that applet. A box with the applet name and view displays in the **CURRENT SCREEN VIEW** box.


Figure 8-13 Select a View



5. Select and drag the **applets** to move them around on the screen.
6. Use the **double arrows** to modify the size of the applet (Figure 8-14).

Figure 8-14 Using Double Arrows to Modify Applet Size




7. Select the **settings** button  to change the applet view or to delete it from the workspace.
8. Repeat **steps 2-7** until the workspace layout is complete.
9. Select the **DONE** button in the upper right hand corner of the screen to open the workspace in the Main Application Window.

Note: User-defined workspaces can be edited after they have been customized.

Launch Workspaces

Once the user-defined workspace has been customized, you are able to open the workspace in

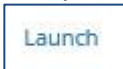the Main Application Window using the Launch link  from the Workspace Manager.

You can also edit the user-defined workspace on the Main Application Window.

To edit the user-defined workspace from the Main Application Window:

1. Launch the **workspace** from the Workspace Manager, or from the Workspace dropdown menu. The workspace opens in the Main Application Window (Figure 8-15).



Figure 8-15 Launch Workspaces

The screenshot shows the eHMP Main Application Window. At the top, there is a patient selection bar with a search icon, a patient profile (Eight, Patient, DOB: 04/07/1935, SSN: 666-00-0008), and a status bar (Inpatient, POSTINGS: C, W, A, D, F). Below this is a navigation bar with a date range (12/08/2014 - 12/08/2016) and a dropdown menu for 'User Defined Workspace 3'. The main content area is divided into two sections: 'NARRATIVE LAB RESULTS' and 'PROBLEMS'. The 'NARRATIVE LAB RESULTS' section contains a table with columns for Date, Description, Type, Author or Verifier, and Facility. The 'PROBLEMS' section contains a table with columns for Problem, Status, and Action. The bottom of the window shows the eHMP version (1.2.6.72715) and a 'Refresh All Data' button.

Date	Description	Type	Author or Verifier	Facility
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1
02/03/2015 - 12:24	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST2
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST1
02/03/2015 - 12:15	LR MICROBIOLOGY REPORT	Laboratory Report	None	TST2

Problem	Status	Action
To Hand	Active	DOD
Bone Pain (Finding)	Chronic	Active
Pain In Limb (Finding)	Chronic	Active
Foot Pain (Finding)	Chronic	Active
Neck Pain (Finding)	Chronic	Active
Sinus Pressure	Chronic	Active

2. Select and drag the **applet** to move it on the screen.

3. Select the **settings**  button on the applet to select a different view, or to delete the applet from the workspace.
4. Hover the mouse on the bottom left corner of the applet, then **select and drag** to re-size the applet.
5. Select the **Workspace Editor** button  to add applets to the workspace.
6. Repeat **steps 2-5** until the customization is complete.

Chapter 9: Other eHMP Features


Search Record Field

The Search Record field allows you to search within a selected patient's record for specific information.

To search the selected patient's record:

1. Enter the **specific item** (i.e., medication, document, immunization) in the *Search Record* field. A list of suggested terms drops down from the field as you enter the text (Figure 9-1).

Figure 9-1 Search Record Field Dropdown

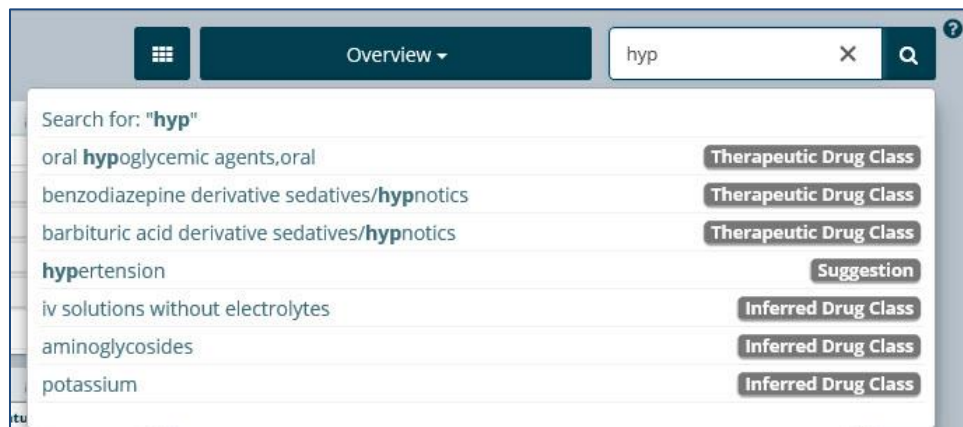

2. Select the **item** from the list, or press **Enter** to run the search.
3. The number of results display and are categorized (Figure 9-2).

Figure 9-2 Search Record Results



4. Use the **preset date ranges** or the **calendar dropdowns** to filter the results for a specific timeframe.
5. Select the **dropdown** arrows to display the results found within each category (Figure 9-3). The selected keyword is highlighted.

Figure 9-3 Selected Keyword

9 results

Search is still under development and may not return all expected results. Please review the patient record manually to ensure accuracy.

DISCHARGE SUMMARY			
03/25/2004 - 19:17	Vehu,Ten ... CHF exacerbation 3) Hyperlipidemia 4) Hypertension 5 for CHF, patient has a scale at home. 3) Hypertension Patient has controlled hypertension. On last admission, pt...	ABILENE (CAA)	
03/25/2004 - 19:17	Vehu,Ten ... CHF exacerbation 3) Hyperlipidemia 4) Hypertension 5 for CHF, patient has a scale at home. 3) Hypertension Patient has controlled hypertension. On last admission, pt...	ABILENE (CAA)	
03/25/2004 - 19:17	Vehu,Ten ... CHF exacerbation 3) Hyperlipidemia 4) Hypertension 5 for CHF, patient has a scale at home. 3) Hypertension Patient has controlled hypertension. On last admission, pt...	ABILENE (CAA)	

PROBLEM

6. Select an **item** from the results list to open a detail dialog box with the selected term highlighted (Figure 9-4).

Figure 9-4 Selected Keyword Detail Dialog Box

Discharge Summary Details

Facility: ABILENE (CAA)
Author: Vehu,Ten
Status: Completed
Attending: Vehu,Ten
Date/Time: 03/25/2004 - 19:17
Expected Cosigner:

Admission Diagnosis:
 1) Chest pain
 2) Shortness of breath

Discharge Diagnosis:
 1) Acute coronary syndromes
 2) CHF exacerbation
 3) Hyperlipidemia
 4) Hypertension
 5) Controlled Diabetes Mellitus II:

Home Outpatient Medications:
 1) Metoprolol 50mg PO BID
 2) Metformin HCL 500 mg PO BID
 3) Simvastatin 40mg PO QPM
 4) Aspirin 81mg PO QDAY

Discharge Medications:
 1) Metoprolol 50 mg PO BID
 2) Metformin HCL 500 mg PO BID
 3) Simvastatin 40mg PO QPM

Close

7. Select **CLOSE** to close the detail dialog box and return to the results list.

Patient Record Refresh and Status Bar

The Refresh All Data button, Status Bar, and the eHMP Data Sources button appear in the bottom right-hand corner of the eHMP application (Figure 9-5). The Refresh All Data button provides the opportunity to refresh the selected patient's data for the most up-to-date information from all sources. The status bar displays the status of each repository that contributes to the aggregated data displayed in eHMP. Green bullets with check marks indicate that data is being synced from that source.

Figure 9-5 eHMP Status Bar



The repositories used to aggregate data:

- My Site: Your current site
- All VA: All VA sites
- DoD: Department of Defense sites
- Communities: Community Health Partners of the VA

The eHMP Data Sources status screen (Figure 9-6) provides detailed information on the refresh status of data by source.


Figure 9-6 eHMP Data Sources

eHMP Data Sources			
Source	My Site		
	Refresh My Site		
<div> <div> <div></div> <div>My Site</div> </div> <div> <div>New data since</div> <div>Current</div> </div> </div>	Domain	Last Refresh	New Data Since
	Allergies	Current	Current
	Appointments	Current	Current
<div> <div> <div></div> <div>All VA</div> </div> <div> <div>Last refresh < 1h ago</div> <div>New data since</div> <div>Current</div> </div> </div>	Consult	Current	Current
	CPT	Current	Current
	Document	Current	Current
<div> <div> <div></div> <div>DoD</div> </div> <div> <div>Last refresh < 1h ago</div> </div> </div>	Health Factor	Current	Current
	Image	Current	Current
	Lab Results	Current	Current
<div> <div> <div></div> <div>Communities</div> </div> <div> <div>Last refresh < 1h ago</div> </div> </div>	Active Medications	Current	Current
	Orders	Current	Current
	Patient	Current	Current
	Purpose of Visit	Current	Current
	Active Problems	Current	Current
	Surgery	Current	Current
	Visit	Current	Current
	Vitals	Current	Current

Refresh All Data

Close

To view the refresh status of patient data:

1. Select the **data sources** button  to open a new window with the sync statuses for all sources. The new window displays the source, when the data was last refreshed, and if there is new data since the last refresh.
2. Select the **Source title** to display the detailed status for that source. You can refresh all data, or refresh each source individually, from this window.
3. Select the **Close** button to return to the current workspace.