Joint Longitudinal Viewer (JLV) 3.6.0.0 User Guide



November 2023 Version 1.0

Department of Veterans Affairs

Office of Information and Technology (OIT)

Revision History

Date	Version	Document	Author
10/12/2023	10/12/2023 1.0 Submitting the document for approval		Booz Allen Hamilton
10/03/2023	0.1	Initial draft of document from last approved	Booz Allen Hamilton

Artifact Rationale

Per the Veteran-focused Integrated Process (VIP) Guide, the user guide is required to be completed prior to Critical Decision Point #2 (CD2), with the expectation that it is updated as needed. A user guide is a technical communication document intended to give assistance to people using a system, such as Veterans Health Information Systems and Technology Architecture (VistA) end users. It is usually written by a technical writer, although it can also be written by programmers, product or project managers, or other technical staff. Most user guides contain both a written guide and the associated images. It is usual to include screenshots of the human-machine interfaces, and hardware manuals often include clear, simplified diagrams. The language used is matched to the intended audience, with jargon kept to a minimum or explained thoroughly. The user guide is a mandatory, build-level document, and should be updated to reflect the contents of the most recently deployed build. The sections documented herein are required if applicable to your product.

Table of Contents

1.	Introduction	1
	1.1. Purpose of the Guide	1
	1.1.1. Guide Conventions	2
	1.1.2. Terminology	2
	1.2. Assumptions	2
	1.3. System Requirements	2
	1.4. Getting Help	2
2.	Logging Into JLV	3
	2.1. VHA User Login	3
	2.1.1. Link Your PIV Card and VistA Account	5
	2.2. CAPRI-Claims User Login	6
	2.2.1. Editing Your User Profile at Login	9
	2.3. Logging Out of the Current JLV Session	10
3.	The Elements and Functionality of JLV	10
	3.1. The Patient Data and Widgets Page	10
	3.1.1. Default View	10
	3.1.2. System Menu	11
	3.1.3. Viewing System Status	12
	3.1.4. Patient Context Management	14
	3.1.4.1. Suspending Context	15
	3.1.4.2. Establishing Context	15
	3.1.4.3. Unable to Communicate with Vault	16
	3.1.4.4. Context Synchronization Failure	16
	3.2. Customizing JLV	17
	3.2.1. Cache Manager	19
	3.2.2. About Data Display Settings	20
	3.2.3. About UI Themes	21
	3.2.4. Using the Accessible Theme's Functionality	23
	3.2.4.1. Keyboard Focus and Navigation	24
	3.2.4.2. Accessible Theme Navigation	26
	3.2.4.3. Table Functionality	27
	3.2.4.4. Widget Move	28
	3.2.4.5. Widget Resize	29
	3.3. User-Restricted Access to Patient Data	30
	3.3.1. VHA Users Viewing DOD-only Patient Data	30
	3.3.2. CAPRI-Claims Users Viewing Patients with DOD Data	31
	3.3.2.1. Patient Registered in MPI	31

	3.3.2.2. Patient Not Registered in MPI	31
	3.3.3. Access to DOD Sensitive Records	32
	3.4. Patient Searches	33
	3.4.1. PDWS Search Rules and Supported Searches	33
	3.4.2. MPI Search Guidelines	34
	3.4.3. Patient Search Using DOD ID	34
	3.4.3.1. Invalid and Maximum Characters	35
	3.4.4. Patient Search Using SSN	37
	3.4.5. Advanced Patient Search Options	39
	3.4.6. Family Member Search	40
	3.4.7. Recently Viewed Patients List	42
4.	Widgets	43
	4.1. Accessing and Opening Widgets	
	4.2. Widget Navigation and Display Options	
	4.3. Widget Toolbars	
	4.4. Minimized vs. Expanded Widget View	46
	4.5. Sort and Filter Options	
	4.5.1. Widget Date Range Filters	49
	4.5.1.1. Quick Date Range Filter	50
	4.5.1.2. Start and End Date Filter	51
	4.5.1.3. Tab Date and Text Filter Tool	52
	4.5.2. Text Filter	53
	4.6. Viewing Connection Status	55
	4.7. Displaying Widgets on Custom Tabs	58
	4.7.1. Exporting Tabs	59
	4.7.2. Importing Tabs	60
	4.8. Refreshing Data in a Widget	62
	4.9. Terminology Normalization	62
	4.10. Image Support	63
	4.11. Using Report Builder	64
	4.11.1. Adding One Record	66
	4.11.2. Removing One Record Directly from a Widget	66
	4.11.3. Adding Multiple Records with EZ Select	67
	4.11.4. Adding Details from Multiple Dates	68
	4.11.5. Adding Sensitive Records	69
	4.11.6. Adding Only Record Details or Notes	69
	4.11.7. Adding Patient Demographics	70
	4.11.8. Generating a Report	71
	4.11.9. Opening a Report	72

	4.11.10. Delete a Report	73
	4.11.11. Downloading and Printing a Report	73
	4.12. Printing	76
5.	Widget Functionality	76
	5.1. Patient Demographics Widget	
	5.1.1. Patient Flags	78
	5.1.2. Advance Directives	79
	5.1.3. Viewing Third-Party Insurance Information	80
	5.1.4. Viewing Military Service Data	81
	5.1.5. Viewing Eligibility and Enrollment Data	82
	5.2. Cardiology StudiesFEHR and VA MUSE Only Widget	83
	5.2.1. Cardiology StudiesFEHR and VA MUSE Only Widget Data	83
	5.3. Community Health Summaries and Documents	83
	5.3.1. Community Health Summaries and Documents Widget Data	84
	5.3.2. Viewing VA Community Health Summary Documents	85
	5.4. Imaging Widget	86
	5.4.1. Imaging Widget Data	87
	5.4.2. Viewing Images	87
	5.5. Documents & Images Widget	90
	5.5.1. Documents & Images Widget Data	90
	5.5.2. Viewing Documents & Images	92
	5.6. Federal EHRH Military Health System (MHS) GENESIS Widget	92
	5.6.1. Federal EHR/MHS GENESIS Widget Data	93
	5.6.2. Viewing Federal EHR/MHS GENESIS Documents	93
	5.7. Health Summaries and Reports—VistA Only Widget	94
	5.7.1. Health Summaries and Reports—VistA Only Widget Data	94
	5.7.2. Viewing Health Summaries	96
	5.8. Lab Results Widget	97
	5.8.1. Lab Results Widget Data	97
	5.8.2. Viewing Lab Results Details	99
	5.8.3. Viewing Lab Panel Results Details	99
	5.8.4. Viewing Abnormal Results in the Lab Results Widget	100
	5.8.5. Viewing Results and Interpretations in the Lab Results Widget	101
	5.8.6. Graph/Table View from the Lab Results Widget	102
6.	Troubleshooting	103
	6.1. Special Instructions for Error Correction	104
	6.1.1. Login Page Error Messages	105
	612 System Error Messages	107

Α.	Acronyms and Abbreviations1	09
B.	Patient Data and Widgets page Widgets, Columns, and Supporte	d:
	Views	
		• •
	Table of Figures	
E:	Table of Figures	4
_	e 1: SSOi Page for VA Userse 2: Home VistA Site Dialog for VHA Users	
_	e 3: JLV Link My Account Page	
_	e 4: IAM Provisioning Services Link VistA User Page	
_	e 5: Home VistA Site Dialog for CAPRI-Claims	
•	e 6: The JLV Login Page	
_	e 7: User Profile Data Source and User Interface Option Fields	
_	e 8: JLV Logout Icon	
_	e 9: Default View	
_	e 10: JLV Tools on the Patient Data and Widgets Page	
•	e 11: More JLV Functionality	
_	e 12: System Status on the JLV Login Page	
•	e 13: System Status on the JLV Toolbar	
_	e 14: Context Indicator	
•	e 15: Suspending Context	
•	e 16: CCOW Warning Message	
•	e 17: CCOW Error Message	
_	e 18: Context Synchronization Failure	
•	e 19: JLV Settings Dialog Box	
_	e 20: VHA JLV User Management Settings Tab	
_	e 21: Admin Cache Manager	
_	e 22: Blue (Default) UI Theme	
_	e 23: Green UI Theme	
_	e 24: Dark Gray UI Theme	
	e 25: Gray UI Theme	
_	e 26: Accessible (Section 508-Compliant) UI Theme	
_	e 27: Navigation Using Elements with Keyboard Focus	
_	e 28: Document Details with Keyboard Focus	
_	e 29:Table Column Header	
_	e 30: Widget Column Select Dropdown	
_	e 31: Widget Move Icon	
_	e 32: Widget Move Dropdown	
_	e 33: Widget Resize Icon	
	e 34: Widget Resize Dropdown	
_	e 35: VA User Restricted Access Dialog	
_	e 36: CAPRI-Claims User Audit Notification	
•	e 37: Restricted Access Notification	32

Figure 38: DOD Sensitive Record and Audit Dialog Box	33
Figure 39: DOD ID (EDIPI) Search	35
Figure 40: Invalid Character Entry	36
Figure 41: Maximum Characters Reached	
Figure 42: SSN Search	38
Figure 43: Additional Identifiers Needed	38
Figure 44: Advanced Search Options	39
Figure 45: Sponsor SSN Search	41
Figure 46: Family Members Dialog Box	42
Figure 47: Recently Viewed Patients List.	43
Figure 48: Widget Tray	44
Figure 49: Dynamically Expanded Filter Dropdown in the Documents & Images Widget	49
Figure 50: Date Range Filters Applied	50
Figure 51: Quick Date Range Filter Options	51
Figure 52: Start and End Date Display Filter	52
Figure 53: Tab Date and Text Filter Tool	53
Figure 54: Documents & Images Widget Text Filter	54
Figure 55: Demographics Widget, PCMM Connection Status Unavailable	56
Figure 56: Connection Status Indicator	56
Figure 57: Connection Status Details	57
Figure 58: Connection Status Window	58
Figure 59: Add Custom Tab Icon	58
Figure 60: Add Tab Dialog Box	59
Figure 61: Export Tabs	60
Figure 62: Import Tabs	61
Figure 63: Import Successful	61
Figure 64: Import Failed	62
Figure 65: Camera Icon in the Image Column	63
Figure 66: Thumbnail and Link in Details View of the Documents & Images Widget	64
Figure 67: Report Builder Icon	64
Figure 68: Report Builder Header Icons	64
Figure 69: New Report Tab	65
Figure 70: Saved Reports Tab.	66
Figure 71: Report Builder EZ Select	67
Figure 72: Records Added to Report Builder	68
Figure 73: Adding Date-Flexible Details to Report Builder, Medicine Full Report	69
Figure 74: Report Builder Options	70
Figure 75: Demographics Widget Add to Report Builder Icon	70
Figure 76: Name Report Dialog Box	71
Figure 77: Saved Reports Tab Report Ready and Delete Option	72
Figure 78: Report Ready Message and Indicator	72
Figure 79: Report Contents Window	73
Figure 80: Delete Reports Dialog Box	73
Figure 81: Report Hyperlink	74

Figure 82: Printing the Downloaded Report	
Figure 83: Download Selected Option	75
Figure 84: Opening the PDF	75
Figure 85: Print Icon on the PDF	76
Figure 86: Patient Demographics Widget in the JLV Toolbar	77
Figure 87: Patient Demographics Widget Site Tabs	77
Figure 88: Patient Flags Indicator	78
Figure 89: Patient Flags Details	79
Figure 90: Advanced Directives Dialog	80
Figure 91: Insurance Information	81
Figure 92: Military Service Window	82
Figure 93: Eligibility and Enrollment Window	82
Figure 94: Cardiology StudiesFEHR and VA MUSE Only Widget	83
Figure 95: Community Health Summaries and Documents, Minimized View	84
Figure 96: Community Health Summaries and Documents, Expanded View	85
Figure 97: Community Health Document	86
Figure 98: Imaging Widget, Minimized View	
Figure 99: Imaging Widget	
Figure 100: View of an Image	
Figure 101: View of a PDF-formatted Document	
Figure 102: Documents & Images Widget, Minimized View	
Figure 103: Documents & Images Widget, Expanded View	
Figure 104: Details View of a Document	
Figure 105: Federal EHR/MHS GENESIS Widget	
Figure 106: C-CDA Document	
Figure 107: Health Summaries and Reports—VistA Only Widget, Minimized View	
Figure 108: Health Summaries and Reports—VistA Only Widget, Expanded View	
Figure 109: Health Summary from a Single, Local Site	
Figure 110: National Health Summary, Multiple Sites	
Figure 111: Lab Results Widget, Expanded View	
Figure 112: Lab Results Details	
Figure 113: Lab Panel Results Detail	
Figure 114: Lab Results Widget, Abnormal Results in Expanded View	
Figure 115: Lab Results Widget, Expanded View, Interpretations Available	
Figure 116: Lab Results Interpretation Details	
Figure 117: Lab Results Graph/Table View	
Figure 118: Lab Results Graph Table View, Select a Lab Test	
Figure 119: Login Error—CAPRI Patient Selection Restrictions	
Tigule 117. Bogin Bilot Cili Iti I utote Solection Resultations	
Table of Tables Table 1: JLV User Profiles	1
Table 2: System Status Notices by Service	
Table 3: Data Sources	
1 au 10 J. Data Dutters	∠U

Table 4: Accessible Theme Keyboard Shortcuts – Patient Data& Widget Navigation	24
Table 5: Accessible Theme Keyboard Shortcuts - Tabs	24
Table 6: Accessible Theme Keyboard Shortcuts – Windows or Dialog Boxes	24
Table 7: Accessible Theme Keyboard Shortcuts – Adding Widgets to the Patient Data and Wi	dgets Page
	25
Table 8: Navigation Icons and Actions	45
Table 9: Widget Toolbar Icons	45
Table 10: Available Widget Filters	48
Table 11: Terminology Specific to Widgets	62
Table 12: Gender ID (SIGI) Code Values	77
Table 13: Error Messages and Resolution Steps	105
Table 14: System Error Messages	107
Table 15: Acronyms and Abbreviations	
Table 16: Widget Details	

1. Introduction

Born from a joint Department of Defense (DOD)—Department of Veterans Affairs (VA) venture called JANUS, Joint Longitudinal Viewer (JLV) was directed by the Secretary of the VA and the Secretary of Defense in early 2013 to further support interoperability between the two departments. JLV is a centrally hosted, Java-based web application managed as a single code baseline and deployed in separate DOD and VA environments. Its browser-based, graphical user interface (GUI) provides an integrated, read-only view of Electronic Health Record (EHR) data from the VA, DOD, and community partners within a single application.

JLV eliminates the need for VA and DOD clinicians to access disparate viewers. The GUI retrieves clinical data from several native data sources and systems, then presents it to the user via widgets, each corresponding to a clinical data domain. Users can create and personalize tabs, drag, and drop widgets onto tabs, sort data within a widget's columns, set date filters, and expand a widget for a detailed view of patient information.

This document is intended for VA JLV users supporting the:

- Veterans Health Administration (VHA)
- Veterans Benefits Administration (VBA)

<u>Table 1</u> describes authorized JLV users and their responsibilities.

User	Description and Responsibilities	
Compensation and Pension Record Interchange (CAPRI)-Claims	VA administrative staff who access patient EHRs to assist in VBA processes	
VHA	VA clinicians and administrative staff who access patient EHRs to assist in providing health care services	
DOD Clinician	DOD clinicians who access patient EHRs	

Table 1: JLV User Profiles

1.1. Purpose of the Guide

The purpose of the user guide is to familiarize VA users with the important features and navigational elements of JLV.

The major features of the JLV GUI include:

- Access to patient data (through widgets)
- Patient context management
- User-restricted access (Break the Glass)
- Use of the Report Builder
- Option to use the accessible (Section 508-compliant) interface

1.1.1. Guide Conventions

This document is designed for both online and hardcopy consumption.

- <u>Cross References</u> are indicated by blue, underlined text and provide a hyperlink to figures, tables, and other sections within this guide.
- Emphasis is expressed by **bold** and *italicized* text.
- The information symbol i calls the reader's attention to additional information.

1.1.2. Terminology

The following standard terms are used throughout this guide:

- Patient Data and Widgets: The page that displays patient data through widgets.
- **Tab:** A component of the JLV interface that displays patient data and user-customized views of widgets.
- **Widget:** A component of the JLV interface that enables a user to view information or perform a function.
- **Widget Tray:** An expandable and collapsible tray on the Patient Data and Widgets page that provides access to all widgets available for placement on the Patient Data and Widgets page.

1.2. Assumptions

The user guide is written from the perspective of VA users, assuming the following:

- You can open, navigate, and use a web browser.
- You can use web-based applications, their menu options, and navigation tools.
- You have the Uniform Resource Locator (URL) for the JLV **Login** page, a system username, and the required Access and Verify codes.
- You are using the functionality of JLV to support VHA and/or VBA workflows.

1.3. System Requirements

JLV is optimized for use with the VA standard browsers, Microsoft Edge, Chrome, Internet Explorer (IE) version 11 and Safari. Other web browsers may function but are not certified for support. It is recommended that you view JLV using:

- A monitor set to a minimum screen resolution of 1024 x 768.
- An approved browser, with the browser's zoom percentage set to 100%.

It is also recommended that Adobe Reader be installed on your system.

1.4. Getting Help

Please see online help for how-to information within JLV. Please visit REDACTED for access to JLV training videos and additional JLV training materials.

Prior to contacting the Enterprise Service Desk (ESD) for support, please refer to <u>Logging Into JLV</u> for detailed information about how to access JLV and to <u>Troubleshooting</u> for suggested resolution steps and troubleshooting information.

If you are an authorized user having trouble logging into JLV or experiencing other application issues, please contact the ESD via telephone or by using the YourIT self-service portal for assistance.



YourIT Self-Service Portal: REDCACTED

If you are unable to retrieve community partner documents for a patient, please contact your local Veterans Health Information Exchange (VHIE) Coordinator. If you need assistance locating your local contact, please e-mail REDACTED.

2. Logging Into JLV

Users who have access to the Computerized Patient Record System (CPRS), CAPRI, and Cerner PowerChart are automatically authorized to use JLV.

0


NOTE: There is no functionality available to change a user identification (ID) and password in JLV. See <u>Troubleshooting</u> for information about login errors, Access code errors, or Verify code errors.

2.1. VHA User Login

To log into JLV through SSOi, you must first link your Personal Identity Verification (PIV) card to your VistA account (see <u>Link Your PIV Card and VistA Account</u>). If you can log into CPRS or JLV without entering your Access and Verify code, you are already linked.

- 1. Enter the JLV URL (REDACTED) into the address bar of your Internet browser, select the **JLV** button CPRS if available, or select the **JLV** button in PowerChart.
 - a. You are redirected to the VA Single Sign-On Internal (SSOi) page (Figure 1).

Figure 1: SSOi Page for VA Users U.S. Department

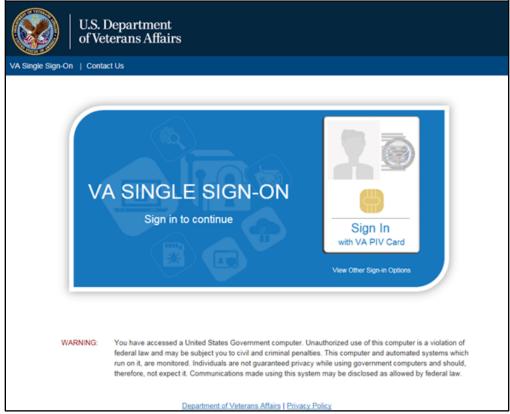

- 2. Click the **Sign in with VA PIV Card** button.
- 3. Select the current authentication certificate Veterans Affairs User CA B1, and click **OK.**
- 4. Enter your Personal Identification Number (PIN), and click **OK.**
 - a. If the **Home VistA Site** dialog appears (Figure 2), select your site, and click **Save.**

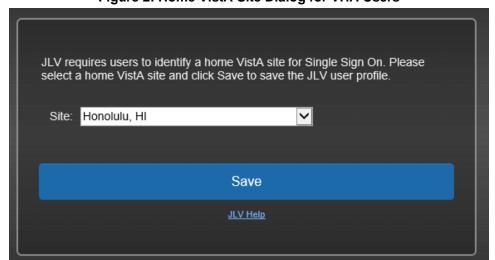


Figure 2: Home VistA Site Dialog for VHA Users

- b. If your PIV card is properly linked to the home VistA site you selected, the **Government Information System Acknowledgement** displays.
- c. If your PIV is not properly linked to the home VistA site you selected, the **Link My Account** page displays (Figure 3, Link Your PIV Card and VistA Account).
- 5. Read the important user consent information regarding access to a Government information system, and click **I Accept.**
 - a. You are redirected to JLV (Figure 9).

2.1.1. Link Your PIV Card and VistA Account

You need your Access and Verify codes for your VistA account and a PIV card PIN to link your PIV credentials to your VistA account. If you do not remember your PIN or experience any other issue related to your PIV card, visit your local PIV office.

- 1. Complete steps 1–4 listed in <u>VHA User Login</u>.
 - b. Link My Account page displays (Figure 3).



Figure 3: JLV Link My Account Page

- 6. Select the **Link My Account** button.
 - a. JLV redirects the user to the Identity Access Management (IAM) Provisioning Services page (Figure 4).

IAM Provisioning Services [VAEC Development] Tuesday May 5th, 2020 VistA instances that are already bound-Name Instance My User Account Information Network Id Not You? User Name VA Email Address @va.gov User Account Request Information-* Link Account Select a VistA Instance * Access Code * Verify Code Submit

Figure 4: IAM Provisioning Services Link VistA User Page

- NOTE: Close all other open applications and browser windows. You may experience a failure to link your account(s) if you are logged into CPRS and your VistA application. If this occurs, you should log out and not log back into CPRS until you receive an email stating your account is linked.
- 7. Click the Link Vista User link (Figure 4)
 - NOTE: Ignore the "No VistA stations linked to your account in provisioning" message that displays.
- 8. Select the appropriate instance of VistA from the **Link Account** dropdown.
 - NOTE: The VistA instance list is sorted by station number. If your station is a divisional site of the main PARENT VA site, select the parent site from the list. Example: if your divisional site is Albany, select the parent site 528 Upstate New York HCS from the list.
- 9. Enter your Access and Verify codes for the selected VistA instance.
- 10. Click Submit.


2.2. CAPRI-Claims User Login

The following instructions apply to CAPRI-Claims users accessing JLV.

- NOTE: Users with any CAPRI patient selection restrictions must use the **Joint Longitudinal Viewer** tab in CAPRI to access JLV.
- 1. Enter the URL into the address bar of your Internet browser (REDACTED), or select the **Joint Longitudinal Viewer** tab in CAPRI.

- a. You are redirected to the **VA SSOi** page (Figure 1).
- 11. Click the **Sign in with VA PIV Card** button.
- 12. Select the authentication certificate *Veterans Affairs User CA B1*, and click **OK**.
- 13. Enter your PIN, and click **OK**.
 - a. If the **Home VistA Site** dialog appears (<u>Figure 5</u>), select *CAPRI-Claims or All VBA*, and click **Save**. All VBA will be the first entry in the drop-down list.

Figure 5: Home VistA Site Dialog for CAPRI-Claims

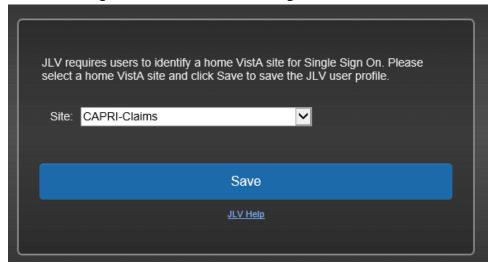

- iNOTE: When CAPRI/Claims is stored as a user's JLV Login site and the user's PIV card is linked to the Claims system in IAM, the user is only prompted for PIV Cert and PIN to authenticate to JLV.
- Note: A new Cerner user, or a user without a login site saved to their profile, will have the option to select EHRM/Oracle Cerner as their site. Users with access to Cerner and Cerner as the login site will have the ability to log into standalone JLV. Users without Cerner access as determined by IAM will be prohibited access to the Cerner site login.
- b. You are redirected to the JLV **Login** page (Figure 6).

Joint Longitudinal Viewer

Accessibility Mode

VistA Access Code:

VistA Verify Code:

Agency: VA Site: CAPRI-Claims

Edit Profile

Login

System Status

JIV data sources available.

JIV Help

Build 2.9.1.0.0.20201123

Figure 6: The JLV Login Page

- 14. Read the important user consent information regarding access to a Government information system, and click **I Accept**.
- 15. If desired, check the **Accessibility Mode** option.
 - a. If selected, JLV opens and displays the application in the Section 508-compliant Accessible UI theme (See <u>Using the Accessible Theme's Functionality</u> for detailed information).
- 16. Enter your login credentials:
 - a. Enter your CAPRI Access code.
 - b. Enter your CAPRI Verify code.
 - c. Select VA from the **Agency** field.
 - d. Select CAPRI-Claims or All VBA from the Site dropdown list.
- 17. Click Login.
 - a. You are redirected to JLV (Figure 9).

2.2.1. Editing Your User Profile at Login

One method of setting application configuration options (i.e., customizing) is by using the **Edit Profile** link (Figure 6) to access the user profile fields through the JLV **Login** page (Figure 7). User profile options can also be set within JLV. See <u>Customizing JLV</u>.

- NOTE: Only CAPRI-Claims users and those who haven't yet linked their VistA accounts to their PIV can load the **Login** page and access the **Edit Profile** link.
- NOTE: The selected options are saved to your user profile by clicking Save and Login and remain set for all subsequent JLV sessions until you change them.

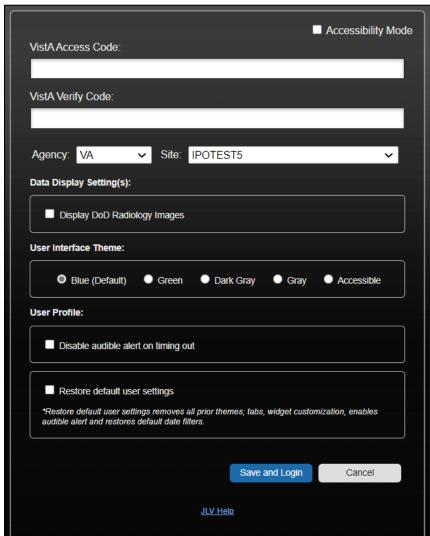


Figure 7: User Profile Data Source and User Interface Option Fields

2.3. Logging Out of the Current JLV Session

Click (highlighted in <u>Figure 8</u>) at any time to end the current session and exit JLV.

Figure 8: JLV Logout Icon



A JLV session ends automatically after 60 minutes of inactivity. Actions in JLV that are considered activity are changing patients, refreshing widgets, opening, or closing widgets, adding, or removing widgets from the Patient Data and Widgets page, and adjusting date ranges. You receive an audible and visual warning that the current session is about to terminate after 55 minutes of inactivity. Automatic termination of a session logs you out of the application, closes all JLV tabs and windows, and displays the **Logout** page. Close and reopen the browser and return to the JLV **Login** page to initiate a new JLV session. If SSOi is enabled, you are redirected to the **VA Single Single-On Logout** page.

NOTE: The audible alert can be disabled by selecting **Disable audible alert on time out** in the **Settings** panel.

3. The Elements and Functionality of JLV

3.1. The Patient Data and Widgets Page


After a successful log in, your username displays on the JLV toolbar, and the Patient Data and Widgets page is displayed.

3.1.1. Default View

Upon log in, you see the **Q Patient Search** dialog box by default.

NOTE: CPRS users who select a patient within CPRS, then launch JLV are taken directly to the Patient Data and Widgets page, where the records for the patient selected in CPRS are displayed.

Figure 9: Default View

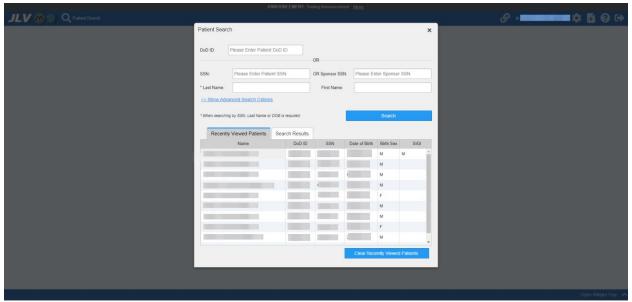


Clicking a patient name in the **Patient Search** window displays the selected patient's clinical records and other patient-centric information.



Note: The Patient Flags window opens and displays by default when the Patient Data and Widgets page opens if there are warnings associated with the selected patient. For more information, see Patient Flags and Alerts.

New JLV users see a preconfigured, default widget layout on the Patient Data and Widgets page when a patient is selected. The default layout includes the **Allergies**, **Problem List**, **Documents & Images**, **Medications - Outpatient**, and **Lab Results** widgets, as well as the **Patient Flags** widget if applicable.

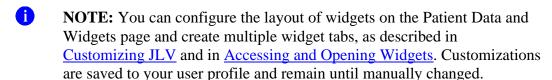

3.1.2. System Menu

Each Patient Data and Widgets page has a toolbar (highlighted in <u>Figure 10</u>) that enables quick access to the following features:

- View Context Management status
- **Settings:** Click to set user preferences, including the UI theme (Refer to <u>Customizing</u> <u>JLV</u>).
- **Report Builder:** Click to create custom reports (Refer to <u>Using Report Builder</u>).

- Online Help: Click **2** to open online help.
- Logout: Click to log out and close the current JLV session (Refer to Logging Out of the Current JLV Session).

Figure 10: JLV Tools on the Patient Data and Widgets Page



The following functions can be completed from the Patient Data and Widgets pages (highlighted in Figure 11):


- Perform a **Q Patient Search**.
- View the system status 1.
- Open the widget tray using the <u>Open Widget Tray ^ link</u> on the Patient Data and Widgets page (as shown in <u>Accessing and Opening Widgets</u>, <u>Figure 48</u>).

Figure 11: More JLV Functionality



3.1.3. Viewing System Status

The system status reports the condition of the JLV application. Hovering over **System Status** provides additional information.

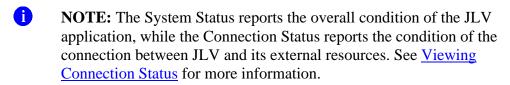

JLV displays the system status on the Login page and on the Patient Data and Widgets toolbar.

- When all monitored systems and services are online and connected, ✓ appears with the status message, "JLV data sources are available".
- When one or more of the monitored systems or services is offline or unavailable, appears with the status message, "Services are not available to JLV".
- When JLV is unavailable, 2 appears with the message, "System status is unavailable".
 - You may not be able to log into JLV or view patient data until the connection is restored.

These warnings may be accompanied by one or more of the notices listed in <u>Table 2</u>, depending on which services are unavailable. <u>Figure 12</u> shows a system status message on the **Login** page, and <u>Figure 13</u> shows the **System Status** hover text on the Patient Data and Widgets toolbar.

Table 2: System Status Notices by Service

Service	Notice		
jMeadows	"The service that gathers patient data and prepares it for display in JLV is unavailable. JLV is available for login but may not function otherwise."		
Master Person Index (MPI)	"VA Patient Identity Service is unavailable. Patient search may not be available, and VA and community partner data may not display."		
Patient Discovery Web Service (PDWS)	"DOD Patient Identity Service is unavailable. Patient search may not be available and DOD data may not display."		
Relay Service	"All DOD, Federal EHR, and Community Partner data in widgets is currently unavailable. The source connection is down and DOD, Federal EHR, and Community Partner data of all types from all sites may not display."		
VistA Data Service (VDS)	"The connection to VA sites is unavailable. VA patient records of all types from all sites may not display."		
Data Exchange Service (DES)	"Some Federal EHR, DOD, and Community Partner data are currently unavailable. The source connection is down and some DOD and Community Partner, and some FEHR data of all types from all sites may not display in widgets.		
Electronic Health Record Modernization (EHRM)	down and some Federal EHR data may not display."		

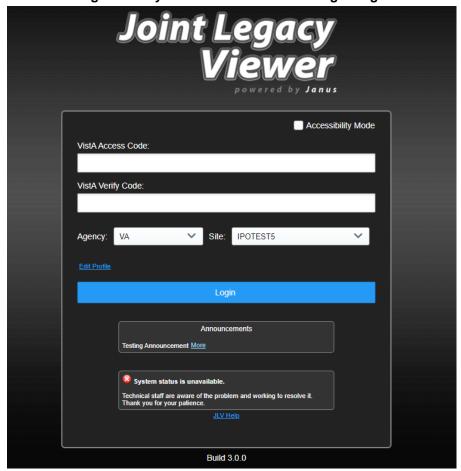


Figure 12: System Status on the JLV Login Page

Figure 13: System Status on the JLV Toolbar



A System Status icon appears in the JLV Navigation Bar to indicate the status.

3.1.4. Patient Context Management

Clinical Context Management (CCOW) is a way for applications to synchronize the display of clinical context, based on the Health Level 7 (HL7) CCOW standard. If CCOW-compliant applications are sharing context and one of the applications changes the data display to a different patient, the other applications switch to display that same patient.

JLV is context management enabled. When context is enabled, a patient change made in one context-enabled application (i.e., CPRS) triggers the same change in JLV. The same effect is shown vice versa; that is, when a patient is selected in JLV, it triggers a patient change in other context-enabled applications.

Context management is enabled by default, and JLV attempts to connect to the context vault upon a valid login. The context status appears on the Patient Data and Widgets toolbar. When context is established, is shown. When context is suspended, is shown.

<u>Figure 14</u> shows the context indicator location on a JLV page. Context is established in this example.

Figure 14: Context Indicator



3.1.4.1. Suspending Context

Context can be manually suspended.

- 1. Click .
- 2. The CCOW warning message, "Are you sure you want to break patient link?" appears (Figure 15).
- 3. Click **Yes** to continue.

Once context is suspended, displays on the Patient Data and Widgets page. When context is suspended, patient changes are not reflected in either JLV or other CCOW-enabled applications.

CCOW Warning

Are you sure you want to break patient link?

OK

Cancel

Figure 15: Suspending Context


3.1.4.2. Establishing Context

When context is suspended, it can be reestablished manually.

- 1. Click %.
- 2. The CCOW warning message, "Are you sure you want to rejoin patient link?" appears (Figure 16).
- 3. Click Yes.

Once context is established, on the Patient Data and Widgets page.


Figure 16: CCOW Warning Message

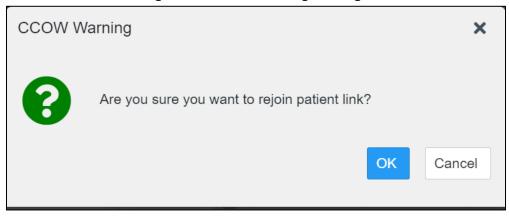

3.1.4.3. Unable to Communicate with Vault

If JLV is unable to initiate context management upon log in, you see the CCOW Error message (Figure 17). Click **OK** to continue. displays on the Patient Data and Widgets page. If the CCOW error persists, contact your local service desk for assistance.

Figure 17: CCOW Error Message

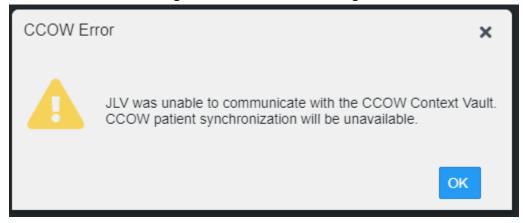

NOTE: This error message automatically displays if you have more than one active JLV session open.

3.1.4.4. Context Synchronization Failure

A CCOW error message is presented if JLV is unable to reestablish context (Figure 18). This indicates context is suspended. Click **OK** to continue. displays on the Patient Data and Widgets page. CAPRI-Claims users may encounter this error message at each log in, as there are usually no CCOW system components installed on their computers.

Figure 18: Context Synchronization Failure



3.2. Customizing JLV

Configuration options can be set in the **JLV Settings** dialog (<u>Figure 19</u>). These settings are saved to your user profile and persist.

- 1. Click the ink on the **Patient Data and Widgets** page (Figure 10).
 - a. The **JLV Settings** dialog box opens (all selections are optional).
- 2. Select the data display setting checkbox for DOD Radiology Images to display within the Radiology Reports Widget.
- 3. Select the radio button to change the color theme.
- 4. To disable audible alert on timing out, select the check box.
- 5. To restore your user profile default settings, click the checkbox.
- 6. You can view your JLV Login Site
- 7. To export or import custom tabs, click the **Export Tabs** or **Import Tabs** buttons.
- 8. Click **Save** to apply the selected data display setting and/or theme **OR**
- 9. Click **Cancel** to exit the dialog box and restore the previous data display and/or theme configuration.

Additional information on each of these selections is contained in the following sections.

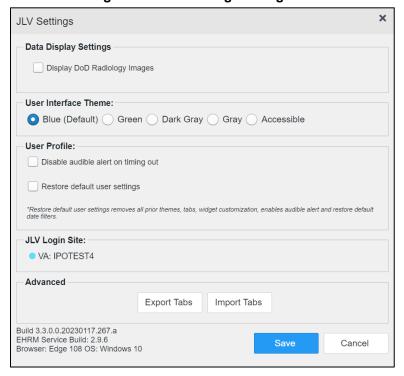


Figure 19: JLV Settings Dialog Box

NOTE: Users logging into JLV with a PIV card will have the option of changing their site location. Users should only select sites linked to their PIV card. The new selection will not take effect until the next login and will remain the user's default site until updated or reset.

You may reset your user profile to JLV defaults by selecting the **Restore Default User Settings** checkbox and then the **Save** button in the **JLV Settings** dialog. This feature resets the theme and data display to default; enables the audible alert; removes any custom tabs; and resets all widgets to their default configurations, which removes any custom sort or filter settings.

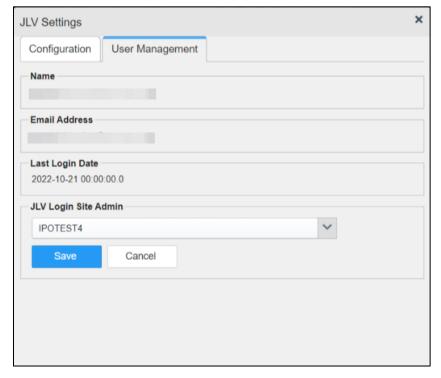


Figure 20: VHA JLV User Management Settings Tab

JLV users that have privileges to edit login sites for end users will have the additional User Management tab. Users will be prompted to enter an email address to search for the associated current entry in the JLV user table. If there is no associated email address, a popup will appear, "User not found".


If the search locates a current entry in the JLV user table, the user's email, username, and last login date will appear in a read only mode as in <u>Figure 20</u>. A login site will appear prepopulated with the selected end user's login site from the user table and allow selection of any site in the JLV site drop-down list. Clicking the save button will update the login site for the user.

3.2.1. Cache Manager

Administrative users will be able to manage cache by selecting the Show Cache Manager checkbox, found under the Configuration Tab on the JLV Settings dialogue. Once selected, the Cache Manager button will be available on the top right JLV toolbar.

The Cache Manager button opens a Cache Manager window displaying the current state of the cache, with the option to manually expire a cache chunk, as in <u>Figure 21</u>.

Figure 21: Admin Cache Manager

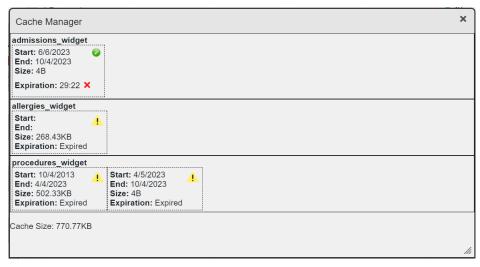

3.2.2. About Data Display Settings

The data display setting is part of your user profile. The data display setting enables you to specify whether to include DOD radiology images. This option is disabled by default.

JLV displays data source indicators in the **Site** columns of individual widgets (<u>Figure 22</u>). A blue circle indicates VA data; an orange square indicates DOD data; a purple hexagon indicates community partner data (if enabled); and a green triangle indicates Cerner Millennium Federal Electronic Health Record (FEHR) data (<u>Table 3</u>). Hovering over individual data source indicators provides more specific information about the data source.

Table 3: Data Sources

Icon	Sources	Hover information	Data Supplied
XXX XXX = the three-letter VistA site abbreviation	CPRS/VistA system data	VA: location Examples: VA: Walla Walla, WA VA: Central Texas HCS ¹	Supplies all records from CPRS/VistA for all sites. These records will remain available even after a site transitions to VA EHRM (Cerner).
DoD	CHCS ² (CDR ³), AHLTA ⁴ , Essentris, TMDS ⁵ , various clinical applications	DOD source system: Full name of DOD Medical Treatment Facility Examples: CDR: Fairchild AFB ⁶ AHLTA: JB ⁷ San Antonio	Supplies all records from all sites from all DOD facilities for all legacy applications. These records will remain available even after a site transitions to MHS GENESIS (Cerner).

¹ Health Care System

² Composite Health Care System

³ Clinical Data Repository

⁴ Armed Forces Health Longitudinal Technology Application

⁵ Theater Medical Data Store

⁶ Air Force Base

⁷ Joint Base

Icon	Sources	Hover information	Data Supplied
FEHR	DOD MHS GENESIS (Cerner) and VA EHRM (Cerner) records stored in Cerner Millennium	FEHR: VA or DOD facility name/location Examples: FEHR: Spokane, WA FEHR: Madigan AMC ⁸	Supplies records from VA and DOD facilities that have migrated to the new Federal EHR, called MHS GENESIS in DOD and VA EHRM (Cerner) in VA.
Name Name = Community Partner name	All DOD and VA Community Partners (VHIE/VLER ⁹), Sequoia Project and DOD Health Information Exchange)	Federal HIE ¹⁰ : Community Partner Organization Name Examples: • Federal HIE: UC ¹¹ San Diego • Federal HIE: Mayo Clinic	Supplies continuity of care documents, care summaries, and notes when provided by partners. Partner capabilities and responses vary.

3.2.3. About UI Themes

The UI theme setting is also part of your user profile. Choosing a UI theme sets the font, foreground, and background colors of the widgets, toolbars, and dialog boxes. Once selected, the UI theme is stored in your user profile and remains set until you change it. The available UI themes are Blue (Default), Green, Dark Gray, Gray, and Accessible (Section 508-compliant).

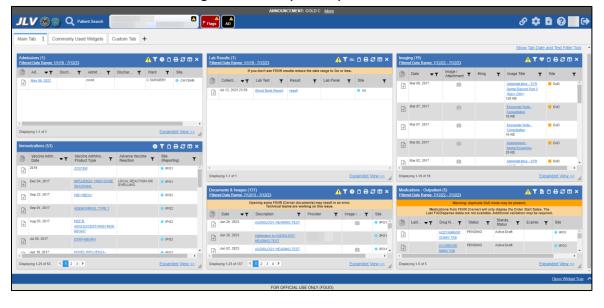


Figure 22: Blue (Default) UI Theme¹²

⁸ Army Medical Center

⁹ Virtual Lifetime Electronic Record

¹⁰ Health Information Exchange

¹¹ University of California

¹² Primary Care Provider (PCP), Patient-Aligned Care Team (PACT), Social Security Number (SSN), Date of Birth (DOB)

Figure 23: Green UI Theme

Figure 24: Dark Gray UI Theme



Figure 25: Gray UI Theme

Figure 26: Accessible (Section 508-Compliant) UI Theme



3.2.4. Using the Accessible Theme's Functionality

The Accessible theme includes larger fonts, enhanced tab and keystroke functionality, and screen reader compatibility in compliance with Federal accessibility requirements and Section 508 standards.

NOTE: The Accessible theme is available for users who have adaptive technologies installed on their computers.

3.2.4.1. Keyboard Focus and Navigation

Keyboard focus is the highlighting of an element, field, or control that enables interaction with, and navigation of, JLV using the keyboard and keystrokes. When navigating interactable elements with the Tab key, the item in focus is outlined in gold as a visual indicator.

JLV uses standard accessibility software key commands for navigation and interaction.

Pressing **Enter** or the **Spacebar** when an element that provides action is in keyboard focus performs the associated action; for example, pressing **Enter** or the **Spacebar** when the is in focus opens the **JLV Settings** dialog box. Using the arrow keys or the **Tab** key moves between keyboard focus items to navigate through screen elements.

The following tables provide a complete list of accessible keystrokes.

Table 4: Accessible Theme Keyboard Shortcuts - Patient Data& Widget Navigation

Keystroke	Description
ALT + W	Press to enter widget body.
ENTER	Press to initiate the action associated with the item in focus.
TAB	Press to transfer keyboard focus to other user interface items.
ESC	Press to return keyboard focus to the panel containing the user interface item with keyboard focus or to exit a window, widget, or tab, and close dropdowns.
SPACEBAR	Press to activate any user interface item (for example, click a button).
Arrow Keys	When focused on a dropdown list, press the down arrow key to view list contents.

Table 5: Accessible Theme Keyboard Shortcuts - Tabs

Keystroke	Description
Arrow Keys	Use the left and right arrow keys to change tab panels.
ТАВ	 Press Tab until + is in focus Press Enter The new tab dialog box opens and prompts you to enter a name for the new tab Tab to the Add button to confirm the new tab name or tab to the Cancel button to discard the new tab
ENTER	Press to initiate data load of tab.
1	On custom tabs, press 1 to rename the tab.
DELETE	On custom tabs, press DELETE to remove the tab.

Table 6: Accessible Theme Keyboard Shortcuts – Windows or Dialog Boxes

Keystroke	Description
ТАВ	Press to transfer keyboard focus to other user interface items within the window or dialog box.

Table 7: Accessible Theme Keyboard Shortcuts – Adding Widgets to the Patient Data and Widgets Page

Keystroke	Description
=	Press to load a dialogue box with options to place a widget.

3.2.4.2. Accessible Theme Navigation

An element with keyboard focus is the starting point for navigation. Pressing the **ALT+W** keys when a widget is in focus, will take you to the first interactable element in the widget.


Figure 27 depicts the **Documents & Images** widget with keyboard focus. Focus can be changed by pressing the **Tab** key. Pressing the **Tab** or Shift Tab key, in this example, shifts the keyboard focus from one interactable element to the next. Press the escape key to return focus to the widget body.



Figure 27: Navigation Using Elements with Keyboard Focus

The **Document** details window has keyboard focus in <u>Figure 28</u>. Move through the toolbar buttons by using the arrow keys, or by pressing **Tab**. Move to **X** and press **Enter** to close the details window. After closing the window, keyboard focus returns to the last element that had keyboard focus prior to opening the window.

Figure 28: Document Details with Keyboard Focus

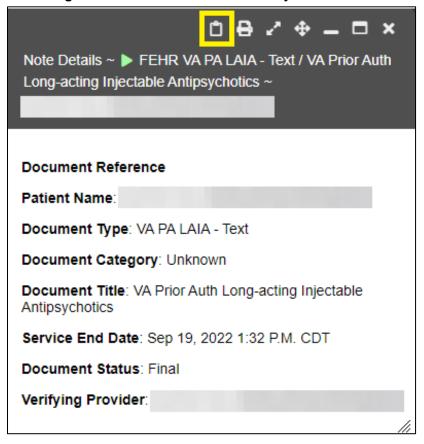

3.2.4.3. Table Functionality

Table Column Headers (<u>Figure 29</u>) can be sorted by selecting the column header and pressing the space bar or enter key to sort the column. Upon a sort, (<u>Figure 30</u>) a user will have to leave the table and reenter for JAWS to receive the updated Table Data.


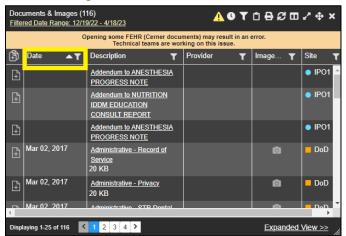


Figure 29:Table Column Header

Figure 30: Widget Column Select Dropdown

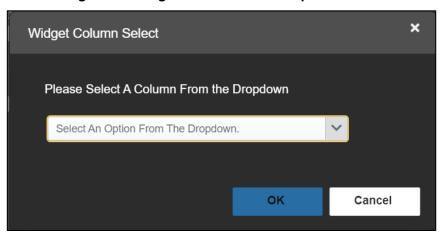

NOTE: Not all table columns are sortable. EZ Select All Column Header and Report Builder Selected Items Column Header are not sortable.


3.2.4.4. Widget Move

In the widget tools menu, pressing the four-arrow icon or move widget button (Figure 31) will open a dialogue with options to move the widget on the screen or inside the Patient Data and Widgets page (Figure 32). This function is available in both the minimized and expanded view.

Documents & Images (116) <u>^ O T Ů ⊖ ♂ □ /</u> Filtered Date Range: 12/19/22 - 4/18/23 Opening some FEHR (Cerner documents) may result in an error. Technical teams are working on this issue (A) Date Description Provider Site Image... Apr 12, 2023 RADRPT / US Abdomen ▶ FEHR Limited **Details** Apr 12, 2023 RADRPT / US **▶** FEHR Encephalogram **Details** RADRPT / CT Abdomen and Apr 12, 2023 ▶ FEHR Pelvis w/ Contrast Details Apr 10, 2023 Unspecified abdominal pain DoD 09:49 Details Apr 05, 2023 IPO1 Addendum to NUTRITION **HYPERLIPIDEMIA EDUCATION CONSULT** < 1 2 3 4 > Displaying 1-25 of 116 Expanded View >>

Figure 31: Widget Move Icon

Figure 32: Widget Move Dropdown

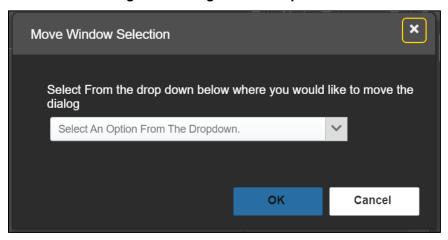


3.2.4.5. Widget Resize

In the widget tools menu, pressing the two-arrow icon or resize widget button (Figure 33) will open a dialogue with options to resize the widget on the screen or inside the Patient Data and Widgets page (Figure 34). This function is available in both the minimized and expanded view.

Documents & Images (116) Filtered Date Range: 12/19/22 - 4/18/23 Opening some FEHR (Cerner documents) may result in an error. Technical teams are working on this issue. Date Provider Description Image.. RADRPT / US Abdomen Apr 12, 2023 ▶ FEHR Limited RADRPT / US Apr 12, 2023 **▶** FEHR Encephalogram **Details** RADRPT / CT Abdomen and Apr 12, 2023 **▶** FEHR Pelvis w/ Contrast **Details** Apr 10, 2023 Unspecified abdominal pain DoD 09:49 Details Apr 05, 2023 IPO1 Addendum to NUTRITION **HYPERLIPIDEMIA EDUCATION CONSULT** < 1 2 3 4 > Displaying 1-25 of 116 Expanded View >>

Figure 33: Widget Resize Icon

Figure 34: Widget Resize Dropdown

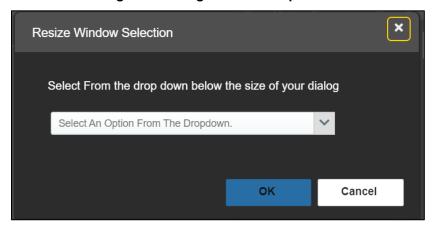

3.3. User-Restricted Access to Patient Data

3.3.1. VHA Users Viewing DOD-only Patient Data

VHA users can view the records of DOD-only patients, meaning they are not registered in MPI, but VA requires that these actions be audited. If an attempt is made to access DOD-only patient records, the VHA user is asked to specify the purpose for access. The purpose options presented to the user are: Emergent Care, Clinical Care, or Authorized Administrative Use (Figure 35).

Once the purpose for access is selected, clicking **Agree to be Audited & Access this Patient Record** (Figure 35), displays the requested data, and audits the action. The access purpose, organization of the VA provider, date, username, IP address, user's Internal Entry Number (IEN), host system's ID, and patient's Electronic Data Interchange Personal Identifier (EDIPI) are captured in the audit record.

This patient is NOT a registered VA Patient.
Select an authorized purpose for accessing this DoD record.

Emergent Care
Clinical Care
Authorized Administrative Use
i.e. claim processing

Your access to this DoD patient record will be recorded and audited by the Department of Defense. Unauthorized access may result in disciplinary actions to include termination of employment.

Organization: VHA
Username:
Date/Time: 2019-07-15 14:11:05
IP Address:

Agree to be Audited & Access this Patient Record

Figure 35: VA User Restricted Access Dialog

3.3.2. CAPRI-Claims Users Viewing Patients with DOD Data

CAPRI-Claims users may see one of two different dialog boxes, depending on whether the patient is registered for care in the VA. If the patient is registered for care in the VA, they are in MPI.

3.3.2.1. Patient Registered in MPI

CAPRI-Claims users are permitted to view the records of a patient who is registered in MPI, but the VA requires that these actions be audited. After performing a patient search and selecting a patient from the search results, a CAPRI-Claims user sees the audit notification (<u>Figure 36</u>) if the patient they selected is registered in MPI. After agreeing to the audit, the CAPRI-Claims user can access the patient's record.


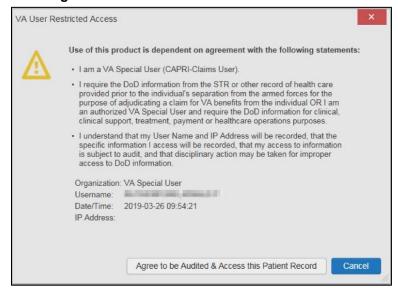


Figure 36: CAPRI-Claims User Audit Notification

3.3.2.2. Patient Not Registered in MPI

CAPRI-Claims users are not permitted to access the records of a patient who is not registered in MPI. After performing a patient search and selecting a patient from the search results, a CAPRI-Claims user sees the restricted access notification (Figure 37) if the selected patient is not registered in MPI. Clicking **Cancel** returns the user to the previous screen.

Figure 37: Restricted Access Notification

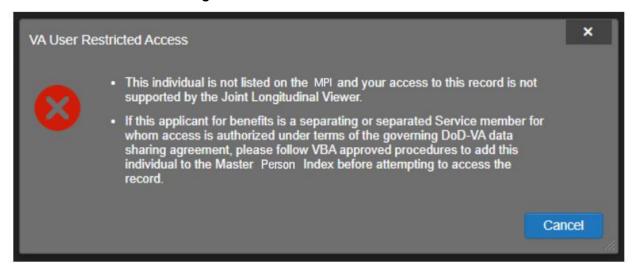

3.3.3. Access to DOD Sensitive Records

All JLV users are audited each time a sensitive DOD record is accessed. Auditing applies to sensitive documents, lab results, outpatient encounters, or progress notes records. JLV records the user's name, organization, SSN/EDIPI (for DOD users), PIV (for VA users), location, patient identifiers (patient last name, first name, middle initial [MI], SSN/EDIPI [DOD only], MPI [VA only], and date of birth [DOB]), data accessed, date and time, and reason for access for each attempt to access sensitive data.

JLV notifies you before you are audited. The message is triggered the first time you click a masked record that is marked ** Sensitive**. Your acknowledgement of the audit notification covers all subsequent DOD sensitive records in any widget during the remainder of the session, and it only expires when you switch patients or log off. An example of a masked record audit notification is depicted in Figure 38.

After you select the purpose for viewing the record and agree to be audited, the sensitive data displays in the **Details** view of the widget.

This record is marked sensitive and should only be viewed if necessary, as per the HIPAA principle of minimum use.

Select an authorized purpose for accessing this record.

Emergent Care
Clinical Care
Authorized Administrative Use
i.e. claim processing

You have requested to view sensitive data. This usage is subject to audit and will apply to all other sensitive data viewed for this patient. Do you wish to proceed?

Organization:
Username:
Date/Time: 2021-07-27 10:42:43
IP Address:

Agree to be Audited & Access this Sensitive Record

Cancel

Figure 38: DOD Sensitive Record and Audit Dialog Box

3.4. Patient Searches

The core function of JLV is to display patient information on the Patient Data and Widgets page. Use the **Q Patient Search** feature to perform a patient search, perform a family search, or use the advanced search options.

JLV determines which enterprise service to use for a patient search, based upon the search criteria:

- If the EDIPI or Sponsor SSN is entered, JLV utilizes PDWS for the search (See <u>PDWS</u> Search Rules and Supported Searches for details).
- If the EDIPI or Sponsor SSN fields are empty, JLV utilizes MPI for the search (See MPI Search Guidelines for details).

3.4.1. PDWS Search Rules and Supported Searches

When JLV utilizes PDWS for patient search, the applicable rules and rule sets required for the PDWS interface must be used:

- **Rule Set 1:** When searching by SSN, the full nine-digit SSN must be entered, and the Last Name or DOB must also be supplied.
 - This is true for search by patient SSN as well as Sponsor SSN.
- **Rule Set 2**: If you do not have the patient's SSN, you must supply ALL the following identifiers: Last Name, First Name, DOB, and Birth Sex.
 - Blanks in any of these four required fields generates an error
- A 10-digit string is required for telephone number.

The following search combinations are supported when PDWS is utilized:

- DOD ID, also referred to as the EDIPI
- Patient SSN and Last Name
- Patient SSN and DOB
- Sponsor SSN and Last Name
- Sponsor SSN and DOB
- Last Name, First Name, DOB, and Birth Sex (ALL identifiers must be entered)

3.4.2. MPI Search Guidelines

JLV utilizes MPI for a patient search when the search parameters include the patient's name and SSN. MPI requires the full last name is entered when the patient's SSN is entered. MPI may also require entering the patient's first name and DOB to increase the efficiency of the search, and the likelihood of locating the patient. An example would be searching for a patient with a common last name, such as Smith or Jones.

3.4.3. Patient Search Using DOD ID

- 1. Click Q Patient Search.
 - a. The **Patient Search** dialog opens.
- 2. Enter the 10-digit patient DOD ID (EDIPI) in the **DOD ID** field (<u>Figure 39</u>).

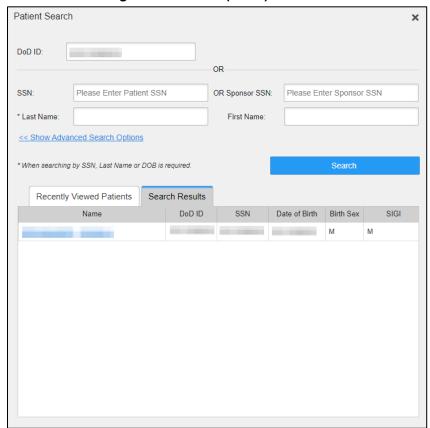


Figure 39: DOD ID (EDIPI) Search¹³

3. Click Search.

- a. The search results display on the **Search Results** tab.
- 4. Click a patient's name in the search results list to open the associated record.
 - a. The **Name**, **SSN**, **DOB**, **Birth Sex**, and **SIGI** columns in the results list provide information to assist with patient identification.
- 5. After selecting a patient, the Patient Data and Widgets page opens and displays the selected patient's data.

3.4.3.1. Invalid and Maximum Characters

Text Fields on Patient Search will alert the user if an invalid key is entered into the text box (<u>Figure 40</u>). Text fields will inform the user when they have reached the maximum character length for entry (<u>Figure 41</u>).

¹³ Self-Identified Gender Identity (SIGI)

× Patient Search DoD ID: Please Enter Patient DoD ID Invalid Character Entered OR-SSN: Please Enter Patient SSN OR Sponsor SSN: Please Enter Sponsor SSN * Last Name: First Name: << Show Advanced Search Options Search * When searching by SSN, Last Name or DOB is required. Recently Viewed Patients Search Results Name SSN DoD ID Date of Birth Birth Sex SIGI No Results

Figure 40: Invalid Character Entry

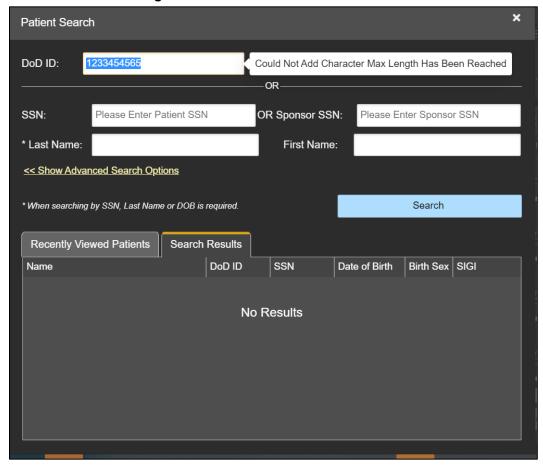


Figure 41: Maximum Characters Reached

3.4.4. Patient Search Using SSN

- 1. Click **Q Patient Search**.
 - a. The **Patient Search** dialog box opens.
- 6. Enter the patient's full, nine-digit SSN in the SSN field (Figure 42).
 - a. When searching by SSN, either the last name or the DOB *must* be entered to perform the search.

Patient Search × Please Enter Patient DoD ID DoD ID: OR OR Sponsor SSN: Please Enter Sponsor SSN * Last Name First Name << Show Advanced Search Options * When searching by SSN, Last Name or DOB is required. Recently Viewed Patients Search Results DoD ID Name SSN Date of Birth Birth Sex SIGI

Figure 42: SSN Search

7. Click Search

- a. If the required patient identifiers are not provided, hover text appears, indicating the information necessary to complete the patient search.
- b. If a match of patient names cannot be made, enter additional identifiers (first name, DOB) to refine the search (Figure 43).

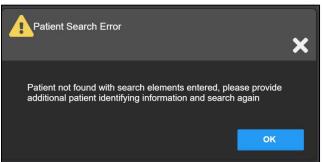


Figure 43: Additional Identifiers Needed

8. The search results display on the **Search Results** tab.

- 9. Click a patient's name in the search results list to open the associated record.
 - a. The **Name**, **SSN**, **DOB**, **Birth Sex**, and **SIGI** columns in the results list provide information to assist with patient identification.
- 10. After selecting a patient, the Patient Data and Widgets page opens and displays the selected patient's data.

3.4.5. Advanced Patient Search Options

Advanced patient search options in the **Patient Search** dialog box allow for the use of additional patient identifiers when performing a patient search.

- 1. Click Q Patient Search.
 - a. The **Patient Search** dialog box opens.
- 2. Click the << Show Advanced Search Options link under the *Last Name field (seen in Figure 42).
- 3. Enter the patient identifiers in the appropriate search fields, as desired.
 - a. Advanced search option fields include *DOB, MI, Birth Sex, Address, State, City, Zip, and Telephone (seen in Figure 44).

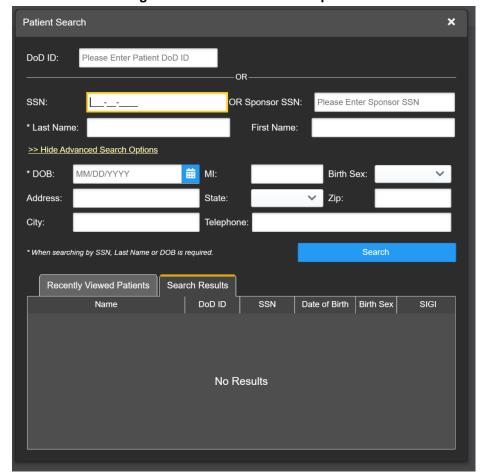


Figure 44: Advanced Search Options

4. Click Search

- a. If the required patient identifiers are not provided, hover text appears, indicating the information necessary to complete the patient search.
- b. The search results display on the **Search Results** tab.
- 5. Click a patient's name in the search results list to open the associated record.
 - a. The **Name**, **SSN**, **DOB**, **Birth Sex**, and **SIGI** columns in the results list provide information to assist with patient identification.
- 6. After selecting a patient, the Patient Data and Widgets page opens and displays the selected patient's data.

3.4.6. Family Member Search

A family member search can be performed using the **Sponsor SSN** field in the **Patient Search** dialog box.

- 1. Click **Q Patient Search**.
 - a. The **Patient Search** dialog box opens.
- 2. Enter the SSN in the **Sponsor SSN** field.
 - a. The **Last Name** or **DOB** of the patient must also be entered in the appropriate fields (Figure 45).
- 3. Click **Search**.
 - a. If the required patient identifiers are not provided, hover text appears indicating the information necessary to complete the patient search.
 - b. The search results display on the **Search Results** tab (Figure 45).

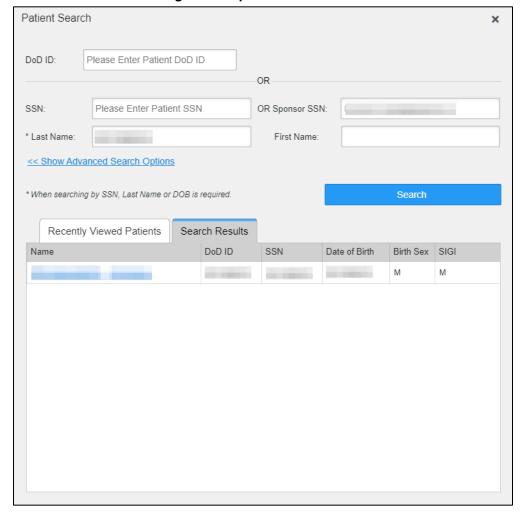


Figure 45: Sponsor SSN Search

- 4. Click a name in the **Search Results** list.
 - a. The **Family Members** dialog opens (<u>Figure 46</u>).
 - b. The sponsor's identifiers and a list of dependents display in the dialog.
- 5. Click a name in the list of family members to open the associated record.
- 6. After selecting a patient, the Patient Data and Widgets page opens and displays the selected patient's data.

Family Members

Sponsor:

Name

SSN

Date of Birth

F

Name

DoD ID

SSN

Date of Birth

Birth Sex

SIGI

M

M

M

Cancel

Figure 46: Family Members Dialog Box

3.4.7. Recently Viewed Patients List

Users who have previously logged into JLV, searched for a patient, and viewed that patient's records can see and access a list of recently viewed patients. A recent patient is defined as a patient whose record has been viewed (opened) by the user.

- NOTE: The Recently Viewed Patients list does not include search history or recent search results, only a list of patients whose data was accessed and viewed. The list is limited to 10 patients. You can clear this list by selecting the Clear Recently Viewed Patients button.
- 1. Click Q Patient Search.
 - a. The **Patient Search** dialog box opens.
- 2. Click the **Recently Viewed Patients** tab in the **Patient Search** dialog box.
 - a. A list of recently viewed patients displays (Figure 47).
 - b. The **Name**, **DOD ID**, **SSN**, **DOB**, **Birth Sex**, and **SIGI** columns in the results list provide information to assist with patient identification.
- 3. Click a patient's name in the **Recently Viewed Patients** list to open the associated record.
- 4. After selecting a patient, the Patient Data and Widgets page opens and displays the selected patient's data.

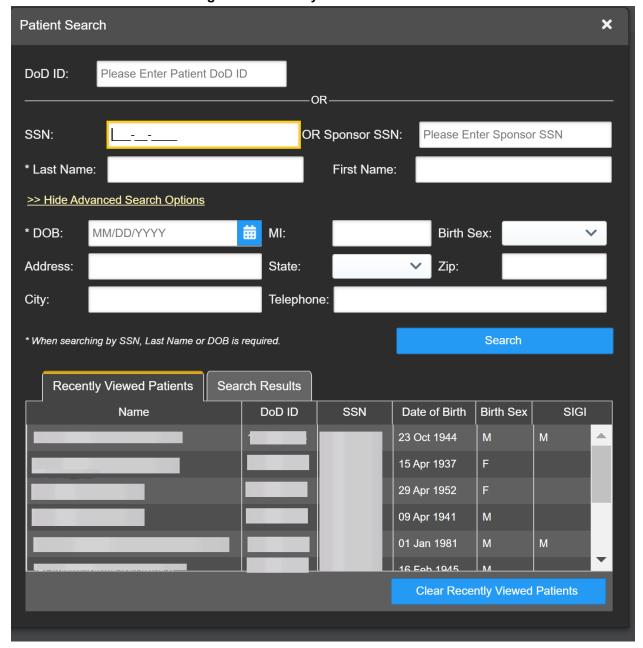


Figure 47: Recently Viewed Patients List

4. Widgets

Widgets are elements on the JLV Patient Data and Widgets page that display data specific to a clinical domain. By default, widgets are displayed in minimized view on the Patient Data and Widgets page (Figure 22), but they can be expanded to view additional details. Available widgets are listed in and launched from the widget tray.

O

NOTE: Widgets may take anywhere between 10 seconds to 90 seconds to load, depending on the date range settings.

4.1. Accessing and Opening Widgets

All available widgets display in the widget tray (<u>Figure 48</u>). Hover over a widget icon to view a description of the widget. Widgets can be opened, rearranged, and closed. JLV can display up to 12 widgets per tab.

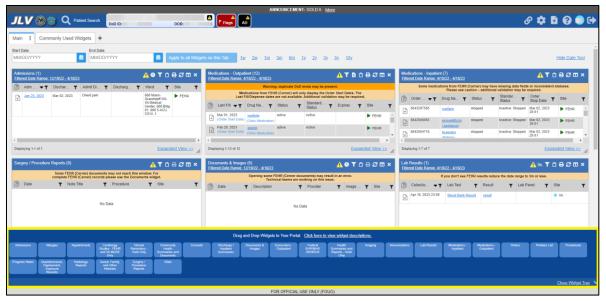


Figure 48: Widget Tray

- 1. Open the widget tray by clicking the **Open Widget Tray** link. You can also click the **:** symbol located on the selected tab, then select the **Open Widget Tray** link.
- 2. Use the scroll bar arrows <> at either end of the tray to see all widget choices, or click the **More Widgets** >> link to scroll through the list of widgets.
- 3. Click and hold a widget icon in the tray, drag it to the Patient Data and Widgets page, and drop it in the desired location.
 - a. The widget is docked on the page and opens in minimized view.
- 4. Close the widget tray by clicking the **Close Widget Tray** link.
 - NOTE: JLV displays a **Duplicate Widget** notice if you attempt to add another instance of a widget to a tab.

Widgets can be closed (removed from the Patient Data and Widgets page) by clicking **X**. Rearrange the widgets displayed on the page by dragging and dropping them.

4.2. Widget Navigation and Display Options

Each widget has tools and display options available to navigate through, and change, the display of data. The vertical scroll bar allows you to move through, and view, the widget's data. Navigation icons and actions are detailed in <u>Table 8</u>.

Table 8: Navigation Icons and Actions

Navigation Icon	Action	Description
<	Changes the focus of the widget to the previous group or page of records within the results display.	Go to Previous Page
1	Changes the focus of the widget to the page number selected.	Jump to Page
>	Changes the focus of the widget to the next group or page of records within the results display.	Go to Next Page
Expanded View >>	Available in minimized views only. Opens the expanded view of the widget in a secondary window.	Go to Expanded View
1-25 of 55	Indicates the number of records displayed in that widget page out of the total number of results for that widget.	Record Display Indicator
Show All/ Show Paged	Click Show All to open all records for a given widget in a scrollable window. Click Show Paged to return to the display of records grouped by pages.	Display Setting

When there are more than 25 records available in a widget, they are grouped in **Show Paged** view. Records 1 through 25 can be viewed by using the widget's vertical scroll bar. Records 26+ can be viewed by using or to view to subsequent pages. The expanded view of a widget contains a **Show All** link, which opens all records for that widget, in a scrollable window. Click **Show Paged** to return to the display of records grouped by pages.

4.3. Widget Toolbars

There is a toolbar on both the minimized and expanded views of each widget and most dialog boxes. Toolbar buttons vary by widget, dialog box, and window. <u>Table 9</u> describes the functionality of each toolbar button.

Table 9: Widget Toolbar Icons

Icon	Name	Function
+	Add to Report Builder	Adds the information displayed in the widget to the Report Builder, including any Details or Notes .
0	Clinical Reminders	Opens the Clinical Reminders widget. NOTE: This icon only appears on the Admissions, Appointments, Encounters-Outpatient, and the Documents & Images widgets.
×	Close	Removes the widget from the Patient Data and Widgets screen or closes a dialog box.

Icon	Name	Function	
	Column Settings	Configures the columns within the widget. Turn the columns on or off by checking the column names that appear in the pop-up box and clicking Apply. NOTE: Column widths and window sizes are auto saved to user profiles for both the minimized and expanded views.	
7	Configure Filter	Filters on specific record types or other elements within the widget, including a date filter. If a date range filter is applied in an open widget, the date range is displayed in the widget header. NOTE: After setting a filter, the Close Filter option is enabled in the widget. Clicking Close Filter restores the full display of records within the widget.	
	Copy to Clipboard	Copies the selected content of the open window to the clipboard for pasting into another application. Copy to Clipboard is disabled in the widget toolbar after clicking Show All within a widget. The clipboard is cleared after closing a patient record, selecting a new patient record, and logging off JLV. Once a user chooses the column selections for Copy to Clipboard, the new selections for that widget will be retained even when switching to a new patient, unless restore default settings is executed.	
1	Connection Status	Both icons provide a status indicator for DOD and VA data sources. indicates all sources are available. indicates one or more data sources are unavailable. Clicking either status icon opens the connection status details in a separate window.	
÷	Print	Prints the selected contents of the open window. Print is disabled in the widget toolbar after clicking Show All within a widget. Once a user chooses the column selections for Print, the new selections for that widget will be retained even when switching to a new patient, unless restore default settings is executed. NOTE: It is recommended that you have the latest Adobe Reader installed for the Report Builder and other printing features.	
	Graph	Opens a dialog to configure the display of multiple data of the same type in a graph.	
C	Refresh	Refreshes the widget or window display. Only the data within that widget's dataset is updated.	
ČĎ	Copy Active Rx to Clipboard	Converts the medication data within the widget to text and copies the text to the user's clipboard. This option is available from the Medications - Outpatient widget. NOTE: Copy Active Rx is disabled when a date range filter is set to less than 120 days.	
+	Widget Move	Opens a dialogue box with options to move the widget on the screen or inside the portal. NOTE: This is only available in the accessible theme.	
~	Widget Resize	Opens a dialogue box with options to resize the widget on the screen or inside the portal. NOTE: This is only available in the accessible theme.	

4.4. Minimized vs. Expanded Widget View

Each widget can be displayed in either a minimized (default) or an expanded view. Minimized view displays a simple list of the available patient data for a clinical domain, arranged in reverse

chronological order by default. Expanded view provides a detailed list of similarly arranged patient data. When the expanded view of a widget is launched, it opens in a new, separate window that displays more attributes of the records in the widget, including additional, sortable columns of data.



NOTE: The Cardiology Studies - FEHR and VA MUSE Only, **Community Health Summaries and Documents**, and **Federal EHR/MHS GENESIS** widgets do not have an expanded view. Instead, a record opened from the minimized view of the widget is displayed in a new, separate window.

Click the **Expanded View >>** link in minimized view to launch expanded view and see the additional display and functionality options.

4.5. Sort and Filter Options

Widgets have a variety of methods to sort and filter data.

The Filter Indicator (+Filter(s)) in the widget header shows that column or text filters have been applied.

In general, you can:

- Click a column title to sort records according to data in that column.
 - If you sort a column that appears in both minimized and expanded views of the widget, the sort saves in your user profile and persists. Data views persist from minimized to expanded views and expanded to minimized views.
 - If you sort a column that is only in the expanded view of the widget, your next session restores the default widget sort, or the last sort saved to your user profile.
 - If you sort a **Date** column that has some partial or missing date values, understand:
 - Dates that have month and year, but lack a day are treated as though the day were "00" (e.g., July 2018 sorts between June 30, 2018 and July 01, 2018).
 - Dates that have a year only are treated as though the month is January and the day is "00" (e.g., 2018 sorts between December 31, 2017 and January 01, 2018.
- Click on the widget toolbar to show or hide columns within the widget.
 - Check/uncheck the column names that appear in the pop-up dialog box, then click **Close**.
- If a widget's column data contains blanks, the column Filter popup will move all blanks to the bottom of the list. If there are no blanks in the widget, the Filter will not be available as a selection in the column.

Sort and filter capabilities are specific to the data presented in each widget and may vary between the minimized and expanded views of an individual widget or between widgets. <u>Table 10</u> lists available filters and their corresponding widgets.

Table 10: Available Widget Filters

	l able 10: Available Widget Filters		
Filter By	Widgets		
Site	Admissions, Allergies, Appointments, Consults, Documents & Images, Health Summaries and Reports, Medications - Inpatient, Discharge/Inpatient Summaries, Lab Results, Orders, Encounters - Outpatient, Medications - Outpatient, Problem List, Procedures, Progress Notes, Questionnaires, Deployment, Exposure Records, Radiology Reports, Social, Family, and Other Histories, Surgery/Procedure Reports, Vitals		
Site (Reporting)	Immunizations		
Provider	Discharge/Inpatient Summaries, Progress Notes, Surgery/Procedure Reports		
Provider Specialty	Appointments, Consults, Documents & Images, Discharge/Inpatient Summaries, Orders, Encounters - Outpatient, Procedures, Progress Notes		
Clinic	Appointments, Encounters - Outpatient		
VA Health Factors	Admissions, Appointments, Documents & Images, Encounters - Outpatient		
Consult Order	Consults		
Description	Documents & Images		
Problem Description	Problem List		
Standardized Description	Documents & Images, Problem List		
Location	Admissions, Allergies, Appointments, Assessments, Consults, Documents & Images, Discharge/Inpatient Summaries, Encounters – Outpatient, Federal HER MHS GENESIS, Imaging, Lab Results, Medications - Inpatient, Medications – Outpatient, Orders, Problem List, Procedures, Progress Notes, Questionnaires, Deployment, Exposure Records, Radiology Reports, Social, Family, and Other Histories, Surgery/Procedure Reports, Vitals		
Source System	Documents & Images		
Vaccine Administered Product Type	Immunizations		
Standardized Vaccine Product Type	Immunizations		
Drug Name	Medications - Inpatient, Medications - Outpatient		
Standardized Drug Name	Medications -Inpatient, Medications -Outpatient		
Ordering Health Care Provider (HCP) Specialty	Medications - Inpatient, Medications - Outpatient		
Note Type	Discharge/Inpatient Summaries		
Lab Test	Lab Results		
Туре	Lab Results; Orders; Problem List ¹⁴ ; Procedures; Social, Family, and Other Histories; Vitals		
Standardized Type	Orders, Vitals		
Status	Problem List		
Document Type/Title	Progress Notes		

-

¹⁴ Filter applies to DOD data only

Filter By	Widgets
Standardized Document Type	Progress Notes
Exam	Radiology Reports
Standardized Radiology Exam	Radiology Reports
Standardized Finding	Social, Family, and Other Histories
Note Title	Surgery/Procedure Reports
Procedure	Surgery/Procedure Reports

Filter dropdowns dynamically expand and sometimes have a horizontal scroll bar to display long lines of data by which records can be filtered (Figure 49).

Documents & Images (54) State Code Training. 17(1927-41927)

Stat date

L2719-2022 Out 1920-2023 Out 1920-2

Figure 49: Dynamically Expanded Filter Dropdown in the Documents & Images Widget

Some widgets also have the capability to configure data to display multiple records of the same type in a graph or table view.

4.5.1. Widget Date Range Filters

Patient data displayed within a widget can be filtered by date range. If a date range filter is applied, the selected range is indicated in the widget header. Click $\widehat{\mathbf{Y}}$ or the **Filtered Date Range** link on the widget toolbar to change the date range of the data displayed (Figure 50).

Documents & Images (919) Filtered Date Range: 5/30/22 - 9/27/22 If the information you expect does not appear, expand the date range End date: Start date: **Apply** 05/30/2022 09/27/2022 1w 2w 1m 3m 6m 1y 2y 3y 5y 10y Remove All Text Q Enter text to filter Remove All Filters **Hide Filters** Opening some FEHR (Cerner documents) may result in an error. Technical teams are working on this issue. Description Provider 2 3 4 > Displaying 1-25 of 919 Expanded View >>

Figure 50: Date Range Filters Applied

The two ways to filter the date range of data displayed in a widget are the **Quick Date Range** filter and the **Start and End Date** filter. If the expected information does not appear in the widget display, use the date filter options to change the date range.

- **NOTE:** Changes made to the date range in a widget persist from minimized to expanded view, expanded to minimized view, patient to patient, and session to session, until the widget is either removed or the date range is changed.
- **NOTE:** When looking for Essentris information associated with a specific admission in the **Discharge/Inpatient Summaries** widget, expand your search date ranges to include additional days before and after the admission. This includes preadmission notes as well as notes signed after the patient was discharged in the filtered widget display.

4.5.1.1. Quick Date Range Filter

Use the preset, **Quick Date Range** filter to set or limit the display of patient records within a widget. This option is available by clicking the **Filtered Date Range** link or by clicking **T** in a widget on the Patient Data and Widgets page.

Clicking a **Quick Date Range** filter refreshes the data displayed to show only the records for the selected range. Selecting **Cancel**, which is the **Apply** button while the data refresh is running, stops the query. The options represent time counting back from the present day (for example, selecting 2y displays only records within the last 2 years). Examples of preset date ranges include 1w (last 7 days), 6m (last 6 months), and 3y (last 3 years).

NOTE: The Quick Date Range filter is contextual. If today is February 2, 2021, the preset filter counts backwards using that date. If today is April 16, 2021, the preset filter counts backwards using that date.

Figure 51 highlights the **Quick Date Range** filter options in the expanded view of the **Documents & Images** widget.

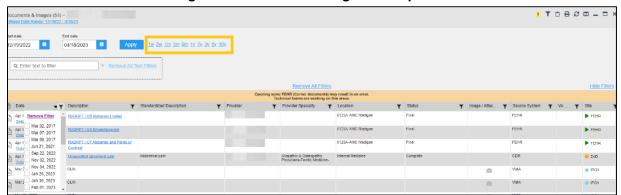


Figure 51: Quick Date Range Filter Options

4.5.1.2. Start and End Date Filter

Use the **Start Date and End Date** filter to display data for a specific time range. Click the **Filtered Date Range** link or **Y**.

NOTE: The Start and End Date filter is static. If the start date is set to June 2, 2016 and the end date is set to June 2, 2017, the data in that date range displays, no matter what today's date is.

Figure 52 provides an example of the start and end date filter fields in the expanded view of the **Admissions** widget.

Start Date: End Date: 04/04/2007 04/12/202 Apply to all Widgets on this Tab Medications - Ou 2007 April A O Y O B S O × Filtered Date Ran Wed Thur nit Di Dischard. 2 3 7 668 Mann-▶ FEHR Grandstaff W∆ VA Medical 10 11 12 14 13 Last Fill Center, 668 Bldg 01, 668 5 ACU, 21 15 18 19 16 20 C514. 1 Mar 01, 2023 (Order Start I 25 27 28 22 26 Feb 28, 2023 5 29 30 10 11 12 Expanded View >> Displaying 1-12 of

Figure 52: Start and End Date Display Filter

- 1. Click mext to Start Date.
- 2. Select a month, day, and year start date for the display filter.
- 3. Click mext to **End Date**.
- 4. Select the month, day, and year end date for the display filter.
- 5. When both a start and end date have been chosen, click **Apply**.
- 6. (Optional) Click Cancel to stop the query.
 - NOTE: Prior day selections persist if you change only the month or year. If you change the month or year and the prior day value does not exist, the calendar will default to the last day of the new month. For example, changing only the month value of an existing August 31, 2018 entry to February would set the date to February 28, 2018 because there are less days in February. Similarly, if you had February 29, 2016 selected, a day which only exists in a leap year, switching the year to 2015 would set the date to February 28, 2015, the last day of the month in a standard year.

The widget refreshes and displays only the records that fall within the custom date range. The date range in use is displayed in the widget header.

4.5.1.3. Tab Date and Text Filter Tool

All tabs in the Patient Data and Widgets page have a **Tab Date and Text Filter** tool that you can use to filter all widgets on a single tab by the same date range or text. The **Tab Date and Text Filter** tool has the same **Start and End Date** and **Quick Date Range** filters available on most individual widgets, but it applies the date filter to all widgets on an individual tab.

- 1. Select the **Show Tab Date and Text Filter Tool** hyperlink, or the **:** button and then select the **Show Tab Date and Text Filter Tool** hyperlink.
- 2. Enter dates in the **Start Date** and **End Date** fields or select a **Quick Date Range** value, select the applicable widgets, then select the **Apply to Selected Widgets on This Tab**

- button. To reset the date ranges of widgets, select **Restore Default Date Ranges**, select the applicable widgets, then select **Restore Default Date Range for Selected Widgets**. **OR**
- 3. Enter text in the **Text Filter** section and select the **Apply** button. Select the widgets you want to apply the text filter to, then select **Apply Text Filter to Selected Widgets on this Tab**. If the widgets are grayed out, check if you have already applied this filter. You can **Lock, Unlock, Remove, or Unlock and Remove** text filters once it is applied.

The **Tab Date and Text Filter** tool (Figure 53) is hidden by default. After you apply a date range or text filter, the tool remains open unless closed manually.

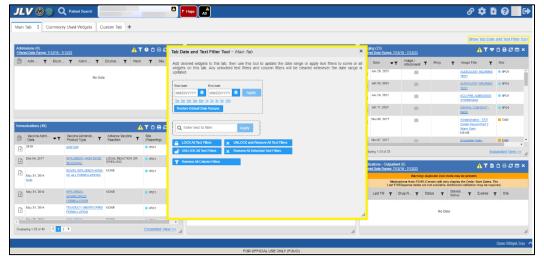


Figure 53: Tab Date and Text Filter Tool

The **Tab Date and Text Filter** tool only filters widgets currently displayed on the active tab that support date range filtering. Widgets that display all data, without a date range, will continue to do so. Filtering only applies to widgets on the tab at the time you apply the **Tab Date and Text Filter** tool and will not affect new widgets added to the tab. You may still apply widget-level date and text filters to individual widgets even after using the **Tab Date and Text Filter** tool. All your widget date and text settings will save to your profile regardless of whether they were set with the **Tab Date and Text Filter** tool or widget-level date or text filters.

4.5.2. Text Filter

Widgets provide a text filter option (Figure 54), allowing you to quickly locate relevant patient data using keywords or characters by applying multiple text filters. Text filters are performed on the information displayed in all columns within a widget, including hidden columns. Filters are not performed on notes, attachments, or details associated with patient records.

Documents & Images (1) ~

Filtered Date Ranger 12/19/22 4/19/2 4/

Figure 54: Documents & Images Widget Text Filter

Click **T** in the minimized view of a widget to display the **Enter text to filter** field. If the filter options are not displayed in the expanded view of the widget, clicking **T** also opens the filter options.

- A text filter box is displayed in the widget filter section. When text is entered, a pop-up box displays with "Press Enter or the + button to apply to filter."
- The applied search term is displayed below the text filter box with a red "X" next to it which allows a user to remove it and updates the search results accordingly.
- The funnel shaped icon and the **+Filters** displays with a yellow highlighted color to indicate that a filter has been applied in the widget.
- If another search term is typed in the Text Filter Box and the user clicks "enter" or "+", a second text filter is applied against the ORIGINAL DATA SET. Filtered results from a second text filter search are displayed in addition to the results from the first text filter search.
- If a third text filter is applied against the ORIGINAL DATA SET, filtered results from the first and second text filter searches are displayed in addition to the results from the third text filter search.
- The text filter box is disabled after a total of 10 text filters search terms are applied. Each search term can be up to 20 characters. The Text Filter Box remains disabled until one (or all) search terms are removed. A message displays with instructions that the maximum number of search terms have been applied, and to remove one or all terms to apply a different search term.
- The word "or" will appear between the search bubbles when multiple text filters are used.
- Text filtering applies only within the widget being searched.
- Text is filtered within the subset of filtered data if a date range filter or other dropdown filters are applied within the widget.
- Text filters are performed on data within all columns. Some widgets have columns that are not shown by default in either the minimized or expanded view.
- Text filters entered in the minimized view of a widget persist when you open the widget in expanded view, except for the Vitals and Lab Results widgets.

- On widgets with site tabs, text filters entered on one site tab persist when you select other site tabs.
- All records across the widget pages are filtered if the widget has multiple page views.
- The **Enter text to filter** field is not case sensitive and allows up to 25 characters, including numeric and special characters.
- JLV begins to filter when you pause while typing the term in the **Enter text to filter** field.
- Text filtering is not available in the **Demographics** widget.
- If a column is hidden, all filters are removed from that column in minimized and expanded views.
- Filtered text is highlighted.
- Users can lock and persist multiple text filters in a widget. This allows automatic application of key search terms to any patient chart.
- If a text filter is locked, the word will be highlighted in yellow with the locked icon closed and the "X" in gray. This text filter will not be able to be removed until unlocked.

4.6. Viewing Connection Status

JLV widgets display the status of their connection to VA and DOD data sources. Connection status information is available for each widget.



NOTE: The Connection Status reports the condition of the connection between JLV and its external resources, while System Status reports the overall condition of the JLV application. See <u>Viewing System Status</u> for more information.

The **Connection Status** indicator icon on a widget toolbar indicates the state of the connection to VA and DOD data sources. There are two status conditions:

- **1** All sources are connected (available).
- • One or more sources are not connected (unavailable).

Figure 55 shows ! in the **Demographics** widget. When the Primary Care Management Module (PCMM) service is unavailable, the primary care fields in the **Demographics** widget display this error message: "Team assignment system is currently unavailable."

<u>A</u>₿_□× Patient Demographics ~ DAYTSHR ▶ FEHR SPOEHR03 DoD CHYSHR VA: CHYSHR + Preferred Name: Pronouns: DOB Asian VA Care Teams: Team assignment system is currently unavailable. Military Service Birth Sex Self-Identified Gender Identity Eligibility and Enrollment Insurance DoD Sharing Opted Out Status: (Jan 26, 2022) Details:

Figure 55: Demographics Widget, PCMM Connection Status Unavailable

<u>Figure 56</u> shows a connection status indicator icon displayed on the **Patient Flags** widget toolbar.


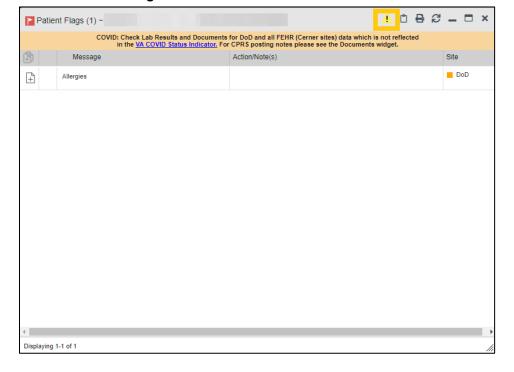


Figure 56: Connection Status Indicator

Clicking a **Connection Status** indicator icon (<u>Figure 56</u>) opens a standalone window with information about the connection to data sources. Clicking the < <u>Hide All Active Interfaces</u> link (<u>Figure 57</u>) in the **Connection Status** window hides the detailed **All Active Connections** view and displays only the connection errors and connection time out errors.

Figure 57: Connection Status Details

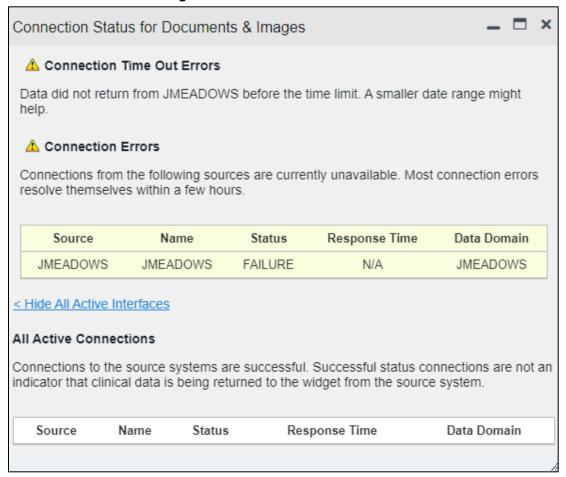

Clicking the > Show All Active Interfaces link in the Connection Status window (Figure 58) opens a Connection Status Details view (Figure 57).

_ 🗆 × Connection Status for Lab Results Connection Errors Connections from the following sources are currently unavailable. Most connection errors resolve themselves within a few hours. Name Data Domain Source Status Response Time VA CHYSHR Connection Unavailable 00:21 Lab Results IPOTEST1 Connection Unavailable 00:21 Lab Results VA VA IPOTEST2 Connection Unavailable 00:21 Lab Results <u>DOD</u> **TMDS** NOT_AVAILABLE 01:06 Lab Chemistries **FEHR FEHR FAILURE** Lab Chemistries 01:06 **FEHR FEHR FAILURE** 01:01 Lab Results > Show All Active Interfaces

Figure 58: Connection Status Window

4.7. Displaying Widgets on Custom Tabs

You can use the *custom tabs* feature to create additional widget configurations for ease of use and quick reference. Changes made to the Patient Data and Widgets page layouts (widget layouts and custom tabs) are saved to your user profile and displayed in future sessions.

- 1. Click + beside the existing tabs on the Patient Data and Widgets page (Figure 59).
 - a. The **Add Tab** dialog box opens (<u>Figure 60</u>).



Figure 59: Add Custom Tab Icon

Figure 60: Add Tab Dialog Box



- 2. Type the name of the new tab in the **Add Tab** dialog box, along with a description, and click **Save**. A message will display if there are any duplicate Tab names, or if max character limits have been reached for the name or description for the custom tab.
 - a. The custom tab opens with a blank workspace, with instructions on how to add widgets and an open widget tray.
 - b. After the tab has been created, hovering on the tab name will display the tab description.
- 3. Click, hold, and drag each of the desired widgets to the screen (max 12).
- 4. Click between tabs at any time, without losing each tab's configuration. Click, hold, and drag the tabs into any order, when there is more than a single tab.
 - a. Tab layouts persist, even when a patient change is made.
- 5. Hovering over the Tab Menu area of the selected tab, the options to Edit Name and Description, Create a Copy, and Delete will become available.
 - a. When creating a copy of a tab, a copy of the tab name will be created with all the settings from the copied widgets, filters, etc.
 - b. Once the copy has been created, the two tabs will be independent of each other.
 - NOTE: Tab configurations remain set until you manually change them or restore your profile to default in the JLV Settings dialog, which will remove all custom tabs.

4.7.1. Exporting Tabs

Users can export one or more tabs using the Export Tabs option, located under Settings. Upon clicking Export Tabs, a scrollable list of tabs will be available to export in the same order as in the user's JLV application, with a tab description available as a hover hint.

Export Tabs

User Tabs

Export All Tabs

Main

Description (Optional)

Commonly Used Widgets

User Roles/Workflow Tags (max 50 characters)

Figure 61: Export Tabs

×

By default, no tabs are preselected within the list, but users may choose one, multiple, or all tabs as seen in (Figure 61). The Export button will become highlighted when one or more Tabs are selected. An optional description may be added with up to 50 characters that will be included in the Export.


Cancel

After clicking Export, a Save As dialogue box will appear with a default export file name, <userLastNameFirstName>_JLV_tab_config. Users will be able to edit the name of the exported tab however, the extension will be fixed to .JSON and is not changeable. Users can then save the exported config to a local or shared drive.

4.7.2. Importing Tabs

JLV users will be able to import one or more tabs using the Import Tabs option (Figure 62), located under Settings. When this option is selected, the user sees a dialogue box where JSON files will be displayed, any available tab descriptions will be displayed as a hover hint on the tab name. After the user clicks Import, the user's current configuration is saved and JLV starts the tab import.


Figure 62: Import Tabs

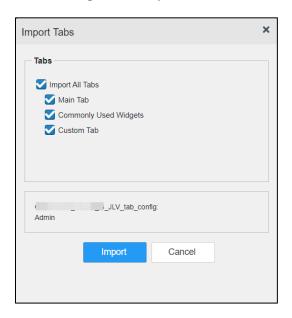

The imported tab will display to the right of the last tab. If the tab is imported and it does not exist only the <tabname> will appear. If a tab with the same name already exists in the user's JLV workspace, the tab being imported will be named <tabname>*, and on subsequent duplicates, the imported tabs will be named with added asterisk, <tabname>**.

An information dialogue popup will be displayed that the import was successful (<u>Figure 63</u>). When the user logs out and logs back in, the imported tabs will still be displayed.

Figure 63: Import Successful

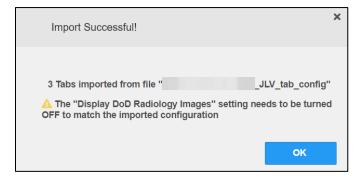

If the import is unsuccessful due to an incorrect file type, an alternate popup will appear (<u>Figure</u> 64).

Figure 64: Import Failed



4.8. Refreshing Data in a Widget

The data displayed in a widget can be refreshed by clicking **2**. This action retrieves data from VA and DOD sources and refreshes the individual widget's display.

4.9. Terminology Normalization

Table 11: Terminology Specific to Widgets

Normalization Standard	Description	Widgets
Centers for Disease Control (CDC) Race Standard	Code set based on the current Federal standards to classify race and ethnicity	Patient Demographics
Current Procedural Terminology (CPT) Standard	Uniform code to describe medical, surgical, and diagnostic services	Procedures, Radiology Reports15
Vaccine Administered (CVX) Standard	Numeric string that identifies the type of vaccine product used	Immunizations
Logical Observation Identifiers Names and Codes (LOINC) Standard	Universal standard for identifying health measurements, observations, and documents	Discharge/Inpatient Summaries, Documents & Images, Lab Results, Orders, Progress Notes, Questionnaires, Deployment, Exposure Records, Assessments, Radiology Reports, Vitals
National Uniform Claim Committee (NUCC) Taxonomy Standard	Unique, 10-character alphanumeric code that identifies a provider grouping, classification, and area of specialization	Admissions, Appointments, Consults, Documents & Images, Medications - Inpatient, Discharge/Inpatient Summaries, Orders, Encounters - Outpatient, MedicationsOutpatient, Procedures, Progress Notes
RxNorm Standard	Standard names given to allergens, clinical drugs, and drug delivery devices in the United States	Allergies, MedicationsInpatient, Orders, MedicationsOutpatient

¹⁵ If LOINC mappings are not available for an entry in the Radiology Reports widget, JLV looks for CPT mappings and displays CPT normalization for standardized type and standardized code, if found.

Normalization Standard	Description	Widgets
SNOMED CT Standard	Core, general terminology used in electronic health records	Admissions, Encounters - – Outpatient, Problem List
X12 Health Insurance Type Standard	Defines electronic data interchange standards for health care insurance	Patient Demographics

4.10. Image Support

JLV integrates access to the VistA Imaging Viewer and MUSE, allowing you access to VA imaging artifacts for supported clinical domains (widgets). Please see the REDACTED for more information regarding the VistA Imaging Viewer.

You can open one or more images associated with a record by selecting the displayed in the Image column (Figure 65) of the Advanced Directives, Cardiology Studies - FEHR and VA MUSE Only, Documents & Images, Imaging, Procedures, Progress Notes, Radiology Reports, and Surgery/Procedure Reports widgets. Clicking from the Image column opens the image(s) in a standalone window. Like a hyperlink, changes color after you launch an image so you can track which images you've viewed in the minimized or expanded widget view of the current session. Changing between minimized and expanded view resets the color.



Figure 65: Camera Icon in the Image Column

If a record has one or more images associated with it, the **Details** view of the record includes a thumbnail with a <u>Click image to open viewer</u> link (<u>Figure 66</u>). Clicking either the thumbnail or the link opens the VistA Imaging Viewer and displays the associated image(s).

Document & Images Details ~ • IPO1 NUTRITION HYPERLIPIDEMIA EDUCATION CONSULT REPORT ~

LOCAL TITLE: NUTRITION HYPERLIPIDEMIA EDUCATION CONSULT REPORT STANDARD TITLE: NUTRITION EDUCATION NOTE DATE OF NOTE: JUL 03, 2022@09:00 ENTRY DATE: JUL 03, 2022@15:46:03 AUTHOR: EXP COSIGNER: URGENCY: STATUS: COMPLETED

VistA Imaging - Manual Closure

Administrative Closure: 07/03/2022 by:

Figure 66: Thumbnail and Link in Details View of the Documents & Images Widget

4.11. Using Report Builder

The **Report Builder** feature is used to create custom PDF reports using patient data and records displayed in JLV widgets. Content for reports can be selected from either the minimized or expanded view of a supported widget. You can launch **Report Builder** from the toolbar on the **Patient Data and Widgets page** (highlighted in Figure 67).

Figure 67: Report Builder Icon



The **Report Builder** pane is comprised of the **New Report** tab (<u>Figure 69</u>) and the **Saved Reports** tab (<u>Figure 70</u>). When open, the **Report Builder** appears docked on the left side of the Patient Data and Widgets page. It can be moved and docked anywhere on the page. JLV will save the location the user selects for the RB window, but will be reset with the option to restore user settings.

Figure 68: Report Builder Header Icons



Click the show () or hide icon () to display or hide the **Report Builder**.

Click the Pin icon () to pin or unpin the Report Builder to a location in the interface. When unpinned, users are able to resize the Report Builder window horizontally from both the left and right edges, and vertically from the lower edge. The unpinned Report Builder will not be able to be positioned where it could hide the tab name or the top JLV banner. The bottom edge will

always be clearly visible so it can be resized without moving it. The unpinned size will be saved for the current and subsequent JLV sessions.

Click the Close icon () to close the Report Builder. Any report configurations in progress are saved when collapsing or closing the Report Builder and restored when the Report Builder is displayed again.

If Report Builder is hidden or out of focus, clicking on the Report Builder icon or on any visible part of the Report Builder popup will bring it back into focus.

The **New Report** tab (Figure 69) is used to build the custom report. The patient data and/or records are listed here as you select and add them. The **print linked artifacts in the report** box will be checked automatically, with all attachments that are unable to be printed denoted by **. The option to **Clear All** records after Build will need to be manually checked if the functionality is desired. The records can be arranged in the desired order for the final report. The total number of New Reports will appear on the tab. Once a record is selected, users may use the to clear a record. If the selected record is no longer visible after scrolling the user will be unable to take any actions on the selected record until it is visible again. Users will be able to unselect the already selected record by clicking on the record again, If Report Builder is closed, the record will be automatically unselected.

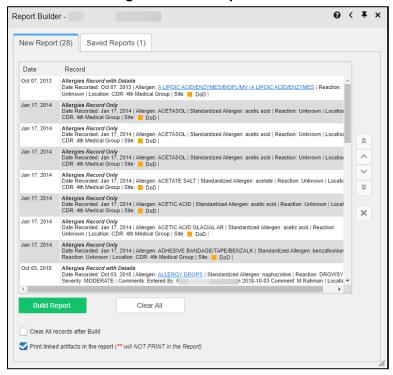


Figure 69: New Report Tab

The **Saved Reports** tab (<u>Figure 70</u>) presents a list of all completed reports with the option to select one or all of the reports for action, the status of the report generation progress, the date the report expires which is seven days, the report size, and a hyperlink to the contents of the

completed report. If one or more is selected, Saved Reports can be deleted or downloaded using the buttons.

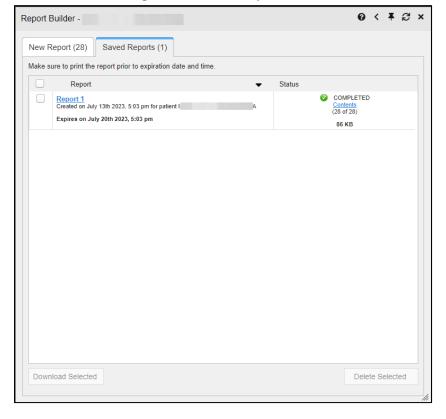


Figure 70: Saved Reports Tab

You can perform actions within widgets to add patient data when the Report Builder is open or closed. JLV enforces a maximum number of 200 records per report, where no more than 100 of these records have attachments. A message appears when this limit is reached.

4.11.1. Adding One Record

Click next to a record in a widget to add data from that record to the Report Builder, including the details and notes associated with the record. in the row changes to confirm the record has been added to the Report Builder (Figure 72). This applies to both minimized and expanded widget views and the **Patient Flags** dialog box.

• **NOTE:** A single site report can now be added to the Report Builder from the Blood Bank Report in the Lab Results widget.

4.11.2. Removing One Record Directly from a Widget

Click to remove a singular record from Report Builder directly from a widget. The record will be removed from Report Builder and will no longer display in the Report Builder selected records. The icon will change back to .

4.11.3. Adding Multiple Records with EZ Select

EZ Select allows you to add all records displayed on a single page of a single widget to the Report Builder with one click. Click (Figure 71) to add data from all records shown in the widget's page display to the Report Builder, including the details and notes associated with the record. EZ Select is available from the minimized and expanded views as well as the **Patient Flags** dialog box.

0

NOTE: Sensitive records and Medication Administration History (MAH) / Medication Administration Log (MAL) reports accessed from the Medications - Inpatient widget cannot be added using EZ Select. You must add them to a report one at a time, and you can only add MAH/MAL to reports from Detail view in the Medications - Inpatient widget. However, you can use EZ Select to add MAH/MAL reports to Report Builder from the Health Summaries and Reports widget.

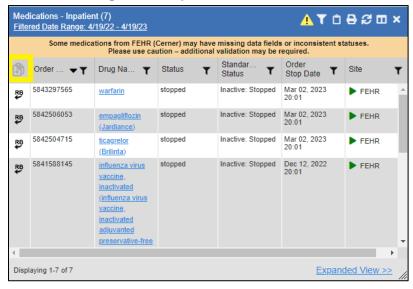


Figure 71: Report Builder EZ Select

in each row changes to be to confirm the record has been added (Figure 72).

Immunizations (161) ATOBSEX Vaccine Admini...
Date Vaccine Administe.. Adverse Vaccine (Reporting) INFLUENZA, INJECTABLE, Mar 22, 2023 CHYSHR QUADRIVALENT, PRESERVATIVE FREE COVID-19 (PFIZER), CHYSHR MRNA, LNP-S, PF, 30 MCG/0.3 ML DOSE, TRIS-SUCROSE (AGES 12+ INFLUENZA, SEASONAL CHYSHR INJECTABLE COVID-19 (MODERNA), CHYSHR MRNA, LNP-S, PF, 100 MCG/0.5ML DOSE OR 50 COVID-19 CHYSHR (ASTRAZENECA) VECTOR-NR, RS-Displaying 1-25 of 161 < 1 2 3 4 > Expanded View >

Figure 72: Records Added to Report Builder

NOTE: Records previously added to a report remain in the Report Builder after a widget refresh. Clicking in the widget refreshes the clinical data, refreshes the widget data display, and resets to , even if the record has been added to the Report Builder. Records previously added to a report prior to the refresh remain in the Report Builder after a widget refresh.

4.11.4. Adding Details from Multiple Dates

When the **Details** view of a record has a date filter (<u>Figure 73</u>)—as in the **Health Summaries** and **Reports, Medications - –Inpatient**, and **Procedures** widgets—resets to each time you change the date range, allowing you to build date-flexible reports using records from multiple date ranges.

T 🖶 📋 🛨 Medicine Full Report ~ National Start date: End date: 鯔 09/12/2019 09/12/2016 Apply 1w 2w 1m 3m 6m 1y 2y 3y 5y 10y Close Filter IPO4 IPO3 IPO5 Printed for data from 09/12/2016 to 09/12/2019 ******************** CONFIDENTIAL SUMMARY pg. 1 ******************** DOB: C MEDICINE ---- MEDF - Med Full Report -JUL 5,2019@13:01 CARDIOLOGY Date/time: JUL 5,2019@13:01 Medical Patient: Indication Comment: JLV TESTING - VISTAWEB MEF REPORT Provider/physician:

Figure 73: Adding Date-Flexible Details to Report Builder, Medicine Full Report

4.11.5. Adding Sensitive Records

Sensitive records can be added to the Report Builder by for a sensitive record in the widget. When prompted, click **Agree to be Audited & Access This Sensitive Record** to continue. in the row changes to to confirm the record has been added (Figure 72). Sensitive records must be added to the Report Builder individually; therefore, repeat as necessary for additional sensitive items.

4.11.6. Adding Only Record Details or Notes

When adding records to a report using or in the widget, JLV by default adds the data displayed in the widget as well as any details and notes associated with the record. This applies to both minimized and expanded widget views and the **Patient Flags** dialog box. If you would like to add just the details and notes or the record row data:

- 2. Select and hold $\stackrel{\square}{=}$ to access these options (Figure 74).
- 3. Select **Add Detail** or **Add Note** to add the details and/or notes associated with the record as well as the record row data to the Report Builder.
- 4. Select the **Add Record Row** option to add only the data displayed in the widget display for that record and not the information provided in the details and/or notes.

Lab Results (62) <u>♠<u>♥</u>▼☆₽♂□×</u> Filtered Date Range: 4/19/22 - 4/19/23 If you don't see FEHR results reduce the date range to 3m or less Collectio. Lab Test Result ▼ Lab Panel Site Add Detail Blood Bank Report result Add Record 16:49 97 mL/min/1.73 m2 eGFR CKD EPI ► FEHR вмР ▶ FEHR 140 mmol/L ВМР ▶ FEHR 31.0 mmol/L **BMP** ▶ FEHR 14 mEq/L BMP ▶ FEHR 10.0 mg/dL ВМР ► FEHR Displaying 1-25 of 62 < 1 2 3 > Expanded View >

Figure 74: Report Builder Options

NOTE: Once you add a record to Report Builder using either the Add Detail/Note, or Add Record Row option, you cannot change the record data included in the report without first clearing the record from the report and adding it again.

4.11.7. Adding Patient Demographics

Clicking in the **Demographics** widget (Figure 75) adds the patient's demographics details for that site to the Report Builder. changes to to confirm the demographics data has been added. Click additional site tabs within the **Demographics** widget and repeat as desired to add demographics details from other sites within the patient's health record.



Figure 75: Demographics Widget Add to Report Builder Icon

4.11.8. Generating a Report

The selected records appear in the **Report Builder** pane (Figure 69).

- 1. Use the navigation arrows to arrange the records in the desired order in the generated report. When a record is selected for the report, the complete record will be visible including all the columns in the widget view, hyperlinks for details and attachments in the **Record** column of the **Report Builder** pane. The **Record** column will be sortable.
- 2. Click **Build Report** to validate the selected records and prepare the report. If a record is added to the report without an error, a **Ready** notation displays next to it in the **Status** column. If a record cannot be added to the report, an **Error** notation appears next to the record. Users will have the option to **Clear All** of the records after the report starts to build. If the user cancels the build, the records will not be cleared.
- 3. Name the report when prompted, and click **OK** (Figure 76)

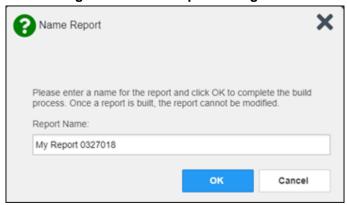


Figure 76: Name Report Dialog Box

- 4. The **Saved Reports** tab opens and displays a report processing indicator in the **Status** column.
- 5. When processing is complete, the **Status** column displays either *COMPLETED* or *ERROR* (Figure 77).
 - An *ERROR* in the **Status** column does not indicate the report failed to build, it is an indicator that one or more records could not be included in the generated report.
 - Both the *COMPLETED* and *ERROR* entries include a <u>Contents</u> link (<u>Figure 77</u>), which provides a list of the records that appear in the generated report.


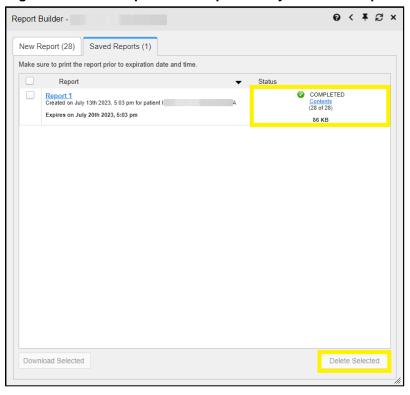


Figure 77: Saved Reports Tab Report Ready and Delete Option

6. A report ready message appears for 6 seconds on the Patient Data and Widgets page when the report has been built and is ready to be printed (<u>Figure 78</u>).

Figure 78: Report Ready Message and Indicator

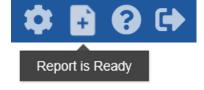

NOTE: Once a report is created, it is available on the Saved Reports tab for seven days. After seven days, the report expires, is removed from the Saved Reports tab, and can no longer be accessed.

4.11.9. Opening a Report

- 1. Click the report name in the **Report** column of the **Saved Reports** tab (Figure 77) **OR** -
- 2. Click the <u>Contents</u> link in the **Status** column of the **Saved Reports** tab to open a list of the records included in the report in the **Report Contents** window (<u>Figure 77</u>) **OR** -
- 3. Click the blue, hyperlinked report name in the **Report Contents** window (Figure 79).

Figure 79: Report Contents Window




NOTE: The Report Builder prints the report content to a file in PDF format. It is recommended that you have the latest Adobe Reader installed on the system from which you access JLV to utilize the Report Builder and other JLV features.

4.11.10. Delete a Report

- 1. Select the checkbox for the report name on the **Saved Reports** tab (<u>Figure 77</u>).
- 2. Click the **Delete Selected** button.
- 3. Click **OK** to confirm deletion (Figure 80).

Figure 80: Delete Reports Dialog Box

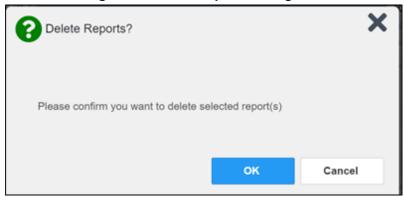

4.11.11. Downloading and Printing a Report

There are two ways to download a report.

Option1: Once a new report has been created, click on the hyperlink to open the report (<u>Figure 81</u>).

•

0 > # 2 × Report Builder -Saved Reports (1) New Report (29) Make sure to print the report prior to expiration date and Report Status COMPLETE <u>report</u> June 27th 2023, 12:06 pm for Contents (29 of 29) Expires on July 4th 2023, 12:06 pm 97 KB **Download Selected** Delete Selected

Figure 81: Report Hyperlink

Click the ... option (Figure 82) within the report window to download an Adobe PDF report created in Report Builder. A printing option will be visible within this view.

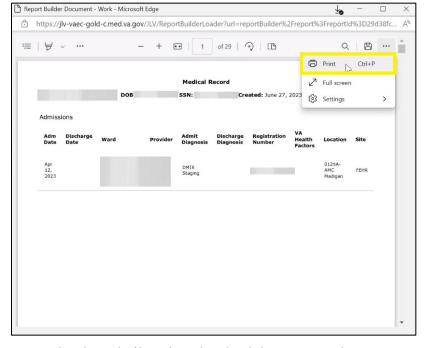


Figure 82: Printing the Downloaded Report

Option 2: After a report has been built, select the check box next to the report name. Then, click on the **Download Selected** button (Figure 83).

Report Builder
New Report (29) Saved Reports (1)

Make sure to print the report prior to expiration date and time.

Report Status


report
Created on June 27th 2023, 12:06 pm for patient
Expires on July 4th 2023, 12:06 pm

Delete Selected

Figure 83: Download Selected Option

Once the file is finished downloading, select Open (Figure 84).

Figure 84: Opening the PDF

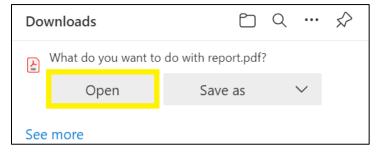

To Print the report, click on the print icon in the toolbar (<u>Figure 85</u>).

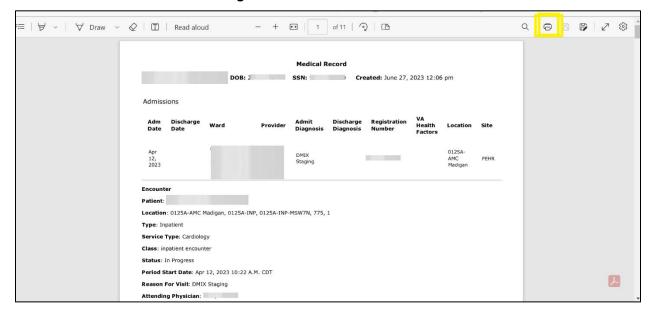


Figure 85: Print Icon on the PDF

4.12. Printing

The data list of a widget in either minimized or expanded view, the details window of a widget, and reports created in **Report Builder** can be printed.

- Click on the desired widget's toolbar to print the data list of a widget in either minimized or expanded view.
- Click **o** on the details window's toolbar to print the details window of a widget.
- Click the (PDF) **Print** icon within the report window to print a report created in **Report Builder**.

The columns currently displayed in the widget will be checked by default in the column popup box. The print selections for the widget will remain, even when switching patients, unless restore default settings is executed.

NOTE: Each report generated using **Report Builder** includes this disclaimer: "The information contained in this transmission may contain privileged and confidential information, including patient information protected by federal and state privacy laws."

5. Widget Functionality

The following subsections detail the functionality in a sampling of the widgets available in JLV.

5.1. Patient Demographics Widget

The **Patient Demographics** widget displays a summary of the patient's nonclinical, personal data on the primary JLV toolbar (<u>Figure 86</u>). Clicking the Patient Demographics widget will open the expanded view of the widget (<u>Figure 87</u>).

Figure 86: Patient Demographics Widget in the JLV Toolbar



The site tabs provide a demographics summary for each site (VA, FEHR, or DOD) the patient has visited.

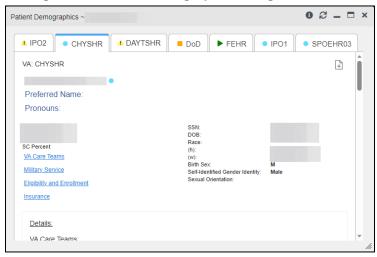


Figure 87: Patient Demographics Widget Site Tabs

Nonclinical, personal data displayed for the patient varies by site and could include site name, patient name, address, SSN, DOB with the age in parenthesis, race, phone numbers, birth sex, and gender ID (SIGI). SIGI code values are listed in Table 12. If the patient is deceased, a DECEASED label will be displayed along with the date of death and age at death in parenthesis next to it in red. The deceased patient's age will not be displayed next to the DOB.

Table 12: Gender ID (SIGI) Code Values

Code	Description
M	Male
F	Female
TM	Transmale/Transman/Female-to-Male
TF	Transfemale/Transwoman/Male-to-Female
0	Other
NB	Nonbinary
Any Other Code	Unknown

The **Patient Demographics** widget provides additional functionality:

- Clicking Patient Demographics opens details in a standalone window, including inpatient providers.
- Clicking or not opens the connection status details in a separate window.

- Clicking Flags opens clinical flag details in a new window.
- Clicking AD opens the **Advance Directives** dialog in a new window.
- Clicking refreshes Patient Demographics and Patient Flags simultaneously.

If <u>A</u> displays on one or more tabs within the **Demographics** widget, it is an indication that the patient is not registered at that VA site or DOD facility.

5.1.1. Patient Flags

If the patient has one or more clinical warnings, alerts, or flags in their record from VA and DOD sites, flags displays next to the **Demographics** widget banner on the Patient Data and Widgets page (Figure 88). A spinner will appear on the Flags icon while data is being retrieved. If the patient does not have flags, the Flags icon will appear greyed out in the toolbar and will not be clickable and Flags. A hover hint for the greyed-out icon will appear, "Patient does not have any Flags".

The **Patient Flags** window opens and displays by default when the patient data appears, if there are warnings associated with the selected patient. If there are any connection errors in Flags, the yellow connection error icon will be visible on the Flags icon Flags.

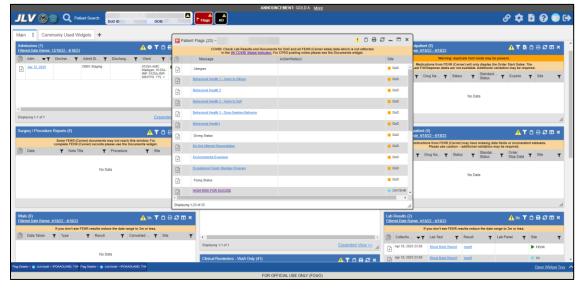


Figure 88: Patient Flags Indicator

Click Flags to open the **Patient Flags** window (Figure 88). The window displays a list of alerts and flags within the patient's record. Selecting the hyperlinked message for each patient alert opens the **Flag Details** window (Figure 89). Selecting the Action/Note(s) link (Figure 88) will open documents pertaining to placement of the flag. Category I Patient flags that are present at multiple sites will be grouped under the owner site.

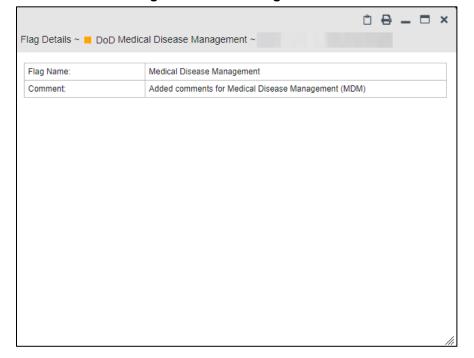


Figure 89: Patient Flags Details

5.1.2. Advance Directives

You can view Advance Directives (Figure 90) by selecting AD on Demographics widget. Advance Directives are pulled from DOD, VA, and Cerner data sources. A spinner will appear on the Advanced Directives icon while data is being retrieved and remains greyed out and unclickable AD if Advance Directives were not found for the patient. A hover hint will appear, "Patient does not have any Advanced Directives". If there are any connection errors, the yellow connection error icon will be visible on the AD icon AD.

<u>1 û ₽ 2 - □ ×</u> Advance Directives (5) ~ (A) Date ▼ ▼ Document Title ▼ Site Dec 15, 2022 Advance Directive Documents / ► FEHR Linux Migration 12/15/22 Details Jun 30, 2020 ADVANCE DIRECTIVE CHYSHR Oct 04, 2017 ADVANCE DIRECTIVE CHYSHR ADVANCE DIRECTIVE CHYSHR ADVANCE DIRECTIVE CHYSHR Displaying 1-5 of 5

Figure 90: Advanced Directives Dialog

5.1.3. Viewing Third-Party Insurance Information

Third-party payers and insurance information is available from the **Patient Demographics** widget. Click the **Insurance** link in the **Demographics** widget.

A new window opens with the following insurance information, some of which is displayed in (Figure 91):

- Health Plan Type
- Health Plan Name
- Standardized Insurance Type (See <u>Terminology Normalization</u> for more information.)
- Plan Effective Date
- Plan Expiration Date
- Site
- Group Number
- Member ID
- Subscriber ID
- RxBIN
- RxPCN
- Notes
- Comments
- Subscriber Date of Birth
- Subscriber's Relationship to Insurer
- Health Plan Mailing Address
- Health Plan Contact Information

<u>. ▼ ů 8 8 ⊞ = □ ×</u> Insurance Information (2) ~ Remove All Filters Hide Filters Health Plan Ty... ▼ Health Plan N... ▲▼ Standardized Insurance Ty... ▼ Plan Effective Date ▼ Plan Expiration Date ▼ Site ▼ Group Number ▼ Member ID ▼ Subsc Health Plan Type was 000 CIVILIAN EMERGENCY RX Mar 24, 2023 ▶ FEHR not supplied Health Plan Type was not supplied CIVILIAN EMERGENCY 000 Mar 24, 2023 ▶ FEHR Displaying 1-2 of 2


Figure 91: Insurance Information

Use the horizontal scroll bar in the window to view the columns not seen within the default window size.

5.1.4. Viewing Military Service Data

Military service data (<u>Figure 92</u>) are available from the **Patient Demographics** widget. Click the <u>Military Service</u> link in the **Demographics** widget.


Figure 92: Military Service Window

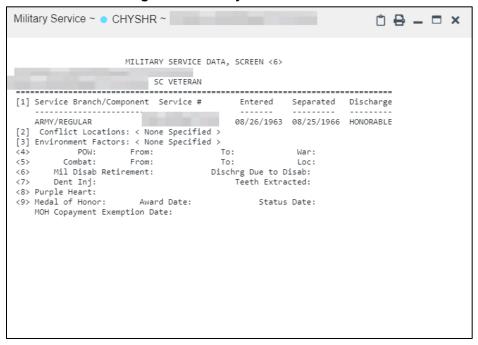

5.1.5. Viewing Eligibility and Enrollment Data

Eligibility and enrollment data (Figure 93) are available from the **Patient Demographics** widget. Click the **Eligibility & Enrollment** link in the **Demographics** widget.

Figure 93: Eligibility and Enrollment Window

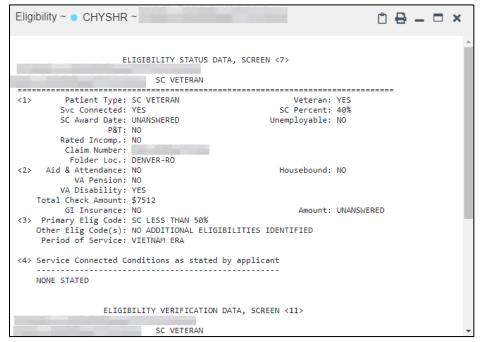

5.2. Cardiology Studies - -FEHR and VA MUSE Only Widget

The Cardiology Studies - FEHR and VA MUSE Only widget (Figure 94) displays electrocardiograms (EKGs) and cardiology studies with related images from all FEHR and VA MUSE sites and systems. The information displays in reverse chronological order by date with one tab for each site where the patient is registered.

Clicking from the **Image** column opens the image(s) in a standalone window.

Report Builder functionality is not available for this widget.

5.2.1. Cardiology Studies - -FEHR and VA MUSE Only Widget Data

The information in the Cardiology Studies - –FEHR and VA MUSE Only widget includes:

- Date
- Description (Test Type)
- Status
- Image
- Site

There is no expanded view for this widget.

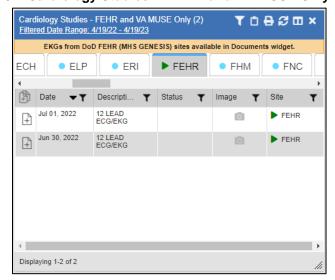


Figure 94: Cardiology Studies - -FEHR and VA MUSE Only Widget

5.3. Community Health Summaries and Documents

The Community Health Summaries and Documents widget (Figure 95) displays the patient's community partner information, including Continuity of Care Documents (CCDs) in C32 and C62 formats and HL7 Consolidated-Clinical Document Architecture (C-CDA) structured documents available for the patient. The information in this widget is sorted by the title of the entry listed in the **Document Type/Title** column.

NOTE: By default, CAPRI-Claims users cannot see the Community Health Summaries and Documents widget. To enable access to this widget, CAPRI-Claims users must be assigned the security key, "DVBA ALLOW JLV COMM DATA".

5.3.1. Community Health Summaries and Documents Widget Data

Data displayed within the **Community Health Summaries and Documents** widget (<u>Figure 95</u>) in the minimized view includes:

- Date
- Document Type/Title
- Site

The purple hexagon beside entries in the **Site** column denotes the source of the data is outside the VA.

Community Health Summaries and Documents (32) TOBSOX To request VHIE <u>training</u> and troubleshooting support for this widget, please contact <u>VHIECoordinators@va.gov</u>. Date ▼ ▼ Document Type/Title Apr 19, 2023 Big Sky Care Connect Longitudinal Patient Big Sky Care Connect Summary Patient Summary East Tennessee HIN Continuity of Care Document InnovAge Continuity of Care Document HealthInfoNet Continuity of Care Document Baseline West Medical Center Feb 07, 2023 attachment Kentucky HIE Office Visit Summary InnovAge < 1 2 > Displaying 1-25 of 32 Expanded View >>

Figure 95: Community Health Summaries and Documents, Minimized View


Click <u>Expanded View >></u> in the minimized view of the **Documents & Images** widget to open the expanded view (<u>Figure 96</u>). Information in expanded view includes:

- Date
- Document Type/Title
- Site

Figure 96: Community Health Summaries and Documents, Expanded View

5.3.2. Viewing VA Community Health Summary Documents

Click on the hyperlinked entries in the **Document Type/Title** column within the **Community Health Summaries and Documents** widget, to open the documents in a new browser window. (Figure 97).

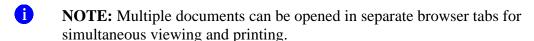

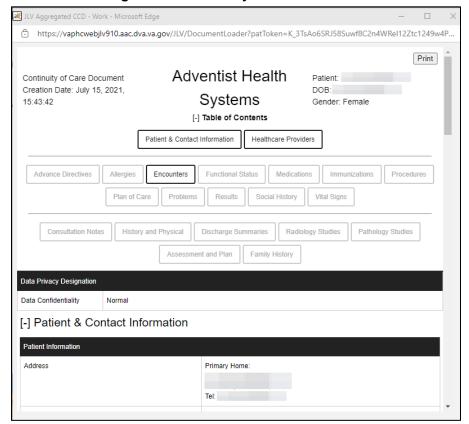


Figure 97: Community Health Document

5.4. Imaging Widget

The **Imaging** widget (<u>Figure 98</u>) displays images from multiple sources, including all Vista Imaging, Cerner imaging, DOD Radiology, and HAIMS artifacts. When the widget is displayed, data is sorted by date.

Record data can be added to Report Builder, although images cannot be added.

1. **NOTE:** The **Imaging** widget is only available to VA JLV users. DoD JLV users do not have access to the widget and do not see it in the Patient Data and Widgets widget tray.

Imaging (19) <u>Filtered Date Range: 4/18/22 - 4/18/23</u> <u>▲♥▼☆₽₽□×</u> Image / Date #Img Image Title Site Mar 08, 2017 0 Administrative - STR DoD **Dental Record Part 3** (Navy Only) 128 KB Mar 07, 2017 0 Encounter Note -DoD Consultation 0 Encounter Note -DoD Consultation 16 KB DoD 0 Assessment -Intake/Screening 20 KB Mar 02, 2017 0 Administrative - STR DoD AHLTA.pdf 20 KB Displaying 1-19 of 19 Expanded View >>

Figure 98: Imaging Widget, Minimized View

5.4.1. Imaging Widget Data

The minimized view of the **Imaging** widget displays the following information:

- Date
- Image/Attachment (See Image Support)
- #Img
- Image Title
- Site

Clicking the con in the toolbar displays the Cardiology Studies - -FEHR and VA MUSE Only Widget which has electrocardiograms (EKGs) and cardiology studies, and related images from all FEHR and VA MUSE sites and systems.


5.4.2. Viewing Images

Click **Expanded View** >> in the minimized view of the **Imaging** widget to open the expanded view (<u>Figure 99</u>). Information in expanded view includes:

- Date
- Image/Attachment (See Image Support)
- #Img
- Image Title
- Standardized Description
- Description

- Procedure
- Origin
- Specialty
- Event
- Type
- Capture Date
- Source System
- Location
- Site



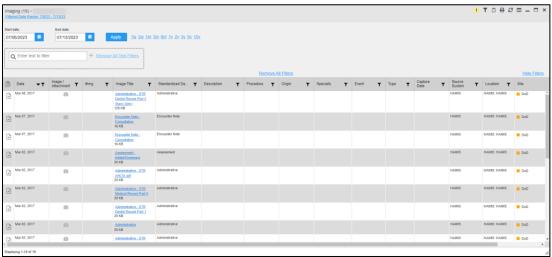

Click a hyperlinked entry in the **Image/Attachment** or the **Image Title** column in either the minimized or expanded view of the **Imaging** widget to access additional information for a listed record in text, PDF, graphic, or document format. The additional information will open in a separate browser window when the image source is VistA Imaging. <u>Figure 100</u> and <u>Figure 101</u> display two types of data accessed through the **Imaging** widget.

CT BRAIN/HEAD W/O CONTRAST

TOWN Summary

To

Figure 100: View of an Image

Audit By Patient 29 Oct 2015 1105 EDT Requested By: Criteria: Patient Id(s) = From Date = 29 Oct 2015, To Date = 29 Oct 2015 FMP/Sponsor SSN: Grade: Captain 03 DOR: Gender: F User Name | User ID | Action Clinical Event ID HAIR ID Sensitive BTG Reason n Date / h01982fb2a310e649b 9af80cb1de10584540 101 29 Oct 2015 1101 h01982fb2a310e649b 9af80cb1de10584540 101 ate Asset(s) 29 Oct 2015 vent 1101 h0134ccf2ffbf9841d9 9cac74ba6383e62201 01 Disassociate Asset 29 Oct 2015 (s) From Event 1059 Event Type: Encounte Event Id: Associate Asset(s) 29 Oct 2015 h0134ccf2ffbf9841d9 h01982fb2a310e649b 9af80cb1de10584540 h0176f18a14a11f4ca 8b2550637d3213f0e0 101 THIS INFORMATION IS PROTECTED BY THE PRIVACY ACT OF 1974 (PL-93-579). UNAUTHORIZED ACCESS TO THIS INFORMATION IS A VIOLATION OF FEDERAL LAW, VIOLATORS WILL BE PROSECUTED. FOR OFFICIAL USE ONLY Page 1 of 2

Figure 101: View of a PDF-formatted Document

5.5. Documents & Images Widget

The Documents & Images widget (Figure 102) includes documents from multiple clinical domains, including radiology reports (exams), progress notes, outpatient encounters, consults, discharge summaries (inpatient notes), questionnaires, and HAIMS records, displayed in reverse chronological order by the document date. HAIMS records retrieved by JLV include scanned paper records, imported paper records, advance directives, scanned non-radiology images, and imported non-radiology images.

5.5.1. Documents & Images Widget Data

The minimized view of the **Documents & Images** widget displays the following information:

- Date
- Description
- Provider
- Image/Attachment (See Image Support)

• Site

Figure 102: Documents & Images Widget, Minimized View



Outpatient encounter records displayed in the **Documents & Images** widget may have a <u>Details</u> link enabled in the **Date** column. Where available, click the link to open a **Details** window for records of this type.

Click **Expanded View >>** in the minimized view of the **Documents & Images** widget to open the expanded view (<u>Figure 103</u>). Information in expanded view includes:

- Date
- Description
- Standardized Description
- Provider
- Provider Specialty
- Location
- Status
- Image/Attachment
- Source System
- VA Health Factors
- Site

Document & Analysis (A1) - Color State Col

Figure 103: Documents & Images Widget, Expanded View

5.5.2. Viewing Documents & Images

Click a hyperlinked entry in the **Description** column in either the minimized or expanded view of the **Documents & Images** widget to access additional information for a listed record. Depending on the record type, the additional information may open in either a **Details** window or a separate browser window. <u>Figure 104</u> displays outpatient encounters details, accessed through the **Documents & Images** widget.

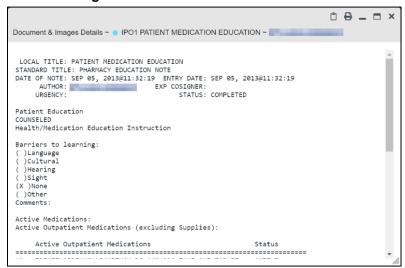


Figure 104: Details View of a Document

5.6. Federal EHRH Military Health System (MHS) GENESIS Widget

The Federal EHR/MHS GENESIS widget displays DOD patient documents and dental summaries from DOD's MHS GENESIS system and VA's Cerner PowerChart. The information is displayed in reverse chronological order by date.

5.6.1. Federal EHR/MHS GENESIS Widget Data

The minimized view of the **Federal EHR/MHS GENESIS** widget includes the following information:

- Date
- Document Title
- Document Type
- Site

<u>Figure 105</u> highlights the **Federal EHR/MHS GENESIS** widget in minimized view. There is no expanded view for the widget; instead, the C-CDA document selected opens in a new browser window (Figure 106).



NOTE: The **Medications** - **Outpatient** widget displays an orange banner with a warning that reads, "Duplicate DOD meds may be present" due to the migration of data into FEHR/MHS GENESIS (Cerner.)

Medical CCDs, Dental CCDs, and aggregate CCDs with blank entries in the **Date** column for the record appear first in the widget by default.

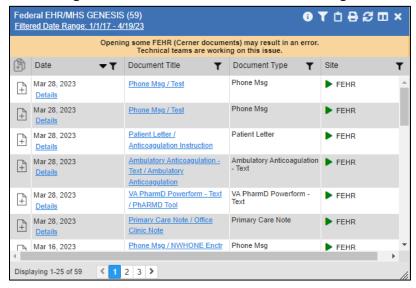


Figure 105: Federal EHR/MHS GENESIS Widget

5.6.2. Viewing Federal EHR/MHS GENESIS Documents

Click a hyperlink in the **Document Title** column in either view of the **Federal EHR/MHS GENESIS** widget to display a document (<u>Figure 106</u>) in a new browser window.

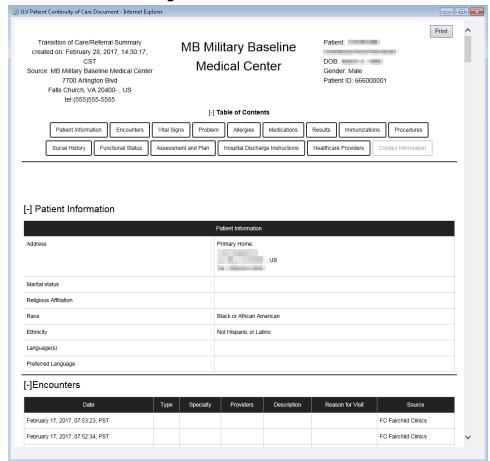


Figure 106: C-CDA Document

5.7. Health Summaries and Reports—VistA Only Widget

5.7.1. Health Summaries and Reports—VistA Only Widget Data

The Health Summaries and Reports—VistA Only widget displays the patient's national and local health summaries from VistA sites where the patient has been registered, including the Autopsy, Blood Bank, Cytopathology, Lab Summary—Cumulative, MAH¹⁶, MAL¹⁷, Medicine Full, Surgical Pathology, Transfers, and Unit Dose reports. This widget is available to VA users only.

JLV sorts records by **Site** by default. If a patient has national health summaries, JLV sorts to display National records first. National records are named in the **Site** column in the minimized view and in the **Type** and **Site** columns in the expanded view.

The **Description** column displays the document title, where available. Click a hyperlinked entry in the **Description** column in either the minimized or expanded views of the widget to open the document.


¹⁶ Limited to any 7-day period

¹⁷ Limited to any 14-day period

Information in the minimized view of the **Health Summaries and Reports—VistA Only** widget (Figure 107) includes:

- Description
- Site

Figure 107: Health Summaries and Reports—VistA Only Widget, Minimized View

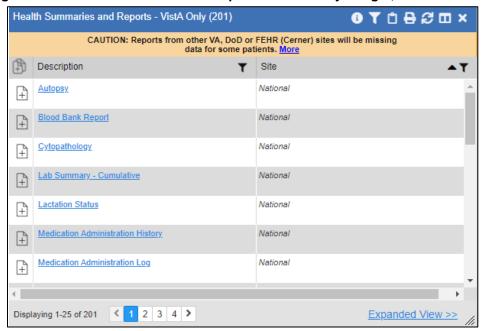

Click <u>Expanded View >></u> from the minimized view of the **Health Summaries and Reports—VistA Only** widget to open an expanded view of the widget. Information in the expanded view (<u>Figure 108</u>) includes:

- Description
- Type
- Site

Figure 108: Health Summaries and Reports—VistA Only Widget, Expanded View

5.7.2. Viewing Health Summaries

Click a hyperlinked entry in the **Description** column of either the minimized or expanded view of the **Health Summaries and Reports—VistA Only** widget to view additional information for the record. This widget is available to VA users only.

Figure 109 displays a health summary record from a single, local site.

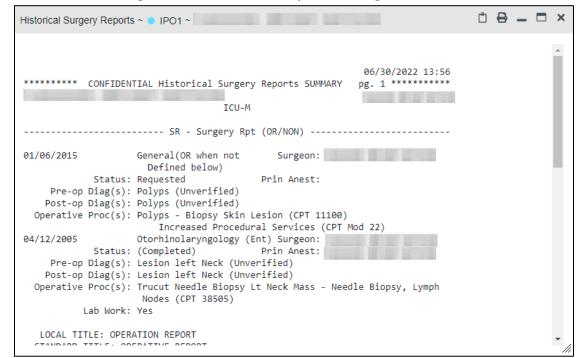


Figure 109: Health Summary from a Single, Local Site

<u>Figure 110</u> is an example of JLV's display of a national health summary. Within the window, tabs represent the VistA sites where the patient is registered. Click each tab to view patient records for that site.

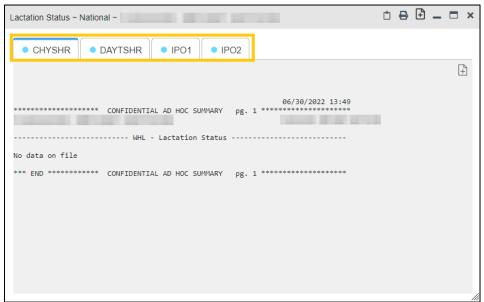


Figure 110: National Health Summary, Multiple Sites

5.8. Lab Results Widget

5.8.1. Lab Results Widget Data

The Lab Results widget displays the patient's lab results information, as well as skin test data and blood transfusion history when available in the patient's record. The information is displayed in reverse chronological order by collection date. Lab data received by DOD sources that include a sensitive flag are masked in the widget views and access to this data is subject to auditing. Additional data may be available in the Community Health Summaries and Documents widget.

The information in the minimized view of the **Lab Results** widget includes:

- Collection Date
- Lab Test
- Result
- Lab Panel
- Site

Click **Expanded View** >> from the minimized view of the **Lab Results** widget to open the expanded view (Figure 111). The lab results information in the expanded view includes:

Collection Date

- Specimen Source
- Lab Test
- Result
- Interpretation
- Units
- Ref Range
- Lab Panel
- Type
- Order Number
- Ordering HCP
- Status
- Location
- Site
 - NOTE: While the Standardized Lab Test and Ordering HCP Specialty columns are not listed by default, you may add them to the expanded view by selecting their corresponding checkboxes in ...

Figure 111: Lab Results Widget, Expanded View

The following lab test abbreviations are used in the **Type** column of the expanded view. Hovering over an entry in the **Type** column displays the full lab type name.


- CHEM for CH, COAG, HEM, HE, TOX, RIA, SER, and SEND
- MICRO for MI, MICROBIOLOGY, and BACT
- AP for ANATOMIC PATHOLOGY
- EM for ELECTRON MICROSCOPY
- SP for SURGICAL PATHOLOGY

- ST for SKIN or SKIN TEST
- CY for CYTOLOGY
- AU for AUTOPSY
- BB for BLOOD BANK
- H.I.V for HIV

5.8.2. Viewing Lab Results Details

Clicking a hyperlinked entry in the **Lab Test** column of the **Lab Results** widget opens a separate window containing a report of the record (Figure 112).

Figure 112: Lab Results Details

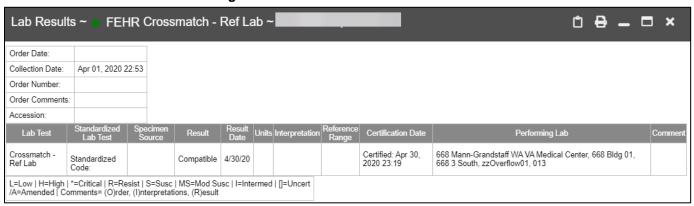

NOTE: A single site report can now be added to the Report Builder from the Blood Bank Report in the Lab Results widget.

5.8.3. Viewing Lab Panel Results Details

Clicking a hyperlinked entry in the **Lab Panel** column of the **Lab Results** widget opens a separate window containing a detailed lab panel report (<u>Figure 113</u>).

Figure 113: Lab Panel Results Detail



5.8.4. Viewing Abnormal Results in the Lab Results Widget

JLV highlights the row for that record in the minimized and expanded views of the widget for records in the **Lab Results** widget where result values fall within a calculated reference range for abnormal.

Pink highlighting represents an abnormal result (i.e., high (H) or low (L)). Red highlighting represents that the abnormal indicator has an asterisk (i.e., H* or L*) or the result is critical.

Examples of records with abnormal results are shown in the expanded view of the **Lab Results** widget in Figure 114.

Remove All Filters

| Remove All Filters | Show | Cardiology | Cardiol

Figure 114: Lab Results Widget, Abnormal Results in Expanded View

5.8.5. Viewing Results and Interpretations in the Lab Results Widget

JLV displays links in the **Result** and **Interpretation** columns in the expanded views of the widget to indicate additional information is available for records in the **Lab Results** widget where result values and interpretation are available in a record. <u>Figure 115</u> highlights multiple result and interpretation links within the widget.



Figure 115: Lab Results Widget, Expanded View, Interpretations Available

Figure 116 displays an example lab interpretation accessed by clicking a link in the **Interpretation** column. The contents of the window may vary depending on the lab type and what is contained in the patient record.

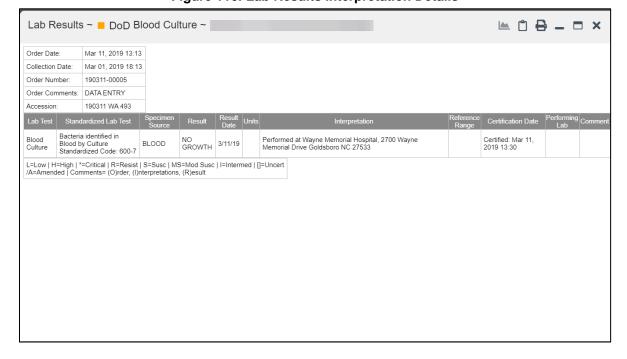


Figure 116: Lab Results Interpretation Details

5.8.6. Graph/Table View from the Lab Results Widget

The **Lab Results** widget provides the option to display multiple results for the same lab test in a graph and table.

- 1. Select in from the widget toolbar in either minimized or expanded view OR -
- 2. Click <u>Expanded View >></u> from the minimized view of the **Lab Results** widget to open an expanded view of the widget, then select the <u>Show Graph/Table</u> link.
- 3. Select a value from the **Select a Lab Test** dropdown in the **Lab Results Graph/Table** dialog.
 - a. The graph displays with corresponding table values beneath.
- 4. Select the **Show Labels** checkbox to toggle graph data labels on or off.
- 5. Select to exit **Graph/Table** view and return to the minimized or expanded view of the **Lab Results** widget.

Figure 117 displays an example **Graph/Table** view of the lab test results.

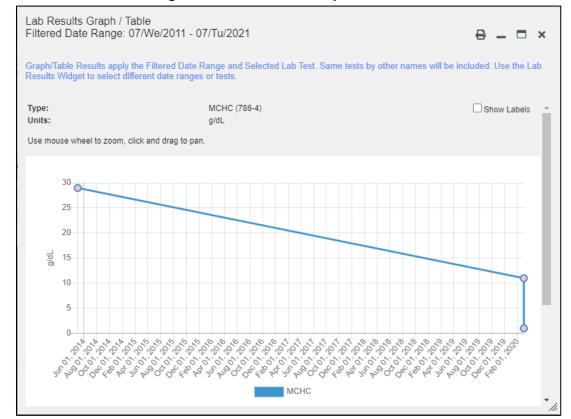


Figure 117: Lab Results Graph/Table View

In the expanded view of the **Lab Results** widget, selecting one or more values in the **Filter by Lab Test** dropdown will disable all other filter dropdowns until the filtered selections are removed. The lab test values in the **Filter by Lab Test** dropdown dynamically increase or decrease depending on the date filter criteria you apply. Tests with the same LOINC code (shown in parentheses following the lab test name) as those selected will be automatically included in the filtered results displayed in the widget.

You may only use the graph/table feature for a single lab test LOINC value. All lab tests with the same LOINC will be graphed even if the lab test names differ. While you can select multiple lab test values in the **Filter by Lab Test** dropdown, you will be forced to select a single lab test value from the **Select a Lab Test** dropdown in the **Lab Results Graph/Table** dialog (**Figure**118) once you select either or the **Show Graph/Table** link. The lab test values in the **Select a Lab Test** dropdown also dynamically increase or decrease depending on the filter criteria applied in the **Lab Results Widget**.



Figure 118: Lab Results Graph Table View, Select a Lab Test

You may also generate a graph for all lab tests with the same LOINC code by selecting the Lab Results Details view (Figure 112).

6. Troubleshooting

The following subsections provide information about troubleshooting common errors in JLV.

Please see online help for how-to information within JLV. Please visit REDACTED to see JLV training videos and access additional JLV training materials.

6.1. Special Instructions for Error Correction

JLV utilizes access control and authentication services to limit access to registered, authorized users. When enabled, JLV validates you against information retrieved from your Smart Card. If you are having trouble logging into JLV, and have used the correct URL for your agency, please review the following before contacting the ESD:

VHA users:

- Have accessed JLV using the JLV URL (REDACTED) or the JLV button in CPRS if available
- Have entered their existing CPRS Access and Verify codes correctly in the fields on the JLV Login page
- Have selected a site from the **Site** dropdown list on the JLV **Login** page (Most users must specify the individual facility name or their parent healthcare system name)

CAPRI-Claims users:

- Have accessed JLV using the JLV URL(REDACTED) or the Joint Longitudinal Viewer tab in CAPRI
- Have entered their existing CAPRI Access and Verify codes correctly in the fields on the JLV **Login** page
- Have selected CAPRI-Claims from the **Site** dropdown list
- Have connected PIV card to the Claims system in IAM.
- Do not have any CAPRI patient selection restrictions

Users with any CAPRI patient selection restrictions (restricted users) must use the **Joint Longitudinal Viewer** tab in CAPRI to access JLV. Restricted users will see the message shown in <u>Figure 119</u> on the JLV **Login** page if attempting to access JLV via the JLV URL.

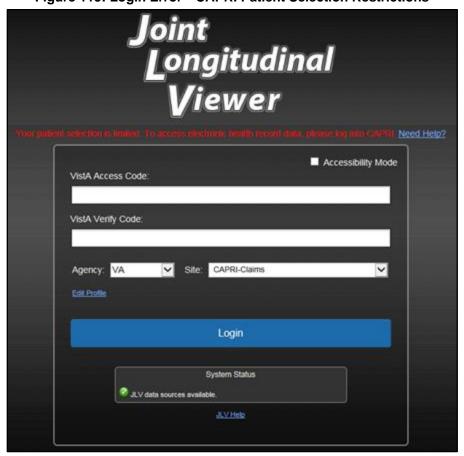


Figure 119: Login Error—CAPRI Patient Selection Restrictions

6.1.1. Login Page Error Messages

Troubleshooting steps for error messages received at the **Login** page are provided in <u>Table 13</u>.

Table 13: Error Messages and Resolution Steps

Error Message/Behavior	Description/Resolution Steps
No access allowed WHY? You have not signed into CAPRI or CPRS for 90 days and account is inactive.	
	FIX IT: Contact the ESD and tell them that either your access to CAPRI or CPRS (whichever you use) has been deactivated.

Error Message/Behavior	Description/Resolution Steps			
Access denied. You are not an authorized user.	WHY? There is a PIV card processing problem, you selected an invalid security certificate, you may not have the OR CPRS GUI Chart VistA option enabled (VHA), or there may be an issue with your single sign-on account. FIX IT:			
	Reinsert PIV and Choose a Valid Authentication Certificate			
	Close your browser window			
	Reinsert your PIV card and relaunch JLV			
	Select a certificate that is not expired and specifies: "Issuer: Veterans Affairs User CA B1"			
	Confirm Your VistA Account is Linked			
	Go to the REDACTED website			
	Select a certificate that is not expired and specifies: "Issuer: Veterans Affairs User CA B1"			
	Select Link VistA User			
	Select the VA Medical Center to link			
	Enter your Access/Verify codes			
	Click Submit			
	If neither of these methods resolve the error, contact the ESD.			
Could not save User Profile.	WHY? The error occurred during PIV authentication. It is either a PIV card processing problem, or you chose the wrong security certificate. FIX IT:			
	Close your browser window			
	Reinsert your PIV card and relaunch JLV			
	Select a certificate that is not expired and specifies: "Issuer: Veterans Affairs User CA B1"			
	If unsuccessful, close all open browser windows/tabs, then open the browser and try JLV again			
Not a valid ACCESS/VERIFY CODE pair.	WHY? JLV could not match your Access and Verify codes to the site selected in the dropdown, or you entered username and password instead of Access/Verify codes. FIX IT:			
	VHA users must select the parent VistA host site for their facility Reenter your Access and Verify codes (CAPRI codes for CAPRI-Claims users, CPRS codes for VHA users)			
Page cannot be	WHY? The JLV URL requires certain IE settings.			
displayed.	FIX IT:			
	Click the Tools menu (press ALT+X) in IE			
	Select Internet Options			
	When the dialog box opens, click the Advanced tab			
	Scroll down in the list until you see the Secure Sockets Layer (SSL) 2.0 setting, and ensure it is NOT checked			
	Ensure that the Transport Layer Security (TLS) 1.0, 1.1 and 1.2 settings ARE checked			
	Click OK to close the dialog box, then relaunch JLV (The page should load)			

Error Message/Behavior	Description/Resolution Steps
Smart Card required.	WHY? Your Smart Card was not read by Windows Security and JLV before opening the JLV URL. FIX IT: 1. Close all browser sessions and browser-based applications Reinsert your PIV card, and relaunch the JLV URL If unsuccessful, close all open IE windows/tabs, then open IE and try JLV again
VERIFY CODE must be changed before continued use.	WHY? Your CPRS or CAPRI Verify code has expired. FIX IT: 1. Open CPRS, VistA, or CAPRI (CAPRI-Claims users) If prompted for a PIV card certificate by CPRS, click Cancel You are prompted to create a new Verify code Once your Verify code has been changed for CPRS or CAPRI, relaunch JLV, which recognizes the new code immediately

6.1.2. System Error Messages

 $\underline{\text{Table 14}}$ summarizes system messages that may be presented to JLV users.

Table 14: System Error Messages

Error Message/Behavior	Description/Resolution Steps
No message displays but behavior occurs: The browser window does not allow the user to scroll to widgets or access the JLV icons on the Patient Data and Widgets page.	 WHY? You may be using an unsupported browser, or your screen resolution settings do not match the recommended configuration. FIX IT: 1. Access JLV using IE 11 or Chrome Use a monitor with a minimum screen resolution of 1024 x 768 Set the browser zoom feature to 100% If browser display issues continue, log out of JLV, close the browser window, open a new browser window, and log into JLV again
Warning: An error occurred while attempting to retrieve VistA Imaging Viewer URL.	WHY? JLV is unable to retrieve the necessary data object to launch a VistA Imaging Viewer instance. FIX IT: Try again later; either VDS or CVIX is offline.
MPI Error: Your query yields too many results. Please modify your search parameters to narrow the search.	WHY? JLV adheres to VA guidelines and blocks search results when a search request yields more than 10 patients. FIX IT: Modify the information in the Patient Search dialog box using additional patient identifiers to reduce the number of results.
MPI Error: There are no patients found using the current parameters.	 WHY? No patient records were found using the information entered in the Patient Search dialog box. FIX IT: 1. Verify the accuracy of the information Reenter patient identifiers in the Patient Search dialog box fields

Error Message/Behavior	Description/Resolution Steps
MPI Error: Application Reject. There was an error attempting to process your query. Please modify your search parameters and try again.	 WHY? There may be one or more errors in the Patient Search dialog fields. FIX IT: 1. Validate the patient identifiers, and try the search again If problems persist, there may be an error between MPI and the Defense Enrollment Eligibility Reporting System (DEERS)
MPI Error: Patient not found with search elements entered, please provide additional patient identifying information and search again.	 WHY? No unique patient is found with the identifiers used in the search. FIX IT: 1. Click OK to return to the Patient Search dialog Provide additional patient identifiers, like the full first name or DOB
You do not have authorization to view this record. Security regulations prohibit computer access to your own medical record.	WHY? VA security regulations prevent you from accessing your own medical records.
You do not have authorization to view this record. Your SSN is missing from the NEW PERSON file. Contact your ADP Coordinator.	WHY? Per VA policy, JLV does not allow access to patient records if the JLV user's SSN is not in their VistA profile. FIX IT: Contact your Automated Data Processing Application Coordinator (ADPAC).
An error occurred while attempting to retrieve the document.	WHY? This error occurs when something goes wrong when JLV tries to retrieve a selected document. FIX IT: 1. Try again If the error persists, contact local support or the ESD
An error occurred while attempting to display the document.	 WHY? This error occurs when JLV tries to retrieve a document but has trouble converting the document to a viewable format. FIX IT: 1. Try again If the error persists, contact local support or the ESD
Your patient selection is limited. To access electronic health record data, please log into CAPRI.	WHY? Some CAPRI-Claims users have patient or site restrictions attached to their profile. FIX IT: Log into CAPRI and select the Joint Longitudinal Viewer tab to view patient records.

A. Acronyms and Abbreviations

 $\underline{\text{Table 15}}$ lists the acronyms and abbreviations used throughout this document and their descriptions.

Table 15: Acronyms and Abbreviations

Acronym	Description
ADPAC	Automated Data Processing Application Coordinator
AFB	Air Force Base
AHLTA	Armed Forces Health Longitudinal Technology Application
AMC	Army Medical Center
ВМІ	Body Mass Index
ВР	Blood Pressure
BSA	Body Surface Area
CAPRI	Compensation and Pension Record Interchange
CCD	Continuity of Care Document
C-CDA	Consolidated Clinical Document Architecture
CCOW	Clinical Context Object Workgroup
CD2	Critical Decision Point #2
CDC	Centers for Disease Control
CDR	Clinical Data Repository
CHCS	Composite Health Care System
CPRS	Computerized Patient Record System
СРТ	Current Procedural Terminology
CVX	Vaccine Administered
DES	Data Exchange Service
DEERS	Defense Enrollment Eligibility Reporting System
DOB	Date of Birth
DOD	Department of Defense
EDIPI	Electronic Data Interchange Pelnal Identifier
EHR	Electronic Health Record
EHRM	Electronic Health Record Modernization
EKG	Electrocardiogram
ESD	Enterprise Service Desk
FEHR	Federal Electronic Health Record
GUI	Graphical User Interface
HAIMS	Healthcare Artifact and Image Management Solution
НСР	Health Care Provider
HCS	Health Care System
HIE	Health Information Exchange
HL7	Health Level Seven International

Acronym	Description
IAM	Identity Access Management
ID	Identification
IE	Internet Explorer
IEN	Internal Entry Number
IP	Internet Protocol
JB	Joint Base
JLV	Joint Longitudinal Viewer
LOINC	Logical Observation Identifiers Names and Codes
МАН	Medication Administration History
MAL	Medication Administration Log
MHS	Military Health System
МІ	Middle Initial
MPI	Master Person Index
NUCC	National Uniform Claim Committee
OIT	Office of Information and Technology
PACT	Patient-Aligned Care Team
PCMM	Primary Care Management Module
PCP	Primary Care Provider
PDWS	Patient Discovery Web Service
PIN	Personal Identification Number
PIV	Personal Identity Verification
SIGI	Self-Identified Gender Identity
SSL	Secure Sockets Layer
SSN	Social Security Number
SSOi	Single Sign-On Internal
TLS	Transport Layer Security
TMDS	Theater Medical Data Store
UC	University of California
UI	User Interface
URL	Uniform Resource Locator
VA	Department of Veterans Affairs
VBA	Veterans Benefits Administration
VDS	VistA Data Service
VHIE	Veterans Health Information Exchange
VIP	Veteran-Focused Integrated Process
VLER	Virtual Lifetime Electronic Record
VHA	Veterans Health Administration
VistA	Veterans Health Information Systems and Technology Architecture

B. Patient Data and Widgets page Widgets, Columns, and Supported Views

<u>Table 16</u> details each data column within the minimized and expanded views of each widget. Data columns available in each view are ordered from left to right.

- An asterisk (*) next to a column title in the table indicates the data for that record is available in "detail view" from within the widget; clicking a link opens the contents of that record in a standalone window.
- A dagger (†) next to a column title in the table indicates the data in that column is normalized (mapped to standards).
- A double dagger (\ddagger) next to a column title in the table indicates the column is hidden by default and must be added from \square .

Table 16: Widget Details

Widget	Description	Minimized View	Expanded View
Admissions	Displays the patient's inpatient admissions information including VA expanded admission discharge transfer history and detailed discharge diagnosis data	Adm Date* Discharge Date Admit Diagnosis† Discharge Diagnosis Ward Site	Adm Date* Discharge Date Ward Provider† Diagnosis† Discharge Diagnosis† Registration Number VA Health Factors Location Site
Allergies	Displays the patient's allergy information	Date Recorded Allergen*† Standardized Allergen† Site	Date Recorded Allergen*† Standardized Allergen† Reaction Severity Comments Location Site
Appointments	Displays the patient's past and future appointments information	Scheduled Date/Time* Clinic Provider† ¹⁸ Site	Scheduled Date/Time* Clinic Provider† ¹⁸ Provider Specialty† ¹⁸ Appointment Status Type Reason VA Health Factors Location Site
Cardiology Studies - FEHR and VA MUSE Only	Displays EKGs and cardiology studies with related images from all FEHR and VA MUSE sites and systems	Date Description (Test Type) Status Image	N/A
Clinical Reminders—VistA Only	Displays the patient's clinical reminders information from all VistA sites	Reminder Site*19	Reminder Site*19

¹⁸ Column hover display includes provider address and phone number.

JLV 3.6.0.0 User Guide 111 October 2023

¹⁹ Lists a sortable, filterable reminder column for each site to which the patient is registered. A maximum of three site columns display in minimized view. All site columns display in expanded view.

Widget	Description	Minimized View	Expanded View
Community Health Summaries	Displays the patient's community health documents from VHIE partners including any CCDs, in C32 and C62 formats, and HL7 C-CDA structured documents	Date Document* Source	N/A
Consults	Displays the patient's outpatient consult information	Date Consult Order* Status Site	Date Consult Order* Provider† Provider Specialty† Status Location Site
Demographics	Displays the patient's non-clinical or personal data including Patient Inquiry, Demographics, Disabilities, Primary Care assignments, Insurance, Military Service, and Enrollment and Eligibility	The site tabs provide a summary from the site (VA or DOD) where the patient is registered	Click the demographics tab to see patient demographic details (VistA patient inquiry for each VA site and DOD DEERS), as well as inpatient and attending providers
Discharge/Inpatient Summaries	Displays the patient's discharge summaries, history, the LOINC Standardized Name, and physical summaries	Date Note Type/Title* Site	Date Note Type/Title* Standardized Note Type / Title Provider† Provider Specialty† Location Visit/Adm Date Status Site
Documents & Images	Displays documents from multiple clinical domains including radiology reports (exams), progress notes, encounters - outpatient, consults, discharge summaries, inpatient notes, questionnaires, and HAIMS and MHS GENESIS records	Date (Note[s]*) Description*† Provider† Image Site	Date (Note[s]*) Description*† Standardized Description† Provider† Provider Specialty† Location Status Image/Attachment Source System VA Health Factors Site
Encounters— Outpatient	Displays records of the patient's outpatient encounters	Encounter Date (Details* Note*) Clinic Provider† Diagnosis† Attachment Site	Encounter Date (Details* Note*) Clinic Status Type Provider† Provider Specialty† Reason Diagnosis† Standardized Diagnosis† Attachment VA Health Factors Location Site
Federal EHR/MHS GENESIS	Displays patient documents and dental summaries from the MHS GENESIS system	Date Document Title* Document Type Location Site	N/A
Health Summaries and Reports— VistA Only	Displays the patient's national and local health summaries and reports from VistA sites where the patient has been registered	Description* Site	Description* Type Site

Widget	Description	Minimized View	Expanded View
Imaging	Displays images from multiple sources, including all Vista Imaging, Cerner imaging, DOD Radiology, and HAIMS artifacts.	Date Image/Attachment #image Image Title Site	Date Image/Attachment #image Image Title Description Class Package Procedure Origin Status Specialty Event Type Capture Date Source System Location Site
Immunizations	Displays the patient's immunization history, including skin test data (where available) in the patient's record	Vaccine Administered Date (Note[s]*) Vaccine Administered Product Type*† Adverse Vaccine Reaction Site (Reporting)	Vaccine Administered Date (Note[s]*) Vaccine Administered Product Type*† Standardized Vaccine Product Type† Series Immunization Result Adverse Vaccine Reaction Exemption/Refusal Reason Exemption/Refusal Date Site of Admin Location Site (Reporting) ²⁰
Lab Results	Displays the patient's lab results information by individual lab test for all test types (i.e., Microbiology, Anatomical Pathology, Cytopathology, and Surgical Pathology) as well as skin test data and blood bank (blood type testing and transfusion) history, when available in the patient's record	Collection Date Lab Test*† Result* Lab Panel* Site	Collection Date Specimen Source Lab Test*† Standardized Lab Test†‡ Result* Interpretation* Units Ref Range Lab Panel* Type Order Number Ordering HCP† Ordering HCP Specialty†‡ Status Location Site
Medications— Inpatient	Displays a patient's inpatient medications information including inpatient infusions	Order Number Drug Name*† Status Standardized Status † Order Stop Date Site	Order Number Order Start Date Order Stop Date Drug Name*† Standardized Drug Name† Status Standardized Status † Ordering HCP† Ordering HCP Specialty† Schedule MAH/MAL Quantity Location Site
Medications— Outpatient	Displays the patient's outpatient medications information, including those discontinued or expired within the past 120 days	Last Fill Drug Name*† Status Standardized Status † Expires Site	Last Fill Drug Name*† Standardized Drug Name† Prescription Number Sig Indication Quantity Days' Supply Refills Left* Status Standardized Status † Ordering HCP† Ordering HCP Specialty† Expires Location Site

²⁰ The **Immunizations** widget differentiates between the reporting site and the site of administration of a vaccine.

Widget	Description	Minimized View	Expanded View
Orders	Displays the patient's medication, consult, radiology, and lab orders	Order Date Description*† Status/Priority Type Site	Order Date Order Number Description*† Standardized Description† Status/Priority Start Date Stop Date Type Provider† Provider Specialty† Location Site
Problem List	Displays the patient's problem list information	Updated Date Problem Description*† Status Linked Items ²¹ Site	Updated Date Onset Problem Description*† Standardized Description† ICD Code ICD Description Severity Type Status Linked Items ²¹ Location Site
Procedures	Displays the patient's procedures for all Current Procedural Terminology (CPT) codes	Procedure Date (Medicine Full Reports* Note[s]*) Location ²² Provider† Procedure Description† Image Site ²³	Procedure Date (Medicine Full Reports* Note[s])* Location ²² Provider† Provider Specialty† CPT Code ²⁴ Procedure Description† Standardized Procedure Description† Image Type Site ²³
Progress Notes	Displays the patient's progress notes information as well as advance directives, clinical warnings, and crisis notes	Date Document Type/Title*† Provider† Image Site	Date Document Type/Title*† Standardized Document Type† Provider† Provider Specialty† Clinic Image Location Site
Questionnaires, Deployment, Exposure Records (AHLTA Only)	Displays the patient's history of questionnaires and DOD pre- and post-deployment assessments	Date Document Title*† Site	Date Entered By Document Title*† Standardized Document Title† Location Site
Radiology Reports	Displays the patient's radiology exams information	Date Exam*† Image Site	Date Exam Number Exam*† Standardized Radiology Exam† CPT Description† Status Results Image Location Site
Social, Family, and Other Histories	Displays the patient's history records from DOD only	Date Reported Type Findings Status Site	Date Reported Type Findings Standardized Findings Comments Status Site

Linked Items contains DOD data only
Location represents a ward or clinic within a given DOD or VA facility where the procedure took place.
Site represents the source of procedure (i.e., VA or DOD). Site is different from Location or Facility.
CPT Code is limited to clinically relevant medical or administrative procedure types.

Widget	Description	Minimized View	Expanded View
Surgery/Procedure Reports	Displays surgical and operative reports for all dates	Date Note Title Procedure Site	Date Note Title Standardized Title Procedure Provider Image Location Site
Vitals	Displays the patient's vital signs information from inpatient or outpatient settings	Date Taken Type*† Result ²⁵ Converted Result ²⁶ Site	Date Taken Type*† Standardized Type† Standardized Group Result ²⁵ Converted Result ²⁶ Location Site

-

²⁵ **Result** displays all results received natively as metric, received as imperial and converted to metric, and those values that are neither metric nor imperial.

²⁶ **Converted Result** displays all results received natively as imperial and all imperial conversions of metric results with the applicable, concatenated unit (height, weight, temperature, etc.), but remains blank for all numeric or text results that are neither metric or imperial and metric results that have no imperial equivalent/conversion (e.g., blood pressure [BP], body mass index [BMI], body surface area [BSA]).