



Vitals / Measurements Installation Guide

GMRV*5.0*23

September 2009

Department of Veterans Affairs
Office of Information & Technology
Office of Enterprise Development

Revision History

Date	Description	Authors
September 2009	Patch GMRV*5.0*23 released. <ul style="list-style-type: none"> - Added Windows XP to list of operating systems, p. 3 - Added names of files for TM and UM changes pages, p. 5 - Removed TEST V1 from installation example, p.6 - Modified Item 6 instructions, p. 10 	Alfred Bustamante, Robert Sutton, Frank Traxler

This page intentionally left blank for double-sided printing.

Table of Contents

Overview	1
Purpose.....	1
Introduction.....	1
Related Manuals.....	1
System Requirements	3
M Server Requirements.....	3
Client Requirements	3
Installation Instructions	5
Pre-Installation Instructions	5
M Server Installation	6
Client Installation	8
Post-Installation Considerations.....	13
Customizing the Client Installation.....	13

This page intentionally left blank for double-sided printing.

Overview

Purpose

The purpose of this installation guide is to provide instructions for use by the technical staff who will install GMRV*5.0*23. Please note that this installation guide does NOT contain instructions for installing the Vitals 5.0 package; this guide only covers Patch 23 installation.

Introduction

Patch GMRV*5.0*23 contains changes to the Gen. Med. Rec. – Vitals package (aka Vitals/Measurements). The package namespaces are GMRV and GMV.

- This patch includes a Dynamic Link Library (DLL) file. This file is used by CPRS.
- This patch modifies the current Vitals.exe file which is the Graphical User Interface (GUI) for entering patient vitals/measurements data via the Vitals/Measurements package.

Related Manuals

The following documents can be downloaded from the VistA Documentation Library (VDL) at <http://www.va.gov/vdl/> under Vitals/Measurements:

- *Vitals/Measurements User Manual 5.0*
- *Vitals/Measurements Release Notes GMRV*5.0*23*
- *Vitals/Measurements Technical Manual and Package Security Guide 5.0*
- *Portable Vital Signs Monitors Interface Specifications (October 2005)*

This page intentionally left blank for double-sided printing.

System Requirements

M Server Requirements

The following packages and patches must be installed and fully patched for the installation environment:

1. VA FileMan V. 22 or greater
2. Kernel V. 8.0 or greater
3. Kernel Toolkit V. 7.3 or greater
4. Kernel RPC Broker V. 1.1 or greater
5. PIMS V. 5.3 or greater
6. Intake and Output V. 4.0
7. Health Summary V. 2.7 or greater
8. Nursing V. 4.0 or greater

Client Requirements

The client (disk) storage requirements are approximately:

Type of Data	Size
Vitals.exe	1900 k
VitalsManager.exe	1200 k
GMV_VitalsViewEnter.dll	1500 k
VITALS.HLP	41 k
VITALSMANAGER.HLP	22 k
GMV_VitalsViewEnter.hlp	36 k

The installation environment on the VistA client workstation requires the following:

1. Workstations must be running Windows NT (V4 or later), Windows 2000, or Windows XP.
2. Workstations must be connected to the local area network (LAN).
3. 12 megabytes of disk space must be available.
4. RPC Broker Workstation may be installed (optional).

This page intentionally left blank for double-sided printing.

Installation Instructions

This section contains instructions for preparing for the installation, installing the M Server, installing the Client, and customizing the Client. Examples and screen shots are included for your reference.

Pre-Installation Instructions

Before attempting to install this patch, complete the following steps:

1. Coordinate the installation with the Clinical Application Coordinator (CAC), package Automated Data Processing Application Coordinator (ADPAC) and Information Resource Management Systems (IRMS).
2. Forward the GMRV*5.0*23 patch message from FORUM to your system. The FORUM message contains the KIDS build.
3. Download the VITL5_P23.ZIP file. The following files are included in this ZIP file:

<u>File Name</u>	<u>Contents</u>	<u>Retrieval Format</u>
GMV_VitalsViewEnter.dll	Dynamic Link Library file	Binary
GMV_VitalsViewEnter.hlp	Help file for DLL	Binary
GMV_VitalsViewEnter.cnt	Help file TOC for DLL	Binary
VITL5_P23_IG.PDF	Patch GMRV*5.0*23 Installation Guide	Binary
VITL5_P23_RN.PDF	Patch GMRV*5.0*23 Release Notes	Binary
VITL5_P23.EXE	Installation Wizard	Binary
VITL5_P23.SRC.ZIP	Source Code	Binary
VITL5_TM.PDF	Technical Manual (all pages)	Binary
VITL5_UM.PDF	User Manual (all pages)	Binary
VITL_5_P23_TM_CHANGE_PAGES.PDF	Technical Manual (change pages)	Binary
VITL_5_P23_UM_CHANGE_PAGES.PDF	User Manual (change pages)	Binary

The preferred method is to FTP the ZIP file from:

- download.vista.med.va.gov

This transmits the ZIP file from the first available FTP server. Sites may also elect to retrieve the ZIP directly from a specific server:

CIO Field Office	FTP Address	Directory
Albany	ftp.fo-albany.med.va.gov	[anonymous.software]
Hines	ftp.fo-hines.med.va.gov	[anonymous.software]
Salt Lake City	ftp.fo-slc.med.va.gov	[anonymous.software]

4. After you have saved the ZIP file, double click on it to “Unzip” it. Highlight all of the files and click the Extract button. Save the files to a directory of your choice.

M Server Installation

Note: The M Server installation must be done before the Client installation.

1. On the VistA system, set the variables `DUZ` and `DUZ(0)` by executing the command `D ^XUP`. Verify `DUZ(0) = @`.
2. Load the `GMRV*5.0*23 KIDS` build from the MailMan message.
3. Use the `KIDS` installation menu option [`XPD MAIN`] and select **Installation** and then **Install Package(s)** and select `GMRV*5.0*23`. See the M Server sample installation below for additional information.

Users may remain on the system at the time of installation, though it should be installed when entry of patient vitals data is low. The software should be installed when use of the Vitals/Measurements package is minimal. Follow your facility's policy regarding the rebuilding of the menu trees upon patch completion.

M Server – sample installation

```
> D ^XUP

Setting up programmer environment
Terminal Type set to: C-VT100

Select OPTION NAME: XPD MAIN      Kernel Installation & Distribution System

      Edits and Distribution ...
      Utilities ...
KIDS  Installation ...
      Patch Monitor Main Menu ...

Select Kernel Installation & Distribution System Option: INSTALLATION

1      Load a Distribution
2      Verify Checksums in Transport Global
3      Print Transport Global
4      Compare Transport Global to Current System
5      Backup a Transport Global
6      Install Package(s)
      Restart Install of Package(s)
      Unload a Distribution

Select Installation Option: 6  Install Package(s)
Select INSTALL NAME: GMRV*5.0*23      Loaded from Distribution 2/20/08@14:20:03

=> GMRV*5*23

This Distribution was loaded on Feb 20, 2008@14:20:03 with header of
GMRV*5*23
It consisted of the following Install(s):
GMRV*5.0*23
```

Checking Install for Package GMRV*5.0*23

Install Questions for GMRV*5.0*23

Want KIDS to INHIBIT LOGONS during the install? NO// NO
Want to DISABLE Scheduled Options, Menu Options, and Protocols? NO// NO

Enter the Device you want to print the Install messages.
You can queue the install by enter a 'Q' at the device prompt.
Enter a '^' to abort the install.

DEVICE: HOME// TELNET

Install Started for GMRV*5.0*23 :
Feb 20, 2008@14:20:33

Build Distribution Date: Feb 20, 2008

Installing Routines:
Feb 20, 2008@14:20:33

Installing PACKAGE COMPONENTS:
GMRV*5.0*23

Installing REMOTE PROCEDURE
Feb 20, 2008@14:20:33

Running Post-Install Routine: ^GMV23PST

Updating system parameters.

Updating DLL parameter.


Updating Routine file...

Updating KIDS files...

GMRV*5.0*23 Installed.
Feb 20, 2008@14:20:33

Install Message sent #68780

100%
Complete

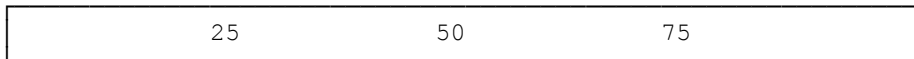

Install Completed

Client Installation

Note: The M Server installation must be done before the VistA Client installation.

1. Double click on the VITL5_P23.exe file in the directory where you placed it. This file is an installation wizard that contains the updated versions for:
 - Vitals.exe
 - VitalsManager.exe
 - Vitals.hlp
 - Vitals.cnt
 - VitalsManager.hlp
 - VitalsManager.cnt
 - Roboex32.dll
2. An InstallShield welcome screen opens. Click **Next** to start the installation.

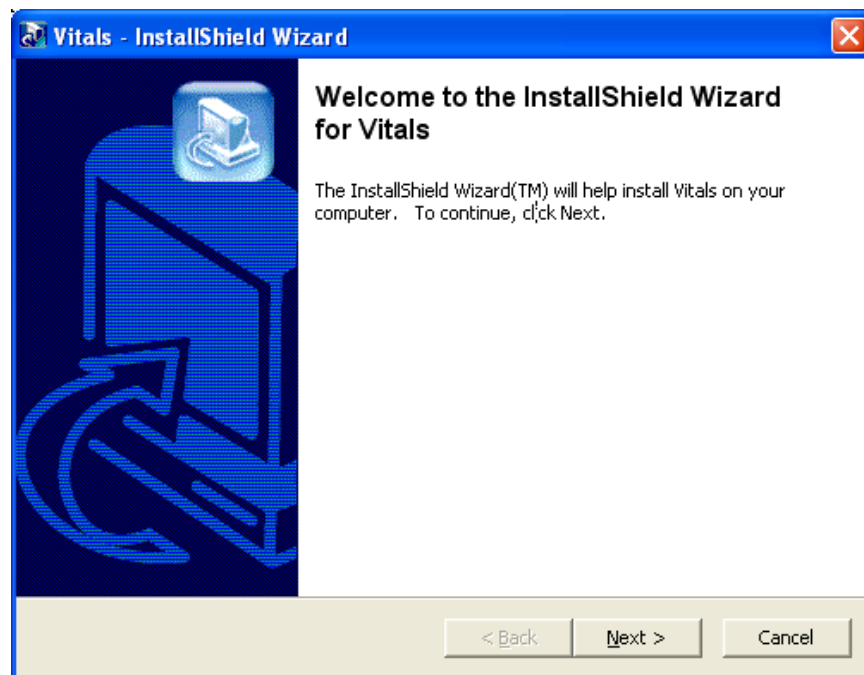


Figure 1

3. Several dialog boxes will then quickly flash across the screen before the “Modify, repair or remove the program” dialog box appears.

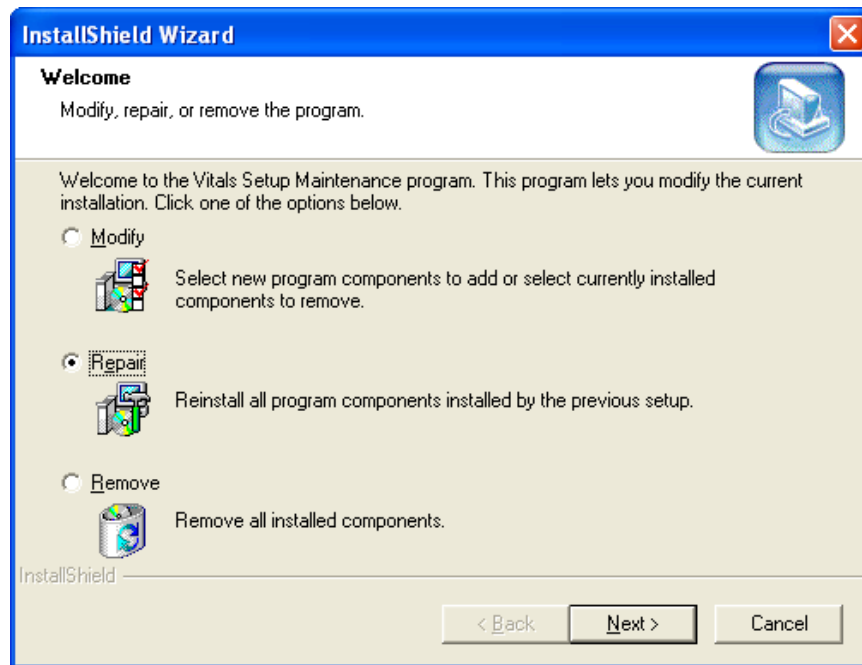


Figure 2

4. Select **Repair**, then click **Next**. The Install wizard will verify settings and replace all existing files on your workstation with the newer versions. If the installation files are located on the client PC, the installation should complete in less than one minute. Installations over the network may be slower because of server traffic or connectivity issues.
5. When all files have been copied, the InstallShield Wizard Complete screen will open. Click **Finish** to finalize the client installation.

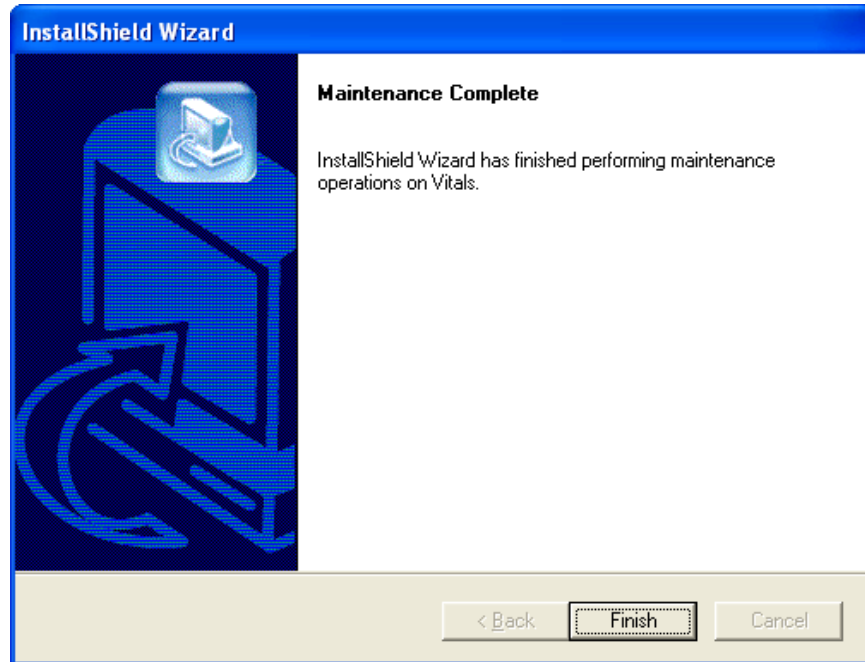


Figure 3

6. The installation wizard creates the following files:
 - \Program Files\Vista\Vitals\Vitals.exe
 - \Program Files\Vista\Vitals\VitalsManager.exe
 - \Program Files\Vista\Vitals\Help\Vitals.hlp
 - \Program Files\Vista\Vitals\Help\Vitals.cnt
 - \Program Files\Vista\Vitals\Help\VitalsManager.hlp
 - \Program Files\Vista\Vitals\Help\VitalsManager.cnt
 - \Program Files\Vista\Vitals\Help\Roboex32.dll

The exe files are placed in the “Program Files\vista\Vitals” directory on the workstation. The other files are placed in the “Program Files\vista\Vitals\Help” directory.

If you are running the software from a server, you should move these files to the server.

If you are running the software from individual workstations and want to push these files, you should include them in your script.

7. Push the GMV_VitalsViewEnter.dll, GMV_VitalsViewEnter.hlp, and GMV_VitalsViewEnter.cnt files to the “Program Files\vista\Common Files” directory of the workstation where CPRS is located. For example, if CPRS is installed on the C:\ drive, the three DLL related files must be copied into the C:\Program Files\vista\Common Files directory. If this directory path does not exist, you must create it.

IMPORTANT: These files must not be installed in the same folder as CPRS. CPRS looks for these three files in the specified directory path. If CPRS cannot find the DLL and its supporting files, the users will not be able to enter patient vitals/measurements data.

Note: The GMV_VitalsViewEnter.dll does not have to be registered in the Windows registry.
The Client installation is complete.

Post-Installation Considerations

Customizing the Client Installation

By default, the client installation creates the icons and “Program Files\vista\Vitals” folders without any command line switches. Vitals/Measurements utilizes the ServerList utility of the RPC Broker for selecting a server to connect to if it is configured on the client workstation. Instructions for configuration and utilization of the ServerList utility can be found in the RPC Broker documentation located on the VDL.

If the ServerList utility has not been configured on the client, both Vitals and Vitals Manager applications will, by default, attempt to connect to the server identified in the users HOSTS file as BROKERSERVER on Listener Port 9200.

To override these default parameters, use the following procedure to add command line parameters to the application shortcuts.

1. On the client desktop, right-click the Vitals icon and select **Properties**. The Vitals Properties window opens. Click the **Shortcut** tab to display the current target settings.

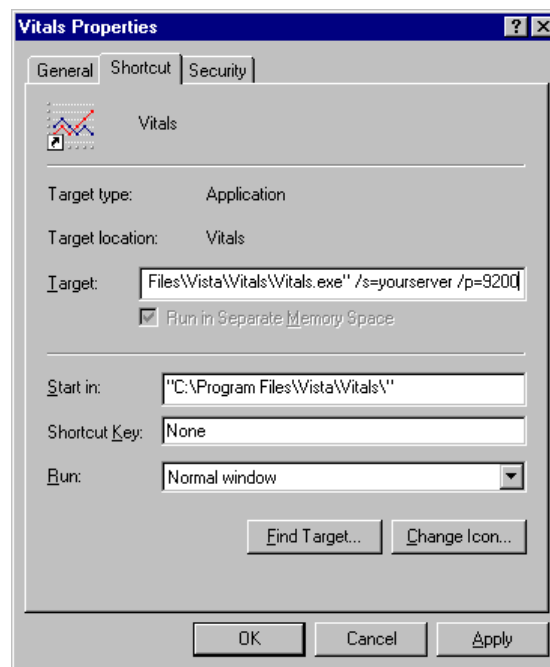


Figure 4

2. In Figure 4, the application will attempt to connect to the server identified in your HOSTS file as *yourserver* and will use listener port 9200. In this example, the complete Target line would read:

```
"C:\Program Files\Vista\Vitals\Vitals.exe" /s=yourserver /p=9200
```

Post-Installation Considerations

3. Enter a different parameter or switch in the Target field. The command line parameters available from the command prompt or within Windows shortcut definitions are:

Vitals.exe **[/server=servername] [/port=listenerport]**
 [/tempdir=temporarydirectory] [/helpdir=helpdirectory]
 [/debug={on|off}] [/noccow] [/ccow=patientonly]

VitalsManager.exe **[/server=servername] [/port=listenerport] [/helpdir=helpdirectory]**
 [/debug={on|off}]

The following table describes each of the available parameters and switches.

Switches	Description	Example
/server	Specifies an alternate server to connect to. The server must be defined in the clients hosts file. Default Hosts. file locations: NT 4.0/W2K = c:\winnt\system32\drivers\etc\hosts. Windows 9x = c:\windows\hosts. Default = BROKERSERVER	/server=vista
/port	Specifies an alternate listener port on the selected server. This is the TCP/IP port that the broker is running on VistA server. Default = 9200	/port=9200
/tempdir	Location accessible to the client workstation and current user for storage of temporary scratch files. Default = <i>application_directory</i> \temp	/tempdir=C:\temp
/helpdir	Location of the Vitals/Measurements windows help files. Default = <i>application_directory</i> \help	/helpdir=C:\help
/debug	Set the debug mode for both the RPC Broker and the Vitals/Measurements application. Default = Off.	/debug=On
/noccow	The application will not check the CCOW context at all. This switch will force the user to sign on and select a patient when invoking the Vitals GUI.	/noccow
/ccow=patientonly	The application will use CCOW, but will be set to check for patient context only. Automatic sign on will be disabled, but the automatic selection of a patient will be enabled. If a patient is already selected in an open application, Vitals will automatically open the patient being used by that application.	/ccow=patientonly